

This Page Is Inserted by IFW Operations
and is not a part of the Official Record

BEST AVAILABLE IMAGES

Defective images within this document are accurate representation of
The original documents submitted by the applicant.

Defects in the images may include (but are not limited to):

- BLACK BORDERS
- TEXT CUT OFF AT TOP, BOTTOM OR SIDES
- FADED TEXT
- ILLEGIBLE TEXT
- SKEWED/SLANTED IMAGES
- COLORED PHOTOS
- BLACK OR VERY BLACK AND WHITE DARK PHOTOS
- GRAY SCALE DOCUMENTS

IMAGES ARE BEST AVAILABLE COPY.

**As rescanning documents *will not* correct images,
please do not report the images to the
Image Problem Mailbox.**

This Page Blank (uspto)

(19) World Intellectual Property Organization
International Bureau



(43) International Publication Date
30 August 2001 (30.08.2001)

PCT

(10) International Publication Number
WO 01/63244 A1

(51) International Patent Classification⁷: G01N 9/00, 33/48

(74) Agent: ARNOLD, Craig, J.; Amster, Rothstein & Ebenstein, 90 Park Avenue, New York, NY 10016 (US).

(21) International Application Number: PCT/US01/05150

(81) Designated States (*national*): AE, AG, AL, AM, AT, AU, AZ, BA, BB, BG, BR, BY, BZ, CA, CH, CN, CR, CU, CZ, DE, DK, DM, DZ, EE, ES, FI, GB, GD, GE, GH, GM, HR, HU, ID, IL, IN, IS, JP, KE, KG, KP, KR, KZ, LC, LK, LR, LS, LT, LU, LV, MA, MD, MG, MK, MN, MW, MX, MZ, NO, NZ, PL, PT, RO, RU, SD, SE, SG, SI, SK, SL, TJ, TM, TR, TT, TZ, UA, UG, UZ, VN, YU, ZA, ZW.

(22) International Filing Date: 16 February 2001 (16.02.2001)

(25) Filing Language: English

(26) Publication Language: English

(30) Priority Data:
09/514,026 25 February 2000 (25.02.2000) US

(84) Designated States (*regional*): ARIPO patent (GH, GM, KE, LS, MW, MZ, SD, SL, SZ, TZ, UG, ZW), Eurasian patent (AM, AZ, BY, KG, KZ, MD, RU, TJ, TM), European patent (AT, BE, CH, CY, DE, DK, ES, FI, FR, GB, GR, IE, IT, LU, MC, NL, PT, SE, TR), OAPI patent (BF, BJ, CF, CG, CI, CM, GA, GN, GW, ML, MR, NE, SN, TD, TG).

(71) Applicant: AMERICAN HOME PRODUCTS CORPORATION [US/US]; Five Giralda Farms, Madison, NJ 07940-0874 (US).

Published:

- with international search report
- before the expiration of the time limit for amending the claims and to be republished in the event of receipt of amendments

For two-letter codes and other abbreviations, refer to the "Guidance Notes on Codes and Abbreviations" appearing at the beginning of each regular issue of the PCT Gazette.



(72) Inventors: CHEN, James, M.; 7 Sgt. David Stoddard Court, Bedminster, NJ 07921 (US). MOBILIO, Dominick; 35 Sneider Road, Warren, NJ 07059 (US). MOY, Franklin, J.; 37 Burch Street, Arlington, MA 02414 (US). PARRIS, Kevin, D.; 112 Woodbine Street, Auburndale, MA 02466 (US). POWERS, Robert; 3 Magnolia Drive, Westford, MA 01866 (US). BAO XU, Zhang; 40 Fieldston Circle, Tewksbury, MA 01876 (US).

WO 01/63244 A1

(54) Title: SOLUTION AND CRYSTAL STRUCTURES OF MMP-13 ACTIVE SITE AND USES THEREOF

(57) Abstract: The present invention relates to the three dimensional structure of human collagenase 3 (MMP-13), as well as to (i) methods of using the MMP-13 structure to rationally design or identify compounds or molecules that inhibit or activate MMP-13 activity, and (ii) compounds identified using said methods.

USSN: 10/075,909 FILED: 02/13/02
A0000517-01-CFP

SOLUTION AND CRYSTAL STRUCTURES OF MMP-13 ACTIVE SITE AND USES THEREOF

5 Field of the Invention

The present invention relates to the three dimensional structure of human collagenase 3 (MMP-13), as well as to (i) methods of using the MMP-13 structure to rationally design or identify compounds or molecules that inhibit or activate MMP-13 activity, and (ii) compounds identified using said methods.

10 Background of the Invention

Human collagenase-3 (MMP-13) is a member of the matrix metalloproteinase (MMP) family which includes the collagenases, stromelysins and gelatinases. The MMPs are involved in the degradation of the extracellular matrix and are associated with normal tissue remodeling processes such as pregnancy, wound healing, and angiogenesis. MMP expression and activity is highly controlled because of the degradative nature of these enzymes, where an apparent loss in MMP regulation results in the pathological destruction of connective tissue and the ensuing disease state. Accordingly, MMPs are a highly active set of targets for the design of therapeutic agents for the disease areas of arthritis and oncology (for reviews, see Woessner, J. F., FASEB 1991; Ries, C., and Petrides, E., Biol. Chem. Hoppe-Seyler 1995; Browner, M. F., Perspect. Drug Discovery Des. 1995; Morphy, *et al.*, Curr. Med. Chem. 1995; and Zask, *et al.*, Curr. Pharm. Des. 1996).

25 MMP-13 was identified on the basis of differential expression in normal breast tissues and in breast carcinoma. In addition, its expression has been reported in squamous cell carcinomas of the larynx, head and neck, in HCS-2/8 human chondrosarcoma cells, during fetal ossification, and in articular cartilage of arthritic patients.

30 There have been a number of X-ray and NMR structures solved for the catalytic domain of MMPs complexed with a variety of inhibitors (see e.g., Bode, *et al.*, EMBO J. 1994; Gooley, *et al.*, Nat. Struct. Biol. 1994; Lovejoy, *et al.*, Science 1994; Lovejoy, *et al.*, Ann. N. Y. Acad. Sci. 1994; Lovejoy, *et al.*,

-2-

Biochemistry 1994; Spurlino, *et al.*, Proteins: Struct., Funct., Genet. 1994; Stams, *et al.*, Nat. Struct. Biol. 1994; Becker, *et al.*, Protein Sci. 1995; Gonnella, *et al.*, Proc. Natl. Acad. Sci. U.S.A. 1995; Van Doren, *et al.*, Protein Sci. 1995; Botos, *et al.*, Proc. Natl. Acad. Sci. USA 1996; Broutin, *et al.*, Acta Crystallogr., Sect. D: Biol. Crystallogr. 1996; Gooley, *et al.*, J. Biomol. NMR 1996; Betz, *et al.*, Eur. J. Biochem. 1997; Gonnella, *et al.*, Bioorg. Med. Chem. 1997; and Moy, *et al.*, Biochemistry 1998). There is a close similarity in the overall three-dimensional fold for these proteins consistent with the relatively high sequence homology (> 40%). Despite this similarity in the MMP structures, there is a distinct substrate specificity between these enzymes indicative of specific biological roles for the various MMPs and a corresponding association with unique disease processes. One example of this potential specificity is the over-expression of MMP-13 in breast carcinoma and MMP-1 in papillary carcinomas. Therefore, the current paradigm in the development of MMP inhibitors is to design specificity into the structures of the small molecule instead of developing a broad spectrum MMP inhibitor (Birkedal-Hansen, *et al.*, Crit. Rev. Oral Biol. Med. 1993; and Rockwell, *et al.*, J. Am. Chem. Soc. 1996). The rationale behind this approach is that an inhibitor specific for the MMP uniquely associated with a disease process may potentially minimize toxic side effects. Therefore, extensive structural information for the various MMPs is critical for a structure-based approach in designing inhibitor selectivity (Birkedal-Hansen, *et al.*, Crit. Rev. Oral Biol. Med. 1993; Rockwell, *et al.*, J. Am. Chem. Soc. 1996; Ghose, *et al.*, J. Am. Chem. Soc. 1995; Hajduk, *et al.*, J. Am. Chem. Soc. 1997; and Olejniczak, *et al.*, J. Am. Chem. Soc. 1997).

This concept has been facilitated by the extensive structural data available for the MMPs where a significant difference in the size and shape of the S1' pocket has been observed (Moy, *et al.*, Biochemistry 1998; Bode, *et al.*, EMBO J. 1994; Gooley, *et al.*, Nat. Struct. Biol. 1994; Lovejoy, *et al.*, Ann. N.Y. Acad. Sci. 1994; Lovejoy, *et al.*, Biochemistry 1994; Lovejoy, *et al.*, Science 1994; Spurlino, *et al.*, Proteins: Struct., Funct., Genet. 1994; Stams, *et al.*, Nat. Struct. Biol. 1994; Becker, *et al.*, Protein Sci. 1995; Gonnella, *et al.*, Proc. Natl.

-3-

Acad. Sci. U.S.A. 1995; Van Doren, *et al.*, Protein Sci. 1995; Botos, *et al.*, Proc. Natl. Acad. Sci. U.S.A. 1996; Broutin, *et al.*, Acta Crystallogr., Sect. D: Biol. Crystallogr. 1996; Gooley, *et al.*, J. Biomol. NMR 1996; Betz, *et al.*, Eur. J. Biochem. 1997; and Gonnella, *et al.*, Bioorg. Med. Chem. 1997). This structural
5 difference across the MMP family provides an obvious approach for designing specificity into potent MMP inhibitors by designing compounds that appropriately fill the available space in the S1' pocket while taking advantage of sequence differences. A number of examples have been previously reported using this approach where some selectivity between MMPs has been achieved by
10 incorporating a biphenyl into the S1' pocket (*see e.g.*, Hajduk, *et al.*, J. Am. Chem. Soc. 1997; and Olejniczak, *et al.*, J. Am. Chem. Soc. 1997).

The inventors have determined both the solution and crystal structures of MMP-13, and, using rational drug design methods, have designed a novel, potent inhibitor that is highly selective for MMP-13.

15

Summary of the Invention

The present invention relates to the three dimensional structure of human collagenase 3 (MMP-13), and more specifically, to the crystal and solution structures of MMP-13 complexed with the inhibitor N-Hydroxy-2-[(4-methoxy-benzenesulfonyl)-pyridin-3-ylmethyl-amino]-3-methyl-benzamide
20 (hereinafter referred to as "Compound A"), as determined using crystallography, spectroscopy and various computer modeling techniques. Particularly, the invention is directed to an MMP-13 active site comprised of the three dimensional structures of various binding pockets located both to the right (S1', S2', S3') and left (S1, S2, S3) of the catalytic zinc of MMP-13, and most
25 particularly is directed to the three dimensional structure of an MMP-13 active site comprising the catalytic zinc and the S1' binding pocket, which is critical to the design and selection of inhibitors with increased potency and specificity for MMP-13, or conversely, for the design and selection of inhibitors of matrix
30 metalloproteinases that are specific against MMP-13.

-4-

Accordingly, the present invention discloses a solution comprising a biologically active catalytic fragment of human collagenase-3 (MMP-13) complexed with Compound A, as well as a crystallized catalytic fragment of MMP-13 complexed with Compound A. The three dimensional structure of the catalytic fragment of MMP-13 is provided by the relative atomic structural coordinates of Figure 4, as obtained from spectroscopy data, and Figure 5, as obtained from crystallography data. Also provided is an active site of MMP-13, characterized by a catalytic zinc, a beta strand, a Ca^{2+} binding loop, an alpha helix and a random coil region, wherein the beta strand of said active site preferably comprises residues N14, L15, T16, Y17, R18, I19, and V20 according to Figure 1, the Ca^{2+} binding loop comprises residues F75, D76, G77, P78, and S79 according to Figure 1, the alpha helix comprises residues N112, L113, F114, L115, V116, A117, A118, H119, E120, F121, G122, and H123 according to Figure 1, and the random coil region comprises residues P139, I140, and Y141 according to Figure 1. Said active site is further characterized by a three dimensional structure comprising the relative structural coordinates of the catalytic zinc and amino acid residues N14, L15, T16, Y17, R18, I19, V20, F75, D76, G77, P78, S79, N112, L113, F114, L115, V116, A117, A118, H119, E120, F121, G122, H123, P139, I140, and Y141 according to the solution or crystal coordinates of Figures 4 or 5, respectively, in each case, \pm a root mean square deviation from the catalytic zinc and conserved backbone atoms of said amino acids of not more than 1.5\AA .

In an alternate embodiment of the invention, an active site of MMP-13 is characterized by a three dimensional structure comprising the relative structural coordinates of the catalytic zinc and amino acid residues L81, L82, L115, V116, H119, L136 and I140 according to the solution or crystal coordinates of Figures 4 or 5, respectively, in each case, \pm a root mean square deviation from the catalytic zinc and conserved backbone atoms of said amino acids of not more than 1.5\AA .

The solution or crystal structural coordinates of MMP-13 or portions thereof as provided by this invention may be stored in a

machine-readable form on a machine-readable storage medium, e.g. a computer hard drive, diskette, DAT tape, etc., for display as a three-dimensional shape or for other uses involving computer-assisted manipulation of, or computation based on, the structural coordinates or the three-dimensional structures they

5 define. By way of example, the data defining the three dimensional structure of MMP-13 or an MMP-13 complex of the present invention, or of a portion of MMP-13 or an MMP-13 complex as disclosed herein, may be stored in a machine-readable storage medium, and may be displayed as a graphical three-dimensional representation of the relevant structural coordinates, typically

10 using a computer capable of reading the data from said storage medium and programmed with instructions for creating the representation from such data.

Accordingly, the present invention provides a machine, such as a computer, programmed in memory with the coordinates of the MMP-13 molecule or molecular complex, or portions thereof (such as, by way of

15 example, the coordinates of the MMP-13 catalytic zinc with adjacent S1', S2' and/or S3' binding pockets), together with a program capable of converting the coordinates into a three dimensional graphical representation of the structural coordinates on a display connected to the machine. A machine having a memory containing such data aids in the rational design or selection of

20 inhibitors or activators of MMP-13 activity, including the evaluation of ability of a particular chemical entity to favorably associate with MMP-13 or an MMP-13 complex as disclosed herein, as well as in the modeling of compounds, proteins, complexes, etc. related by structural or sequence homology to MMP-13.

The present invention is additionally directed to a method of

25 determining the three dimensional structure of a molecule or molecular complex whose structure is unknown, comprising the steps of first obtaining crystals or a solution of the molecule or molecular complex whose structure is unknown, and then generating X-ray diffraction data from the crystallized molecule or molecular complex and/or generating NMR data from the solution of the

30 molecule or molecular complex. The generated diffraction or spectroscopy data from the molecule or molecular complex can then be compared with the known

-6-

three dimensional structure of MMP-13 as disclosed herein, and the three dimensional structure of the unknown molecule or molecular complex conformed to the known MMP-13 structure using standard techniques such as molecular replacement analysis, 2D, 3D and 4D isotope filtering, editing and
5 triple resonance NMR techniques, and computer homology modeling. Alternatively, a three dimensional model of the unknown molecule may be generated by generating a sequence alignment between MMP-13 and the unknown molecule, based on any or all of amino acid sequence identity, secondary structure elements or tertiary folds, and then generating by computer
10 modeling a three dimensional structure for the molecule using the three dimensional structure of, and sequence alignment with, MMP-13.

The present invention further provides a method for identifying a potential inhibitor or activator of MMP-13, comprising the steps of using a three dimensional structure of MMP-13 as defined by the relative structural
15 coordinates of amino acids encoding MMP-13 to design or select a potential inhibitor or activator, and synthesizing or obtaining said potential inhibitor or activator. The inhibitor or activator may be selected by screening an appropriate database, may designed *de novo* by analyzing the steric configurations and charge potentials of an empty MMP-13 active site in
20 conjunction with the appropriate software programs, or may be designed using characteristics of known inhibitors or activators to MMP-13 or other collagenases in order to create "hybrid" activators or inhibitors. The method of the present invention is preferably used to design or select inhibitors of MMP-13 activity.

25 Alternatively, the present invention provides a method for identifying a potential inhibitor or activator that is selective for one or more members of the matrix metalloproteinase family except MMP-13, comprising the steps of (i) using the three dimensional structures of MMP-13 and the desired target matrix metalloproteinase(s) as defined by the relative structural
30 coordinates of amino acids encoding MMP-13 and the target matrix metalloproteinase(s) in order to design or select such a potential inhibitor or

-7-

activator, and (ii) synthesizing or obtaining said potential inhibitor or activator. In this case, the potential inhibitor or activator is designed to incorporate chemical or steric features favorable for association with an active site of the desired matrix metalloproteinase(s) and unfavorable for association with an

5 MMP-13 active site, preferably where said active site comprises the MMP-13 S1' pocket. The inhibitor or activator may be selected by screening an appropriate database, may designed *de novo* by analyzing the steric configurations and charge potentials of empty MMP-13/matrix metalloproteinase active sites in conjunction with the appropriate software programs, or may be designed using

10 characteristics of known inhibitors or activators to MMP-13 or other collagenases in order to create "hybrid" activators or inhibitors.

Also provided by the present invention are the inhibitors and activators designed or selected using the methods disclosed herein.

15 Brief Description of the Figures

Figure 1 depicts the amino acid sequence encoding the catalytic fragment of human MMP-13.

Figure 2 depicts the sequence based alignment between (A) MMP-13 and MMP-8 and (B) MMP-13 and MMP-1 used for the MMP-13 homology

20 model.

Figure 3 is an illustration of the sulfonamide derivative of the hydroxamic inhibitor N-Hydroxy-2-[(4-methoxy-benzenesulfonyl)-pyridin-3-ylmethyl-amino]-3-methyl-benzamide (Compound A), with the corresponding proton labels.

25 Figure 4 lists the atomic structure coordinates for the restrained minimized mean structure of MMP-13 complexed with Compound A as derived by NMR spectroscopy. "Atom type" refers to the atom whose coordinates are being measured. "Residue" refers to the type of residue of which each measured atom is a part - i.e., amino acid, cofactor, ligand or solvent. The "x, y and z"

30 coordinates indicate the Cartesian coordinates of each measured atom's location

-8-

(Å). All non-protein atoms (Compound A, zinc and calcium) are listed as HETATM instead of atoms using PDB conventions.

Figure 5 lists the atomic structure coordinates for MMP-13 as derived by X-ray diffraction of a crystallized MMP-13:Compound A complex.

- 5 Figure headings are as noted above, except "Occ" indicates the occupancy factor, and "B" indicates the "B-value", which is a measure of how mobile the atom is in the atomic structure (\AA^2). "MOL" indicates the segment identification used to uniquely identify each molecule in the crystal.

- Figure 6 is an illustration of the Compound B inhibitor, with the
10 corresponding proton labels.

- Figure 7 is a design scheme dividing 2-[Benzyl-(4-phenethyloxy-benzenesulfonyl)-amino]-N-hydroxy-3,5-dimethyl-benzamide (hereinafter referred to as "Compound C") into two components corresponding to its potency component (2-[Benzyl-(4-methoxy-benzenesulfonyl)-amino]-N-hydroxy-3,5-
15 dimethyl-benzamide, hereinafter referred to as "Compound D") and its selectivity component, thereby providing the basis for the design of a hybrid inhibitor with Compound B.

- Figure 8A is a design scheme showing the flow from Compound B and Compound C to the hybrid inhibitor benzofuran-2-carboxylic acid (2-{4-
20 [benzyl-(2-hydroxycarbamoyl-4,6-dimethyl-phenyl)-sulfamoyl]-phenoxy}-ethyl)-amide (hereinafter referred to as "Compound E"). Figure 8B illustrates an expanded view of the NMR MMP-13:Compound B complex overlaid with the MMP-13:Compound D model, demonstrating the approach to forming the hybrid inhibitor Compound E. The MMP-13 active site is shown as a grid
25 surface with Compound B and Compound D shown as liquorice bonds. The view is looking at the S1' pocket.

Detailed Description of the Invention

As used herein, the following terms and phrases shall have the meanings set forth below:

“Compound A” is N-Hydroxy-2-[(4-methoxy-benzenesulfonyl)-pyridin-3-ylmethyl-amino]-3-methyl-benzamide, as shown in Figure 3.

“Compound B” is the compound having the chemical structure shown in Figure 6. “Compound C” is 2-[Benzyl-(4-phenethyloxy-benzenesulfonyl)-amino]-N-hydroxy-3,5-dimethyl-benzamide, as shown in Figure 7. “Compound D” is 2-[Benzyl-(4-methoxy-benzenesulfonyl)-amino]-N-hydroxy-3,5-dimethyl-benzamide, also shown in Figure 7. “Compound E” is Benzofuran-2-carboxylic acid (2-{4-[benzyl-(2-hydroxycarbamoyl-4,6-dimethyl-phenyl)-sulfamoyl]-phenoxy}-ethyl)-amide, as shown in Figure 8A. “Compound F” is 2-(Benzyl-4-(3-phenyl-propoxy)-benzenesulfonyl]-amino)-N-hydroxy-3,5-dimethyl-benzamide.

Unless otherwise noted, “MMP-13” includes both human collagenase 3 as encoded by the amino acid sequence of Figure 1 (including conservative substitutions thereof), as well as “MMP-13 analogues”, defined herein as proteins comprising an MMP-13 like active site as defined by the present invention, including, but not limited to, an active site characterized by a three dimensional structure comprising the relative structural coordinates of the catalytic zinc and amino acid residues L81, L82, L115, V116, H119, L136 and I140 according to the solution or crystal coordinates of Figures 4 or 5, respectively, in each case, \pm a root mean square deviation from the catalytic zinc and conserved backbone atoms of said amino acids of not more than 1.5Å, or more preferably, not more than 1.0Å, or most preferably, not more than 0.5Å. Alternatively, an MMP-13 analogue of the present invention is a protein which comprises an MMP-13 like active site characterized by a catalytic zinc, a beta strand, a Ca^{2+} binding loop, an alpha helix and a random coil region, or, more particularly, comprising an active site characterized by a three dimensional structure comprising the relative structural coordinates of the catalytic zinc and of amino acid residues N14, L15, T16, Y17, R18, I19, V20,

-10-

F75, D76, G77, P78, S79, N112, L113, F114, L115, V116, A117, A118, H119, E120, F121, G122, H123, P139, I140, and Y141 according to Figures 4 or 5, or more preferably, where said three dimensional structure further comprises the relative structural coordinates of amino acid residues G80, L81, L82, A83, H84, 5 A85, K109, G110, Y111, S124, L125, G126, L127, D128, H129, S130, K131, D132, P133, G134, A135, L136, M137, F138, T142, Y143, T144, and G145 according to Figures 4 or 5, or most preferably, where said three dimensional structure still further comprises the relative structural coordinates of F149 and P152 according to Figures 4 or 5, in each case, \pm a root mean square deviation 10 from the catalytic zinc and the conserved backbone atoms (N, C α , C, and O) of said amino acids of not more than 1.5Å (or more preferably, not more than 1.0Å, or most preferably, not more than 0.5Å).

Unless otherwise indicated, "protein" or "molecule" shall include a protein, protein domain, polypeptide or peptide.

15 "Structural coordinates" are the Cartesian coordinates corresponding to an atom's spatial relationship to other atoms in a molecule or molecular complex. Structural coordinates may be obtained using x-ray crystallography techniques or NMR techniques, or may be derived using molecular replacement analysis or homology modeling. Various software 20 programs allow for the graphical representation of a set of structural coordinates to obtain a three dimensional representation of a molecule or molecular complex. The structural coordinates of the present invention may be modified from the original sets provided in Figures 4 or 5 by mathematical manipulation, such as by inversion or integer additions or subtractions. As 25 such, it is recognized that the structural coordinates of the present invention are relative, and are in no way specifically limited by the actual x, y, z coordinates of Figures 4 and 5. Further, it is recognized that the structural coordinates taken from Figure 5 may be from either molecule of MMP-13 catalytic fragment in the MMP-13:Compound A crystal (i.e., from A-13 or B-13).

30 An "agent" shall include a protein, polypeptide, peptide, nucleic acid, including DNA or RNA, molecule, compound, antibiotic or drug.

-11-

“Root mean square deviation” is the square root of the arithmetic mean of the squares of the deviations from the mean, and is a way of expressing deviation or variation from the structural coordinates described herein.

It will be obvious to the skilled practitioner that the numbering of the amino acid residues in the various isoforms of MMP-13 or in MMP-13 analogues covered by the present invention may be different than that set forth herein, or may contain certain conservative amino acid substitutions that yield the same three dimensional structures as those defined by Figures 4 or 5 herein. Corresponding amino acids and conservative substitutions in other isoforms or analogues are easily identified by visual inspection of the relevant amino acid sequences or by using commercially available homology software programs. “Conservative substitutions” are those amino acid substitutions which are functionally equivalent to the substituted amino acid residue, either by way of having similar polarity, steric arrangement, or by belonging to the same class as the substituted residue (e.g., hydrophobic, acidic or basic), and includes substitutions having an inconsequential effect on the three dimensional structure of MMP-13 with respect to the use of said structure for the identification and design of MMP-13 activators or inhibitors, for molecular replacement analyses and/or for homology modeling.

An “active site” refers to a region of a molecule or molecular complex that, as a result of its shape and charge potential, favorably interacts or associates with another agent (including, without limitation, a protein, polypeptide, peptide, nucleic acid, including DNA or RNA, molecule, compound, antibiotic or drug). As such, the active site may include both the actual site of substrate cleavage or collagenase activity, as well as certain or all binding sites or pockets adjacent to the site of substrate cleavage that nonetheless may affect MMP-13 activity upon interaction or association with an agent, either by direct interference with the site of substrate cleavage or by indirectly affecting the steric conformation or charge potential of the MMP-13 molecule. The catalytic center of the MMP-13 molecule is characterized by a zinc atom chelated by H119, H123 and H129. MMP-13 binding sites or pockets located to the right of

-12-

the catalytic zinc include S1', S2' and S3'. Binding sites or pockets to the left of the catalytic zinc include S1, S2 and S3.

The present invention relates to the three dimensional structure of human collagenase 3 (MMP-13) or an MMP-13 analogue, and more specifically, to the crystal and solution structures of MMP-13 complexed with an inhibitor, referred to herein as "Compound A", as determined using crystallography, spectroscopy and various computer modeling techniques. The three dimensional solution and crystal structures of the MMP-13:Compound A complex (as disclosed herein at Figures 4 or 5, respectively) and the uncomplexed MMP-13 catalytic fragment (which may be computationally derived from the structural coordinates of Figures 4 or 5) are useful for a number of applications, including, but not limited to, the visualization, identification and characterization of MMP-13 active sites, including the MMP-13 catalytic zinc chelated by H119, H123 and H129, as well as the various MMP-13 binding pockets adjacent to the catalytic zinc of the MMP-13 molecule. The active site structures may then be used to predict the orientation and binding affinity of a designed or selected activator or inhibitor of the MMP-13 protein. Accordingly, the invention is particularly directed to the three dimensional structure of an MMP-13 active site, including but not limited to the S1', S2', S3', S1, S2 and/or S3 binding pockets, taken separately or together with the catalytic zinc of the MMP-13 molecule.

The present invention provides a solution comprising a biologically active catalytic fragment of human collagenase-3 (MMP-13) complexed with Compound A. In a particular embodiment, the catalytic fragment of MMP-13 comprises the amino acid residues of Figure 1, or conservative substitutions thereof. Preferably, the solution provided for herein comprises MMP-13 complexed with Compound A in a 1:1 molar ratio, and more preferably comprises 1 mM MMP-13 in an equimolar complex with Compound A, in a buffer comprising 10mM deuterated Tris-Base, 100mM NaCl, 5mM CaCl₂, 0.1mM ZnCl₂, 2mM NaN₃, and 10 mM deuterated DTT in either 90% H₂O/10% D₂O or 100% D₂O, at a preferred pH of 6.5. The concentration of

-13-

MMP-13:Compound A in the solution should be high enough to yield a good signal-to-noise ratio in the NMR spectrum, but not so high as to result in precipitation or aggregation of the protein. Further, the MMP-13 of the solution may be either ^{15}N enriched or $^{15}\text{N},^{13}\text{C}$ enriched. As exemplified below, NMR spectra from the solution of the present invention are preferably obtained at a temperature of 35°C.

The secondary structure of the catalytic fragment used in the solution of the present invention comprises three alpha helices and a mixed parallel and anti-parallel beta sheet comprising five beta strands, configured in the order $\beta_I, \alpha_A, \beta_{II}, \beta_{III}, \beta_{IV}, \beta_V, \alpha_B$, and α_C . The three alpha helices correspond to residues 28-44 (α_A), 112-123 (α_B) and 153-163 (α_C) of Figure 1, and the five beta strands correspond to residues 83-86 (β_I), 95-100 (β_{II}), 59-66 (β_{III}), 14-20 (β_{IV}), and 49-53 (β_V) of Figure 1, respectively. While the solution of the present invention comprises MMP-13 in a 1:1 molar ratio with Compound A, it is understood that one of ordinary skill in the art may devise additional solutions using alternate inhibitors or ligands in the appropriate molar concentrations, thereby preventing the auto-degradation of MMP-13 and creating a solution of sufficient stability and concentration to obtain a usable NMR spectrum.

The protein used in the solution of the present invention includes MMP-13, as well as MMP-13 analogues, where said protein comprises an active site characterized by the three dimensional structure comprising the relative structural coordinates of the catalytic zinc and amino acid residues L81, L82, L115, V116, H119, L136 and I140 (or conservative substitutions thereof) according to the solution coordinates of Figure 4, \pm a root mean square deviation from the catalytic zinc and the conserved backbone atoms of said amino acids of not more than 1.5Å, or more preferably, not more than 1.0Å, or most preferably, not more than 0.5Å. These residues comprise the residues most closely associated with Compound A in the MMP-13:Compound A complex, as determined from the observed NOEs between MMP-13 and Compound A (Table 1).

-14-

Alternatively, a protein used in the solution of the present invention comprises an active site characterized by a catalytic zinc, a beta strand (comprising amino acid residues N14, L15, T16, Y17, R18, I19, and V20 or conservative substitutions thereof), a Ca^{2+} binding loop (comprising amino acid residues F75, D76, G77, P78, and S79 or conservative substitutions thereof), an alpha helix (comprising amino acid residues N112, L113, F114, L115, V116, A117, A118, H119, E120, F121, G122, and H123 or conservative substitutions thereof) and a random coil region (comprising amino acid residues P139, I140, and Y141 or conservative substitutions thereof), or, more particularly,

characterized by a three dimensional structure comprising the relative structural coordinates of the catalytic zinc and the amino acid residues N14, L15, T16, Y17, R18, I19, V20, F75, D76, G77, P78, S79, N112, L113, F114, L115, V116, A117, A118, H119, E120, F121, G122, H123, P139, I140, and Y141 according to Figure 4, or more preferably, where said three dimensional structure further comprises the relative structural coordinates of amino acid residues G80, L81, L82, A83, H84, A85, K109, G110, Y111, S124, L125, G126, L127, D128, H129, S130, K131, D132, P133, G134, A135, L136, M137, F138, T142, Y143, T144, and G145 according to Figure 4 (incorporating an S1' pocket in the active site), or most preferably, where said three dimensional structure still further comprises the relative structural coordinates of F149 and P152 according to Figure 4 (further defining a hydrophobic area at the bottom of the S1' pocket), including, in each case, conservative substitutions of said amino acids and, in each case, \pm a root mean square deviation from the catalytic zinc and the conserved backbone atoms (N, $\text{C}\alpha$, C, and O) of said amino acids of not more than 1.5\AA (or more preferably, not more than 1.0\AA , or most preferably, not more than 0.5\AA). Finally, in the most preferred embodiment, the protein used in the solution of the present invention comprises the complete structural coordinates according to Figure 4, \pm a root mean square deviation from the conserved backbone atoms of said amino acids (or conservative substitutions thereof) of not more than 1.5\AA (or more preferably, not more than 1.0\AA , and most preferably, not more than 0.5\AA).

-15-

Also provided by the present invention is a crystallized catalytic fragment of MMP-13 complexed with Compound A. The crystal of the present invention effectively diffracts X-rays for the determination of the structural coordinates of the MMP-13:Compound A complex, and is characterized as being
5 in orthorhombic form with space group P21212, and having unit cell parameters of $a=108.3\text{\AA}$, $b=79.8\text{\AA}$, and $c=36.1\text{\AA}$. Further, the crystal complex of the present invention consists of two molecules of MMP-13:Compound A complex in the asymmetric crystal unit.

In a preferred embodiment, the MMP-13 of the crystal complex of
10 the present invention comprises the amino acid residues of Figure 1 (or conservative substitutions thereof), and is characterized by a secondary structure comprising three alpha helices and a mixed parallel and anti-parallel beta sheet comprising five beta strands, configured in the order β_I , α_A , β_{II} , β_{III} , β_{IV} , β_V , α_B , and α_C . Further, the three alpha helices preferably correspond to
15 residues 28-44 (α_A), 112-123 (α_B) and 153-163 (α_C) of Figure 1, and the five beta strands correspond to residues 83-86 (β_I), 95-100 (β_{II}), 59-66 (β_{III}), 14-20 (β_{IV}), and 49-53 (β_V) of Figure 1, respectively.

The protein used in the crystal or crystal complex of the present invention includes MMP-13, as well as MMP-13 analogues, where said protein
20 comprises an active site characterized by the three dimensional structure comprising the relative structural coordinates of the catalytic zinc and amino acid residues L81, L82, L115, V116, H119, L136 and I140 (or conservative substitutions thereof) according to the crystal coordinates of Figure 5, \pm a root mean square deviation from the catalytic zinc and the conserved backbone
25 atoms of said amino acids of not more than 1.5\AA , or more preferably, not more than 1.0\AA , or most preferably, not more than 0.5\AA .

Alternatively, a protein used in the crystal or crystal complex of the present invention comprises an active site characterized by a catalytic zinc, a beta strand (comprising amino acid residues N14, L15, T16, Y17, R18, I19, and
30 V20 or conservative substitutions thereof), a Ca^{2+} binding loop (comprising amino acid residues F75, D76, G77, P78, and S79 or conservative substitutions

-16-

thereof), an alpha helix (comprising amino acid residues N112, L113, F114, L115, V116, A117, A118, H119, E120, F121, G122, and H123 or conservative substitutions thereof) and a random coil region (comprising amino acid residues P139, I140, and Y141 or conservative substitutions thereof), or, more

5 particularly, characterized by a three dimensional structure comprising the relative structural coordinates of the catalytic zinc and amino acid residues N14, L15, T16, Y17, R18, I19, V20, F75, D76, G77, P78, S79, N112, L113, F114, L115, V116, A117, A118, H119, E120, F121, G122, H123, P139, I140, and Y141 according to Figure 5, or more preferably, where said three dimensional

10 structure further comprises the relative structural coordinates of amino acid residues G80, L81, L82, A83, H84, A85, K109, G110, Y111, S124, L125, G126, L127, D128, H129, S130, K131, D132, P133, G134, A135, L136, M137, F138, T142, Y143, T144, and G145 according to Figure 5 (incorporating an S1' pocket in the active site), or most preferably, where said three dimensional structure

15 still further comprises the relative structural coordinates of F149 and P152 according to Figure 5 (further defining a hydrophobic area at the bottom of the S1' pocket), in each case, including conservative substitutions of the said amino acids and, in each case, \pm a root mean square deviation from the catalytic zinc and the conserved backbone atoms of said amino acids of not more than 1.5Å

20 (or more preferably, not more than 1.0Å, or most preferably, not more than 0.5Å).

Finally, in the most preferred embodiment, the protein used in the crystal of the present invention comprises the complete structural coordinates according to Figure 5, \pm a root mean square deviation from the conserved

25 backbone atoms of said amino acids (or conservative substitutions thereof) of not more than 1.5Å (or more preferably, not more than 1.0Å, and most preferably, not more than 0.5Å).

Molecular modeling methods known in the art may be used to identify an active site or binding pocket of the MMP-13 molecule, MMP-13

30 molecular complex, or an MMP-13 analogue. Specifically, the structural coordinates provided by the present invention may be used to characterize a

-17-

three dimensional model of the MMP-13 molecule, molecular complex or MMP-13 analogue. From such a model, putative active sites may be computationally visualized, identified and characterized based on the surface structure of the molecule, surface charge, steric arrangement, the presence of reactive amino acids, regions of hydrophobicity or hydrophilicity, etc. Such putative active sites may be further refined using chemical shift perturbations of spectra generated from various and distinct MMP-13 complexes, competitive and non-competitive inhibition experiments, and/or by the generation and characterization of MMP-13 mutants to identify critical residues or characteristics of the active site.

10 The identification of putative active sites of a molecule or molecular complex is of great importance, as most often the biological activity of a molecule or molecular complex results from the interaction between an agent and one or more active sites of the molecule or molecular complex. Accordingly, the active sites of a molecule or molecular complex are the best targets to use in the design or selection of activators or inhibitors that affect the activity of the molecule or molecular complex.

 The present invention is directed to an active site of MMP-13 or an MMP-13 analogue, that, as a result of its shape, reactivity, charge potential, etc., favorably interacts or associates with another agent (including, without limitation, a protein, polypeptide, peptide, nucleic acid, including DNA or RNA, molecule, compound, antibiotic or drug). As such, the active site of the present invention includes both the actual site of substrate cleavage or collagenase activity (the catalytic zinc chelated by H119, H123, and H129), as well as binding sites or pockets adjacent to the site of substrate cleavage (i.e., S1', S2', S3', S1, S2, and/or S3) that may nonetheless affect MMP-13 activity upon interaction or association with an agent, either by direct interference with the site of substrate cleavage or by indirectly affecting the steric conformation or charge potential of the MMP-13 molecule. Accordingly, the present invention is directed to an active site of the MMP-13 molecule characterized by a zinc atom chelated by H119, H123 and H129, and preferably the S1' binding pocket to the right of the catalytic zinc.

-18-

In an alternate embodiment, the active site of the present invention is characterized by the three dimensional structure comprising the relative structural coordinates of the catalytic zinc and amino acid residues L81, L82, L115, V116, H119, L136 and I140 (or conservative substitutions thereof) according to the solution or crystal coordinates of Figures 4 or 5, respectively, in each case, \pm a root mean square deviation from the catalytic zinc and the conserved backbone atoms of said amino acids of not more than 1.5Å, or more preferably, not more than 1.0Å, or most preferably, not more than 0.5Å.

Alternatively, the active site of the present invention is characterized by a catalytic zinc, a beta strand (comprising amino acid residues N14, L15, T16, Y17, R18, I19, and V20 or conservative substitutions thereof), a Ca²⁺ binding loop (comprising amino acid residues F75, D76, G77, P78, and S79 or conservative substitutions thereof), an alpha helix (comprising amino acid residues N112, L113, F114, L115, V116, A117, A118, H119, E120, F121, G122, and H123 or conservative substitutions thereof) and a random coil region (comprising amino acid residues P139, I140, and Y141 or conservative substitutions thereof), or, more particularly, is characterized by a three dimensional structure comprising the relative solution or crystal structural coordinates of the catalytic zinc and amino acid residues N14, L15, T16, Y17, R18, I19, V20, F75, D76, G77, P78, S79, N112, L113, F114, L115, V116, A117, A118, H119, E120, F121, G122, H123, P139, I140, and Y141 according to Figures 4 or 5, respectively, or more preferably, where said three dimensional structure further comprises the relative solution or crystal structural coordinates of amino acid residues G80, L81, L82, A83, H84, A85, K109, G110, Y111, S124, L125, G126, L127, D128, H129, S130, K131, D132, P133, G134, A135, L136, M137, F138, T142, Y143, T144, and G145 according to Figures 4 or 5, or most preferably, where said three dimensional structure still further comprises the relative solution or crystal structural coordinates of F149 and P152 according to Figures 4 or 5, in each case, including conservative substitutions of said amino acids, and in each case, \pm a root mean square deviation from the catalytic zinc and the conserved backbone atoms of said amino acids of not more than 1.5Å

-19-

(or more preferably, not more than 1.0Å, or most preferably, not more than 0.5Å).

In order to use the structural coordinates generated for a crystal or solution structure of the present invention as set forth in Figures 4 and 5, respectively, it is often necessary to display the relevant coordinates as, or convert them to, a three dimensional shape or graphical representation, or to otherwise manipulate them. For example, a three dimensional representation of the structural coordinates is often used in rational drug design, molecular replacement analysis, homology modeling, and mutation analysis. This is typically accomplished using any of a wide variety of commercially available software programs capable of generating three dimensional graphical representations of molecules or portions thereof from a set of structural coordinates. Examples of said commercially available software programs include, without limitation, the following: GRID (Oxford University, Oxford, UK); MCSS (Molecular Simulations, San Diego, CA); AUTODOCK (Scripps Research Institute, La Jolla, CA); DOCK (University of California, San Francisco, CA); Flo99 (Thistlesoft, Morris Township, NJ); Ludi (Molecular Simulations, San Diego, CA); QUANTA (Molecular Simulations, San Diego, CA); Insight (Molecular Simulations, San Diego, CA); SYBYL (TRIPOS, Inc., St. Louis, MO); and LEAPFROG (TRIPOS, Inc., St. Louis, MO).

For storage, transfer and use with such programs, a machine, such as a computer, is provided for that produces a three dimensional representation of the MMP-13 molecule, a portion thereof (such as an active site or a binding site), a MMP-13 molecular complex, or an MMP-13 analogue. The machine of the present invention comprises a machine-readable data storage medium comprising a data storage material encoded with machine-readable data. Machine-readable storage media comprising data storage material include conventional computer hard drives, floppy disks, DAT tape, CD-ROM, and other magnetic, magneto-optical, optical, floptical and other media which may be adapted for use with a computer. The machine of the present invention also comprises a working memory for storing instructions for processing the

-20-

machine-readable data, as well as a central processing unit (CPU) coupled to the working memory and to the machine-readable data storage medium for the purpose of processing the machine-readable data into the desired three dimensional representation. Finally, the machine of the present invention

5 further comprises a display connected to the CPU so that the three dimensional representation may be visualized by the user. Accordingly, when used with a machine programmed with instructions for using said data, *e.g.*, a computer loaded with one or more programs of the sort identified above, the machine provided for herein is capable of displaying a graphical three-dimensional

10 representation of any of the molecules or molecular complexes, or portions of molecules of molecular complexes, described herein.

In one embodiment of the invention, the machine-readable data comprises the relative structural coordinates of the catalytic zinc and amino acid residues L81, L82, L115, V116, H119, L136 and I140 according to Figures 4 or

15 5, in each case, including conservative substitutions thereof, and in each case, \pm a root mean square deviation from the catalytic zinc and the conserved backbone atoms of said amino acids of not more than 1.5Å (or more preferably, not more than 1.0Å, and most preferably, not more than 0.5Å), wherein said structural coordinates characterize an active site of MMP-13 or an MMP-13

20 analogue.

In an alternate preferred embodiment, the machine-readable data comprises the structural coordinates of the catalytic zinc and amino acid residues N14, L15, T16, Y17, R18, I19, V20, F75, D76, G77, P78, S79, N112, L113, F114, L115, V116, A117, A118, H119, E120, F121, G122, H123, P139,

25 I140, and Y141 according to Figures 4 or 5, in each case, including conservative substitutions thereof, and in each case, \pm a root mean square deviation from the catalytic zinc and the conserved backbone atoms of said amino acids of not more than 1.5Å (or more preferably, not more than 1.0Å, and most preferably, not more than 0.5Å). In an even more preferred embodiment, the machine-

30 readable data further comprises the relative structural coordinates of amino acid residues G80, L81, L82, A83, H84, A85, K109, G110, Y111, S124, L125, G126,

-21-

L127, D128, H129, S130, K131, D132, P133, G134, A135, L136, M137, F138, T142, Y143, T144, and G145 according to Figures 4 or 5, or most preferably, still further comprises the relative structural coordinates of F149 and P152 according to Figures 4 or 5, in each case, including conservative substitutions of
5 said amino acids, and in each case, \pm a root mean square deviation from the catalytic zinc and the conserved backbone atoms of said amino acids of not more than 1.5Å (or more preferably, not more than 1.0Å, or most preferably, not more than 0.5Å).

Finally, it is most preferred that the machine-readable data
10 comprise the relative structural coordinates of all residues constituting the MMP-13 catalytic fragment according to Figures 4 or 5, in each case, \pm a root mean square deviation from the conserved backbone atoms of said amino acids of not more than 1.5Å. In each case, the noted embodiments comprise conservative substitutions of the noted residues resulting in same structural
15 coordinates within the stated root mean square deviation.

The structural coordinates of the present invention permit the use of various molecular design and analysis techniques in order to (i) solve the three dimensional structures of related molecules, molecular complexes or MMP-13 analogues, and (ii) to design, select, and synthesize chemical agents
20 capable of favorably associating or interacting with an active site of an MMP-13 molecule or MMP-13 analogue, wherein said chemical agents potentially act as activators or inhibitors of MMP-13 or of an MMP-13 analogue.

More specifically, the present invention provides a method for determining the molecular structure of a molecule or molecular complex whose
25 structure is unknown, comprising the steps of obtaining crystals or a solution of the molecule or molecular complex whose structure is unknown, and then generating x-ray diffraction data from the crystallized molecule or molecular complex, and/or generating NMR data from the solution of the molecule or molecular complex. The x-ray diffraction data from the molecule or molecular
30 complex whose structure is unknown is then compared to the x-ray diffraction data obtained from the MMP-13:Compound A crystal of the present invention.

-22-

Alternatively, the NMR data from the molecule or molecular structure whose structure is unknown is then compared with the NMR data obtained from the MMP-13:Compound A solution of the present invention. Then, molecular replacement analysis is used to conform the three dimensional structure
5 determined from the MMP-13:Compound A crystal of solution of the present invention to the x-ray diffraction data from the unknown molecule or molecular complex, or, alternatively, 2D, 3D and 4D isotope filtering, editing and triple resonance NMR techniques are used to conform the three dimensional structure determined from the MMP-13:Compound A solution of the present invention to
10 the NMR data from the solution molecule or molecular complex.

Molecular replacement analysis uses a molecule having a known structure as a starting point to model the structure of an unknown crystalline sample. This technique is based on the principle that two molecules which have similar structures, orientations and positions will diffract x-rays similarly. A
15 corresponding approach to molecular replacement is applicable to modeling an unknown solution structure using NMR technology. The NMR spectra and resulting analysis of the NMR data for two similar structures will be essentially identical for regions of the proteins that are structurally conserved, where the NMR analysis consists of obtaining the NMR resonance assignments and the
20 structural constraint assignments, which may contain hydrogen bond, distance, dihedral angle, coupling constant, chemical shift and dipolar coupling constant constraints. The observed differences in the NMR spectra of the two structures will highlight the differences between the two structures and identify the corresponding differences in the structural constraints. The structure
25 determination process for the unknown structure is then based on modifying the NMR constraints from the known structure to be consistent with the observed spectral differences between the NMR spectra.

Accordingly, in one non-limiting embodiment of the invention, the resonance assignments for the MMP-13:Compound A complex provide the
30 starting point for resonance assignments of MMP-13 in a new MMP-13:"unsolved agent" complex. Chemical shift perturbances in two dimensional

-23-

$^{15}\text{N}/^1\text{H}$ spectra can be observed and compared between the MMP-13:Compound A complex and the new MMP-13:agent complex. In this way, the affected residues may be correlated with the three dimensional structure of MMP-13 as provided by the relevant residues of Figure 4. This effectively identifies the region of the MMP-13:agent complex that has incurred a structural change relative to the MMP-13:Compound A complex. The ^1H , ^{15}N , ^{13}C and ^{13}CO NMR resonance assignments corresponding to both the sequential backbone and side-chain amino acid assignments of MMP-13 may then be obtained and the three dimensional structure of the new MMP-13:agent complex may be generated using standard 2D, 3D and 4D triple resonance NMR techniques and NMR assignment methodology, using the MMP-13:Compound A structure, resonance assignments and structural constraints as a reference. Various computer fitting analyses of the new agent with the three dimensional model of MMP-13 may be performed in order to generate an initial three dimensional model of the new agent complexed with MMP-13, and the resulting three dimensional model may be refined using standard experimental constraints and energy minimization techniques in order to position and orient the new agent in association with the three dimensional structure of MMP-13.

The present invention further provides that the structural coordinates of the present invention may be used with standard homology modeling techniques in order to determine the unknown three-dimensional structure of a molecule or molecular complex. Homology modeling involves constructing a model of an unknown structure using structural coordinates of one or more related protein molecules, molecular complexes or parts thereof (i.e., active sites). Homology modeling may be conducted by fitting common or homologous portions of the protein whose three dimensional structure is to be solved to the three dimensional structure of homologous structural elements in the known molecule, specifically using the relevant (i.e., homologous) structural coordinates provided by Figures 4 and/or 5 herein. Homology may be determined using amino acid sequence identity, homologous secondary structure elements, and/or homologous tertiary folds. Homology modeling can

include rebuilding part or all of a three dimensional structure with replacement of amino acids (or other components) by those of the related structure to be solved.

Accordingly, a three dimensional structure for the unknown
5 molecule or molecular complex may be generated using the three dimensional structure of the MMP-13:Compound A complex of the present invention, refined using a number of techniques well known in the art, and then used in the same fashion as the structural coordinates of the present invention, for instance, in applications involving molecular replacement analysis, homology modeling, and
10 rational drug design.

Determination of the three dimensional structure of MMP-13 and its catalytic active site as disclosed herein is critical to the rational identification and/or design of therapeutic agents that may act as inhibitors or activators of MMP-13 enzymatic activity. Alternatively, using conventional drug assay
15 techniques, the only way to identify such an agent is to screen thousands of test compounds, either in culture or by administration to suitable animal models in a laboratory setting, until an agent having the desired inhibitory or activating effect on a target compound is identified. Necessarily, such conventional screening methods are expensive, time consuming, and do not elucidate the
20 method of action of the identified agent on the target compound.

However, advancing X-ray, spectroscopic and computer modeling technologies allow researchers to visualize the three dimensional structure of a targeted compound. Using such a three dimensional structure, researchers identify putative binding sites and then identify or design agents to interact with
25 these binding sites. These agents are then screened for an activating or inhibitory effect upon the target molecule. In this manner, not only are the number of agents to be screened for the desired activity greatly reduced, but the mechanism of action on the target compound is better understood.

Accordingly, the present invention further provides a method for
30 identifying a potential inhibitor or activator of MMP-13, comprising the steps of using a three dimensional structure of MMP-13 as defined by the relative

-25-

structural coordinates of amino acids encoding MMP-13 to design or select a potential inhibitor or activator, and synthesizing or obtaining said potential inhibitor or activator. The inhibitor or activator may be selected by screening an appropriate database, may be designed *de novo* by analyzing the steric configurations and charge potentials of an empty MMP-13 active site in conjunction with the appropriate software programs, or may be designed using characteristics of known inhibitors or activators to MMP-13 or other collagenases in order to create "hybrid" activators or inhibitors. The method of the present invention is preferably used to design or select inhibitors of MMP-13 activity.

An agent that interacts or associates with an active site of MMP-13 or an MMP-13 analogue may be identified by determining an active site of MMP-13 or of the MMP-13 analogue from a three dimensional model of the MMP-13 or MMP-13 analogue, and performing computer fitting analyses to identify an agent which interacts or associates with said active site. Computer fitting analyses utilize various computer software programs that evaluate the "fit" between the putative active site and the identified agent, by (a) generating a three dimensional model of the putative active site of a molecule or molecular complex using homology modeling or the atomic structural coordinates of the active site, and (b) determining the degree of association between the putative active site and the identified agent. The degree of association may be determined computationally by any number of commercially available software programs, or may be determined experimentally using standard binding assays.

Three dimensional models of the putative active site may be generated using any one of a number of methods known in the art, and include, but are not limited to, homology modeling as well as computer analysis of raw structural coordinate data generated using crystallographic or spectroscopy techniques. Computer programs used to generate such three dimensional models and/or perform the necessary fitting analyses include, but are not limited to: GRID (Oxford University, Oxford, UK), MCSS (Molecular Simulations, San Diego, CA), AUTODOCK (Scripps Research Institute, La Jolla,

CA), DOCK (University of California, San Francisco, CA), Flo99 (Thistlesoft, Morris Township, NJ), Ludi (Molecular Simulations, San Diego, CA), QUANTA (Molecular Simulations, San Diego, CA), Insight (Molecular Simulations, San Diego, CA), SYBYL (TRIPOS, Inc., St. Louis, MO) and LEAPFROG (TRIPOS, Inc.,
5 St. Louis, MO).

In a preferred method of the present invention, the identified active site of MMP-13 or the MMP-13 analogue comprises a catalytic zinc, a beta strand, a Ca^{2+} binding loop, an alpha helix and a random coil region. More preferably, the identified active site comprises a catalytic zinc, a beta strand
10 comprising residues N14, L15, T16, Y17, R18, I19, and V20 according to Figure 1 (or conservative substitutions thereof), a Ca^{2+} binding loop comprising residues F75, D76, G77, P78, and S79 according to Figure 1 (or conservative substitutions thereof), an alpha helix comprising residues N112, L113, F114, L115, V116, A117, A118, H119, E120, F121, G122, and H123 according to
15 Figure 1 (or conservative substitutions thereof), and a random coil region comprising residues P139, I140, and Y141 according to Figure 1 (or conservative substitutions thereof).

More specifically, the identified active site of the present method comprises the relative structural coordinates of the catalytic zinc and amino acid
20 residues N14, L15, T16, Y17, R18, I19, V20, F75, D76, G77, P78, S79, N112, L113, F114, L115, V116, A117, A118, H119, E120, F121, G122, H123, P139, I140, and Y141 according to Figures 4 or 5, in each case, including conservative substitutions of said amino acids, and in each case, \pm a root mean square deviation from the catalytic zinc and the conserved backbone atoms of said
25 amino acids of not more than 1.5\AA (or more preferably, not more than 1.0\AA , or most preferably, not more than 0.5\AA). In an alternate preferred embodiment, the identified active site further comprises the relative structural coordinates of amino acid residues G80, L81, L82, A83, H84, A85, K109, G110, Y111, S124, L125, G126, L127, D128, H129, S130, K131, D132, P133, G134, A135, L136,
30 M137, F138, T142, Y143, T144, and G145 according to Figures 4 or 5, in each case, including conservative substitutions of said amino acids, and in each case,

-27-

± a root mean square deviation from the conserved backbone atoms of said amino acids of not more than 1.5Å (or more preferably, not more than 1.0Å, or most preferably, not more than 0.5Å). In yet a third preferred embodiment, the identified active site of the present method further comprises the relative structural coordinates of amino acid residues F149 and P152 according to Figures 4 or 5, in each case, including conservative substitutions of said amino acids, and in each case, ± a root mean square deviation from the conserved backbone atoms of said amino acids of not more than 1.5Å (or more preferably, not more than 1.0Å, or most preferably, not more than 0.5Å). Embodiments comprising conservative substitutions of the noted amino acids result in the same structural coordinates of the corresponding residues in Figures 4 or 5 within the stated root mean square deviation.

The effect of such an agent identified by computer fitting analyses on MMP-13 (or MMP-13 analogue) activity may be further evaluated computationally, or experimentally by contacting the identified agent with MMP-13 (or an MMP-13 analogue) and measuring the effect of the agent on the enzyme's activity. Depending upon the action of the agent on the active site of MMP-13, the agent may act either as an inhibitor or activator of MMP-13 activity. Standard enzymatic assays may be performed and the results analyzed to determine whether the agent is an inhibitor of MMP-13 activity (i.e., the agent may reduce or prevent binding affinity between MMP-13 and the relevant substrate, and thereby reduce the level or rate of MMP-13 activity compared to baseline), or an activator of MMP-13 activity (i.e., the agent may increase binding affinity between MMP-13 and the relevant substrate, and thereby increase the level or rate of MMP-13 activity compared to baseline). Further tests may be performed to evaluate the selectivity of the identified agent to MMP-13 with regard to the other metalloproteinases.

Agents designed or selected to interact with MMP-13 must be capable of both physically and structurally associating with MMP-13 via various covalent and/or non-covalent molecular interactions, and of assuming a three

-28-

dimensional configuration and orientation that complements the relevant active site of the MMP-13 molecule.

Accordingly, using these criteria, the structural coordinates of the MMP-13:Compound A complex as disclosed herein, and/or structural

5 coordinates derived therefrom using molecular replacement analysis or homology modeling, agents may be designed to increase either or both of the potency and selectivity of known inhibitors or activators, either by modifying the structure of known inhibitors or activators or by designing new agents *de novo* via computational inspection of the three dimensional configuration and

10 electrostatic potential of an MMP-13 active site.

Accordingly, in one embodiment of the invention, the structural coordinates of Figures 4 or 5 of the present invention, or structural coordinates derived therefrom using molecular replacement or homology modeling techniques as discussed above, are used to screen a database for agents that

15 may act as potential inhibitors or activators of MMP-13 activity (or the activity of MMP-13 analogues). Specifically, the obtained structural coordinates of the present invention are read into a software package and the three dimensional structure is analyzed graphically. A number of computational software packages may be used for the analysis of structural coordinates, including, but

20 not limited to, Sybyl (Tripos Associates), QUANTA and XPLOR (Brunger, A.T., (1993) XPLOR Version 3.1 Manual, Yale University, New Haven, CT). Additional software programs check for the correctness of the coordinates with regard to features such as bond and atom types. If necessary, the three dimensional structure is modified and then energy minimized using the

25 appropriate software until all of the structural parameters are at their equilibrium/optimal values. The energy minimized structure is superimposed against the original structure to make sure there are no significant deviations between the original and the energy minimized coordinates.

The energy minimized coordinates of MMP-13 complexed with a

30 "solved" inhibitor or activator are then analyzed and the interactions between the solved ligand and MMP-13 are identified. The final MMP-13 structure is

modified by graphically removing the solved inhibitor or activator so that only MMP-13 and a few residues of the solved agent are left for analysis of the binding site cavity. QSAR and SAR analysis and/or conformational analysis may be carried out to determine how other inhibitors or activators compare to the

5 solved inhibitor or activator. The solved agent may be docked into the uncomplexed structure's binding site to be used as a template for data base searching, using software to create excluded volume and distance restrained queries for the searches. Structures qualifying as hits are then screened for activity using standard assays and other methods known in the art.

10 Further, once the specific interaction is determined between the solved inhibitor or activator, docking studies with different inhibitors or activators allow for the generation of initial models of new inhibitors or activators in complex with MMP-13. The integrity of these new models may be evaluated a number of ways, including constrained conformational analysis

15 using molecular dynamics methods (*i.e.*, where both MMP-13 and the complexed activator or inhibitor are allowed to sample different three dimensional conformational states until the most favorable state is reached or found to exist between the protein and the complexed agent). The final structure as proposed by the molecular dynamics analysis is analyzed visually to

20 make sure that the model is in accord with known experimental SAR based on measured binding affinities. Once models are obtained of the original solved agent bound to MMP-13 and computer models of other molecules bound to MMP-13, strategies are determined for designing modifications into the activators or inhibitors to improve their activity and/or enhance their selectivity.

25 Once an MMP-13 binding agent has been optimally selected or designed, as described above, substitutions may then be made in some of its atoms or side groups in order to improve or modify its selectivity and binding properties. Generally, initial substitutions are conservative, *i.e.*, the replacement group will have approximately the same size, shape, hydrophobicity and charge

30 as the

-30-

original group. Such substituted chemical compounds may then be analyzed for efficiency of fit to MMP-13 by the same computer methods described in detail above.

Alternatively, the present invention provides a method for

5 identifying a potential inhibitor or activator that is selective for one or more members of the matrix metalloproteinase family except MMP-13, comprising the steps of (i) using the three dimensional structures of MMP-13 and the desired target matrix metalloproteinase(s) as defined by the relative structural coordinates of amino acids encoding MMP-13 and the target matrix

10 metalloproteinase(s) in order to design or select such a potential inhibitor or activator, and (ii) synthesizing or obtaining said potential inhibitor or activator. In this case, the potential inhibitor or activator is designed to incorporate chemical or steric features favorable for association with an active site of the desired matrix metalloproteinase(s) and unfavorable for association with an

15 MMP-13 active site, preferably where said active site comprises the MMP-13 S1' pocket. The inhibitor or activator may be selected by screening an appropriate database, may designed *de novo* by analyzing the steric configurations and charge potentials of empty MMP-13/matrix metalloproteinase active sites in conjunction with the appropriate software programs, or may be designed using

20 characteristics of known inhibitors or activators to MMP-13 or other collagenases in order to create "hybrid" activators or inhibitors.

Various molecular analysis and rational drug design techniques are further disclosed in U.S. Patent Nos. 5,834,228, 5,939,528 and 5,865,116, as well as in PCT Application No. PCT/US98/16879, published as WO

25 99/09148, the contents of which are hereby incorporated by reference.

The present invention may be better understood by reference to the following non-limiting Examples. The following Examples are presented in order to more fully illustrate the preferred embodiments of the invention, and should in no way be construed as limiting the scope of the present invention.

Example 1

¹H, ¹⁵N and ¹³CO Assignments and Secondary Structure Determination of MMP-13 Complexed with Compound A

5

Methods and Results: The uniform ¹⁵N and ¹³C- labeled 165 amino-acid catalytic fragment of human collagenase-3 (MMP-13) was expressed in *E. coli* strain BL21(DE3) containing the plasmid pProMMP-13 according to a published method (Freije *et al.*, *J. Biol. Chem.* 1994). MMP-13 was purified as previously described (Moy *et al.*, *J. Biomol.* 1997) with minor modifications. N-terminal amino acid sequencing was performed to confirm the protein's identity while the uniform ¹⁵N and ¹³C labeling of MMP-13 was confirmed by MALDI-TOF mass spectrometry (PerSeptive Biosystems). The sulfonamide derivative of the hydroxamic acid compound, N-Hydroxy-2-[(4-methoxy-benzenesulfonyl)-pyridin-3-ylmethyl-amino]-3-methyl-benzamide, was prepared from 2-amino-3-methyl-benzoic acid methyl ester and p-methoxybenzenesulfonyl chloride followed by alkylation with 3-picolyl chloride, hydrolysis (LiOH/THF) to afford the carboxylic acid and conversion to the hydroxamic acid (oxalyl chloride/DMF/NH₂OH). Formation of the HCl salt yielded Compound A as shown in Figure 3.

The NMR samples contained 1 mM of MMP-13 determined spectrophotometrically in a equimolar complex with Compound A in a buffer containing 10 mM deuterated Tris-Base, 100 mM NaCl, 5 mM CaCl₂, 0.1 mM ZnCl₂, 2 mM NaN₃, 10 mM deuterated DTT, in either 90% H₂O/ 10% D₂O or 100% D₂O at pH 6.5. All NMR spectra were recorded at 35°C on a Bruker AMX-2 600 spectrometer equipped with a triple-resonance gradient probe.

Spectra were processed using the NMRPipe software package (Delaglio *et al.*, *J. Biomol. NMR* 1995) and analyzed with PIPP (Garrett *et al.*, *J. Magn. Reson.* 1991), NMRPipe and PEAK-SORT, an in-house software package. The assignments of the ¹H, ¹⁵N, ¹³CO, and ¹³C resonances were based on the following experiments: CBCA(CO)NH, CBCANH, C(CO)NH, HC(CO)NH,

30

-32-

HBHA(CO)NH, HNCO, HCACO, HNHA, HNCA, HCCH-COSY and HCCH-TOCSY (for reviews, see Bax *et al.*, Methods Enzymol. 1994; and Clore & Gronenborn, Methods Enzymol. 1994). The accuracy of the MMP-13 NMR assignments was further confirmed by sequential NOEs in the ^{15}N -edited NOESY-HSQC spectra.

- 5 Prior to analysis of the MMP-13 NMR structure, the structure determination of the inhibitor-free catalytic fragment of MMP-1 has been reported (Moy *et al.*, Biochemistry 1998; Moy *et al.*, J. Biomol. NMR 1997) (30 simulated annealing structures deposited with Protein Data Bank, Accession No. 1AYK; restrained minimized mean structure deposited with Protein Data Bank, 10 Accession No. 2AYK). Because the MMPs are highly autocatalytic, the NMR analysis of the inhibitor-free MMP-1 was accomplished by establishing buffer conditions where the enzyme was still active but the rate of self-cleavage of the enzyme had been diminished. This was achieved by the addition of DTT which significantly diminished self-aggregation of the enzyme and by lowering the pH 15 of the sample to 6.5, just above the pH where the enzyme was known to be inactivated because of the loss of the catalytic zinc. Under these conditions, an MMP-1 NMR sample was typically stable for 1-2 months. Unfortunately this was not the case for MMP-13, the protein rapidly degraded within a few hours which required the use of an inhibitor to assign the MMP-13 NMR resonances.
- 20 The secondary structure of the MMP-13:Compound A complex is based on characteristic NOE data involving the NH, H α and H β protons from ^{15}N -edited NOESY-HSQC and ^{13}C -edited NOESY-HMQC spectra, $^3\text{JHN}\alpha$ coupling constants from HNHA, slowly exchanging NH protons and $^{13}\text{C}\alpha$ and $^{13}\text{C}\beta$ secondary chemical shifts (for reviews, see Wishart & Sykes, Methods Enzymol. 1994; and Wuthrich, NMR of Proteins and Nucleic Acids, John Wiley & Sons, 25 New York 1986). It was determined that the MMP-13 NMR structure in the complex is composed of three α -helices corresponding to residues 28-44 (α_a), 112-123 (α_b) and 153-163 (α_c) and a mixed parallel and anti-parallel β -sheet consisting of 5 strands corresponding to residues 83-86 (β_1), 95-100 (β_2), 59-66 30 (β_3), 14-20 (β_4) and 49-53 (β_5). This is essentially identical to the secondary structure observed for other MMP structures.

-33-

There were three distinct regions in the MMP-13:Compound A spectra where the resonance assignments are incomplete. These correspond to residues G70-Y73, P87-N91 and T144-H148. Residues T144-H148 correspond to part of the dynamic loop region previously seen in the MMP-1 structure (Moy
5 *et al.*, *J. Biomol. NMR* 1997). This suggests a similar dynamic profile for this region in the MMP-13 structure even in the presence of a high-affinity inhibitor ($IC_{50} = 33$ nM). Residues P87 to N91 contain a cluster of prolines which disrupt the sequential assignment process because of the missing NH. Residues G70 to Y73 correspond to a loop region in the vicinity of the structural zinc which was
10 readily assigned in the MMP-1 structure. The backbone and side-chain 1H , ^{15}N , ^{13}C , and ^{13}CO assignments are essentially complete for the remainder of the protein.

Example 2

15 High Resolution Solution Structure of the Catalytic Fragment of MMP-13 Complexed with Compound A

Materials and Methods:

Preparation of Compound A: The sulfonamide derivative of the hydroxamic acid
20 compound, Compound A, was prepared according to the procedure noted in Example 1 to yield the compound of Figure 3.

Expression of recombinant ^{15}N and ^{13}C / ^{15}N -labeled MMP-13: A 169-residue C-terminally truncated human collagenase-3 (MMP-13) was expressed in *E. coli*.
25 The coding sequence of a C-terminally truncated procollagenase was amplified by PCR from the plasmid pNot3a, that contains the entire coding sequence of MMP-13 (Frieje, *et al.*, *J. Biol. Chem.* 1994). The PCR primers contained the appropriate restriction sites for ease of cloning. The construct codes for a truncated proMMP-13 with an N-terminal methionine added and a C-terminal
30 proline at residue 169 of the native proMMP-13 sequence. The PCR amplified DNA fragment was the cloned into pET-21a (+) at the Nde I/Sal I sites ,

-34-

resulting in a recombinant plasmid designated as pProMMP-13. *E. coli* bacteria, BL21(DE3), containing the plasmid pProMMP-13, were grown in LB broth supplemented with 100 µg/ml ampicillin. An overnight culture was diluted 1:20 and grown at 37°C to an A₆₀₀ of 0.6-0.8 with vigorous shaking. Isopropyl β-D-
5 galactoside (IPTG) was added to a final concentration of 1 mM and cultures were shaken for 3 h at 37°C. The cells were harvested by centrifugation (7000 Xg for 15 min) at 4°C, washed with PBS, and frozen at -70°C until further use.

Uniform ¹⁵N and ¹³C- labeled ProMMP-13 was obtained by growing BL21(DE3) *E. coli* in defined media containing 2.0 g/l [¹³C₆, 98%+]D-glucose and 1.0 g/l [¹⁵N, 98%+] ammonium chloride as the sole carbon and
10 nitrogen sources, respectively. In addition, the defined media contained M9 salts (Sambrook, *et al.*, Molecular Cloning: A Laboratory Manual, Cold Spring Harbor Laboratory Press, New York, NY 1989), trace elements, vitamins and 100 µg/ml ampicillin. Conditions for induction and growth were the same as
15 above.

Purification of recombinant ¹⁵N and ¹³C MMP-13: MMP-13 was purified according to Moy *et al.*, *J. Biomol. NMR* 1997, with modifications as follows. Frozen cell pellets were thawed on ice. Cells were resuspended by homogenization in lysis
20 buffer (0.1 M Tricine, pH 8.0, 10 mM EDTA, 2mM DTT, 0.5 mM PMSF). Cells were lysed by French Press (2X) followed by treatment with lysozyme (1 mg/ml; final) at room temperature for 30 min. The lysate was centrifuged at 45,000 x g for 30 minutes. The pellet was washed twice with 50 mM Tricine pH 7.5, 0.2 M NaCl₂, 0.5% Triton X-100, resuspended in fresh urea buffer (20 mM Tricine,
25 pH 7.5, 8 M urea, 0.2% NaN₃, 2 mM DTT) and incubated at room temperature for 1 hour. The urea solubilized protein was centrifuged at 45,000 x g for 30 min and the resultant supernatant was filtered and applied to a Hitrap-Q Sepharose (Pharmacia Biotech) anion exchange column equilibrated in 6 M urea buffer. The column was washed with urea buffer and eluted with a 0-0.25
30 M NaCl linear gradient. Fractions containing proMMP-13 were detected by SDS-PAGE, pooled and quickly diluted into 5-fold excess of renaturing buffer

-35-

(50 mM Tricine, pH 7.5, 0.4 M NaCl, 10 mM CaCl₂, 0.1 mM ZnOAc₂, 0.02% NaN₃). After 2 days of dialysis against 25 volumes of renaturing buffer (with three changes), refolded proMMP-13 was concentrated to about 4-10 mg/ml in a Millipore Biomax 5 concentrator. ProMMP-13 was activated to MMP-13CAT
5 (catalytic domain) by an overnight incubation at 37 °C in the presence of 1 mM *p*-aminophenylmercuric acetate (APMA).

The activated protein is then applied onto a Superdex-75 16/60 gel filtration column equilibrated in 2.5 mM Tris-HCl, pH 7.5, 5 mM CaCl₂, 0.4 M NaCl, 2 mM DTT, 0.02% NaN₃ and 0.05 mM ZnOAc₂. The protein is eluted
10 and fractions containing MMP-13CAT were identified by SDS-PAGE. Peak fractions were pooled and the protein was concentrated in a Millipore Biomax concentrator to about 5 mg/ml and stored at -70 °C. N-terminal amino acid sequencing was performed to confirm the protein's identity. The uniform ¹⁵N and ¹³C labeling of MMP-13-CAT was confirmed by MALDI-TOF mass
15 spectrometry (PerSeptive Biosystems).

NMR Sample Preparation: The MMP-13:Compound A NMR sample contained 1mM ¹⁵N-or ¹⁵N/¹³C-labeled MMP-13 with Compound A in a 1:1 ratio. The sample was prepared by repeated buffer exchange using 20-30ml solution
20 containing 10mM deuterated Tris-Base, 100mM NaCl, 5mM CaCl₂, 0.1mM ZnCl₂, 2mM NaN₃, 10mM deuterated DTT, and 0.2mM Compound A in either 90% H₂O/10 % D₂O or 100% D₂O. Buffer exchange was carried out on a Millipore Ultrafree-15 Centrifugal Filter Unit. Excess Compound A was removed by additional buffer exchanges where Compound A was removed from the
25 buffer.

NMR Data Collection: All spectra were recorded at 35°C on a Bruker AMX-2 600 spectrometer using a gradient enhanced triple-resonance ¹H/¹³C/¹⁵N probe. For spectra recorded in H₂O, water suppression was achieved with the WATERGATE
30 sequence and water-flip back pulses (Piotto, *et al.*, *J. Biomol. NMR* 1992; Grzesiek and Bax, *J. Am. Chem. Soc.* 1993). Quadrature detection in the

indirectly detected dimensions were recorded with States-TPPI hypercomplex phase increment (Marion, *et al.*, J. Magn. Reson. 1989). Spectra were collected with appropriate refocusing delays to allow for 0,0 or -90,180 phase correction.

The resonance assignments and bound conformation of

- 5 Compound A in the MMP-1: Compound A complex were based on the 2D $^{12}\text{C}/^{12}\text{C}$ -filtered NOESY (Petros, *et al.*, FEBS Lett. 1992; Gemmecker, *et al.*, J. Magn. Reson. 1992), 2D $^{12}\text{C}/^{12}\text{C}$ -filtered TOCSY (Petros, *et al.*, FEBS Lett. 1992; Gemmecker, *et al.*, J. Magn. Reson. 1992) and $^{12}\text{C}/^{12}\text{C}$ -filtered COSY experiments (Ikura and Bax, J. Magn. Reson. 1992).
- 10 The MMP-13:Compound A structure is based on the following series of spectra: HNHA (Vuister and Bax, J. Am. Chem. Soc. 1993), HNHB (Archer, *et al.*, J. Magn. Reson. 1992), 3D long-range ^{13}C - ^{13}C correlation (Bax and Popchapsky, J. Magn. Reson. 1992), coupled CT-HCACO (Powers, *et al.*, J. Magn. Reson. 1991; Vuister, *et al.*, J. Am. Chem. Soc. 1992), HACAHB-COSY
- 15 (Grzesiek, *et al.*, J. Amer. Chem. Soc. 1995), 3D ^{15}N - (Mario, *et al.*, Biochemistry 1989; Zuiderweg and Fesik, Biochemistry 1989) and ^{13}C -edited NOESY (Zuiderweg, *et al.*, J. Magn. Reson. 1990; Ikura, *et al.*, J. Magn. Reson. 1990), and 3D ^{13}C -edited/ ^{12}C -filtered NOESY (Lee, *et al.*, FEBS Lett. 1994).
- 20 experiments. The ^{15}N -edited NOESY, ^{13}C -edited NOESY and 3D ^{13}C -edited/ ^{12}C -filtered NOESY experiments were collected with 100 msec, 120 msec and 110 msec mixing times, respectively. The acquisition parameters for each of the experiments used in determining the solution structure of MMP-13 complexed with Compound A were as reported previously (Moy, *et al.*, Biochemistry, 1998).
- 25 Spectra were processed using the NMRPipe software package (Delaglio, *et al.*, J. Biomol. NMR, 1995) and analyzed with PIPP (Garrett, *et al.*, J. Magn. Reson., 1991) on a Sun Sparc Workstation. When appropriate, data processing included a solvent filter, zero-padding data to a power of two, linear predicting back one data point of indirectly acquired data to obtain zero phase
- 30 corrections, linear prediction of additional points for the indirectly acquired dimensions to increase resolution. Linear prediction by the means of the mirror

-37-

image technique was used only for constant-time experiments (Zhu and Bax, J. Magn. Reson., 1992). In all cases data was processed with a skewed sine-bell apodization function and one zero-filling was used in all dimensions.

- 5 *Interproton Distance Restraints:* The NOEs assigned from 3D ^{13}C -edited/ ^{12}C -filtered NOESY and 3D ^{15}N -edited NOESY experiments were classified into strong, medium, and weak corresponding to interproton distance restraints of 1.8-2.7 Å (1.8-2.9 Å for NOEs involving NH protons), 1.8-3.3 Å (1.8-3.5 Å for NOEs involving NH protons), and 1.8-5.0 Å, respectively (Williamson, *et al.*, J. Mol. Biol., 1985; Clore, *et al.*, EMBO J., 1986). Upper distance limits for
- 10 distances involving methyl protons and non-stereospecifically assigned methylene protons were corrected appropriately for center averaging (Wuthrich, *et al.*, J. Mol. Biol., 1983).
- 15 *Torsion Angle Restraints and Stereospecific Assignments.* The β -methylene stereospecific assignments and χ_1 torsion angle restraints were obtained primarily from a qualitative estimate of the magnitude of $^3J_{\alpha\beta}$ coupling constants from the HACAHB-COSY experiment (Grzesiek, *et al.*, J. Am. Chem. Soc., 1992) and $^3J_{N\beta}$ coupling constants from the HNHB experiment (Archer, *et al.*, J. Magn.
- 20 Reson., 1991). Further support for the assignments was obtained from approximate distance restraints for intraresidue NOEs involving NH, C α H, and C β H protons (Powers, *et al.*, Biochemistry, 1993).

- The ϕ and ψ torsion angle restraints were obtained from $^3J_{\text{NH}\alpha}$ coupling constants measured from the relative intensity of H α crosspeaks to the
- 25 NH diagonal in the HNHA experiment (Vuister and Bax, J. Am. Chem. Soc., 1993), from a qualitative estimate of the magnitude of $^3J_{\alpha\beta}$ coupling constants from the HACAHB-COSY experiment (Grzesiek, *et al.*, J. Am. Chem. Soc., 1992) and from approximate distance restraints for intraresidue and sequential NOEs involving NH, C α H, and C β H protons by means of the conformational grid
- 30 search program STEREOSEARCH (Nilges, *et al.*, Biopolymers 1990), as described previously (Kraulis, *et al.*, Biochemistry 1989). $^1J_{\text{C}\alpha\text{H}\alpha}$ coupling

constants obtained from a coupled 3D CT-HCACO spectrum were used to ascertain the presence of non-glycine residues with positive ϕ backbone torsion angles (Vuister, *et al.*, J. Am. Chem. Soc. 1992). The presence of a $^1J_{\alpha H\alpha}$ coupling constant greater than 130 Hz allowed for a minimum ϕ restraint of -2° to -178° .

The Ile and Leu χ_2 torsion angle restraints and the stereospecific assignments for leucine methyl groups were determined from $^3J_{\alpha C\delta}$ coupling constants obtained from the relative intensity of $C\alpha$ and $C\delta$ cross peaks in a 3D long-range ^{13}C - ^{13}C NMR correlation spectrum (Bax, *et al.*, J. Am. Chem. Soc. 1992), in conjunction with the relative intensities of intraresidue NOEs (Powers, *et al.*, Biochemistry 1993). Stereospecific assignments for valine methyl groups were determined based on the relative intensity of intraresidue $NH-C\gamma H$ and $C\alpha H-C\gamma H$ NOEs as described by Zuiderweg *et al.* (1985) (Zuiderweg, *et al.*, Biopolymers 1985). The minimum ranges employed for the ϕ , ψ , and χ torsion angle restraints were $\pm 30^\circ$, $\pm 50^\circ$, and $\pm 20^\circ$, respectively (Kraulis, *et al.*, Biochemistry 1989).

Structure Calculations: The structures were calculated using the hybrid distance geometry-dynamical simulated annealing method of Nilges *et al.* (1988) (Protein Eng.) with minor modifications (Clore, *et al.*, Biochemistry 1990) using the program XPLOR (Brunger, X-Plor Version 3.1 Manual, Yale University, New Haven, CT), adapted to incorporate pseudopotentials for $^3J_{NH\alpha}$ coupling constants (Garrett, *et al.*, J. Magn. Reson. Ser. B 1994), secondary $^{13}C\alpha/^{13}C\beta$ chemical shift restraints (Kuszewski, *et al.*, J. Magn. Reson. Ser. B 1995) and a conformational database potential (Kuszewski, *et al.*, Protein Sci. 1996; Kuszewski, *et al.*, J. Magn. Reson. 1997). The target function that is minimized during restrained minimization and simulated annealing comprises only quadratic harmonic terms for covalent geometry, $^3J_{NH\alpha}$ coupling constants and secondary $^{13}C\alpha/^{13}C\beta$ chemical shift restraints, square-well quadratic potentials for the experimental distance and torsion angle restraints, and a quartic van der Waals term for non-bonded contacts. All peptide bonds were constrained to be

-39-

planar and trans. There were no hydrogen-bonding, electrostatic, or 6-12 Lennard-Jones empirical potential energy terms in the target function.

To prevent the Zn and Ca ions from being expelled during the high-temperature simulated annealing stages of the refinement protocol, a
5 minimal number of distance restraints between the His sidechain and Zn and between backbone atoms and C α were included in the XPLOR distance restraint file based on the observed coordination in the X-ray structures (Lovejoy, *et al.*, Science 1994; Lovejoy, *et al.*, Biochemistry 1994; Spurlino, *et al.*, Proteins: Struct., Funct., Genet. 1994; Borkakoti, *et al.*, Nat. Struct. Biol. 1994).

10 The starting MMP-13:Compound A complex structure for the simulated-annealing protocol was obtained by manually docking Compound A into a homology model for MMP-13. The initial orientation of Compound A in the MMP-13 active site was based on the previously reported MMP-1:CGS-27023A structure (Moy, *et al.*, Biochemistry 1999).

15 Homology modeling methods were utilized to generate a three dimensional model of MMP-13. The linear amino acid sequence corresponding to the catalytic domain of MMP-13 was aligned (SYBYL) with the catalytic domains of MMP-1, MMP-7 and MMP-8 based on the availability of their x-ray crystallographic structures (Bode, *et al.*, EMBO J 1994; Spurlino, *et al.*, Proteins: Struct., Funct., Genet. 1994; Betz, *et al.*, Eur. J. Biochem. 1997; Lovejoy, *et al.*,
20 Nat. Struct. Biol. 1999; Borkakoti, *et al.*, Nat. Struct. Biol. 1994; Browner, *et al.*, Biochemistry 1995). The alignments of MMP-13 with MMP-1 and MMP-8 demonstrated the highest homology where the computed identities are 58.9% and 61.4%, respectively (Figure 2).

25 The X-ray structure of MMP-8 was selected to be used as the template for homology modeling the structure of MMP-13. This decision was based mainly on the sequence alignment shown in Figure 2B where no insertions (labeled “###”) are found in the critical specificity loop (Labeled Underlined and Boldface). In Figure 2A, the region labeled “##” in the
30 specificity loop shows that there is an “insertion” of 2 additional amino acid residues compared to the sequence length of MMP-1. Based on our analysis of

-40-

the alignments, MMP-8 would allow for a more accurate modeling of the inhibitor binding pockets since no predictions have to be made within this loop region.

COMPOSER (SYBYL) was used to construct the initial homology
5 model of MMP-13. The only insertion was a serine (labeled '**' in Figure 2B) at position 32 of MMP-13. The insertion of S32 occurs within a coiled region which is at the entrance of a long alpha helix and about 17 angstroms from the S' specificity loop. The model of MMP-13 was then energy minimized utilizing a set of nested refinement procedures (Chen, *et al.*, J. Biomol. Struct. Dyn. 1995),
10 but where the protein backbone heavy atoms were constrained as close as possible to their original positions.

The MMP-13:Compound A model was then subjected to a 1000 steps of CHARMM minimization with the 5 intramolecular NOE restraints and the 47 distance restraints observed between MMP-13 and Compound A where
15 the coordinates for MMP-13 were kept fixed. This approach approximated the positioning of Compound A in the active site of MMP-13 without distorting the MMP-13 structure. The final structure was exported as a PDB file and used as the starting point for XPLOR simulated annealing protocol where all the residues in the structure were free to move. Since the initial stage of the
20 simulated annealing protocol corresponds to high-temperature dynamics (1500 K) with a relatively weak XPLOR NOE force constant (Ries and Petrides, Biol. Chem. Hoppe-Seyler 1995), the initial MMP-13:Compound A structure does not bias the structure determination process since the structure is effectively free to explore the available conformational space. Additionally, each iteration of the
25 simulated annealing process begins with a random trajectory for the molecular dynamics. The fact that these trajectories differ by upwards of 10 Å assures a distinct exploration of conformational space for the ensemble of MMP-13:Compound A structures determined from the simulated annealing protocol.

Results and Discussion

Compound A Resonance Assignments and Bound Conformation: The primary structure of Compound A along with the proton naming convention is shown in Figure 3. The NMR assignments for Compound A in the MMP-13 complex followed established protocols using 2D ^{12}C -filtering experiments (Petros, *et al.*, FEBS Lett. 1992; Gemmecker, *et al.*, J. Magn. Reson. 1992; Ikura and Bax, J. Am. Chem. Soc. 1992) since the NMR sample was composed of $^{13}\text{C}/^{15}\text{N}$ labeled MMP-13 and unlabeled Compound A. Thus, traditional 2D-NOESY, COSY and TOCSY spectra of Compound A in the presence of MMP-13 yielded straight-forward assignments for Compound A along with assignments for free Compound A (data not shown). The only notable difference in the assignments for free and bound Compound A is the observation of two distinct resonances for 2HB1/2 in the complex (4.91 ppm; 4.67 ppm). The missing resonance in the free Compound A may simply be obscured by water. Also, an observation that the protons on the p-methoxyphenyl ring are degenerate suggests rapid ring flips when complexed to MMP-13. This was also seen with CGS-27023A complexed with both MMP-1 and stromelysin (Gonnella, *et al.*, Bioorg. Med. Chem. 1997; Moy, *et al.*, Biochemistry 1998; Moy, *et al.*, Biochemistry 1999).

Compound A does not adopt a preferred conformation in the absence of MMP-13 as evident by the lack of structural NOEs. Only a minimal number of intramolecular NOEs were observed for Compound A in the MMP-13 complex which were relevant to the bound conformation of Compound A (data not shown). The minimal number of structural NOEs is a result of the Compound A conformation, structure and chemical shift degeneracy. A number of the observed NOEs correspond to a sequential interaction which have no effect on the overall conformation of the inhibitor and were not used in the refinement of Compound A or the complex. The structural intramolecular NOEs observed are primarily between 1HH* and the pyridine ring and between 2HB1/2 and both the p-methoxyphenyl and aryl ring. These NOEs are consistent with the "splayed" conformation previously observed for CGS-27023A bound to both MMP-1 and stromelysin, but the bound conformation of

-42-

Compound A is predominately determined from the intermolecular NOEs between Compound A and MMP-13 (Table 1).

Structure Determination: The NMR structure determination methodology is an iterative procedure where the current state of the structure is used to analyze the ambiguous NOE data. In essence, the structure is used as a distance filter to sort through the ambiguous NOE list where the first structure is determined from unambiguous data. For the refinement of MMP-13, the initial structure was a homology model based on the MMP-8 X-ray structure. This was justified by the overall similarity in previously reported MMP structures and from the secondary structure assignments by NMR for MMP-13. The regular secondary structure elements of MMP-13 were identified from a qualitative analysis of sequential and inter-strand NOEs, NH exchange rates, $^3\text{J}_{\text{HN}\alpha}$ coupling constants (Clare, *et al.*, Crit. Rev. Biochem. Mol. Biol. 1989) and the $^{13}\text{C}\alpha$ and $^{13}\text{C}\beta$ secondary chemical shifts (Spera and Bax, J. Am. Chem. Soc. 1991). The deduced secondary structure is essentially identical to the inhibitor-free MMP-1 NMR structures previously reported.

The final 30 simulated annealing structures calculated for residues 7-164 were based on 3279 experimental NMR restraints, consisting of 2561 approximate interproton distance restraints, 51 distance restraints between MMP-13 and Compound A, 88 distance restraints for 44 backbone hydrogen bonds, 391 torsion angle restraints, 103 $^3\text{J}_{\text{NH}\alpha}$ restraints 123 $\text{C}\alpha$ restraints and 108 $\text{C}\beta$ restraints. Stereospecific assignments were obtained for 81 of the 100 residues with β -methylene protons, for the methyl groups of 5 of the 6 Val residues, and for the methyl groups of 12 of the 13 Leu residues. In addition, 12 out of the 12 Phe residues and 7 out of the 8 Tyr residues were well defined making it possible to assign NOE restraints to only one of the pair of $\text{C}\delta\text{H}$ and $\text{C}\epsilon\text{H}$ protons and to assign a χ^2 torsion angle restraint. Similarly, χ^2 torsion angle restraints were assigned for the three Trp residues. The atomic rms distribution of the 30 simulated annealing structures about the mean coordinate positions for residues 7-164 is $0.43 \pm 0.06 \text{ \AA}$ for the backbone atoms, $0.81 \pm$

-43-

0.09 Å for all atoms, and 0.47 ± 0.04 Å for all atoms excluding disordered surface side chains. The mean standard deviation for the ϕ and ψ backbone torsion angles of residues 7-164 are $6.2 \pm 11.3^\circ$ and $7.1 \pm 11.8^\circ$, respectively. The high quality of the MMP-13 NMR structure is also evident by the results of PROCHECK analysis and by a calculated, large negative value for the Lennard-Jones-van der Waals energy (-695 ± 11 kcal mol⁻¹). For the PROCHECK statistics, an overall G-factor of 0.16 ± 0.16 , a hydrogen bond energy of 0.82 ± 0.05 and only 7.8 ± 1.0 bad contacts per 100 residues are consistent with a good quality structure comparable to ~ 1 Å X-ray structure.

10 The high quality of the MMP-13 NMR structure is also evident by the very small deviations from idealized covalent geometry, by the absence of interproton distance and torsion angle violations greater than 0.1 Å and 1° , respectively and by the fact that most of the backbone torsion angles for non-glycine residues lie within expected regions of the Ramachandran plot (not shown). 91.5% of the residues lie within the most favored region of the Ramachandran ϕ , ψ plot and 7.8% in the additionally allowed regions. ¹JCαHα coupling constants from the coupled CT-HCACO experiment indicated that all non-glycine residues have negative ϕ torsion angles.

20 The quality of the NMR data to properly define the complex is also supported by the well-defined coordinates for Compound A and the active site residues, where the atomic rms distribution is 0.47 ± 0.08 Å and 0.18 ± 0.03 Å for the heavy atoms of Compound A and MMP-13 backbone atoms, respectively.

Description of the MMP-13:Compound A Structure: The overall fold of MMP-13 is essentially identical to previously reported MMP structures (Bode, *et al.*, EMBO J. 1994; Gooley, *et al.*, Nat. Struct. Biol. 1994; Lovejoy, *et al.*, Science 1994; Lovejoy, *et al.*, Ann. N. Y. Acad. Sci. 1994; Lovejoy, *et al.*, Biochemistry 1994; Spurlino, *et al.*, Proteins: Struct., Funct., Genet. 1994; Stams, *et al.*, Nat. Struct. Biol. 1994; Becker, *et al.*, Protein Sci. 1995; Gonnella, *et al.*, Proc. Natl. Acad. Sci. U. S. A. 1995; Van Doren, *et al.*, Protein Sci. 1995; Botos, *et al.*, Proc. Natl. Acad. Sci. USA 1996; Broutin, *et al.*, Acta Crystallogr., Sect. D: Biol. Crystallogr.

-44-

1996; Gooley, *et al.*, J. Biomol. NMR 1996; Betz, *et al.*, Eur. J. Biochem. 1997; Gonnella, *et al.*, Bioorg. Med. Chem. 1997; Moy, *et al.*, Biochemistry 1998 and Moy, *et al.*, Biochemistry 1999). The MMP-13 NMR structure is composed of three α -helices corresponding to residues 28-44 (α_A), 112-123 (α_B) and 153-163 (α_C) and a mixed parallel and anti-parallel β -sheet consisting of 5 strands corresponding to residues 83-86 (β_1), 95-100 (β_2), 59-66 (β_3), 14-20 (β_4) and 49-53 (β_5). The active site of MMP-13 is bordered by β -strand IV, the Ca^{+2} binding loop, helix B and a random coil region from residues P139-Y141. The catalytic zinc is chelated by H119, H123, and H129 while the structural zinc is chelated by H69, H84 and H97. The calcium ion is chelated in a loop region consisting of residues D75 to G79. An interesting feature of the MMP active-site structure is an apparent kink in the backbone that occurs between the Ca^{+2} binding loop and β -strand IV. In the case of MMP-13, this results in the NHs of both L82 and A83 facing toward the active site of the enzyme. An important feature of substrate and inhibitor binding to the MMPs are hydrogen bonding interactions with β -strand IV which is facilitated by this unusual kink conformation (Lovejoy, *et al.*, Science 1994; Lovejoy, *et al.*, Biochemistry 1994; Spurlino, *et al.*, Proteins: Struct., Funct., Genet. 1994; and Borkakoti, *et al.*, Nat. Struct. Biol. 1994).

20 The interaction of Compound A in the active site of MMP-13 was determined by 5 intramolecular NOEs for Compound A and by a total of 47 intermolecular distance restraints between MMP-13 and Compound A. The key MMP-13 residues involved in the interaction with the inhibitor correspond to three distinct MMP-13 regions: residues L81, L82 and A83 from β -strand IV; residues L115, V116, and H119 from α -helix II; and L136, I140 and Y141 from the active site loop which comprise the S1' and S2' pockets of MMP-13. A unique feature of the MMP-13 structure is the large S1' pocket which nearly reaches the surface of the protein.

30 Compound A binds to the right-side of the catalytic Zn where the p-methoxyphenyl of Compound A sits in the S1' pocket of the MMP-13 active site. This positioning is evident from the observed NOEs from 3HH*, 3HE1/2

-45-

and 3HD1/2 to L115, V116, H119, L136, and Y141. The aryl group primarily interacts with the side-chain of L81 as evident by the strong NOEs between 1HH*, 1HE2 and 1HZ and the L81 spin-system. Finally, the pyridine ring is essentially solvent exposed but interacts with the side-chain of I140. These interactions position Compound A such that the hydroxamic acid moiety of Compound A chelates to the "right" of the catalytic zinc and the sulfonyl oxygens are in hydrogen-bonding distance to the backbone NH of L82.

It is interesting to note that the active site loop is highly dynamic in both the inhibitor-free and CGS-27023A structures based on S^2 order-parameters (Moy, *et al.*, J. Biomol. NMR 1997). This region in the MMP-13:Compound A structure appears to be significantly less mobile by the observation that most of the residues in this loop region were easily observable in the ^1H - ^{15}N HSQC spectra and readily assigned. One possible explanation for this difference is the hydrophobic interaction between the pyridine ring of Compound A and the side-chain for Ile-140. In MMP-1, I140 is replaced by a serine which essentially eliminates this beneficial interaction.

Another unique feature of the MMP-13 NMR structure is the apparent dynamic nature of residues H69 to Y73. These residues are completely disordered due to the lack of any assignment information and the resulting absence of any constraint information presumably a result of the flexible nature of these residues. Residues H69 to Y73 occur between the Ca^{+2} binding loop and the structural zinc where the corresponding region in the previously solved MMP-1 NMR structures is well defined. There is no apparent explanation for this change in mobility between the two NMR structures but it may contribute to the observed difference in the physical behavior of MMP-1 and MMP-13. Under identical conditions, inhibitor-free MMP-1 is stable for upwards of two months whereas inhibitor-free MMP-13 degrades immediately.

Comparison of the MMP-13:Compound A and MMP-1:CGS-27023A Structures:

The high-resolution NMR structure for the MMP-13:Compound A complex was effectively and efficiently determined by using a homology model based on the

-46-

MMP-1 NMR structure as an initial structure to analyze ambiguous NOESY data. This result is evident of the high structural and sequence similarity between members of the MMP family and consistent with the previously observed best-fit superposition of the backbone atoms for MMP-1, stromelysin, matrilysin and
5 neutrophil collagenase (Moy, *et al.*, Biochemistry 1998; Moy, *et al.*, Biochemistry 1999).

The strong similarity between the various MMP structures creates an initial difficulty in designing specific MMP inhibitors. This is exemplified by the high sequence similarity among the MMPs in the active site. Comparison of
10 the sequence similarity between MMP-13 and MMP-1 illustrates this difficulty. There are only a few significant residue differences between the two enzymes where these modifications results in a significant change in the local environment of the active site. The R114 to V115 modification results in a conversion from a hydrophilic to a hydrophobic environment at the base of the
15 S1' pocket between MMP-1 and MMP-13, respectively. Similarly, the N80 to L81 substitution places a bulkier hydrophobic residue in the S2' pocket for MMP-13 compared to a more hydrophilic environment for MMP-1. Similarly in the active loop region, I140 a bulky hydrophobic residue in MMP-13 replaces the smaller hydrophilic S139 residue in MMP-1. Clearly, it is feasible to
20 incorporate substituents into a small molecule to take advantage of these spatial distinct environmental changes between MMP-1 and MMP-13. Nevertheless, when these sequence and environmental differences are averaged across the MMP family it becomes less discriminating and extremely difficult to design an inhibitor to a specific MMP subtype based strictly on the small sequence
25 differences.

Conversely, the most distinct structural difference between the MMPs and readily amenable to incorporating specificity in drug design is the relative size and shape of the S1' pocket. This is clearly evident by comparison of the defined S1' pockets for MMP-13 and MMP-1. The large difference in size
30 in the S1' pockets between the MMP-13 and MMP-1 NMR structures is striking. The S1' pocket for MMP-13 nearly reaches the outer surface of the protein and

-47-

is greater than twice the size of MMP-1. The additional size of the MMP-13 S1' pocket relative to MMP-1 is best illustrated by the filling capacity of the two inhibitors. In the MMP-1:CGS-27023A NMR structure, the p-methoxyphenyl effectively fills the available S1' pocket for MMP-1. Conversely, in the MMP-13:Compound A complex the p-methoxyphenyl only partially fills the available space within the MMP-13 S1' pocket. The size of the MMP-13 pocket is actually similar in size to stromelysin where the design of stromelysin inhibitors has taken advantage of this deeper S1' pocket by using a biphenyl substituent in another series instead of the p-methoxyphenyl in Compound A to bind into the S1' pocket (Hajduk, *et al.*, J. Am. Chem. Soc. 1997; Olejniczak, *et al.*, J. Am. Chem. Soc. 1997). Thus, the NMR structures for MMP-13 and MMP-1 suggest that a ready approach to designing specificity between these MMPs is to take advantage of the significantly different sized S1' pockets. The high mobility of the MMP-1 active site presents a potential caveat to this analysis of the static images of the MMP-1 and MMP-13 structures. It is probable that the MMP-1 active site is capable of accommodating a S1' substituent larger than implied from its current structure due to its increased mobility in both free and inhibited structures.

Examination of the binding mode of Compound A in the MMP-13:Compound A complex suggests a conformation generally similar to CGS-27023A in the MMP-1:CGS-27023A NMR structure previously reported (30 simulated annealing structures deposited with Protein Data Bank, Accession No. 4AYK; restrained minimized mean structure deposited with Protein Data Bank, Accession No. 3AYK). Compound A and CGS-27023A are structurally very similar with the only difference being the nature of the substituent binding in the S2' pocket where an aryl group in Compound A replaces the isopropyl group in CGS-27023A. The strong resemblance between the binding mode of Compound A and CGS-27023A is apparent from the nearly identical intermolecular NOE patterns observed between the inhibitors and the proteins. The key MMP-13 residues involved in the interaction with Compound A correspond to L81, L82 and A83 from β -strand IV; residues L115, V116, and

-48-

H119 from α -helix II; and L136, I140 and Y141 from the active site loop. Similarly, the MMP-1 residues involved in the interaction with CGS-27023A correspond to residues N80, L81, A82 and H83 from β -strand IV; residues R114, V115, H118 and E119 from α -helix II; and L135, P138, Y137, S139 and
5 Y140 from the dynamic flexible loop.

As stated previously, there are three distinct residue changes between MMP-13 and MMP-1 in the active site. The R114 to L115 change between MMP-1 and MMP-13, respectively, has a significant impact on the environment at the base of the S1' pocket but since Compound A only partially
10 fills the MMP-13 S1' pocket this change should not effect the binding conformation of Compound A relative to CGS-27023A. Conversely, the N80 to L81 substitution directly interacts with the inhibitors in the S2' pocket and may result in an effective change in the binding mode of the inhibitors. To complicate the analysis, the only change in the inhibitors are the substituents
15 that bind the S2' pocket. For the MMP-1:CGS-27023A complex, the isopropyl group interacts with both the sidechains of N80 and H83 where the aryl group from Compound A only interacts with L81 in MMP-13. Additionally, CGS-27023A is in hydrogen-bonding distance to both L81 and A82, whereas Compound A appears to form a bifurcated hydrogen bond with L82. This
20 analysis suggests that CGS-27023A binds closer to β -strand IV since the S2' pocket is more accessible in MMP-1 due to the absence of the bulky L81 side-chain and the presence of the aryl group in Compound A. A direct comparison of the bound conformations suggest only a subtle difference in the relative orientation of the inhibitors. The S139 to I140 difference between MMP-1 and
25 MMP-13, respectively, appears to be related to a mobility change as opposed to a structural change. In the MMP-1:CGS-27023A structure the pyridine ring position is essentially undefined and solvent exposed this compares to the MMP-13:Compound A structure where the pyridine ring binds with the side-chain of I140. Clearly, Ile is a bulkier more hydrophobic group relative to Ser which
30 would provide a beneficial hydrophobic interactions with the pyridine ring. The more interesting observation is the apparent decrease in mobility for the active

-49-

loop in the MMP-13 structure which may be related the pyridine ring I140 interaction. This appears to be consistent with previously inhibited MMP X-ray structures (Spurlino, *et al.*, Proteins: Struct., Funct., Genet. 1994) where the inhibitor may extend the formation of a β -sheet between b-strand IV and the active loop region which results in low B-factors in the X-ray structure. This may suggest that the mobility of the active loop region is easily removed with any positive interaction with the inhibitor.

There are apparently some interesting differences between the mode of binding for the two inhibitors in the MMP-13:Compound A and MMP-1:CGS-27023A NMR structures. The more striking observation is the overall similarity between the two structures. Despite some significant sequence differences and a large difference in the size and shape of the S1' pocket either inhibitor structure would accurately predict the other structure. This observation seems to indicate that the major contributing factors to inhibitors binding the MMPs is the fit in the S1' pocket and the binding of the hydroxamic acid to the catalytic zinc. The interaction in the S2' pocket appears to have a more subtle impact on inhibitor binding and selectivity since both Compound A and CGS-27023A are low nanomolar inhibitors of MMP-13 and MMP-1, respectively. Therefore, the high-resolution solution structure of the MMP-13:Compound A in conjunction with the previously reported MMP-1 NMR structures suggest that taking advantage of the significant differences in the size and shape of the S1' pocket is a reasonable approach for developing specific MMP inhibitors.

The studies described herein present the high-resolution solution structure of MMP-13 complexed with a sulfonamide derivative of a hydroxamic acid compound (Compound A). The overall fold of MMP-13 is similar to previously reported MMPs structures. The major difference is the large S1' pocket which nearly reaches the surface of the protein. The structure was based on a total of 3279 constraints including 47 distance restraints between MMP-13 and Compound A from X-filtered NOESY experiments. The inhibitor was found to bind to the "right" side of the catalytic Zn such that the p-methoxyphenyl ring

-50-

- sits in the S1' pocket, the aryl moiety interacts with L81 of β IV, the pyridine ring interacts with I140 of the active site loop, hydrogen bond interactions exist between the sulfonamide oxygens with residue L82 and the hydroxamic acid chelates the catalytic Zn. This inhibitor binds MMP-13 similarly to the MMP-1:
- 5 CGS-27023A complex suggesting that appropriately filling the S1' pocket may play a key role in developing selective MMP inhibitors.

Table 1. Observed NOEs Between Compound A and MMP-13

Compound A	MMP-13	NOE Class	Compound A	MMP-1	NOE Class
1HH*	L81 H γ	W	3HH*	Y141 H α	M
1HH*	L81 H δ 1#	W	3HH*	Y141 H β 1	W
1HH*	L81 H δ 2#	M	3HH*	Y141 H β 2	W
1HH*	L81 H α	S	3HH*	Y141 H δ 2	W
1HE2	L81 H δ 1#	W	3HE2	L82 H δ 1#	W
1HE2	L81 H δ 2#	M	3HE1	A83 H β #	W
1HZ	L81 H δ 1#	W	3HE1	H116 H α	W
1HZ	L81 H δ 2#	M	3HE1	H116 H γ 1#	M
2HZ	I140 H γ 2#	W	3HE2	H116 H γ 2#	W
2HE1	I140 H δ 1#	W	3HE2	I140 H γ 2#	W
3HH*	L82 H δ 1#	W	3HE2	Y141 H α	W
3HH*	L115 H β #	W	3HE2	Y141 H β 1	W
3HH*	L115 H γ	W	3HE2	Y141 H β 2	W
3HH*	L115 H δ 1#	W	3HD2	L82 H δ 1#	W
3HH*	L115 H δ 2#	W	3HD1	A83 H β #	W
3HH*	V116 H α	W	3HD1	V116 H γ 1#	W
3HH*	V116 H γ 1#	W	3HD2	V116 H γ 2#	W
3HH*	V116 H γ 2#	M	3HD2	I140 H α	W
3HH*	H119 H α	W	3HD2	I140 H γ 2#	W
3HH*	H119 H δ 2	W	3HD2	Y141 H α	W
3HH*	H119 H β 1	W	3HD2	Y141 H β 1	W
3HH*	H119 H β 2	W	3HD2	Y141 H β 2	W
3HH*	L136 H δ 1#	W	3HD2	Y141 HN	W
3HH*	L136 H δ 2#	W			

Example 3

Structure Based Design of a Novel, Potent, and Selective Inhibitor for MMP-13

The matrix metalloproteinases (MMPs) comprise a family of zinc
5 containing enzymes that cleave a broad range of substrates including collagens,
fibronectin and gelatins where the substrate preference varies for individual
MMPs. The design of MMP inhibitors has been initially based upon imitation of
the binding interaction of natural protein substrates to MMPs where structural
information of MMPs complexed with peptide substrates has been determined
10 by x-ray crystallography and NMR spectroscopy. This structural information has
provided a general description of the MMPs active site.

The active site for the MMPs is composed of a catalytic zinc
chelated by three histidines where three substrate binding pockets are located to
both the right (S1', S2', S3') and left (S1, S2, S3) of the catalytic zinc. The
15 substrate binding pockets were identified by the interactions of side chains from
the peptide substrate with the MMPs. The primary effort in MMP inhibitor
design has focused on compounds that chelate the catalytic zinc while primarily
binding in the S1' and S2' pockets. This has evolved from the observation that
the structural characteristics of the S1' pocket (size, shape, amino acid
20 composition) incurs the greatest variability between the individual MMPs and
this provides an obvious approach in designing selective and specific MMP
inhibitors. Nevertheless, there has also been success in utilizing the binding
pockets to the left of the catalytic zinc in addition to or in combination with the
right handed binding pockets in the design of inhibitors.

25 The underlying challenge in designing MMP inhibitors is the
reasonably high sequence and structural homology observed between the
individual members of the MMP family making it intrinsically difficult to design
an inhibitor that will function against a single MMP in the absence of structural
information. The problem with a non-specific MMP inhibitor as a drug is the
30 high likelihood of serious side-effects because of the large number of enzymes in
the MMP family and their corresponding diversity in targets and function.

-53-

Accordingly, the detailed structural information provided herein is a critical component of an inhibitor design program targeting a particular MMP enzyme.

Materials and Methods:

- 5 *Synthesis of Compound D and Compound E:* The sulfonamide derived from 2-amino-3,5-dimethyl-benzoic acid methyl ester and p-methoxybenzenesulfonyl chloride was N-alkylated with benzyl bromide and the ester group of the resulting intermediate was hydrolyzed (LiOH/THF) to afford the carboxylic acid. The corresponding hydroxamic acid was formed by preparation of the
10 acid chloride (oxalyl chloride/DMF) followed by reaction with hydroxylamine. Compound E was synthesized by reaction of 2-amino-3,5-dimethyl-benzoic acid methyl ester and p-fluorobenzenesulfonyl chloride followed by N-alkylation with benzyl bromide. Hydrolysis of the methyl ester (LiOH/THF) followed by displacement of fluorine with the alkoxide of benzofuran-2-carboxylic acid (2-
15 hydroxy-ethyl)-amide gave, after conversion to the hydroxamic acid and formation of the HCl salt as described above, Compound E.

- NMR Sample Preparation:* Uniformly (>95%) ^{15}N - and $^{15}\text{N}/^{13}\text{C}$ -labeled human recombinant MMP-13 was expressed in *E. coli* and purified as described
20 previously. 1mM $^{13}\text{C}/^{15}\text{N}$ - and ^{15}N - MMP-13 NMR samples were prepared by concentration and buffer exchange using Millipore Ultrafree -10 centrifugal filters into a buffer containing 10mM deuterated Tris-base, 100mM NaCl, 5mM CaCl_2 , 0.1 mM ZnCl_2 , 2 mM NaN_3 , 10mM deuterated DTT in 90% $\text{H}_2\text{O}/10\%$ D_2O or 100% D_2O . The 10:1 Compound B:MMP-13 samples were prepared by
25 addition of Compound B into either a 1mM $^{13}\text{C}/^{15}\text{N}$ - or ^{15}N -MMP-13 sample followed by pH readjustment. The sample to explore the potential of competitive inhibition between Compound B and Compound A was prepared by first adding 1mM of Compound A to a 1 mM ^{15}N - MMP-13 sample followed by the addition of 10mM Compound B. The initial MMP-13:Compound A sample
30 was made by buffer exchange of ^{15}N - MMP-13 into the buffer containing 0.1 mM Compound A followed by additional buffer exchanges to remove excess

Compound A. Finally, 10mM of Compound B was added to the 1mM ^{15}N -MMP-13:Compound A sample followed by pH readjustment.

NMR Data Collection: All spectra were recorded at 35°C on a Bruker AMX-2 600
5 spectrometer using a gradient enhanced triple-resonance $^1\text{H}/^{13}\text{C}/^{15}\text{N}$ probe. For spectra recorded in H_2O , water suppression was achieved with the WATERGATE sequence and water-flip back pulses (Piotto, *et al.*, J. Biomol. NMR 1992; Grzesiek and Bax, J. Am. Chem. Soc. 1993). Quadrature detection in the indirectly detected dimensions were recorded with States-TPPI hypercomplex
10 phase increment (Marion, *et al.*, J. Magn. Reson. 1989). Spectra were collected with appropriate refocusing delays to allow for 0,0 or -90,180 phase correction.

The resonance assignments and bound conformation of Compound A in the MMP-1: Compound A complex were based on the 2D $^{12}\text{C}/^{12}\text{C}$ -filtered NOESY (Petros, *et al.*, FEBS Lett. 1992; Gemmecker, *et al.*, J. Magn. Reson. 1992), 2D $^{12}\text{C}/^{12}\text{C}$ -filtered TOCSY (Petros, *et al.*, FEBS Lett. 1992; Gemmecker, *et al.*, J. Magn. Reson. 1992) and $^{12}\text{C}/^{12}\text{C}$ -filtered COSY
15 experiments (Ikura and Bax, J. Am. Chem. Soc. 1992).

The assignments of the ^1H , ^{15}N , and ^{13}C resonances of MMP-13 in the MMP-13:Compound B complex were based on the previous assignments for
20 the MMP-13:Compound A complex in combination with a minimal set of experiments: 2D ^1H - ^{15}N HSQC, 3D ^{15}N -edited NOESY (Marion, *et al.* Biochemistry 1989; Zuiderweg and Fesik, Biochemistry 1989), CBCA(CO)NH (Grzesiek and Bax, J. Am. Chem. Soc. 1992), C(CO)NH (Grzesiek, *et al.*, J. Magn. Reson., Ser. B 1993), HNHA (Vuister and Bax, J. Am. Chem. Soc. 1993)
25 and HNCA (Kay, *et al.*, J. Magn. Reson. 1990). The acquisition parameters for each of the experiments used in determining the solution structure of the MMP-13:Compound B complex were as reported previously (Moy, *et al.*, Biochemistry 1996).

The MMP-13:Compound B structure is based on observed NOEs
30 from the 3D ^{15}N -edited NOESY (Marion, *et al.* Biochemistry 1989; Zuiderweg and Fesik, Biochemistry 1989) and 3D ^{13}C -edited/ ^{12}C -filtered NOESY (Vuister

-55-

and Bax, J. Am. Chem. Soc. 1993; Lee, *et al.*, FEBS Lett. 1994). The 3D ^{15}N -edited NOESY and 3D ^{13}C -edited/ ^{12}C -filtered NOESY experiments were collected with 100 msec and 110 msec mixing times, respectively.

5 *Molecular Analysis and Design:* The minimized models of Compound B and Compound D complexed to MMP-13 were prepared as previously described (Chen, *et al.*, J. Biomol. Struct. Dyn. 1995; Chen, *et al.*, Biochemistry (in press) 1998). Using molecular dynamics methods (Sybyl v6.4 from Tripos Inc), protein regions within 5 Å from Compound B were sampled along with the inhibitor,
10 whereas everything else remained rigid during the simulations. Upon energy convergence, the last 50 frames from the final 100 picoseconds run was averaged and this averaged structure underwent a final minimization. The final protein-Compound B model appeared to have optimized possible polar and van der waals interactions. The identical procedure was applied to the complex of
15 MMP-13 and Compound D. Since the two complexes used identical MMP-13 structures, the proteins were overlapped to depict the positions of the two inhibitors within the active site. Graphics analysis of the inhibitors showed that the methylene carbon of Compound B containing the 2HB1/2 protons (Figure 6) overlapped identically with the methoxy carbon from Compound D. This
20 analysis indicated the optimal or minimal linkage length of connecting the benzofuran moiety to the methoxy region of Compound D. The final design scheme is shown in Figure 8A for the hybrid inhibitor. The homology model of MMP-9 was constructed using the COMPOSER program (Tripos INC, Sybyl v.6.4)

25

High-throughput Screening Analysis: Compound B was identified as an initial lead from the analysis of the MMP-13 high-throughput screen (HTS). A total of 58079 compounds were screened for their ability to inhibit MMP-13 enzymatic activity where 385 compounds were shown to have $\geq 40\%$ inhibition at 10
30 $\mu\text{g/ml}$ dosage. Compound B was shown to exhibit weak inhibition of MMP-13 (89% at the 10 $\mu\text{g/ml}$), but more intriguing was the observation of a complete

lack of activity against other MMPs (MMP-1, MMP-9 and TACE). The primary structure of Compound B along with the proton naming convention is shown in Figure 6.

The resulting HTS hits were further examined by cluster analysis.

5 The hits were clustered based on structural similarities where the properties of these compounds were compared against the properties of the set of orally available drugs. The properties used to profile the HTS hits consists of: total number of non-hydrogen atoms, number of heteroatoms, number of hydrogen-bond donors and acceptors, calculated logP and molecular weight. This profile
10 analysis provides an initial means to predict the likelihood that an HTS hit may have drug-like characteristics such as bioavailability and in-vivo stability. The profile of Compound B indicates that the compound has properties similar to orally available drugs suggesting that it would be an ideal candidate for optimization of its enzyme potency and selectivity.

15 A common feature of known MMP inhibitor structures is the presence of a Zn-chelator that plays a fundamental role in its activity. In most cases Zn chelation occurs from the presence of a hydroxamic acid in the structure of the small molecule. As apparent from the structure of Compound B, the compound does not contain an obvious substituent that would chelate Zn.
20 Thus, the unique structure of Compound B suggested a potential novel mechanism for inhibition of MMP-13 further strengthening the choice of Compound B as an initial lead candidate. Therefore, the identification of Compound B as a candidate to optimize its activity and selectivity was based on three unique observations: its intrinsic MMP-13 selectivity, its structural profile
25 similar to known bioavailable drugs and finally its apparent novel structure.

NMR Structure of the MMP-13 - Compound B Complex: The NMR binding studies provided critical information pertaining to the mechanism of Compound B inhibition of MMP-13 and the method for designing increase potency. The
30 major question presented when Compound B was identified from HTS was its unknown MMP-13 binding site and its method for inducing MMP-13 inhibition.

Previous work on the NMR structure of MMP-13 complexed with Compound A and MMP-1 complexed with CGS-27023A provided the framework and methodology to analysis Compound B bound to MMP-13 (Moy, *et al.*, Biochemistry Submitted 1999; Moy, *et al.*, Biochemistry 1999).

- 5 The Compound B MMP-13 binding site was initially identified from chemical shift perturbation in the ^1H - ^{15}N HSQC spectra. The observed perturbations were mapped onto a GRASP surface (not shown). It is apparent that the major effect of Compound B on the chemical shifts of MMP-13 occurs in the proximity of the S1' pocket suggesting that Compound B sits in this pocket.
- 10 From the NMR and X-ray structures of MMP-13, it was determined that the S1' pocket for MMP-13 is very deep and linear in shape while nearly reaching the surface of the protein. In fact, a number of residues at the surface of MMP-13 near the base of the S1' pocket show significant chemical shift perturbation in the presence of Compound B. Since Compound B is a linear molecule, docking
- 15 studies would place the inhibitor stretched throughout the linear S1' pocket of MMP-13. The only question remaining was whether to place the morpholine or the benzofuran moiety of Compound B at one end of the pocket, adjacent to the catalytic zinc or the opposite end, distant from the zinc atom. Property analysis of the enzymes S1' pocket depicts that the end adjacent to the zinc is relatively
- 20 polar whereas the opposite end is hydrophobic. This analysis lead us to dock Compound B with the morpholine ring adjacent to the catalytic zinc atom with the benzofuran moiety siting in a hydrophobic pocket formed by L115, L136, F149 and P152 at the base of the S1' pocket. To further verify the proposed binding of Compound B in the S1' pocket of MMP-13, a simple competition
- 25 experiment with Compound A was conducted. The ^1H - ^{15}N HSQC experiment for the MMP-13:Compound B complex was collected in the presence of Compound A. The presence of Compound A displaced all of Compound B as evident by the distinct differences in the ^1H - ^{15}N HSQC spectra which further suggests that both compounds bind in the S1' pocket.
- 30 The relative orientation and binding of Compound B with MMP-13 was further confirmed by the observation of intermolecular NOEs between

-58-

Compound B and MMP-13 from the 3D ^{13}C -edited/ ^{12}C -filtered NOESY experiment. The NOESY spectra was collected in the presence of a ten-fold excess of Compound B because of the weak affinity of Compound B with MMP-13. Nevertheless, a total of 16 NOEs were observed between Compound B and

5 L81, L115, V116, Y141, T142 and Y143 which support the initial positioning of Compound B in the MMP-13 S1' pocket. An expanded 2D plane from the 3D ^{13}C -edited/ ^{12}C -filtered NOESY experiment (not shown) demonstrated examples of some key intermolecular NOEs between Compound B benzofuran group resonances and L115 δ and Compound B resonances proximal to the

10 morpholine ring and L82 δ . The complex of Compound B with MMP-13 was subjected to energy refinement using the NMR results as constraints (Moy, *et al.*, Biochemistry 1999; Chen, *et al.*, J. Biomol. Struct. Dyn. 1995). The modeling results depict the morpholine oxygen forming a hydrogen bond with the backbone amide group of Leu-82 and the benzofuran group packs deep in

15 the S1' pocket with the peptide bond linker portion forming hydrogen bonds with protein backbone groups. The complex shows no apparent interactions between the inhibitor and the catalytic zinc justifying the ligands micromolar potency.

20 *Structures of MMP-1, MMP-9 and MMP-13:* The recent NMR solution structures of MMP-1 and MMP-13 were used as starting points for molecular modeling and analysis (Moy, *et al.*, Biochemistry Submitted 1999; Moy, *et al.*, Biochemistry 1998; Moy, *et al.*, Biochemistry 1999). A homology model for MMP-9 was developed based on its strong homology to MMP-1 (54% identity around the

25 catalytic domain). Based on the homology model, the catalytic site of MMP-9 is similar to the corresponding sites in MMP-1 and MMP-13. All three structures were used as starting points for analysis and synthetic design.

Comparative analysis of the MMP structures shows that residue positions 115 and 144, in addition to the length of the specificity loop,

30 determines the size and shape of the S1' pockets. Alignment of the NMR structures for MMP-1 and MMP-13 shows that MMP-13 contains two additional

-59-

insertions in the specificity loop. The homology model of MMP-9 indicates no additional insertions so its length is identical to MMP-1.

Residue positions 115 and 144 are important in establishing the relative length of the S1' pockets for the MMPs where the larger the side chain at these positions results in a smaller S1' pocket. Since residue 115 is spatially closer to the catalytic zinc than residue 144, a larger side chain for residue 115 will have a greater impact on defining a smaller S1' pocket compared to residue 144. MMP-1 has the largest side chain at position 115, thus its S1' pocket is the smallest. MMP-9 has an Arg at position 144 resulting in its S1' pocket being longer compared to MMP-1. Conversely, MMP-13 has short side chains at both positions 115 and 144. The short side chains combined with an increased length of its specificity loop result in MMP-13 having the largest S1' pocket. To summarize, the size of the MMP S1' pockets are as follows: MMP-13 > MMP-9 > MMP-1 where this structural feature plays a critical role in the design strategy for developing a potent and specific MMP-13 inhibitor.

Design Strategy: A strategy utilizing NMR and molecular modeling was applied towards the design and synthesis of an MMP-13 selective inhibitor lead. The basic approach behind the design strategy is to optimize the affinity of the chemical lead Compound B while maintaining its inherent MMP-13 selectivity. This can be achieved by taking advantage of the distinct structural feature of MMP-13, its deep linear S1' pocket, while combining overlapping structural features of Compound B with other potent inhibitors. Compound C is an example of a potent and selective inhibitor for MMP-9 and MMP-13 (See Table 2). Based on the NMR solution structure of MMP-13 complexed with Compound A (Figure 4), structurally similar inhibitors were positioned into the active site of MMP-13.

Figure 7 shows the critical regions of Compound C, which can be broken down into two components, Compound D which represents the zinc chelating portion of the compound that contributes to the binding potency and the toluene group (1A) which contributes to enhanced ligand selectivity against

MMP-1. The strategy was to design a new inhibitor based on replacing the toluene group (1A) with a component of Compound B critical for binding within the extended S1' pocket of MMP-13. The overlay of the NMR solution structure for Compound B with the model for Compound D is shown in Figure 8B. The close similarity between the positioning of the two structures made it readily apparent that it would be possible to generate a hybrid of the two structures combining the potent Compound D with the selective component of Compound B (Figure 8A). These results were then used to design the proposed hybrid inhibitor Compound E. The assay data in Table 2 clearly shows that the new inhibitor, Compound E, has better potency compared to Compound C in addition to improved selectivity towards MMP-13. Thus, the combination of NMR spectroscopy with molecular modeling techniques resulted in the design of a novel, potent and selective MMP-13 inhibitor (Compound E) which has an IC₅₀ of 17 nM for MMP-13 and showed >5800, 56 and >500 fold selectivity against MMP-1, MMP-9 and TACE, respectively. To the best of our knowledge, this represents the first example of a potent MMP-13 inhibitor that has been shown to be selective against MMP-9.

Table 2 - IC50 and Selectivity Data

Compound	MMP-1	MMP-9	MMP-13	TACE	S-1 ^a	S-9 ^a	S-TACE ^a
C	750nM	46nM	75nM	470nM	10.0x	0.6x	6.3x
D	82nM	21nM	15nM	240nM	5.5x	1.4x	16x
E	NA	945nM	17nM	19%	>5800x	56x	>500x
F	1025nM	71nM	301nM	664nM	3.4x	0.2x	2.2x

^a Selectivity data presented as a ratio of the MMP or TACE IC₅₀ with MMP-13

-61-

Example 4

The X-ray crystal structure of the MMP-13:Compound A complex was determined using the following procedure:

- 5 *Gene/expression system/production:* The cDNA coding for human MMP-13 proenzyme had 85 residues of the PRO domain, followed by 165 residues of the catalytic domain (CAT). The gene was carried on a pET-21a expression plasmid, under the control of a bacteriophage T7 promoter. The expression host was *Escherichia coli* BL21 (DE3), which had a chromosomal copy of T7 RNA
10 polymerase under *lac* control. Cells were grown in nutrient broth, and synthesis of PRO-CAT was induced by isopropyl- β -thiogalactoside. The protein accumulated to 5-10% of total cellular protein, essentially all of which was aggregated into inclusion bodies.

- For potential MAD experiments, the plasmid was transferred into a
15 methionine auxotroph host. PRO-CAT with selenomethionine substitution was produced by induction in a defined medium, with methionine replaced by selenomethionine.

- Purification and refolding of PRO-CAT:* Frozen cells were disrupted
20 mechanically, and inclusion bodies were isolated by centrifugation. PRO-CAT was solubilized with urea containing dithiothreitol to disrupt any disulfide bridges. PRO-CAT was partially purified by anion-exchange chromatography, in urea, on Q Sepharose. The protein was diluted to about 400 μ g/ml in a solution of sodium chloride, calcium chloride, and zinc acetate, buffered with
25 tricine-HCl. Refolding proceeded over 3-4 days, during dialysis, with multiple buffer changes. PRO-CAT was then concentrated for activation and release of CAT.

- Activation of PRO-CAT:* The presently-accepted view of MMPs holds that the
30 proenzyme form is maintained in an inactive state through the coordination of one cysteine from the PRO domain into the active-site zinc. If this cysteine is

-62-

displaced, the enzyme becomes active. In our protocol, aminophenyl mercuric acetate was added to the protein solution to form a mercurial adduct with the cysteine. Progress of activation was monitored by SDS polyacrylamide gel electrophoresis. Results indicated that the CAT domain accumulated and the
5 PRO domain was degraded to small peptides.

Purification of MMP-13 (CAT) – Size Exclusion: Following activation and PRO cleavage, MMP-13 was isolated by size-exclusion chromatography through SuperDex 75 in a solution of sodium chloride, calcium chloride, and zinc
10 acetate, buffered with tris-HCl.

Purification of MMP-13 – Affinity: MMP-13 was further purified by affinity chromatography on an immobilized hydroxamate inhibitor. The affinity matrix was prepared by coupling an hydroxamate inhibitor to Sepharose through the
15 amino group of the piperazine ring. MMP-13 can be absorbed to the matrix and desorbed by displacement using another inhibitor of choice.

Characterization of MMP-13: Protein preparations for crystallization trials were validated by several techniques. Routinely, SDS-PAGE showed a predominant
20 species whose migration was consistent with a molecular weight of around 19,000. MALDITOF mass spectroscopy demonstrates a single species consistent with the expected size of 18,588 amu. (MMP-13 prepared with selenomethionine showed essentially complete replacement). N-terminal sequencing demonstrated that the protein begins with YNVF, as expected for
25 correct cleavage between PRO and CAT. Retention volume in analytical size-exclusion chromatography was consistent with a monomeric protein: no detectable aggregation was observed. The final protein was enzymatically active on a fluorogenic peptide substrate, and degraded denatured collagen.

30 *Crystallization of MMP-13 complex with Compound A:* The MMP-13 protein solution was buffered with 10 mM tris-HCL buffer, pH 7.5, and 0.25 M NaCl.

-63-

The concentration of protein used for crystallization was 20.0 mg/ml. The inhibitor solution was added to a protein solution with a mole ratio (protein:inhibitor) of 1:2, and was mixed for more than 1 hour.

Crystallization conditions were screened by the hanging-drop vapor diffusion method (McPherson, A., Methods Biochem. Anal. 1976). A successful procedure for growing crystals of this complex at room temperature was identified, and crystals were obtained. Specifically, a solution was prepared from 3 μ l of protein solution and 3 μ l of precipitant solution, which consisted of 26% PEG4000, 0.1 M ammonium sulfate, and 0.1 M sodium chloride. A drop of this solution was suspended on a microscope coverslip glass which had been coated with silicone to prevent drop spreading. The reservoir solutions consisted of 0.6 ml precipitant solution. Equilibration was performed at room temperature by vapor diffusion. Crystals began appearing after three days. After two weeks, these crystals stopped growing. The X-ray data which have been processed show that the MMP-13 complex was crystallized in two forms. One crystal form is C-centered orthorhombic; it belonged to space group C2221, and had a cell dimension of $a=36.3 \text{ \AA}$, $b=134.4 \text{ \AA}$, and $c=134.8 \text{ \AA}$. This crystal had high mosaicity; therefore, it would be of little use when working on the structure of the complex. The second crystal form is primitive orthorhombic, from space group P21212, with a cell constant of $a=108.3 \text{ \AA}$, $b=79.8 \text{ \AA}$, and $c=36.1 \text{ \AA}$. This crystal had low mosaicity, but it was very small in most cases.

In order to obtain a big single crystal for X-ray data collections, the seeding technique (Thaller, C., et al., J. Mol. Biol. 1981) was applied. This was accomplished by using both the microseeding and the macroseeding methods. Small seed crystals were transferred to a 20% PEG4000 precipitant solution on a depression slide. A single washed crystal was injected into a hanging-drop solution, which was composed of 3 μ l of MMP-13 complex solution and 3 μ l of precipitant solution. The reservoir solutions consisted of 0.6 ml precipitant solution at pH 8.0. This procedure successfully produced bigger crystals with a maximum edge dimension of up to $0.35 \times 0.1 \times 0.1 \text{ mm}^3$. These crystals diffracted X-ray at a resolution of 2.0 \AA .

-64-

X-Ray Data Collection: X-ray diffraction data from 30.0-2.0 Å resolution for the MMP-13:Compound A complex crystal (P21212 form) was collected by using an RAXIS IIC Image Plate area detector which used graphite monochromatic CuKα radiation from a Rigaku RU200 rotating anode generator (operating at 50 kV, 100 mA) at a low temperature of 100 K. The oscillation angle for each plate was 1 degree, and exposure time was 20 minutes per 'image'. The processing of X-ray diffraction data was accomplished using the HKL programs (Otwinowski, Z. and Minor, W., Methods in Enzymology 276:307-26). The R-merges for full and partial reflections were 4.0% and 6.04% respectively. 18,782 unique reflections (81% complete at 2.0 Å resolutions) were collected.

Structure Determination and Refinement: The MMP-13 complex crystal structure has been determined by a combination of crystallographic modeling and the Molecular Replacement method using models of MMP-13 derived from the MMP-1 and MMP-8 structures. The homology between MMP-13 and MMP-8 is 56% by sequence, and at least 70% by structure. Crystals of the MMP-13 complex have two molecules in the asymmetric unit, i.e., the unit is a dimer. Conventional molecular replacement was not effective for determination of this dimer structure by using a monomer model. There are two reasons for this: (1) the high symmetry of the crystal structure; and (2) the conformations and the configurations of the side chain and the main chain in flexible loops of MMP-13 and MMP-8.

Firstly, the crystal structure of the MMP-13 complex is highly symmetrical. The P21212 crystal has four symmetry operations, and there are eight molecules in a unit cell. A second crystal form, belonging to space group C222, and having eight symmetry operations in a unit cell, has been identified. In this crystal, there are 16 monomers per cell in the dimer structure, and 32 monomers per cell in the tetramer structure. Therefore, the rotation search and especial translation search become more difficult. Secondly, even though the MMP family's catalytic domain structure is highly conserved, the conformations and the configurations of the side chain and the main chain in flexible loops of

MMP-13 and MMP-8 may not be the same. In particular, the similarity between the two structures may not be sufficient to permit the determination of the dimer structure using a monomer as the searching model.

Many attempts at a rotation and translation search were made by
5 using the X-ray data and models of either a monomer of MMP-8 or a dimer of MMP-1. Some rotation solutions were obtained, but no final translation solution has been found by using the monomer model. Accordingly, to determine this structure, it was proposed that a dimer model be constructed first; the molecular replacement method was then applied to solve the structure.
10 The key idea of this proposal was crystal packing. To construct a dimer, the orientations of each monomer were determined on the basis of a rotation search. The positions of each monomer were located on the basis of the molecular packing in unit cell. Many dimer models have been constructed and applied as the 'model' for searching the rotation and translation using
15 program AMORE (Collaborative Computational Project, Number 4 (CCP4) (1994), Acta Cryst. D50:760-763). One dimer model was found to be correct, and finally resulted in the MMP-13 3-D crystal structure using the molecular replacement method. The MMP-13 complex structure was confirmed by observing the most important and significant fact that the positions of the two
20 zinc ions and the two calcium ions could be identified from the difference ($F_o - F_c$) maps with five-sigma cut, where F_o was observed structure factor and F_c was the calculated structure factor of the dimer model without zinc and calcium atoms.

These ions were located in the exact positions where they were
25 observed in other MMP family members. The molecule fits the ($2F_o - F_c$) electron densities very well, both in main chain and in side chain. The molecule fits the $2F_o - F_c$ electron density quite well. All of these MMP molecules are conserved in the core structure region, especially the position of the central helix and the catalytic zinc. The MMP-13 dimer structure was further confirmed
30 by applying the molecular replacement programs XPLOR (Brünger, A.T., XPLOR Version 3.1 Manual, Yale University, New Haven CT) and MERLOT (Fitzgerald,

-66-

P., MERLOT, version 2.4 (Nov. 10, 1991). All of them worked very well, and produced results which were in agreement with the MMP-13 structure.

Structure Refinement: The structure refinement was carried out by the program
5 XPLOR. The initial dimer model included 320 amino acid residues without zinc
and calcium ions. The dimer model was refined against 2.0 Å X-ray data,
collected on an RAXIS IIC area detector at a temperature of 100 K. The progress
of the refinement was evaluated from the quality of the protein molecular
conformations and the electron density maps, and the values of the
10 crystallographic R-factor. The initial R-factor was 52%. After rigid-body
minimization, conjugated-gradient minimization, a heating stage, a slow-cooling
stage in the range from 4000K to 300K, energy minimization, B-factor
refinement, and positional refinement, the R-factor lowered to 0.32. Electron-
density maps with coefficients of (2Fo-Fc) and (Fo-Fc), as well as the phases,
15 were calculated. The difference map shows four zinc ions and four calcium ions
in the dimer structure with five-sigma cut. Some side chain loops and a few
main loops were rebuilt on the interactive graphics system. The rebuilt dimer
plus the zinc and calcium ions, as the new model, was refined. The R-factor
was down to 26.6%. At this stage, a model of inhibitor Compound A was
20 positioned in the active-site region based on the difference electron density.

The complex structure was refined by repeating the above steps,
with the R-factor down to 20%. The water molecules were modeled as oxygen
atoms. Their initial positions were located by searching the peaks in the (Fo-Fc)
difference maps. These positions were then checked by calculating the distance
25 between 'water' and the oxygen and nitrogen of the protein. Together with the
protein (complex) atoms, these 'water' molecules were refined against the X-ray
data. Once the temperature factor of water was higher than 50, this water was
omitted. 120 water molecules near the protein were found, and five water
molecules were identified in the active site of each monomer. The (2Fo-Fc)
30 maps were used to adjust the solvent model and to aid in the placement of new
solvent molecules, as well as to check and correct the whole model. The r.m.s.

-67-

deviations of C α atoms for bond angles and bond distances from ideal geometry were 1.6° and 0.012 Å. The final crystallographic R-factor was 22%, at a resolution of 2.0 Å.

All publications mentioned herein above, whether to issued

5 patents, pending applications, published articles, protein structure deposits, or otherwise, are hereby incorporated by reference in their entirety. While the foregoing invention has been described in some detail for purposes of clarity and understanding, it will be appreciated by one skilled in the art from a reading of the disclosure that various changes in form and detail can be made

10 without departing from the true scope of the invention in the appended claims.

What is claimed is:

1. A solution comprising a biologically active catalytic fragment of human collagenase-3 (MMP-13) complexed with N-Hydroxy-2-[(4-methoxy-benzenesulfonyl)-pyridin-3-ylmethyl-amino]-3-methyl-benzamide ("Compound A").
2. The solution of Claim 1, wherein the catalytic fragment of MMP-13 comprises the amino acid residues of Figure 1.
3. The solution of Claim 2, comprising 1 mM MMP-13 complexed with Compound A in a 1:1 molar ratio, in a buffer comprising 10mM deuterated Tris-Base, 100mM NaCl, 5mM CaCl₂, 0.1mM ZnCl₂, 2mM NaN₃, and 10 mM deuterated DTT in either 90% H₂O/10% D₂O or 100% D₂O.
4. The solution of Claim 3, wherein the MMP-13 is either ¹⁵N enriched or ¹⁵N, ¹³C enriched.
5. The solution of Claim 1, wherein the secondary structure of the catalytic fragment of MMP-13 comprises three alpha helices and a mixed parallel and anti-parallel beta sheet comprising five beta strands.
6. The solution of Claim 5, wherein the alpha helices and beta strands are configured in the order β_I , α_A , β_{II} , β'_{III} , β_{IV} , β_V , α_B , and α_C .
7. The solution of Claim 6, wherein the three alpha helices correspond to residues 28-44 (α_A), 112-123 (α_B) and 153-163 (α_C) of Figure 1, and the five beta strands correspond to residues 83-86 (β_I), 95-100 (β_{II}), 59-66 (β_{III}), 14-20 (β_{IV}), and 49-53 (β_V) of Figure 1.
8. A crystallized catalytic fragment of MMP-13 complexed with N-Hydroxy-2-[(4-methoxy-benzenesulfonyl)-pyridin-3-ylmethyl-amino]-3-

-69-

methyl-benzamide ("Compound A").

9. The crystallized complex of Claim 8, wherein the catalytic fragment of MMP-13 comprises the amino acid residues of Figure 1.

10. The crystallized complex of Claim 9, characterized as being in orthorhombic form with space group P21212, and having unit cell parameters of $a=108.3\text{\AA}$, $b=79.8\text{\AA}$, and $c=36.1\text{\AA}$.

11. The crystallized complex of Claim 10, further characterized as consisting of two molecules of MMP-13:Compound A complex in the asymmetric unit.

12. The crystallized complex of Claim 11, wherein the secondary structure of the catalytic fragment of MMP-13 comprises three alpha helices and a mixed parallel and anti-parallel beta sheet comprising five beta strands.

13. The crystallized complex of Claim 12, wherein the alpha helices and beta strands are configured in the order β_I , α_A , β_{II} , β_{III} , β_{IV} , β_V , α_B , and α_C .

14. The crystallized complex of Claim 13, wherein the three alpha helices correspond to residues 28-44 (α_A), 112-123 (α_B) and 153-163 (α_C) of Figure 1, and the five beta strands correspond to residues 83-86 (β_I), 95-100 (β_{II}), 59-66 (β_{III}), 14-20 (β_{IV}), and 49-53 (β_V) of Figure 1.

15. An active site of MMP-13, characterized by a catalytic zinc, a beta strand, a Ca^{2+} binding loop, an alpha helix, and a random coil region.

-70-

16. The active site of Claim 15, wherein the beta strand comprises residues N14, L15, T16, Y17, R18, I19, and V20 according to Figure 1, the Ca^{2+} binding loop comprises residues F75, D76, G77, P78, and S79 according to Figure 1, the alpha helix comprises residues N112, L113, F114, L115, V116, A117, A118, H119, E120, F121, G122, and H123 according to Figure 1, and the random coil region comprises residues P139, I140, and Y141 according to Figure 1.

17. The active site of Claim 16, wherein said active site comprises the relative structural coordinates of the catalytic zinc and amino acid residues N14, L15, T16, Y17, R18, I19, V20, F75, D76, G77, P78, S79, N112, L113, F114, L115, V116, A117, A118, H119, E120, F121, G122, H123, P139, I140, and Y141 according to the solution or crystal coordinates of Figures 4 or 5, respectively, in each case, \pm a root mean square deviation from the catalytic zinc and the conserved backbone atoms of said amino acids of not more than 1.5Å.

18. The active site of Claim 17, further comprising the relative structural coordinates of amino acid residues G80, L81, L82, A83, H84, A85, K109, G110, Y111, S124, L125, G126, L127, D128, H129, S130, K131, D132, P133, G134, A135, L136, M137, F138, T142, Y143, T144, and G145 according to the solution or crystal coordinates of Figures 4 or 5, respectively, in each case, \pm a root mean square deviation from the conserved backbone atoms of said amino acids of not more than 1.5Å.

19. The active site of Claim 18, further comprising the relative structural coordinates of amino acid residues F149 and P152 according to the solution or crystal coordinates of Figures 4 or 5, respectively, in each case, \pm a root mean square deviation from the conserved backbone atoms of said amino acids of not more than 1.5Å.

-71-

20. An active site of MMP-13 comprising the relative structural coordinates of a catalytic zinc and amino acid residues L81, L82, L115, V116, H119, L136 and I140 according to the solution or crystal coordinates of Figures 4 or 5, respectively, in each case, \pm a root mean square deviation from the catalytic zinc and the conserved backbone atoms of said amino acids of not more than 1.5Å.

21. A method for identifying a potential inhibitor or activator of MMP-13, comprising the steps of:

(a) using a three dimensional structure of MMP-13 as defined by the relative structural coordinates of amino acids encoding MMP-13 according to Figures 4 or 5, \pm a root mean square deviation from the conserved backbone atoms of said amino acids of not more than 1.5Å;

(b) employing said three-dimensional structure to design or select a potential inhibitor or activator; and

(c) synthesizing or obtaining said potential inhibitor or activator.

22. The method according to Claim 21, wherein the potential inhibitor is designed de novo.

23. The method according to Claim 21, wherein the potential inhibitor is designed from a known inhibitor.

24. The method of Claim 22, further comprising the step of contacting the potential inhibitor with MMP-13 in the presence of a substrate to determine the ability of the potential inhibitor to inhibit MMP-13.

25. The method of Claim 23, further comprising the step of contacting the potential inhibitor with MMP-13 in the presence of a substrate to determine the ability of the potential inhibitor to inhibit MMP-13.

-72-

26. The method according to Claim 21, wherein the step of employing the three dimensional structure to design or select the potential inhibitor comprises the steps of:

- (a) identifying chemical entities or fragments capable of associating with MMP-13; and
- (b) assembling the identified chemical entities or fragments into a single molecule to provide the structure of the potential inhibitor.

27. The method according to Claim 26, wherein the potential inhibitor is designed de novo.

28. The method according to Claim 26, wherein the potential inhibitor is designed from a known inhibitor.

29. The method of Claim 27, further comprising the step of contacting the potential inhibitor with MMP-13 in the presence of a substrate to determine the ability of the potential inhibitor to inhibit MMP-13.

30. The method of Claim 28, further comprising the step of contacting the potential inhibitor with MMP-13 in the presence of a substrate to determine the ability of the potential inhibitor to inhibit MMP-13.

31. An inhibitor identified or designed by the method of Claim 21.

32. An inhibitor identified or designed by the method of Claim 26.

1 / 75

YNVFP 5	RTLKW 10	SKMNL 15	TYRIV 20	NYTPD 25
MTHSE 30	VEKAF 35	KKAFK 40	VWSDV 45	TPLNF 50
TRLHD 55	GIADI 60	MISFG 65	IKEHG 70	DFYPF 75
DGPSG 80	LLAHA 85	FPPGP 90	NYGGD 95	AHFDD 100
DETWT 105	SSSKG 110	YNLFL 115	VAAHE 120	FGHSL 125
GLDHS 130	KDPGA 135	LMFPI 140	YTYTG 145	KSHFM 150
LPDDD 155	VQGIQ 160	SLYG 164		

FIG. 1

2 / 75

Sequence 1: MMP-13

Sequence 2: MMP-1

Identity score: 58.9 %

VGEYNVFPRTLKWSKMNLTYRIVNYTPDMTHSEVEKAFKKAFKVWSDVTPLNFTRLHDGIADIMISFGIKEHGDFYFPDG
 LTEGN PR WEQTHLTURIENYTPDLPRADVDAIEKAFQLWSNVTPLTFTKVSEGQADIMISFVRGDHRDNSPFDG

PSGLLAHAFPPGPNYGGDAHFDDDETWTS

SSKGYNLF

LVAAHEFGHSLGLDHSKDPGALMF

PIYTYTGKSHFMLPDDDVQ

PGNLAHAFQPGPGIGGDAHFDEDERWTNNFREYNLHRVAAHELGHSLGLS HST DIGALMYPSYTFESGDVQ LAODD
 ID

###

GIQSLYGPGEDEPN
 GIQAIYGRSQ

FIG. 2A

Sequence 1: MMP-13

Sequence 2: MMP-8

Identity score: 61.4 %

VGEYNVFPRTLKWSKMNLTYRIVNYT PDMTH S EVEKAFKKAFKVWSDVTPLNFTRLHDGIADIMISFGIKEHGDFYFPDG
 NPKWER T NLTYRIRNYTP QLSEA EVERAI KDAFEL WSVASPLI
 FTRISQGEADINIAFYQRDHGDNSPFDG

..

PSGLLAHAFPPGPNYGGDAHFDDDETWTSSSKGYNLFVAAHEFGHSLGLDHSKDPGALMF PIYTYTGKSHFMLPDDDVQ
 PNGILAHAFQPGQGIGGDAHFDAEETWTNTSANYNLFVAA HEFGHSLGLAHSSDPGALMYPNYAF RETSNYSLPODD ID

###

GIQSLYGPGEDEPN
 GIQAIYG

FIG. 2B

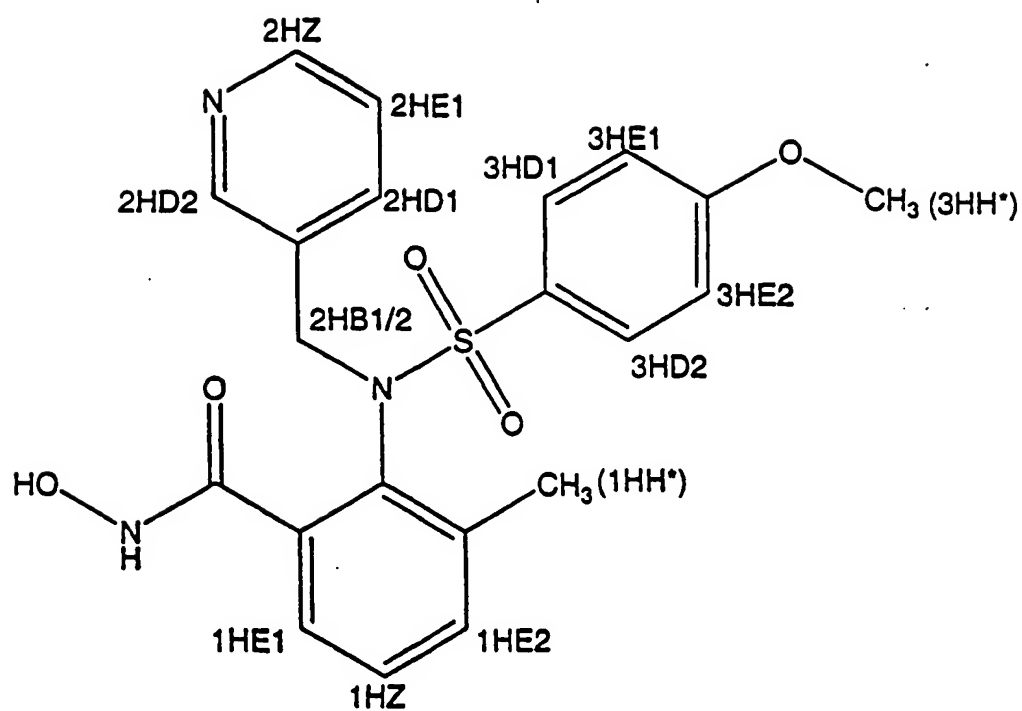


FIG. 3

4/75

	Atom Type	Res.		X	Y	Z		
ATOM	1 N	THR	7	-12.675	-13.911	-8.815	1.00	0.83
ATOM	2 HN	THR	7	-12.001	-14.254	-8.192	1.00	1.22
ATOM	3 CA	THR	7	-14.063	-13.649	-8.340	1.00	0.63
ATOM	4 HA	THR	7	-14.744	-14.330	-8.830	1.00	0.73
ATOM	5 CB	THR	7	-14.132	-13.858	-6.825	1.00	0.61
ATOM	6 HB	THR	7	-13.473	-13.158	-6.335	1.00	0.66
ATOM	7 OG1	THR	7	-13.730	-15.185	-6.514	1.00	0.71
ATOM	8 HG1	THR	7	-13.721	-15.690	-7.330	1.00	1.07
ATOM	9 CG2	THR	7	-15.564	-13.628	-6.336	1.00	0.67
ATOM	10 HG21	THR	7	-15.712	-12.577	-6.139	1.00	1.14
ATOM	11 HG22	THR	7	-15.728	-14.191	-5.429	1.00	1.32
ATOM	12 HG23	THR	7	-16.261	-13.955	-7.093	1.00	1.23
ATOM	13 C	THR	7	-14.451	-12.208	-8.678	1.00	0.52
ATOM	14 O	THR	7	-15.416	-11.962	-9.374	1.00	0.65
ATOM	15 N	LEU	8	-13.704	-11.254	-8.195	1.00	0.47
ATOM	16 HN	LEU	8	-12.927	-11.473	-7.639	1.00	0.61
ATOM	17 CA	LEU	8	-14.027	-9.831	-8.495	1.00	0.42
ATOM	18 HA	LEU	8	-15.098	-9.715	-8.575	1.00	0.43
ATOM	19 CB	LEU	8	-13.495	-8.937	-7.370	1.00	0.52
ATOM	20 HB1	LEU	8	-13.721	-7.905	-7.591	1.00	0.54
ATOM	21 HB2	LEU	8	-12.424	-9.060	-7.292	1.00	0.58
ATOM	22 CG	LEU	8	-14.151	-9.331	-6.042	1.00	0.60
ATOM	23 HG	LEU	8	-13.958	-10.376	-5.844	1.00	0.60
ATOM	24 CD1	LEU	8	-13.566	-8.484	-4.910	1.00	0.74
ATOM	25 HD11	LEU	8	-13.899	-8.875	-3.960	1.00	1.22
ATOM	26 HD12	LEU	8	-13.900	-7.462	-5.016	1.00	1.26
ATOM	27 HD13	LEU	8	-12.488	-8.518	-4.956	1.00	1.31
ATOM	28 CD2	LEU	8	-15.664	-9.096	-6.117	1.00	0.61
ATOM	29 HD21	LEU	8	-15.871	-8.278	-6.791	1.00	1.13
ATOM	30 HD22	LEU	8	-16.040	-8.856	-5.134	1.00	1.18
ATOM	31 HD23	LEU	8	-16.149	-9.991	-6.478	1.00	1.26
ATOM	32 C	LEU	8	-13.374	-9.438	-9.822	1.00	0.40
ATOM	33 O	LEU	8	-12.218	-9.722	-10.064	1.00	0.45
ATOM	34 N	LYS	9	-14.109	-8.795	-10.687	1.00	0.36
ATOM	35 HN	LYS	9	-15.042	-8.581	-10.474	1.00	0.36
ATOM	36 CA	LYS	9	-13.536	-8.393	-12.002	1.00	0.37
ATOM	37 HA	LYS	9	-12.521	-8.050	-11.862	1.00	0.39
ATOM	38 CB	LYS	9	-13.539	-9.599	-12.944	1.00	0.50
ATOM	39 HB1	LYS	9	-12.851	-10.344	-12.573	1.00	0.60

FIG. 4

5/75

ATOM	40	HB2	LYS	9	-13.233	-9.286	-13.932	1.00	0.48
ATOM	41	CG	LYS	9	-14.948	-10.193	-13.007	1.00	0.60
ATOM	42	HG1	LYS	9	-15.632	-9.455	-13.398	1.00	0.66
ATOM	43	HG2	LYS	9	-15.260	-10.482	-12.014	1.00	0.78
ATOM	44	CD	LYS	9	-14.951	-11.421	-13.921	1.00	0.94
ATOM	45	HD1	LYS	9	-13.944	-11.794	-14.033	1.00	1.57
ATOM	46	HD2	LYS	9	-15.344	-11.147	-14.889	1.00	1.62
ATOM	47	CE	LYS	9	-15.829	-12.511	-13.303	1.00	0.57
ATOM	48	HE1	LYS	9	-16.776	-12.086	-13.007	1.00	1.15
ATOM	49	HE2	LYS	9	-15.333	-12.924	-12.437	1.00	1.10
ATOM	50	NZ	LYS	9	-16.060	-13.591	-14.304	1.00	1.61
ATOM	51	HZ1	LYS	9	-15.181	-14.127	-14.445	1.00	2.14
ATOM	52	HZ2	LYS	9	-16.358	-13.168	-15.207	1.00	2.13
ATOM	53	HZ3	LYS	9	-16.802	-14.231	-13.959	1.00	2.14
ATOM	54	C	LYS	9	-14.377	-7.265	-12.605	1.00	0.32
ATOM	55	O	LYS	9	-15.493	-7.021	-12.191	1.00	0.34
ATOM	56	N	TRP	10	-13.850	-6.571	-13.577	1.00	0.31
ATOM	57	HN	TRP	10	-12.947	-6.781	-13.895	1.00	0.33
ATOM	58	CA	TRP	10	-14.618	-5.456	-14.201	1.00	0.30
ATOM	59	HA	TRP	10	-15.030	-4.826	-13.427	1.00	0.29
ATOM	60	CB	TRP	10	-13.684	-4.630	-15.088	1.00	0.29
ATOM	61	HB1	TRP	10	-14.264	-3.917	-15.655	1.00	0.32
ATOM	62	HB2	TRP	10	-13.157	-5.286	-15.765	1.00	0.33
ATOM	63	CG	TRP	10	-12.699	-3.901	-14.230	1.00	0.25
ATOM	64	CD1	TRP	10	-11.516	-4.405	-13.812	1.00	0.30
ATOM	65	HD1	TRP	10	-11.137	-5.390	-14.040	1.00	0.37
ATOM	66	CD2	TRP	10	-12.786	-2.553	-13.683	1.00	0.21
ATOM	67	NE1	TRP	10	-10.872	-3.454	-13.042	1.00	0.30
ATOM	68	HE1	TRP	10	-9.996	-3.569	-12.617	1.00	0.36
ATOM	69	CE2	TRP	10	-11.614	-2.295	-12.934	1.00	0.23
ATOM	70	CE3	TRP	10	-13.758	-1.538	-13.763	1.00	0.24
ATOM	71	HE3	TRP	10	-14.663	-1.706	-14.328	1.00	0.29
ATOM	72	CZ2	TRP	10	-11.412	-1.075	-12.287	1.00	0.22
ATOM	73	HZ2	TRP	10	-10.509	-0.903	-11.720	1.00	0.27
ATOM	74	CZ3	TRP	10	-13.558	-0.309	-13.113	1.00	0.25
ATOM	75	HZ3	TRP	10	-14.310	0.463	-13.181	1.00	0.32
ATOM	76	CH2	TRP	10	-12.387	-0.078	-12.376	1.00	0.23
ATOM	77	HH2	TRP	10	-12.238	0.870	-11.879	1.00	0.26
ATOM	78	C	TRP	10	-15.755	-6.031	-15.050	1.00	0.39
ATOM	79	O	TRP	10	-15.641	-7.098	-15.620	1.00	0.48
ATOM	80	N	SER	11	-16.855	-5.332	-15.132	1.00	0.43
ATOM	81	HN	SER	11	-16.927	-4.476	-14.660	1.00	0.44
ATOM	82	CA	SER	11	-18.006	-5.835	-15.936	1.00	0.52
ATOM	83	HA	SER	11	-18.003	-6.915	-15.930	1.00	0.59
ATOM	84	CB	SER	11	-19.313	-5.330	-15.325	1.00	0.64
ATOM	85	HB1	SER	11	-19.120	-4.425	-14.763	1.00	1.16
ATOM	86	HB2	SER	11	-19.718	-6.079	-14.666	1.00	1.20
ATOM	87	OG	SER	11	-20.246	-5.067	-16.365	1.00	1.39
ATOM	88	HG	SER	11	-19.821	-4.495	-17.008	1.00	1.92
ATOM	89	C	SER	11	-17.893	-5.335	-17.379	1.00	0.47
ATOM	90	O	SER	11	-18.785	-5.528	-18.181	1.00	0.60
ATOM	91	N	LYS	12	-16.808	-4.692	-17.715	1.00	0.42
ATOM	92	HN	LYS	12	-16.101	-4.543	-17.053	1.00	0.51
ATOM	93	CA	LYS	12	-16.646	-4.178	-19.107	1.00	0.41
ATOM	94	HA	LYS	12	-17.243	-4.775	-19.781	1.00	0.47
ATOM	95	CB	LYS	12	-17.116	-2.722	-19.167	1.00	0.43
ATOM	96	HB1	LYS	12	-18.168	-2.674	-18.926	1.00	0.50
ATOM	97	HB2	LYS	12	-16.957	-2.334	-20.163	1.00	0.46
ATOM	98	CG	LYS	12	-16.327	-1.882	-18.160	1.00	0.41
ATOM	99	HG1	LYS	12	-15.275	-1.922	-18.401	1.00	0.37
ATOM	100	HG2	LYS	12	-16.484	-2.272	-17.164	1.00	0.42
ATOM	101	CD	LYS	12	-16.805	-0.430	-18.223	1.00	0.50
ATOM	102	HD1	LYS	12	-17.856	-0.386	-17.981	1.00	0.56
ATOM	103	HD2	LYS	12	-16.648	-0.044	-19.220	1.00	0.65
ATOM	104	CE	LYS	12	-16.018	0.412	-17.218	1.00	0.61
ATOM	105	HE1	LYS	12	-15.054	0.665	-17.636	1.00	1.15
ATOM	106	HE2	LYS	12	-15.879	-0.151	-16.307	1.00	1.16
ATOM	107	NZ	LYS	12	-16.773	1.661	-16.920	1.00	1.39
ATOM	108	HZ1	LYS	12	-16.498	2.018	-15.983	1.00	1.90
ATOM	109	HZ2	LYS	12	-17.794	1.458	-16.927	1.00	1.87
ATOM	110	HZ3	LYS	12	-16.556	2.379	-17.640	1.00	1.97
ATOM	111	C	LYS	12	-15.175	-4.269	-19.521	1.00	0.36
ATOM	112	O	LYS	12	-14.284	-4.250	-18.695	1.00	0.34
ATOM	113	N	MET	13	-14.917	-4.380	-20.796	1.00	0.37
ATOM	114	HN	MET	13	-15.652	-4.402	-21.443	1.00	0.40
ATOM	115	CA	MET	13	-13.506	-4.487	-21.269	1.00	0.38
ATOM	116	HA	MET	13	-12.910	-4.964	-20.506	1.00	0.39

6/75

ATOM	117	CB	MET	13	-13.469	-5.332	-22.543	1.00	0.46
ATOM	118	HB1	MET	13	-12.523	-5.189	-23.043	1.00	0.53
ATOM	119	HB2	MET	13	-14.273	-5.031	-23.199	1.00	0.42
ATOM	120	CG	MET	13	-13.632	-6.809	-22.178	1.00	0.64
ATOM	121	HG1	MET	13	-12.857	-7.097	-21.483	1.00	1.26
ATOM	122	HG2	MET	13	-13.556	-7.411	-23.071	1.00	1.37
ATOM	123	SD	MET	13	-15.252	-7.067	-21.414	1.00	1.22
ATOM	124	CE	MET	13	-14.663	-7.870	-19.903	1.00	0.57
ATOM	125	HE1	MET	13	-14.020	-7.189	-19.362	1.00	1.16
ATOM	126	HE2	MET	13	-14.107	-8.758	-20.158	1.00	1.09
ATOM	127	HE3	MET	13	-15.508	-8.141	-19.286	1.00	1.20
ATOM	128	C	MET	13	-12.936	-3.095	-21.560	1.00	0.32
ATOM	129	O	MET	13	-11.793	-2.957	-21.948	1.00	0.35
ATOM	130	N	ASN	14	-13.718	-2.064	-21.371	1.00	0.28
ATOM	131	HN	ASN	14	-14.635	-2.199	-21.052	1.00	0.29
ATOM	132	CA	ASN	14	-13.217	-0.681	-21.631	1.00	0.26
ATOM	133	HA	ASN	14	-12.359	-0.725	-22.286	1.00	0.29
ATOM	134	CB	ASN	14	-14.319	0.148	-22.297	1.00	0.30
ATOM	135	HB1	ASN	14	-14.025	1.186	-22.318	1.00	0.31
ATOM	136	HB2	ASN	14	-15.235	0.043	-21.735	1.00	0.31
ATOM	137	CG	ASN	14	-14.539	-0.346	-23.729	1.00	0.37
ATOM	138	OD1	ASN	14	-13.677	-0.981	-24.304	1.00	1.16
ATOM	139	ND2	ASN	14	-15.664	-0.077	-24.334	1.00	1.05
ATOM	140	HD21	ASN	14	-16.359	0.435	-23.871	1.00	1.81
ATOM	141	HD22	ASN	14	-15.812	-0.386	-25.252	1.00	1.06
ATOM	142	C	ASN	14	-12.813	-0.024	-20.309	1.00	0.22
ATOM	143	O	ASN	14	-13.566	-0.019	-19.357	1.00	0.23
ATOM	144	N	LEU	15	-11.630	0.533	-20.247	1.00	0.21
ATOM	145	HN	LEU	15	-11.042	0.517	-21.031	1.00	0.24
ATOM	146	CA	LEU	15	-11.171	1.194	-18.987	1.00	0.18
ATOM	147	HA	LEU	15	-12.025	1.447	-18.379	1.00	0.19
ATOM	148	CB	LEU	15	-10.250	0.243	-18.210	1.00	0.18
ATOM	149	HB1	LEU	15	-9.812	0.769	-17.375	1.00	0.19
ATOM	150	HB2	LEU	15	-9.463	-0.102	-18.865	1.00	0.21
ATOM	151	CG	LEU	15	-11.046	-0.964	-17.696	1.00	0.19
ATOM	152	HG	LEU	15	-11.547	-1.442	-18.525	1.00	0.20
ATOM	153	CD1	LEU	15	-10.086	-1.961	-17.044	1.00	0.20
ATOM	154	HD11	LEU	15	-9.726	-1.556	-16.110	1.00	0.98
ATOM	155	HD12	LEU	15	-9.251	-2.141	-17.704	1.00	1.04
ATOM	156	HD13	LEU	15	-10.604	-2.890	-16.857	1.00	1.07
ATOM	157	CD2	LEU	15	-12.083	-0.513	-16.658	1.00	0.21
ATOM	158	HD21	LEU	15	-12.114	-1.228	-15.850	1.00	1.07
ATOM	159	HD22	LEU	15	-13.055	-0.456	-17.122	1.00	1.00
ATOM	160	HD23	LEU	15	-11.814	0.457	-16.268	1.00	1.04
ATOM	161	C	LEU	15	-10.397	2.471	-19.334	1.00	0.18
ATOM	162	O	LEU	15	-9.785	2.570	-20.380	1.00	0.20
ATOM	163	N	THR	16	-10.425	3.447	-18.460	1.00	0.18
ATOM	164	HN	THR	16	-10.929	3.338	-17.627	1.00	0.18
ATOM	165	CA	THR	16	-9.699	4.729	-18.722	1.00	0.19
ATOM	166	HA	THR	16	-9.051	4.617	-19.574	1.00	0.20
ATOM	167	CB	THR	16	-10.716	5.839	-18.996	1.00	0.22
ATOM	168	HB	THR	16	-10.198	6.729	-19.315	1.00	0.24
ATOM	169	OG1	THR	16	-11.445	6.112	-17.808	1.00	0.23
ATOM	170	HG1	THR	16	-11.821	5.286	-17.495	1.00	0.98
ATOM	171	CG2	THR	16	-11.680	5.393	-20.096	1.00	0.26
ATOM	172	HG21	THR	16	-12.200	6.254	-20.489	1.00	1.05
ATOM	173	HG22	THR	16	-12.396	4.696	-19.686	1.00	1.02
ATOM	174	HG23	THR	16	-11.125	4.914	-20.889	1.00	1.05
ATOM	175	C	THR	16	-8.864	5.100	-17.495	1.00	0.17
ATOM	176	O	THR	16	-9.157	4.687	-16.391	1.00	0.16
ATOM	177	N	TYR	17	-7.826	5.878	-17.675	1.00	0.18
ATOM	178	HN	TYR	17	-7.603	6.202	-18.574	1.00	0.19
ATOM	179	CA	TYR	17	-6.981	6.268	-16.507	1.00	0.17
ATOM	180	HA	TYR	17	-7.585	6.233	-15.615	1.00	0.17
ATOM	181	CB	TYR	17	-5.814	5.288	-16.362	1.00	0.19
ATOM	182	HB1	TYR	17	-6.194	4.278	-16.347	1.00	0.19
ATOM	183	HB2	TYR	17	-5.292	5.488	-15.438	1.00	0.20
ATOM	184	CG	TYR	17	-4.857	5.445	-17.520	1.00	0.22
ATOM	185	CD1	TYR	17	-5.037	4.685	-18.682	1.00	0.26
ATOM	186	HD1	TYR	17	-5.867	3.998	-18.755	1.00	0.27
ATOM	187	CD2	TYR	17	-3.782	6.336	-17.426	1.00	0.25
ATOM	188	HD2	TYR	17	-3.643	6.923	-16.530	1.00	0.26
ATOM	189	CE1	TYR	17	-4.143	4.817	-19.751	1.00	0.31
ATOM	190	HE1	TYR	17	-4.282	4.231	-20.647	1.00	0.36
ATOM	191	CE2	TYR	17	-2.888	6.470	-18.496	1.00	0.30
ATOM	192	HE2	TYR	17	-2.059	7.158	-18.424	1.00	0.35
ATOM	193	CZ	TYR	17	-3.068	5.710	-19.658	1.00	0.32

7/75

ATOM	194	OH	TYR	17	-2.186	5.839	-20.711	1.00	0.39
ATOM	195	HH	TYR	17	-1.696	5.016	-20.790	1.00	0.85
ATOM	196	C	TYR	17	-6.448	7.692	-16.690	1.00	0.19
ATOM	197	O	TYR	17	-6.414	8.220	-17.784	1.00	0.21
ATOM	198	N	ARG	18	-6.044	8.320	-15.616	1.00	0.19
ATOM	199	HN	ARG	18	-6.089	7.874	-14.747	1.00	0.19
ATOM	200	CA	ARG	18	-5.523	9.714	-15.712	1.00	0.22
ATOM	201	HA	ARG	18	-5.131	9.877	-16.704	1.00	0.24
ATOM	202	CB	ARG	18	-6.674	10.691	-15.447	1.00	0.27
ATOM	203	HB1	ARG	18	-6.978	10.613	-14.412	1.00	0.31
ATOM	204	HB2	ARG	18	-7.507	10.442	-16.083	1.00	0.30
ATOM	205	CG	ARG	18	-6.229	12.127	-15.733	1.00	0.35
ATOM	206	HG1	ARG	18	-5.504	12.137	-16.531	1.00	0.93
ATOM	207	HG2	ARG	18	-5.790	12.549	-14.843	1.00	0.85
ATOM	208	CD	ARG	18	-7.447	12.946	-16.149	1.00	0.81
ATOM	209	HD1	ARG	18	-8.216	12.867	-15.378	1.00	1.29
ATOM	210	HD2	ARG	18	-7.838	12.561	-17.068	1.00	1.63
ATOM	211	NE	ARG	18	-7.030	14.362	-16.406	1.00	1.52
ATOM	212	HE	ARG	18	-7.071	14.711	-17.318	1.00	2.11
ATOM	213	CZ	ARG	18	-6.561	15.119	-15.456	1.00	2.24
ATOM	214	NH1	ARG	18	-6.119	16.314	-15.736	1.00	3.18
ATOM	215	HH11	ARG	18	-6.142	16.647	-16.679	1.00	3.48
ATOM	216	HH12	ARG	18	-5.760	16.898	-15.009	1.00	3.84
ATOM	217	NH2	ARG	18	-6.564	14.700	-14.220	1.00	2.63
ATOM	218	HH21	ARG	18	-6.928	13.795	-14.000	1.00	2.44
ATOM	219	HH22	ARG	18	-6.205	15.285	-13.493	1.00	3.49
ATOM	220	C	ARG	18	-4.413	9.931	-14.676	1.00	0.21
ATOM	221	O	ARG	18	-4.550	9.576	-13.522	1.00	0.23
ATOM	222	N	ILE	19	-3.314	10.514	-15.079	1.00	0.21
ATOM	223	HN	ILE	19	-3.223	10.794	-16.014	1.00	0.22
ATOM	224	CA	ILE	19	-2.196	10.755	-14.118	1.00	0.23
ATOM	225	HA	ILE	19	-2.200	9.985	-13.360	1.00	0.25
ATOM	226	CB	ILE	19	-0.864	10.721	-14.875	1.00	0.25
ATOM	227	HB	ILE	19	-0.862	11.491	-15.633	1.00	0.25
ATOM	228	CG1	ILE	19	-0.702	9.341	-15.531	1.00	0.29
ATOM	229	HG11	ILE	19	-1.607	9.092	-16.065	1.00	0.82
ATOM	230	HG12	ILE	19	-0.525	8.601	-14.765	1.00	0.97
ATOM	231	CG2	ILE	19	0.291	10.962	-13.893	1.00	0.29
ATOM	232	HG21	ILE	19	1.231	10.914	-14.420	1.00	1.08
ATOM	233	HG22	ILE	19	0.272	10.206	-13.123	1.00	1.09
ATOM	234	HG23	ILE	19	0.187	11.937	-13.440	1.00	1.00
ATOM	235	CD1	ILE	19	0.477	9.345	-16.512	1.00	0.93
ATOM	236	HD11	ILE	19	1.402	9.216	-15.970	1.00	1.59
ATOM	237	HD12	ILE	19	0.501	10.280	-17.050	1.00	1.50
ATOM	238	HD13	ILE	19	0.360	8.533	-17.214	1.00	1.55
ATOM	239	C	ILE	19	-2.381	12.126	-13.454	1.00	0.23
ATOM	240	O	ILE	19	-2.355	13.150	-14.108	1.00	0.23
ATOM	241	N	VAL	20	-2.563	12.152	-12.161	1.00	0.25
ATOM	242	HN	VAL	20	-2.578	11.314	-11.653	1.00	0.27
ATOM	243	CA	VAL	20	-2.746	13.454	-11.454	1.00	0.27
ATOM	244	HA	VAL	20	-3.496	14.035	-11.970	1.00	0.27
ATOM	245	CB	VAL	20	-3.202	13.205	-10.015	1.00	0.31
ATOM	246	HB	VAL	20	-2.522	12.517	-9.534	1.00	0.32
ATOM	247	CG1	VAL	20	-3.216	14.529	-9.247	1.00	0.33
ATOM	248	HG11	VAL	20	-3.607	15.310	-9.883	1.00	0.97
ATOM	249	HG12	VAL	20	-2.211	14.782	-8.944	1.00	1.08
ATOM	250	HG13	VAL	20	-3.842	14.432	-8.372	1.00	1.10
ATOM	251	CG2	VAL	20	-4.612	12.611	-10.028	1.00	0.33
ATOM	252	HG21	VAL	20	-5.296	13.317	-10.476	1.00	1.05
ATOM	253	HG22	VAL	20	-4.924	12.401	-9.016	1.00	1.03
ATOM	254	HG23	VAL	20	-4.612	11.697	-10.602	1.00	1.11
ATOM	255	C	VAL	20	-1.424	14.231	-11.451	1.00	0.27
ATOM	256	O	VAL	20	-1.403	15.435	-11.611	1.00	0.26
ATOM	257	N	ASN	21	-0.321	13.555	-11.259	1.00	0.28
ATOM	258	HN	ASN	21	-0.357	12.585	-11.124	1.00	0.30
ATOM	259	CA	ASN	21	0.992	14.265	-11.235	1.00	0.29
ATOM	260	HA	ASN	21	0.973	15.076	-11.949	1.00	0.26
ATOM	261	CB	ASN	21	1.235	14.829	-9.834	1.00	0.33
ATOM	262	HB1	ASN	21	0.544	15.637	-9.646	1.00	0.33
ATOM	263	HB2	ASN	21	2.249	15.199	-9.766	1.00	0.35
ATOM	264	CG	ASN	21	1.022	13.727	-8.795	1.00	0.40
ATOM	265	OD1	ASN	21	0.459	12.694	-9.097	1.00	1.01
ATOM	266	ND2	ASN	21	1.445	13.908	-7.574	1.00	0.88
ATOM	267	HD21	ASN	21	1.895	14.743	-7.330	1.00	1.50
ATOM	268	HD22	ASN	21	1.312	13.208	-6.901	1.00	0.88
ATOM	269	C	ASN	21	2.116	13.291	-11.606	1.00	0.34
ATOM	270	O	ASN	21	1.929	12.090	-11.619	1.00	0.37

8/75

ATOM	271	N	TYR	22	3.274	13.810	-11.933	1.00	0.38
ATOM	272	HN	TYR	22	3.387	14.783	-11.932	1.00	0.38
ATOM	273	CA	TYR	22	4.417	12.935	-12.340	1.00	0.46
ATOM	274	HA	TYR	22	4.067	11.929	-12.509	1.00	0.45
ATOM	275	CB	TYR	22	5.028	13.481	-13.630	1.00	0.49
ATOM	276	HB1	TYR	22	5.845	12.846	-13.938	1.00	0.56
ATOM	277	HB2	TYR	22	5.397	14.482	-13.457	1.00	0.53
ATOM	278	CG	TYR	22	3.981	13.513	-14.714	1.00	0.43
ATOM	279	CD1	TYR	22	3.684	12.352	-15.436	1.00	0.38
ATOM	280	HD1	TYR	22	4.199	11.430	-15.212	1.00	0.39
ATOM	281	CD2	TYR	22	3.313	14.708	-15.003	1.00	0.46
ATOM	282	HD2	TYR	22	3.543	15.603	-14.445	1.00	0.51
ATOM	283	CE1	TYR	22	2.718	12.386	-16.447	1.00	0.36
ATOM	284	HE1	TYR	22	2.490	11.491	-17.004	1.00	0.36
ATOM	285	CE2	TYR	22	2.345	14.742	-16.013	1.00	0.44
ATOM	286	HE2	TYR	22	1.828	15.663	-16.235	1.00	0.49
ATOM	287	CZ	TYR	22	2.048	13.581	-16.735	1.00	0.39
ATOM	288	OH	TYR	22	1.095	13.615	-17.733	1.00	0.43
ATOM	289	HH	TYR	22	1.173	14.457	-18.187	1.00	0.92
ATOM	290	C	TYR	22	5.499	12.923	-11.258	1.00	0.56
ATOM	291	O	TYR	22	6.554	12.378	-11.470	1.00	1.38
ATOM	292	N	THR	23	5.240	13.544	-10.130	1.00	0.47
ATOM	293	HN	THR	23	4.372	13.987	-10.023	1.00	1.08
ATOM	294	CA	THR	23	6.237	13.623	-9.004	1.00	0.46
ATOM	295	HA	THR	23	5.848	14.338	-8.304	1.00	0.48
ATOM	296	CB	THR	23	6.361	12.265	-8.273	1.00	0.62
ATOM	297	HB	THR	23	5.383	11.969	-7.921	1.00	0.68
ATOM	298	OG1	THR	23	7.223	12.420	-7.156	1.00	0.86
ATOM	299	HG1	THR	23	7.941	11.788	-7.244	1.00	1.28
ATOM	300	CG2	THR	23	6.916	11.159	-9.181	1.00	0.59
ATOM	301	HG21	THR	23	7.753	11.533	-9.748	1.00	1.08
ATOM	302	HG22	THR	23	6.141	10.816	-9.850	1.00	1.16
ATOM	303	HG23	THR	23	7.245	10.332	-8.570	1.00	1.22
ATOM	304	C	THR	23	7.623	14.115	-9.523	1.00	0.40
ATOM	305	O	THR	23	8.077	13.699	-10.565	1.00	0.45
ATOM	306	N	PRO	24	8.302	15.016	-8.823	1.00	0.42
ATOM	307	CA	PRO	24	9.625	15.520	-9.311	1.00	0.42
ATOM	308	HA	PRO	24	9.534	15.918	-10.307	1.00	0.46
ATOM	309	CB	PRO	24	9.924	16.655	-8.335	1.00	0.50
ATOM	310	HB1	PRO	24	9.743	17.605	-8.815	1.00	0.57
ATOM	311	HB2	PRO	24	10.955	16.598	-8.014	1.00	0.49
ATOM	312	CG	PRO	24	8.995	16.507	-7.129	1.00	0.66
ATOM	313	HG1	PRO	24	8.613	17.475	-6.842	1.00	0.84
ATOM	314	HG2	PRO	24	9.537	16.069	-6.303	1.00	0.76
ATOM	315	CD	PRO	24	7.832	15.598	-7.529	1.00	0.56
ATOM	316	HD2	PRO	24	7.675	14.826	-6.786	1.00	0.62
ATOM	317	HD1	PRO	24	6.940	16.183	-7.680	1.00	0.61
ATOM	318	C	PRO	24	10.743	14.470	-9.253	1.00	0.40
ATOM	319	O	PRO	24	11.835	14.692	-9.737	1.00	0.40
ATOM	320	N	ASP	25	10.490	13.337	-8.662	1.00	0.44
ATOM	321	HN	ASP	25	9.608	13.172	-8.270	1.00	0.48
ATOM	322	CA	ASP	25	11.554	12.295	-8.577	1.00	0.48
ATOM	323	HA	ASP	25	12.393	12.695	-8.025	1.00	0.51
ATOM	324	CB	ASP	25	11.016	11.062	-7.847	1.00	0.57
ATOM	325	HB1	ASP	25	11.719	10.249	-7.945	1.00	0.61
ATOM	326	HB2	ASP	25	10.068	10.773	-8.276	1.00	0.56
ATOM	327	CG	ASP	25	10.827	11.394	-6.364	1.00	0.67
ATOM	328	OD1	ASP	25	10.079	10.689	-5.709	1.00	1.23
ATOM	329	OD2	ASP	25	11.437	12.348	-5.908	1.00	1.34
ATOM	330	C	ASP	25	12.025	11.916	-9.985	1.00	0.45
ATOM	331	O	ASP	25	13.179	11.597	-10.191	1.00	0.55
ATOM	332	N	MET	26	11.146	11.948	-10.955	1.00	0.40
ATOM	333	HN	MET	26	10.220	12.209	-10.767	1.00	0.41
ATOM	334	CA	MET	26	11.553	11.590	-12.348	1.00	0.42
ATOM	335	HA	MET	26	12.624	11.686	-12.447	1.00	0.49
ATOM	336	CB	MET	26	11.144	10.149	-12.656	1.00	0.53
ATOM	337	HB1	MET	26	11.282	9.954	-13.709	1.00	0.55
ATOM	338	HB2	MET	26	10.105	10.006	-12.397	1.00	0.51
ATOM	339	CG	MET	26	12.011	9.186	-11.846	1.00	0.71
ATOM	340	HG1	MET	26	11.783	9.288	-10.796	1.00	0.73
ATOM	341	HG2	MET	26	13.053	9.419	-12.009	1.00	0.77
ATOM	342	SD	MET	26	11.683	7.485	-12.380	1.00	0.89
ATOM	343	CE	MET	26	10.000	7.330	-11.728	1.00	0.59
ATOM	344	HE1	MET	26	9.292	7.456	-12.534	1.00	1.25
ATOM	345	HE2	MET	26	9.825	8.084	-10.979	1.00	1.23
ATOM	346	HE3	MET	26	9.877	6.352	-11.285	1.00	1.23
ATOM	347	C	MET	26	10.872	12.530	-13.344	1.00	0.34

9/75

ATOM	348	O	MET	26	9.897	13.184	-13.031	1.00	0.32
ATOM	349	N	THR	27	11.385	12.604	-14.544	1.00	0.33
ATOM	350	HN	THR	27	12.174	12.070	-14.773	1.00	0.38
ATOM	351	CA	THR	27	10.775	13.504	-15.562	1.00	0.32
ATOM	352	HA	THR	27	10.618	14.483	-15.133	1.00	0.35
ATOM	353	CB	THR	27	11.711	13.616	-16.768	1.00	0.39
ATOM	354	HB	THR	27	11.295	14.308	-17.484	1.00	0.42
ATOM	355	OG1	THR	27	11.852	12.338	-17.371	1.00	0.37
ATOM	356	HG1	THR	27	12.765	12.242	-17.653	1.00	0.94
ATOM	357	CG2	THR	27	13.080	14.121	-16.313	1.00	0.51
ATOM	358	HG21	THR	27	13.602	14.553	-17.154	1.00	1.14
ATOM	359	HG22	THR	27	13.655	13.297	-15.918	1.00	1.11
ATOM	360	HG23	THR	27	12.951	14.871	-15.546	1.00	1.12
ATOM	361	C	THR	27	9.436	12.921	-16.013	1.00	0.27
ATOM	362	O	THR	27	9.177	11.743	-15.864	1.00	0.24
ATOM	363	N	HIS	28	8.580	13.740	-16.554	1.00	0.32
ATOM	364	HN	HIS	28	8.807	14.688	-16.657	1.00	0.37
ATOM	365	CA	HIS	28	7.253	13.241	-17.004	1.00	0.34
ATOM	366	HA	HIS	28	6.715	12.833	-16.161	1.00	0.36
ATOM	367	CB	HIS	28	6.457	14.403	-17.601	1.00	0.46
ATOM	368	HB1	HIS	28	5.428	14.104	-17.736	1.00	0.71
ATOM	369	HB2	HIS	28	6.880	14.676	-18.557	1.00	0.88
ATOM	370	CG	HIS	28	6.516	15.583	-16.669	1.00	0.73
ATOM	371	ND1	HIS	28	6.056	16.838	-17.036	1.00	1.66
ATOM	372	HD1	HIS	28	5.659	17.080	-17.898	1.00	2.30
ATOM	373	CD2	HIS	28	6.987	15.716	-15.387	1.00	1.33
ATOM	374	HD2	HIS	28	7.423	14.922	-14.798	1.00	2.01
ATOM	375	CE1	HIS	28	6.258	17.664	-15.993	1.00	1.95
ATOM	376	HE1	HIS	28	5.993	18.711	-15.990	1.00	2.70
ATOM	377	NE2	HIS	28	6.823	17.031	-14.962	1.00	1.71
ATOM	378	C	HIS	28	7.436	12.156	-18.069	1.00	0.30
ATOM	379	O	HIS	28	6.737	11.164	-18.082	1.00	0.30
ATOM	380	N	SER	29	8.362	12.338	-18.970	1.00	0.31
ATOM	381	HN	SER	29	8.912	13.149	-18.952	1.00	0.34
ATOM	382	CA	SER	29	8.567	11.319	-20.039	1.00	0.32
ATOM	383	HA	SER	29	7.660	11.217	-20.615	1.00	0.35
ATOM	384	CB	SER	29	9.699	11.775	-20.959	1.00	0.38
ATOM	385	HB1	SER	29	9.973	10.963	-21.621	1.00	0.39
ATOM	386	HB2	SER	29	10.555	12.056	-20.368	1.00	0.37
ATOM	387	OG	SER	29	9.265	12.896	-21.717	1.00	0.45
ATOM	388	HG	SER	29	9.157	12.614	-22.628	1.00	0.96
ATOM	389	C	SER	29	8.931	9.964	-19.424	1.00	0.26
ATOM	390	O	SER	29	8.479	8.930	-19.876	1.00	0.26
ATOM	391	N	GLU	30	9.747	9.954	-18.405	1.00	0.24
ATOM	392	HN	GLU	30	10.107	10.796	-18.056	1.00	0.25
ATOM	393	CA	GLU	30	10.137	8.657	-17.779	1.00	0.22
ATOM	394	HA	GLU	30	10.484	7.978	-18.542	1.00	0.25
ATOM	395	CB	GLU	30	11.260	8.899	-16.769	1.00	0.23
ATOM	396	HB1	GLU	30	11.424	8.002	-16.191	1.00	0.24
ATOM	397	HB2	GLU	30	10.980	9.707	-16.108	1.00	0.22
ATOM	398	CG	GLU	30	12.547	9.268	-17.510	1.00	0.29
ATOM	399	HG1	GLU	30	12.386	10.165	-18.086	1.00	0.67
ATOM	400	HG2	GLU	30	12.826	8.460	-18.171	1.00	0.68
ATOM	401	CD	GLU	30	13.666	9.509	-16.495	1.00	0.84
ATOM	402	OE1	GLU	30	13.436	9.266	-15.321	1.00	1.49
ATOM	403	OE2	GLU	30	14.731	9.936	-16.908	1.00	1.59
ATOM	404	C	GLU	30	8.935	8.046	-17.051	1.00	0.17
ATOM	405	O	GLU	30	8.715	6.849	-17.082	1.00	0.19
ATOM	406	N	VAL	31	8.163	8.861	-16.387	1.00	0.16
ATOM	407	HN	VAL	31	8.366	9.819	-16.371	1.00	0.17
ATOM	408	CA	VAL	31	6.983	8.341	-15.640	1.00	0.16
ATOM	409	HA	VAL	31	7.292	7.527	-14.999	1.00	0.17
ATOM	410	CB	VAL	31	6.402	9.464	-14.782	1.00	0.20
ATOM	411	HB	VAL	31	6.261	10.344	-15.392	1.00	0.22
ATOM	412	CG1	VAL	31	5.058	9.021	-14.208	1.00	0.23
ATOM	413	HG11	VAL	31	5.135	8.000	-13.867	1.00	0.97
ATOM	414	HG12	VAL	31	4.298	9.090	-14.973	1.00	1.07
ATOM	415	HG13	VAL	31	4.793	9.659	-13.378	1.00	1.07
ATOM	416	CG2	VAL	31	7.364	9.785	-13.636	1.00	0.24
ATOM	417	HG21	VAL	31	7.528	8.897	-13.045	1.00	1.05
ATOM	418	HG22	VAL	31	6.936	10.557	-13.013	1.00	1.03
ATOM	419	HG23	VAL	31	8.304	10.129	-14.040	1.00	0.99
ATOM	420	C	VAL	31	5.911	7.844	-16.617	1.00	0.16
ATOM	421	O	VAL	31	5.293	6.817	-16.406	1.00	0.17
ATOM	422	N	GLU	32	5.672	8.571	-17.677	1.00	0.18
ATOM	423	HN	GLU	32	6.172	9.401	-17.824	1.00	0.19
ATOM	424	CA	GLU	32	4.626	8.146	-18.652	1.00	0.21

10/75

ATOM	425	HA	GLU	32	3.673	8.092	-18.147	1.00	0.24
ATOM	426	CB	GLU	32	4.533	9.170	-19.787	1.00	0.27
ATOM	427	HB1	GLU	32	3.922	8.772	-20.582	1.00	0.31
ATOM	428	HB2	GLU	32	5.524	9.379	-20.164	1.00	0.28
ATOM	429	CG	GLU	32	3.904	10.463	-19.262	1.00	0.29
ATOM	430	HG1	GLU	32	4.456	10.812	-18.405	1.00	0.48
ATOM	431	HG2	GLU	32	2.879	10.272	-18.977	1.00	0.52
ATOM	432	CD	GLU	32	3.937	11.529	-20.359	1.00	0.70
ATOM	433	OE1	GLU	32	4.969	12.161	-20.513	1.00	1.37
ATOM	434	OE2	GLU	32	2.929	11.696	-21.026	1.00	1.45
ATOM	435	C	GLU	32	4.962	6.773	-19.235	1.00	0.20
ATOM	436	O	GLU	32	4.126	5.893	-19.280	1.00	0.20
ATOM	437	N	LYS	33	6.168	6.575	-19.689	1.00	0.20
ATOM	438	HN	LYS	33	6.835	7.293	-19.654	1.00	0.21
ATOM	439	CA	LYS	33	6.518	5.249	-20.269	1.00	0.21
ATOM	440	HA	LYS	33	5.825	5.029	-21.068	1.00	0.24
ATOM	441	CB	LYS	33	7.940	5.281	-20.843	1.00	0.26
ATOM	442	HB1	LYS	33	7.987	6.024	-21.624	1.00	0.31
ATOM	443	HB2	LYS	33	8.179	4.312	-21.257	1.00	0.31
ATOM	444	CG	LYS	33	8.954	5.631	-19.748	1.00	0.26
ATOM	445	HG1	LYS	33	8.823	4.970	-18.906	1.00	0.40
ATOM	446	HG2	LYS	33	8.799	6.648	-19.430	1.00	0.42
ATOM	447	CD	LYS	33	10.380	5.469	-20.291	1.00	0.48
ATOM	448	HD1	LYS	33	10.466	4.517	-20.793	1.00	0.74
ATOM	449	HD2	LYS	33	11.080	5.505	-19.469	1.00	1.11
ATOM	450	CE	LYS	33	10.705	6.593	-21.282	1.00	0.92
ATOM	451	HE1	LYS	33	10.398	7.543	-20.868	1.00	1.52
ATOM	452	HE2	LYS	33	10.184	6.419	-22.211	1.00	1.19
ATOM	453	NZ	LYS	33	12.172	6.614	-21.538	1.00	1.60
ATOM	454	HZ1	LYS	33	12.668	6.957	-20.692	1.00	1.99
ATOM	455	HZ2	LYS	33	12.374	7.247	-22.340	1.00	2.14
ATOM	456	HZ3	LYS	33	12.498	5.653	-21.763	1.00	2.03
ATOM	457	C	LYS	33	6.399	4.158	-19.202	1.00	0.19
ATOM	458	O	LYS	33	6.054	3.035	-19.495	1.00	0.20
ATOM	459	N	ALA	34	6.682	4.471	-17.966	1.00	0.17
ATOM	460	HN	ALA	34	6.965	5.383	-17.740	1.00	0.18
ATOM	461	CA	ALA	34	6.589	3.428	-16.904	1.00	0.16
ATOM	462	HA	ALA	34	7.276	2.625	-17.128	1.00	0.18
ATOM	463	CB	ALA	34	6.952	4.043	-15.551	1.00	0.16
ATOM	464	HB1	ALA	34	6.483	3.476	-14.761	1.00	1.02
ATOM	465	HB2	ALA	34	6.604	5.065	-15.516	1.00	0.98
ATOM	466	HB3	ALA	34	8.024	4.022	-15.423	1.00	1.02
ATOM	467	C	ALA	34	5.164	2.875	-16.844	1.00	0.16
ATOM	468	O	ALA	34	4.954	1.677	-16.847	1.00	0.17
ATOM	469	N	PHE	35	4.182	3.729	-16.792	1.00	0.16
ATOM	470	HN	PHE	35	4.364	4.694	-16.792	1.00	0.16
ATOM	471	CA	PHE	35	2.781	3.230	-16.736	1.00	0.17
ATOM	472	HA	PHE	35	2.690	2.525	-15.924	1.00	0.17
ATOM	473	CB	PHE	35	1.815	4.396	-16.508	1.00	0.18
ATOM	474	HB1	PHE	35	0.802	4.060	-16.672	1.00	0.19
ATOM	475	HB2	PHE	35	2.045	5.192	-17.200	1.00	0.19
ATOM	476	CG	PHE	35	1.953	4.902	-15.089	1.00	0.18
ATOM	477	CD1	PHE	35	1.616	4.071	-14.011	1.00	0.19
ATOM	478	HD1	PHE	35	1.258	3.069	-14.191	1.00	0.19
ATOM	479	CD2	PHE	35	2.415	6.203	-14.849	1.00	0.20
ATOM	480	HD2	PHE	35	2.674	6.847	-15.677	1.00	0.21
ATOM	481	CE1	PHE	35	1.743	4.539	-12.699	1.00	0.21
ATOM	482	HE1	PHE	35	1.484	3.897	-11.870	1.00	0.23
ATOM	483	CE2	PHE	35	2.540	6.670	-13.535	1.00	0.22
ATOM	484	HE2	PHE	35	2.893	7.672	-13.349	1.00	0.24
ATOM	485	CZ	PHE	35	2.205	5.838	-12.460	1.00	0.22
ATOM	486	HZ	PHE	35	2.303	6.198	-11.447	1.00	0.24
ATOM	487	C	PHE	35	2.432	2.524	-18.048	1.00	0.18
ATOM	488	O	PHE	35	1.770	1.507	-18.055	1.00	0.19
ATOM	489	N	LYS	36	2.864	3.053	-19.162	1.00	0.19
ATOM	490	HN	LYS	36	3.394	3.878	-19.144	1.00	0.19
ATOM	491	CA	LYS	36	2.535	2.399	-20.460	1.00	0.22
ATOM	492	HA	LYS	36	1.462	2.358	-20.574	1.00	0.23
ATOM	493	CB	LYS	36	3.135	3.205	-21.614	1.00	0.24
ATOM	494	HB1	LYS	36	3.045	2.641	-22.530	1.00	0.27
ATOM	495	HB2	LYS	36	4.178	3.400	-21.412	1.00	0.24
ATOM	496	CG	LYS	36	2.384	4.530	-21.758	1.00	0.27
ATOM	497	HG1	LYS	36	2.471	5.097	-20.844	1.00	0.69
ATOM	498	HG2	LYS	36	1.341	4.332	-21.963	1.00	0.68
ATOM	499	CD	LYS	36	2.988	5.332	-22.913	1.00	0.75
ATOM	500	HD1	LYS	36	2.898	4.766	-23.828	1.00	1.39
ATOM	501	HD2	LYS	36	4.032	5.525	-22.710	1.00	1.34

11/75

ATOM	502	CE	LYS	36	2.243	6.659	-23.065	1.00	1.15
ATOM	503	HE1	LYS	36	2.728	7.415	-22.464	1.00	1.64
ATOM	504	HE2	LYS	36	1.221	6.540	-22.736	1.00	1.61
ATOM	505	NZ	LYS	36	2.260	7.076	-24.496	1.00	1.99
ATOM	506	HZ1	LYS	36	2.628	6.298	-25.079	1.00	2.51
ATOM	507	HZ2	LYS	36	2.871	7.911	-24.605	1.00	2.40
ATOM	508	HZ3	LYS	36	1.295	7.309	-24.801	1.00	2.38
ATOM	509	C	LYS	36	3.098	0.976	-20.481	1.00	0.21
ATOM	510	O	LYS	36	2.446	0.053	-20.927	1.00	0.23
ATOM	511	N	LYS	37	4.295	0.778	-19.995	1.00	0.21
ATOM	512	HN	LYS	37	4.810	1.527	-19.629	1.00	0.20
ATOM	513	CA	LYS	37	4.864	-0.600	-19.988	1.00	0.22
ATOM	514	HA	LYS	37	4.926	-0.974	-21.000	1.00	0.24
ATOM	515	CB	LYS	37	6.257	-0.581	-19.358	1.00	0.22
ATOM	516	HB1	LYS	37	6.589	-1.596	-19.195	1.00	0.24
ATOM	517	HB2	LYS	37	6.216	-0.061	-18.412	1.00	0.21
ATOM	518	CG	LYS	37	7.244	0.130	-20.285	1.00	0.26
ATOM	519	HG1	LYS	37	6.921	1.140	-20.459	1.00	0.25
ATOM	520	HG2	LYS	37	7.296	-0.398	-21.227	1.00	0.28
ATOM	521	CD	LYS	37	8.625	0.139	-19.628	1.00	0.30
ATOM	522	HD1	LYS	37	8.994	-0.873	-19.551	1.00	0.77
ATOM	523	HD2	LYS	37	8.549	0.570	-18.640	1.00	0.84
ATOM	524	CE	LYS	37	9.594	0.968	-20.473	1.00	0.90
ATOM	525	HE1	LYS	37	10.530	1.076	-19.943	1.00	1.47
ATOM	526	HE2	LYS	37	9.169	1.945	-20.652	1.00	1.59
ATOM	527	NZ	LYS	37	9.836	0.286	-21.774	1.00	1.77
ATOM	528	HZ1	LYS	37	9.798	0.984	-22.543	1.00	2.22
ATOM	529	HZ2	LYS	37	9.106	-0.439	-21.926	1.00	2.28
ATOM	530	HZ3	LYS	37	10.774	-0.161	-21.762	1.00	2.33
ATOM	531	C	LYS	37	3.955	-1.506	-19.158	1.00	0.20
ATOM	532	O	LYS	37	3.689	-2.636	-19.516	1.00	0.21
ATOM	533	N	ALA	38	3.479	-1.013	-18.046	1.00	0.19
ATOM	534	HN	ALA	38	3.711	-0.098	-17.777	1.00	0.19
ATOM	535	CA	ALA	38	2.589	-1.838	-17.182	1.00	0.18
ATOM	536	HA	ALA	38	3.116	-2.727	-16.870	1.00	0.19
ATOM	537	CB	ALA	38	2.183	-1.030	-15.949	1.00	0.19
ATOM	538	HB1	ALA	38	2.831	-0.172	-15.851	1.00	1.05
ATOM	539	HB2	ALA	38	2.270	-1.649	-15.068	1.00	1.00
ATOM	540	HB3	ALA	38	1.161	-0.698	-16.057	1.00	1.06
ATOM	541	C	ALA	38	1.338	-2.238	-17.965	1.00	0.18
ATOM	542	O	ALA	38	0.967	-3.392	-18.012	1.00	0.19
ATOM	543	N	PHE	39	0.688	-1.295	-18.589	1.00	0.18
ATOM	544	HN	PHE	39	1.005	-0.368	-18.547	1.00	0.18
ATOM	545	CA	PHE	39	-0.535	-1.632	-19.367	1.00	0.19
ATOM	546	HA	PHE	39	-1.248	-2.122	-18.720	1.00	0.19
ATOM	547	CB	PHE	39	-1.156	-0.354	-19.937	1.00	0.21
ATOM	548	HB1	PHE	39	-1.883	-0.614	-20.692	1.00	0.24
ATOM	549	HB2	PHE	39	-0.381	0.256	-20.378	1.00	0.21
ATOM	550	CG	PHE	39	-1.836	0.416	-18.829	1.00	0.20
ATOM	551	CD1	PHE	39	-3.010	-0.080	-18.250	1.00	0.25
ATOM	552	HD1	PHE	39	-3.429	-1.014	-18.595	1.00	0.30
ATOM	553	CD2	PHE	39	-1.294	1.627	-18.380	1.00	0.17
ATOM	554	HD2	PHE	39	-0.389	2.012	-18.827	1.00	0.18
ATOM	555	CE1	PHE	39	-3.642	0.633	-17.224	1.00	0.28
ATOM	556	HE1	PHE	39	-4.548	0.250	-16.779	1.00	0.34
ATOM	557	CE2	PHE	39	-1.926	2.341	-17.354	1.00	0.18
ATOM	558	HE2	PHE	39	-1.507	3.275	-17.007	1.00	0.17
ATOM	559	CZ	PHE	39	-3.099	1.843	-16.776	1.00	0.23
ATOM	560	HZ	PHE	39	-3.587	2.394	-15.985	1.00	0.26
ATOM	561	C	PHE	39	-0.154	-2.571	-20.508	1.00	0.18
ATOM	562	O	PHE	39	-0.862	-3.509	-20.817	1.00	0.18
ATOM	563	N	LYS	40	0.963	-2.330	-21.136	1.00	0.19
ATOM	564	HN	LYS	40	1.522	-1.570	-20.870	1.00	0.19
ATOM	565	CA	LYS	40	1.388	-3.214	-22.254	1.00	0.19
ATOM	566	HA	LYS	40	0.642	-3.186	-23.031	1.00	0.20
ATOM	567	CB	LYS	40	2.730	-2.707	-22.804	1.00	0.21
ATOM	568	HB1	LYS	40	3.466	-2.723	-22.014	1.00	0.21
ATOM	569	HB2	LYS	40	2.610	-1.692	-23.155	1.00	0.25
ATOM	570	CG	LYS	40	3.218	-3.588	-23.966	1.00	0.25
ATOM	571	HG1	LYS	40	3.337	-4.604	-23.621	1.00	0.46
ATOM	572	HG2	LYS	40	4.171	-3.218	-24.314	1.00	0.46
ATOM	573	CD	LYS	40	2.213	-3.560	-25.121	1.00	0.38
ATOM	574	HD1	LYS	40	1.840	-2.555	-25.253	1.00	0.54
ATOM	575	HD2	LYS	40	1.392	-4.227	-24.905	1.00	0.56
ATOM	576	CE	LYS	40	2.903	-4.019	-26.407	1.00	0.40
ATOM	577	HE1	LYS	40	3.776	-4.604	-26.158	1.00	1.07
ATOM	578	HE2	LYS	40	3.188	-3.157	-26.005	1.00	1.07

12/75

ATOM	579	NZ	LYS	40	1.958	-4.852	-27.203	1.00	1.40
ATOM	580	HZ1	LYS	40	1.571	-5.607	-26.602	1.00	1.95
ATOM	581	HZ2	LYS	40	2.464	-5.274	-28.009	1.00	1.92
ATOM	582	HZ3	LYS	40	1.181	-4.258	-27.552	1.00	2.02
ATOM	583	C	LYS	40	1.553	-4.648	-21.740	1.00	0.17
ATOM	584	O	LYS	40	1.034	-5.583	-22.314	1.00	0.17
ATOM	585	N	VAL	41	2.271	-4.828	-20.663	1.00	0.17
ATOM	586	HN	VAL	41	2.681	-4.060	-20.214	1.00	0.18
ATOM	587	CA	VAL	41	2.468	-6.204	-20.116	1.00	0.16
ATOM	588	HA	VAL	41	2.953	-6.816	-20.862	1.00	0.17
ATOM	589	CB	VAL	41	3.350	-6.143	-18.868	1.00	0.18
ATOM	590	HB	VAL	41	2.966	-5.393	-18.192	1.00	0.41
ATOM	591	CG1	VAL	41	3.343	-7.508	-18.175	1.00	0.44
ATOM	592	HG11	VAL	41	2.420	-7.631	-17.629	1.00	1.16
ATOM	593	HG12	VAL	41	4.176	-7.571	-17.490	1.00	1.18
ATOM	594	HG13	VAL	41	3.429	-8.289	-18.916	1.00	1.11
ATOM	595	CG2	VAL	41	4.781	-5.785	-19.277	1.00	0.43
ATOM	596	HG21	VAL	41	5.132	-6.492	-20.013	1.00	1.12
ATOM	597	HG22	VAL	41	5.423	-5.820	-18.411	1.00	1.11
ATOM	598	HG23	VAL	41	4.797	-4.790	-19.697	1.00	1.19
ATOM	599	C	VAL	41	1.122	-6.833	-19.751	1.00	0.16
ATOM	600	O	VAL	41	0.887	-7.999	-19.996	1.00	0.17
ATOM	601	N	TRP	42	0.240	-6.080	-19.152	1.00	0.16
ATOM	602	HN	TRP	42	0.448	-5.143	-18.950	1.00	0.17
ATOM	603	CA	TRP	42	-1.079	-6.655	-18.761	1.00	0.17
ATOM	604	HA	TRP	42	-0.927	-7.642	-18.352	1.00	0.17
ATOM	605	CB	TRP	42	-1.739	-5.767	-17.699	1.00	0.18
ATOM	606	HB1	TRP	42	-2.787	-6.018	-17.621	1.00	0.19
ATOM	607	HB2	TRP	42	-1.638	-4.730	-17.983	1.00	0.20
ATOM	608	CG	TRP	42	-1.073	-5.990	-16.377	1.00	0.18
ATOM	609	CD1	TRP	42	-0.311	-5.082	-15.724	1.00	0.22
ATOM	610	HD1	TRP	42	-0.092	-4.084	-16.066	1.00	0.28
ATOM	611	CD2	TRP	42	-1.095	-7.182	-15.539	1.00	0.19
ATOM	612	NE1	TRP	42	0.140	-5.643	-14.543	1.00	0.22
ATOM	613	HE1	TRP	42	0.714	-5.194	-13.887	1.00	0.25
ATOM	614	CE2	TRP	42	-0.315	-6.935	-14.384	1.00	0.20
ATOM	615	CE3	TRP	42	-1.707	-8.441	-15.669	1.00	0.25
ATOM	616	HE3	TRP	42	-2.309	-8.658	-16.539	1.00	0.27
ATOM	617	CZ2	TRP	42	-0.149	-7.903	-13.393	1.00	0.24
ATOM	618	HZ2	TRP	42	0.454	-7.691	-12.521	1.00	0.25
ATOM	619	CZ3	TRP	42	-1.543	-9.418	-14.673	1.00	0.31
ATOM	620	HZ3	TRP	42	-2.018	-10.381	-14.782	1.00	0.39
ATOM	621	CH2	TRP	42	-0.764	-9.149	-13.538	1.00	0.30
ATOM	622	HH2	TRP	42	-0.642	-9.904	-12.775	1.00	0.35
ATOM	623	C	TRP	42	-1.991	-6.754	-19.985	1.00	0.17
ATOM	624	O	TRP	42	-2.726	-7.706	-20.138	1.00	0.18
ATOM	625	N	SER	43	-1.952	-5.782	-20.855	1.00	0.17
ATOM	626	HN	SER	43	-1.352	-5.021	-20.713	1.00	0.17
ATOM	627	CA	SER	43	-2.831	-5.825	-22.062	1.00	0.18
ATOM	628	HA	SER	43	-3.846	-6.028	-21.759	1.00	0.19
ATOM	629	CB	SER	43	-2.779	-4.474	-22.775	1.00	0.20
ATOM	630	HB1	SER	43	-2.965	-3.683	-22.059	1.00	0.21
ATOM	631	HB2	SER	43	-3.533	-4.442	-23.543	1.00	0.23
ATOM	632	OG	SER	43	-1.499	-4.304	-23.368	1.00	0.21
ATOM	633	HG	SER	43	-1.031	-5.140	-23.309	1.00	0.97
ATOM	634	C	SER	43	-2.358	-6.922	-23.019	1.00	0.18
ATOM	635	O	SER	43	-3.085	-7.350	-23.893	1.00	0.21
ATOM	636	N	ASP	44	-1.148	-7.379	-22.866	1.00	0.17
ATOM	637	HN	ASP	44	-0.575	-7.019	-22.156	1.00	0.18
ATOM	638	CA	ASP	44	-0.632	-8.445	-23.770	1.00	0.18
ATOM	639	HA	ASP	44	-0.650	-8.086	-24.788	1.00	0.19
ATOM	640	CB	ASP	44	0.809	-8.793	-23.386	1.00	0.20
ATOM	641	HB1	ASP	44	1.117	-9.683	-23.915	1.00	0.21
ATOM	642	HB2	ASP	44	0.864	-8.969	-22.322	1.00	0.22
ATOM	643	CG	ASP	44	1.734	-7.635	-23.760	1.00	0.24
ATOM	644	OD1	ASP	44	1.340	-6.833	-24.591	1.00	0.85
ATOM	645	OD2	ASP	44	2.820	-7.568	-23.209	1.00	0.84
ATOM	646	C	ASP	44	-1.499	-9.705	-23.665	1.00	0.19
ATOM	647	O	ASP	44	-1.753	-10.366	-24.653	1.00	0.21
ATOM	648	N	VAL	45	-1.927	-10.058	-22.475	1.00	0.21
ATOM	649	HN	VAL	45	-1.689	-9.519	-21.693	1.00	0.21
ATOM	650	CA	VAL	45	-2.749	-11.299	-22.302	1.00	0.26
ATOM	651	HA	VAL	45	-2.833	-11.811	-23.247	1.00	0.28
ATOM	652	CB	VAL	45	-2.045	-12.222	-21.303	1.00	0.30
ATOM	653	HB	VAL	45	-2.645	-13.107	-21.146	1.00	0.37
ATOM	654	CG1	VAL	45	-0.678	-12.626	-21.866	1.00	0.36
ATOM	655	HG11	VAL	45	-0.210	-11.766	-22.322	1.00	0.27

13/75

ATOM	656	HG12	VAL	45	-0.810	-13.400	-22.607	1.00	1.02
ATOM	657	HG13	VAL	45	-0.051	-12.995	-21.068	1.00	1.13
ATOM	658	CG2	VAL	45	-1.855	-11.486	-19.973	1.00	0.32
ATOM	659	HG21	VAL	45	-2.819	-11.303	-19.524	1.00	0.96
ATOM	660	HG22	VAL	45	-1.356	-10.545	-20.149	1.00	1.09
ATOM	661	HG23	VAL	45	-1.258	-12.091	-19.305	1.00	1.11
ATOM	662	C	VAL	45	-4.160	-10.966	-21.790	1.00	0.29
ATOM	663	O	VAL	45	-4.837	-11.819	-21.249	1.00	0.64
ATOM	664	N	THR	46	-4.619	-9.748	-21.963	1.00	0.36
ATOM	665	HN	THR	46	-4.062	-9.076	-22.409	1.00	0.65
ATOM	666	CA	THR	46	-5.998	-9.382	-21.491	1.00	0.38
ATOM	667	HA	THR	46	-6.567	-10.277	-21.320	1.00	0.44
ATOM	668	CB	THR	46	-5.912	-8.577	-20.186	1.00	0.39
ATOM	669	HB	THR	46	-6.889	-8.193	-19.943	1.00	0.46
ATOM	670	OG1	THR	46	-5.018	-7.491	-20.358	1.00	0.36
ATOM	671	HG1	THR	46	-5.532	-6.719	-20.608	1.00	0.94
ATOM	672	CG2	THR	46	-5.430	-9.461	-19.036	1.00	0.43
ATOM	673	HG21	THR	46	-4.929	-10.327	-19.429	1.00	1.08
ATOM	674	HG22	THR	46	-6.277	-9.775	-18.445	1.00	1.15
ATOM	675	HG23	THR	46	-4.746	-8.901	-18.415	1.00	1.05
ATOM	676	C	THR	46	-6.668	-8.482	-22.553	1.00	0.32
ATOM	677	O	THR	46	-6.124	-7.450	-22.892	1.00	0.32
ATOM	678	N	PRO	47	-7.833	-8.829	-23.084	1.00	0.30
ATOM	679	CA	PRO	47	-8.479	-7.951	-24.100	1.00	0.30
ATOM	680	HA	PRO	47	-7.820	-7.790	-24.936	1.00	0.33
ATOM	681	CB	PRO	47	-9.687	-8.773	-24.546	1.00	0.35
ATOM	682	HB1	PRO	47	-9.541	-9.110	-25.561	1.00	0.40
ATOM	683	HB2	PRO	47	-10.579	-8.166	-24.489	1.00	0.37
ATOM	684	CG	PRO	47	-9.825	-9.986	-23.621	1.00	0.35
ATOM	685	HG1	PRO	47	-9.916	-10.885	-24.212	1.00	0.42
ATOM	686	HG2	PRO	47	-10.703	-9.869	-23.001	1.00	0.34
ATOM	687	CD	PRO	47	-8.576	-10.077	-22.739	1.00	0.33
ATOM	688	HD2	PRO	47	-8.853	-10.091	-21.692	1.00	0.31
ATOM	689	HD1	PRO	47	-7.993	-10.946	-22.999	1.00	0.39
ATOM	690	C	PRO	47	-8.933	-6.614	-23.506	1.00	0.25
ATOM	691	O	PRO	47	-9.744	-5.914	-24.080	1.00	0.26
ATOM	692	N	LEU	48	-8.418	-6.252	-22.362	1.00	0.26
ATOM	693	HN	LEU	48	-7.766	-6.828	-21.912	1.00	0.29
ATOM	694	CA	LEU	48	-8.827	-4.960	-21.742	1.00	0.26
ATOM	695	HA	LEU	48	-9.904	-4.905	-21.696	1.00	0.27
ATOM	696	CB	LEU	48	-8.241	-4.858	-20.329	1.00	0.31
ATOM	697	HB1	LEU	48	-8.476	-3.892	-19.909	1.00	0.34
ATOM	698	HB2	LEU	48	-7.167	-4.968	-20.385	1.00	0.33
ATOM	699	CG	LEU	48	-8.816	-5.964	-19.434	1.00	0.34
ATOM	700	HG	LEU	48	-8.808	-6.900	-19.972	1.00	0.32
ATOM	701	CD1	LEU	48	-7.952	-6.091	-18.177	1.00	0.41
ATOM	702	HD11	LEU	48	-8.002	-5.171	-17.613	1.00	1.11
ATOM	703	HD12	LEU	48	-6.928	-6.283	-18.462	1.00	1.05
ATOM	704	HD13	LEU	48	-8.315	-6.906	-17.570	1.00	1.15
ATOM	705	CD2	LEU	48	-10.255	-5.628	-19.016	1.00	0.36
ATOM	706	HD21	LEU	48	-10.569	-4.707	-19.478	1.00	1.10
ATOM	707	HD22	LEU	48	-10.299	-5.524	-17.942	1.00	1.09
ATOM	708	HD23	LEU	48	-10.912	-6.428	-19.325	1.00	1.04
ATOM	709	C	LEU	48	-8.289	-3.806	-22.589	1.00	0.25
ATOM	710	O	LEU	48	-7.174	-3.849	-23.071	1.00	0.26
ATOM	711	N	ASN	49	-9.073	-2.775	-22.762	1.00	0.25
ATOM	712	HN	ASN	49	-9.964	-2.770	-22.355	1.00	0.26
ATOM	713	CA	ASN	49	-8.622	-1.604	-23.568	1.00	0.25
ATOM	714	HA	ASN	49	-7.703	-1.842	-24.082	1.00	0.27
ATOM	715	CB	ASN	49	-9.700	-1.245	-24.593	1.00	0.28
ATOM	716	HB1	ASN	49	-9.390	-0.375	-25.153	1.00	0.30
ATOM	717	HB2	ASN	49	-10.628	-1.033	-24.081	1.00	0.28
ATOM	718	CG	ASN	49	-9.902	-2.419	-25.553	1.00	0.32
ATOM	719	OD1	ASN	49	-9.798	-3.564	-25.161	1.00	1.10
ATOM	720	ND2	ASN	49	-10.186	-2.182	-26.804	1.00	1.14
ATOM	721	HD21	ASN	49	-10.268	-1.258	-27.121	1.00	1.94
ATOM	722	HD22	ASN	49	-10.317	-2.927	-27.427	1.00	1.14
ATOM	723	C	ASN	49	-8.391	-0.417	-22.633	1.00	0.24
ATOM	724	O	ASN	49	-9.290	0.016	-21.939	1.00	0.23
ATOM	725	N	PHE	50	-7.192	0.107	-22.606	1.00	0.24
ATOM	726	HN	PHE	50	-6.485	-0.264	-23.173	1.00	0.26
ATOM	727	CA	PHE	50	-6.896	1.263	-21.710	1.00	0.23
ATOM	728	HA	PHE	50	-7.688	1.380	-20.985	1.00	0.21
ATOM	729	CB	PHE	50	-5.574	1.016	-20.981	1.00	0.24
ATOM	730	HB1	PHE	50	-5.357	1.853	-20.334	1.00	0.25
ATOM	731	HB2	PHE	50	-4.780	0.907	-21.705	1.00	0.27
ATOM	732	CG	PHE	50	-5.676	-0.243	-20.154	1.00	0.22

14/75

ATOM	733	CD1	PHE	50	-6.266	-0.201	-18.886	1.00	0.25
ATOM	734	HD1	PHE	50	-6.652	0.731	-18.500	1.00	0.28
ATOM	735	CD2	PHE	50	-5.176	-1.451	-20.654	1.00	0.22
ATOM	736	HD2	PHE	50	-4.720	-1.483	-21.633	1.00	0.23
ATOM	737	CE1	PHE	50	-6.358	-1.368	-18.117	1.00	0.25
ATOM	738	HE1	PHE	50	-6.813	-1.336	-17.139	1.00	0.28
ATOM	739	CE2	PHE	50	-5.267	-2.618	-19.886	1.00	0.23
ATOM	740	HE2	PHE	50	-4.881	-3.550	-20.272	1.00	0.25
ATOM	741	CZ	PHE	50	-5.858	-2.576	-18.618	1.00	0.24
ATOM	742	HZ	PHE	50	-5.928	-3.476	-18.025	1.00	0.25
ATOM	743	C	PHE	50	-6.777	2.538	-22.545	1.00	0.26
ATOM	744	O	PHE	50	-6.028	2.596	-23.501	1.00	0.31
ATOM	745	N	THR	51	-7.517	3.555	-22.184	1.00	0.24
ATOM	746	HN	THR	51	-8.109	3.468	-21.413	1.00	0.22
ATOM	747	CA	THR	51	-7.470	4.842	-22.940	1.00	0.27
ATOM	748	HA	THR	51	-6.775	4.762	-23.762	1.00	0.31
ATOM	749	CB	THR	51	-8.868	5.153	-23.483	1.00	0.30
ATOM	750	HB	THR	51	-9.562	5.248	-22.663	1.00	0.29
ATOM	751	OG1	THR	51	-9.283	4.100	-24.341	1.00	0.35
ATOM	752	HG1	THR	51	-9.638	4.491	-25.142	1.00	0.84
ATOM	753	CG2	THR	51	-8.835	6.464	-24.273	1.00	0.34
ATOM	754	HG21	THR	51	-9.805	6.640	-24.716	1.00	1.02
ATOM	755	HG22	THR	51	-8.092	6.394	-25.053	1.00	1.07
ATOM	756	HG23	THR	51	-8.588	7.280	-23.611	1.00	1.13
ATOM	757	C	THR	51	-7.024	5.969	-22.001	1.00	0.25
ATOM	758	O	THR	51	-7.553	6.139	-20.920	1.00	0.22
ATOM	759	N	ARG	52	-6.054	6.740	-22.411	1.00	0.29
ATOM	760	HN	ARG	52	-5.645	6.583	-23.287	1.00	0.32
ATOM	761	CA	ARG	52	-5.566	7.861	-21.556	1.00	0.29
ATOM	762	HA	ARG	52	-5.591	7.563	-20.518	1.00	0.27
ATOM	763	CB	ARG	52	-4.128	8.201	-21.955	1.00	0.35
ATOM	764	HB1	ARG	52	-4.125	8.654	-22.935	1.00	0.39
ATOM	765	HB2	ARG	52	-3.539	7.295	-21.977	1.00	0.38
ATOM	766	CG	ARG	52	-3.521	9.177	-20.945	1.00	0.39
ATOM	767	HG1	ARG	52	-3.645	8.787	-19.946	1.00	0.71
ATOM	768	HG2	ARG	52	-4.017	10.134	-21.025	1.00	0.57
ATOM	769	CD	ARG	52	-2.030	9.345	-21.244	1.00	0.79
ATOM	770	HD1	ARG	52	-1.825	9.001	-22.248	1.00	1.45
ATOM	771	HD2	ARG	52	-1.453	8.763	-20.543	1.00	1.39
ATOM	772	NE	ARG	52	-1.656	10.782	-21.120	1.00	1.47
ATOM	773	HE	ARG	52	-2.354	11.468	-21.073	1.00	2.06
ATOM	774	CZ	ARG	52	-0.398	11.127	-21.071	1.00	2.09
ATOM	775	NH1	ARG	52	-0.070	12.385	-20.960	1.00	3.05
ATOM	776	HH11	ARG	52	-0.782	13.084	-20.911	1.00	3.45
ATOM	777	HH12	ARG	52	0.894	12.649	-20.923	1.00	3.60
ATOM	778	NH2	ARG	52	0.532	10.213	-21.138	1.00	2.31
ATOM	779	HH21	ARG	52	0.281	9.249	-21.226	1.00	2.16
ATOM	780	HH22	ARG	52	1.496	10.477	-21.102	1.00	3.05
ATOM	781	C	ARG	52	-6.460	9.090	-21.758	1.00	0.29
ATOM	782	O	ARG	52	-6.719	9.495	-22.875	1.00	0.33
ATOM	783	N	LEU	53	-6.928	9.689	-20.689	1.00	0.26
ATOM	784	HN	LEU	53	-6.702	9.345	-19.798	1.00	0.25
ATOM	785	CA	LEU	53	-7.803	10.896	-20.822	1.00	0.29
ATOM	786	HA	LEU	53	-8.167	10.972	-21.835	1.00	0.32
ATOM	787	CB	LEU	53	-8.992	10.784	-19.862	1.00	0.28
ATOM	788	HB1	LEU	53	-9.579	11.688	-19.908	1.00	0.31
ATOM	789	HB2	LEU	53	-8.624	10.648	-18.855	1.00	0.28
ATOM	790	CG	LEU	53	-9.866	9.587	-20.249	1.00	0.28
ATOM	791	HG	LEU	53	-9.264	8.690	-20.246	1.00	0.29
ATOM	792	CD1	LEU	53	-10.999	9.440	-19.232	1.00	0.29
ATOM	793	HD11	LEU	53	-11.606	8.585	-19.487	1.00	0.95
ATOM	794	HD12	LEU	53	-11.610	10.331	-19.243	1.00	1.05
ATOM	795	HD13	LEU	53	-10.581	9.303	-18.247	1.00	1.07
ATOM	796	CD2	LEU	53	-10.463	9.799	-21.646	1.00	0.36
ATOM	797	HD21	LEU	53	-10.523	10.856	-21.860	1.00	1.01
ATOM	798	HD22	LEU	53	-11.453	9.370	-21.685	1.00	1.09
ATOM	799	HD23	LEU	53	-9.835	9.319	-22.382	1.00	1.14
ATOM	800	C	LEU	53	-7.000	12.154	-20.483	1.00	0.33
ATOM	801	O	LEU	53	-6.315	12.218	-19.482	1.00	0.34
ATOM	802	N	HIS	54	-7.080	13.154	-21.319	1.00	0.41
ATOM	803	HN	HIS	54	-7.637	13.075	-22.121	1.00	0.45
ATOM	804	CA	HIS	54	-6.324	14.413	-21.062	1.00	0.47
ATOM	805	HA	HIS	54	-5.292	14.183	-20.851	1.00	0.54
ATOM	806	CB	HIS	54	-6.407	15.314	-22.297	1.00	0.60
ATOM	807	HB1	HIS	54	-6.018	16.291	-22.052	1.00	0.64
ATOM	808	HB2	HIS	54	-7.438	15.407	-22.603	1.00	0.61
ATOM	809	CG	HIS	54	-5.602	14.726	-23.426	1.00	0.74

15/75

ATOM	810	ND1	HIS	54	-5.645	15.254	-24.707	1.00	1.35
ATOM	811	HD1	HIS	54	-6.172	16.028	-24.996	1.00	1.86
ATOM	812	CD2	HIS	54	-4.740	13.656	-23.493	1.00	0.86
ATOM	813	HD2	HIS	54	-4.480	13.010	-22.668	1.00	1.34
ATOM	814	CE1	HIS	54	-4.834	14.512	-25.481	1.00	1.33
ATOM	815	HE1	HIS	54	-4.670	14.692	-26.533	1.00	1.83
ATOM	816	NE2	HIS	54	-4.257	13.525	-24.792	1.00	0.92
ATOM	817	C	HIS	54	-6.933	15.154	-19.867	1.00	0.43
ATOM	818	O	HIS	54	-6.230	15.714	-19.051	1.00	0.49
ATOM	819	N	ASP	55	-8.236	15.172	-19.767	1.00	0.42
ATOM	820	HN	ASP	55	-8.784	14.719	-20.442	1.00	0.45
ATOM	821	CA	ASP	55	-8.892	15.892	-18.635	1.00	0.49
ATOM	822	HA	ASP	55	-8.217	15.938	-17.796	1.00	0.54
ATOM	823	CB	ASP	55	-9.251	17.314	-19.073	1.00	0.65
ATOM	824	HB1	ASP	55	-9.876	17.774	-18.323	1.00	0.75
ATOM	825	HB2	ASP	55	-9.783	17.277	-20.013	1.00	0.68
ATOM	826	CG	ASP	55	-7.974	18.140	-19.244	1.00	0.71
ATOM	827	OD1	ASP	55	-7.978	19.037	-20.071	1.00	1.19
ATOM	828	OD2	ASP	55	-7.018	17.870	-18.536	1.00	1.28
ATOM	829	C	ASP	55	-10.167	15.156	-18.223	1.00	0.45
ATOM	830	O	ASP	55	-10.638	14.273	-18.912	1.00	0.44
ATOM	831	N	GLY	56	-10.728	15.518	-17.100	1.00	0.46
ATOM	832	HN	GLY	56	-10.328	16.233	-16.563	1.00	0.50
ATOM	833	CA	GLY	56	-11.975	14.848	-16.632	1.00	0.44
ATOM	834	HA1	GLY	56	-12.482	14.399	-17.472	1.00	0.44
ATOM	835	HA2	GLY	56	-12.622	15.579	-16.169	1.00	0.48
ATOM	836	C	GLY	56	-11.624	13.760	-15.614	1.00	0.40
ATOM	837	O	GLY	56	-10.473	13.543	-15.294	1.00	0.42
ATOM	838	N	ILE	57	-12.613	13.078	-15.105	1.00	0.37
ATOM	839	HN	ILE	57	-13.533	13.275	-15.380	1.00	0.39
ATOM	840	CA	ILE	57	-12.352	12.002	-14.106	1.00	0.35
ATOM	841	HA	ILE	57	-11.406	12.184	-13.616	1.00	0.38
ATOM	842	CB	ILE	57	-13.473	12.000	-13.064	1.00	0.41
ATOM	843	HB	ILE	57	-14.415	11.820	-13.561	1.00	0.42
ATOM	844	CG1	ILE	57	-13.508	13.363	-12.360	1.00	0.48
ATOM	845	HG11	ILE	57	-13.512	14.148	-13.101	1.00	0.48
ATOM	846	HG12	ILE	57	-12.631	13.465	-11.737	1.00	0.51
ATOM	847	CG2	ILE	57	-13.216	10.896	-12.037	1.00	0.44
ATOM	848	HG21	ILE	57	-13.315	9.932	-12.513	1.00	1.19
ATOM	849	HG22	ILE	57	-13.934	10.977	-11.235	1.00	1.09
ATOM	850	HG23	ILE	57	-12.218	11.000	-11.639	1.00	1.04
ATOM	851	CD1	ILE	57	-14.765	13.484	-11.488	1.00	0.56
ATOM	852	HD11	ILE	57	-15.459	12.693	-11.728	1.00	1.08
ATOM	853	HD12	ILE	57	-15.235	14.439	-11.668	1.00	1.24
ATOM	854	HD13	ILE	57	-14.487	13.413	-10.447	1.00	1.14
ATOM	855	C	ILE	57	-12.307	10.647	-14.817	1.00	0.30
ATOM	856	O	ILE	57	-13.139	10.353	-15.653	1.00	0.31
ATOM	857	N	ALA	58	-11.337	9.828	-14.493	1.00	0.26
ATOM	858	HN	ALA	58	-10.679	10.096	-13.817	1.00	0.27
ATOM	859	CA	ALA	58	-11.221	8.489	-15.148	1.00	0.23
ATOM	860	HA	ALA	58	-11.957	8.398	-15.932	1.00	0.25
ATOM	861	CB	ALA	58	-9.824	8.339	-15.749	1.00	0.23
ATOM	862	HB1	ALA	58	-9.843	7.585	-16.522	1.00	0.97
ATOM	863	HB2	ALA	58	-9.129	8.044	-14.976	1.00	1.11
ATOM	864	HB3	ALA	58	-9.513	9.280	-16.172	1.00	1.03
ATOM	865	C	ALA	58	-11.443	7.387	-14.114	1.00	0.23
ATOM	866	O	ALA	58	-11.389	7.617	-12.922	1.00	0.27
ATOM	867	N	ASP	59	-11.701	6.189	-14.564	1.00	0.25
ATOM	868	HN	ASP	59	-11.744	6.028	-15.530	1.00	0.28
ATOM	869	CA	ASP	59	-11.934	5.069	-13.613	1.00	0.27
ATOM	870	HA	ASP	59	-12.788	5.296	-12.991	1.00	0.34
ATOM	871	CB	ASP	59	-12.207	3.785	-14.400	1.00	0.33
ATOM	872	HB1	ASP	59	-12.203	2.942	-13.725	1.00	0.34
ATOM	873	HB2	ASP	59	-11.438	3.651	-15.147	1.00	0.32
ATOM	874	CG	ASP	59	-13.572	3.880	-15.084	1.00	0.44
ATOM	875	OD1	ASP	59	-13.791	3.139	-16.028	1.00	1.20
ATOM	876	OD2	ASP	59	-14.374	4.691	-14.653	1.00	1.14
ATOM	877	C	ASP	59	-10.700	4.863	-12.731	1.00	0.22
ATOM	878	O	ASP	59	-10.806	4.767	-11.524	1.00	0.27
ATOM	879	N	ILE	60	-9.534	4.780	-13.326	1.00	0.18
ATOM	880	HN	ILE	60	-9.478	4.850	-14.302	1.00	0.20
ATOM	881	CA	ILE	60	-8.291	4.561	-12.523	1.00	0.22
ATOM	882	HA	ILE	60	-8.554	4.303	-11.512	1.00	0.28
ATOM	883	CB	ILE	60	-7.502	3.404	-13.155	1.00	0.27
ATOM	884	HB	ILE	60	-7.255	3.655	-14.175	1.00	0.28
ATOM	885	CG1	ILE	60	-8.377	2.146	-13.136	1.00	0.30

16/75

ATOM	887	HG12	ILE	60	-8.541	1.839	-12.113	1.00	0.36
ATOM	888	CG2	ILE	60	-6.210	3.127	-12.369	1.00	0.39
ATOM	889	HG21	ILE	60	-6.456	2.704	-11.409	1.00	1.05
ATOM	890	HG22	ILE	60	-5.658	4.043	-12.228	1.00	1.10
ATOM	891	HG23	ILE	60	-5.600	2.428	-12.921	1.00	1.12
ATOM	892	CD1	ILE	60	-7.688	1.015	-13.904	1.00	0.38
ATOM	893	HD11	ILE	60	-7.209	1.413	-14.786	1.00	1.07
ATOM	894	HD12	ILE	60	-8.424	0.280	-14.196	1.00	1.14
ATOM	895	HD13	ILE	60	-6.948	0.549	-13.270	1.00	1.04
ATOM	896	C	ILE	60	-7.438	5.834	-12.518	1.00	0.20
ATOM	897	O	ILE	60	-6.731	6.115	-13.464	1.00	0.25
ATOM	898	N	MET	61	-7.473	6.585	-11.448	1.00	0.20
ATOM	899	HN	MET	61	-8.033	6.326	-10.687	1.00	0.25
ATOM	900	CA	MET	61	-6.641	7.822	-11.373	1.00	0.20
ATOM	901	HA	MET	61	-6.327	8.102	-12.366	1.00	0.19
ATOM	902	CB	MET	61	-7.464	8.963	-10.773	1.00	0.24
ATOM	903	HB1	MET	61	-8.331	9.137	-11.392	1.00	0.35
ATOM	904	HB2	MET	61	-6.860	9.856	-10.743	1.00	0.33
ATOM	905	CG	MET	61	-7.918	8.604	-9.358	1.00	0.31
ATOM	906	HG1	MET	61	-7.146	8.870	-8.653	1.00	0.66
ATOM	907	HG2	MET	61	-8.112	7.544	-9.300	1.00	0.67
ATOM	908	SD	MET	61	-9.433	9.519	-8.967	1.00	0.54
ATOM	909	CE	MET	61	-8.878	11.154	-9.516	1.00	0.40
ATOM	910	HE1	MET	61	-9.492	11.914	-9.056	1.00	1.06
ATOM	911	HE2	MET	61	-8.968	11.227	-10.589	1.00	1.16
ATOM	912	HE3	MET	61	-7.846	11.298	-9.232	1.00	1.12
ATOM	913	C	MET	61	-5.396	7.540	-10.524	1.00	0.20
ATOM	914	O	MET	61	-5.478	6.951	-9.463	1.00	0.22
ATOM	915	N	ILE	62	-4.241	7.937	-11.001	1.00	0.20
ATOM	916	HN	ILE	62	-4.207	8.393	-11.868	1.00	0.21
ATOM	917	CA	ILE	62	-2.971	7.678	-10.252	1.00	0.21
ATOM	918	HA	ILE	62	-3.156	6.982	-9.448	1.00	0.20
ATOM	919	CB	ILE	62	-1.938	7.080	-11.211	1.00	0.24
ATOM	920	HB	ILE	62	-1.753	7.781	-12.012	1.00	0.26
ATOM	921	CG1	ILE	62	-2.480	5.762	-11.785	1.00	0.23
ATOM	922	HG11	ILE	62	-3.479	5.922	-12.162	1.00	0.20
ATOM	923	HG12	ILE	62	-2.508	5.018	-11.003	1.00	0.24
ATOM	924	CG2	ILE	62	-0.635	6.812	-10.455	1.00	0.30
ATOM	925	HG21	ILE	62	-0.863	6.443	-9.466	1.00	1.08
ATOM	926	HG22	ILE	62	-0.070	7.729	-10.375	1.00	1.12
ATOM	927	HG23	ILE	62	-0.052	6.076	-10.988	1.00	0.99
ATOM	928	CD1	ILE	62	-1.584	5.262	-12.927	1.00	0.29
ATOM	929	HD11	ILE	62	-0.979	6.073	-13.305	1.00	1.02
ATOM	930	HD12	ILE	62	-2.201	4.876	-13.724	1.00	1.09
ATOM	931	HD13	ILE	62	-0.941	4.476	-12.559	1.00	1.07
ATOM	932	C	ILE	62	-2.423	8.988	-9.677	1.00	0.22
ATOM	933	O	ILE	62	-2.393	10.004	-10.343	1.00	0.27
ATOM	934	N	SER	63	-1.993	8.976	-8.441	1.00	0.20
ATOM	935	HN	SER	63	-2.028	8.147	-7.916	1.00	0.18
ATOM	936	CA	SER	63	-1.452	10.226	-7.829	1.00	0.22
ATOM	937	HA	SER	63	-0.998	10.836	-8.597	1.00	0.26
ATOM	938	CB	SER	63	-2.597	11.000	-7.176	1.00	0.24
ATOM	939	HB1	SER	63	-3.448	11.012	-7.845	1.00	0.25
ATOM	940	HB2	SER	63	-2.286	12.012	-6.978	1.00	0.29
ATOM	941	OG	SER	63	-2.951	10.369	-5.952	1.00	0.25
ATOM	942	HG	SER	63	-3.682	9.772	-6.127	1.00	0.85
ATOM	943	C	SER	63	-0.404	9.879	-6.764	1.00	0.21
ATOM	944	O	SER	63	-0.364	8.775	-6.259	1.00	0.20
ATOM	945	N	PHE	64	0.440	10.823	-6.419	1.00	0.24
ATOM	946	HN	PHE	64	0.380	11.705	-6.841	1.00	0.27
ATOM	947	CA	PHE	64	1.490	10.569	-5.382	1.00	0.24
ATOM	948	HA	PHE	64	1.560	9.511	-5.179	1.00	0.22
ATOM	949	CB	PHE	64	2.840	11.084	-5.895	1.00	0.28
ATOM	950	HB1	PHE	64	3.564	11.047	-5.097	1.00	0.32
ATOM	951	HB2	PHE	64	2.730	12.103	-6.235	1.00	0.32
ATOM	952	CG	PHE	64	3.316	10.220	-7.040	1.00	0.28
ATOM	953	CD1	PHE	64	4.112	9.096	-6.788	1.00	0.30
ATOM	954	HD1	PHE	64	4.385	8.844	-5.774	1.00	0.32
ATOM	955	CD2	PHE	64	2.963	10.545	-8.355	1.00	0.33
ATOM	956	HD2	PHE	64	2.350	11.412	-8.550	1.00	0.37
ATOM	957	CE1	PHE	64	4.553	8.297	-7.850	1.00	0.36
ATOM	958	HE1	PHE	64	5.166	7.430	-7.656	1.00	0.40
ATOM	959	CE2	PHE	64	3.403	9.747	-9.417	1.00	0.40
ATOM	960	HE2	PHE	64	3.130	9.998	-10.431	1.00	0.47
ATOM	961	CZ	PHE	64	4.198	8.623	-9.165	1.00	0.40
ATOM	962	HZ	PHE	64	4.538	8.007	-9.984	1.00	0.47
ATOM	963	C	PHE	64	1.115	11.318	-4.097	1.00	0.27

17/75

ATOM	964	O	PHE	64	0.924	12.518	-4.108	1.00	0.36
ATOM	965	N	GLY	65	0.996	10.617	-2.996	1.00	0.30
ATOM	966	HN	GLY	65	1.146	9.649	-3.017	1.00	0.33
ATOM	967	CA	GLY	65	0.615	11.282	-1.709	1.00	0.38
ATOM	968	HA1	GLY	65	-0.152	10.697	-1.224	1.00	0.46
ATOM	969	HA2	GLY	65	0.230	12.270	-1.913	1.00	0.45
ATOM	970	C	GLY	65	1.823	11.397	-0.770	1.00	0.32
ATOM	971	O	GLY	65	2.926	11.007	-1.098	1.00	0.40
ATOM	972	N	ILE	66	1.598	11.926	0.408	1.00	0.30
ATOM	973	HN	ILE	66	0.691	12.220	0.635	1.00	0.36
ATOM	974	CA	ILE	66	2.691	12.081	1.417	1.00	0.36
ATOM	975	HA	ILE	66	3.564	11.534	1.093	1.00	0.40
ATOM	976	CB	ILE	66	3.040	13.564	1.571	1.00	0.41
ATOM	977	HB	ILE	66	2.127	14.134	1.656	1.00	0.64
ATOM	978	CG1	ILE	66	3.829	14.026	0.337	1.00	0.68
ATOM	979	HG11	ILE	66	3.301	13.729	-0.557	1.00	0.95
ATOM	980	HG12	ILE	66	4.804	13.561	0.346	1.00	1.01
ATOM	981	CG2	ILE	66	3.886	13.764	2.831	1.00	0.93
ATOM	982	HG21	ILE	66	4.372	14.727	2.790	1.00	1.50
ATOM	983	HG22	ILE	66	4.632	12.986	2.891	1.00	1.41
ATOM	984	HG23	ILE	66	3.249	13.720	3.702	1.00	1.54
ATOM	985	CD1	ILE	66	3.997	15.551	0.343	1.00	0.70
ATOM	986	HD11	ILE	66	4.944	15.806	0.797	1.00	1.22
ATOM	987	HD12	ILE	66	3.196	16.009	0.902	1.00	1.28
ATOM	988	HD13	ILE	66	3.979	15.917	-0.673	1.00	1.23
ATOM	989	C	ILE	66	2.207	11.519	2.760	1.00	0.46
ATOM	990	O	ILE	66	1.021	11.363	2.958	1.00	0.54
ATOM	991	N	LYS	67	3.129	11.205	3.659	1.00	0.59
ATOM	992	HN	LYS	67	4.073	11.343	3.434	1.00	0.64
ATOM	993	CA	LYS	67	2.780	10.630	5.014	1.00	0.74
ATOM	994	HA	LYS	67	3.072	9.594	5.038	1.00	0.83
ATOM	995	CB	LYS	67	3.550	11.404	6.102	1.00	0.90
ATOM	996	HB1	LYS	67	3.237	12.438	6.089	1.00	0.89
ATOM	997	HB2	LYS	67	4.608	11.352	5.891	1.00	0.96
ATOM	998	CG	LYS	67	3.287	10.815	7.504	1.00	1.08
ATOM	999	HG1	LYS	67	2.254	10.524	7.598	1.00	1.31
ATOM	1000	HG2	LYS	67	3.510	11.565	8.249	1.00	1.33
ATOM	1001	CD	LYS	67	4.179	9.590	7.746	1.00	0.98
ATOM	1002	HD1	LYS	67	5.216	9.885	7.694	1.00	1.07
ATOM	1003	HD2	LYS	67	3.979	8.839	6.999	1.00	1.07
ATOM	1004	CE	LYS	67	3.885	9.016	9.135	1.00	1.17
ATOM	1005	HE1	LYS	67	4.331	8.036	9.220	1.00	1.64
ATOM	1006	HE2	LYS	67	2.817	8.938	9.272	1.00	1.50
ATOM	1007	NZ	LYS	67	4.453	9.913	10.180	1.00	1.93
ATOM	1008	HZ1	LYS	67	4.569	10.870	9.792	1.00	2.38
ATOM	1009	HZ2	LYS	67	5.378	9.547	10.485	1.00	2.43
ATOM	1010	HZ3	LYS	67	3.808	9.948	10.995	1.00	2.40
ATOM	1011	C	LYS	67	1.274	10.732	5.280	1.00	0.72
ATOM	1012	O	LYS	67	0.530	9.804	5.035	1.00	0.79
ATOM	1013	N	GLU	68	0.815	11.855	5.760	1.00	0.77
ATOM	1014	HN	GLU	68	1.425	12.601	5.939	1.00	0.84
ATOM	1015	CA	GLU	68	-0.645	12.004	6.011	1.00	0.84
ATOM	1016	HA	GLU	68	-1.014	11.130	6.530	1.00	0.99
ATOM	1017	CB	GLU	68	-0.895	13.254	6.860	1.00	1.05
ATOM	1018	HB1	GLU	68	-0.393	13.149	7.810	1.00	1.23
ATOM	1019	HB2	GLU	68	-1.956	13.370	7.024	1.00	1.10
ATOM	1020	CG	GLU	68	-0.353	14.487	6.134	1.00	1.15
ATOM	1021	HG1	GLU	68	-1.000	14.730	5.304	1.00	1.32
ATOM	1022	HG2	GLU	68	0.642	14.281	5.768	1.00	1.28
ATOM	1023	CD	GLU	68	-0.308	15.669	7.104	1.00	1.75
ATOM	1024	OE1	GLU	68	0.246	16.692	6.736	1.00	2.45
ATOM	1025	OE2	GLU	68	-0.823	15.530	8.202	1.00	2.16
ATOM	1026	C	GLU	68	-1.346	12.132	4.660	1.00	0.76
ATOM	1027	O	GLU	68	-0.899	12.859	3.795	1.00	1.11
ATOM	1028	N	HIS	69	-2.420	11.414	4.454	1.00	0.94
ATOM	1029	HN	HIS	69	-2.755	10.815	5.155	1.00	1.32
ATOM	1030	CA	HIS	69	-3.114	11.487	3.136	1.00	1.04
ATOM	1031	HA	HIS	69	-2.877	12.437	2.679	1.00	1.25
ATOM	1032	CB	HIS	69	-2.545	10.358	2.243	1.00	1.49
ATOM	1033	HB1	HIS	69	-1.750	9.862	2.783	1.00	2.12
ATOM	1034	HB2	HIS	69	-2.131	10.798	1.351	1.00	2.27
ATOM	1035	CG	HIS	69	-3.570	9.333	1.837	1.00	0.95
ATOM	1036	ND1	HIS	69	-3.818	8.195	2.588	1.00	1.43
ATOM	1037	HD1	HIS	69	-3.415	7.972	3.453	1.00	1.83
ATOM	1038	CD2	HIS	69	-4.355	9.223	0.717	1.00	1.04
ATOM	1039	HD2	HIS	69	-4.403	9.946	-0.082	1.00	1.41
ATOM	1040	CE1	HIS	69	-4.715	7.452	1.912	1.00	1.81

18/75

ATOM	1041	HE1	HIS	69	-5.097	6.502	2.257	1.00	2.54
ATOM	1042	NE2	HIS	69	-5.075	8.032	0.765	1.00	1.53
ATOM	1043	C	HIS	69	-4.643	11.435	3.341	1.00	1.14
ATOM	1044	O	HIS	69	-5.392	10.889	2.556	1.00	1.76
ATOM	1045	N	GLY	70	-5.108	12.065	4.393	1.00	1.49
ATOM	1046	HN	GLY	70	-4.487	12.532	4.990	1.00	1.98
ATOM	1047	CA	GLY	70	-6.576	12.123	4.665	1.00	1.86
ATOM	1048	HA1	GLY	70	-7.071	12.633	3.852	1.00	2.28
ATOM	1049	HA2	GLY	70	-6.746	12.667	5.583	1.00	2.09
ATOM	1050	C	GLY	70	-7.155	10.716	4.801	1.00	1.81
ATOM	1051	O	GLY	70	-8.182	10.404	4.232	1.00	2.53
ATOM	1052	N	ASP	71	-6.513	9.863	5.545	1.00	1.55
ATOM	1053	HN	ASP	71	-5.686	10.127	5.999	1.00	1.66
ATOM	1054	CA	ASP	71	-7.047	8.484	5.701	1.00	1.91
ATOM	1055	HA	ASP	71	-8.126	8.513	5.684	1.00	2.42
ATOM	1056	CB	ASP	71	-6.546	7.620	4.546	1.00	2.67
ATOM	1057	HB1	ASP	71	-6.623	6.578	4.813	1.00	3.03
ATOM	1058	HB2	ASP	71	-5.514	7.865	4.341	1.00	2.88
ATOM	1059	CG	ASP	71	-7.397	7.892	3.303	1.00	3.56
ATOM	1060	OD1	ASP	71	-8.476	7.330	3.215	1.00	4.08
ATOM	1061	OD2	ASP	71	-6.960	8.664	2.465	1.00	4.16
ATOM	1062	C	ASP	71	-6.577	7.889	7.028	1.00	1.46
ATOM	1063	O	ASP	71	-5.600	8.323	7.605	1.00	1.78
ATOM	1064	N	PHE	72	-7.260	6.886	7.507	1.00	1.36
ATOM	1065	HN	PHE	72	-8.038	6.546	7.018	1.00	1.67
ATOM	1066	CA	PHE	72	-6.849	6.248	8.786	1.00	1.48
ATOM	1067	HA	PHE	72	-6.504	7.007	9.473	1.00	1.75
ATOM	1068	CB	PHE	72	-8.037	5.503	9.399	1.00	2.01
ATOM	1069	HB1	PHE	72	-8.374	6.028	10.281	1.00	2.58
ATOM	1070	HB2	PHE	72	-7.733	4.503	9.669	1.00	2.43
ATOM	1071	CG	PHE	72	-9.161	5.434	8.395	1.00	2.30
ATOM	1072	CD1	PHE	72	-9.414	4.243	7.704	1.00	2.86
ATOM	1073	HD1	PHE	72	-8.802	3.372	7.887	1.00	3.09
ATOM	1074	CD2	PHE	72	-9.954	6.563	8.158	1.00	2.97
ATOM	1075	HD2	PHE	72	-9.758	7.482	8.691	1.00	3.28
ATOM	1076	CE1	PHE	72	-10.459	4.182	6.775	1.00	3.73
ATOM	1077	HE1	PHE	72	-10.655	3.264	6.242	1.00	4.46
ATOM	1078	CE2	PHE	72	-10.999	6.502	7.229	1.00	3.80
ATOM	1079	HE2	PHE	72	-11.610	7.374	7.045	1.00	4.54
ATOM	1080	CZ	PHE	72	-11.252	5.312	6.537	1.00	4.08
ATOM	1081	HZ	PHE	72	-12.058	5.264	5.821	1.00	4.92
ATOM	1082	C	PHE	72	-5.716	5.266	8.500	1.00	1.41
ATOM	1083	O	PHE	72	-5.384	4.430	9.318	1.00	2.20
ATOM	1084	N	TYR	73	-5.120	5.371	7.338	1.00	1.12
ATOM	1085	HN	TYR	73	-5.412	6.059	6.703	1.00	1.48
ATOM	1086	CA	TYR	73	-3.999	4.457	6.972	1.00	1.25
ATOM	1087	HA	TYR	73	-3.774	3.793	7.790	1.00	1.46
ATOM	1088	CB	TYR	73	-4.391	3.635	5.742	1.00	1.86
ATOM	1089	HB1	TYR	73	-3.531	3.082	5.395	1.00	2.35
ATOM	1090	HB2	TYR	73	-4.726	4.300	4.961	1.00	2.46
ATOM	1091	CG	TYR	73	-5.498	2.670	6.089	1.00	2.08
ATOM	1092	CD1	TYR	73	-5.241	1.585	6.934	1.00	2.58
ATOM	1093	HD1	TYR	73	-4.252	1.444	7.347	1.00	2.82
ATOM	1094	CD2	TYR	73	-6.779	2.853	5.553	1.00	2.85
ATOM	1095	HD2	TYR	73	-6.978	3.691	4.901	1.00	3.24
ATOM	1096	CE1	TYR	73	-6.264	0.683	7.244	1.00	3.48
ATOM	1097	HE1	TYR	73	-6.066	-0.155	7.896	1.00	4.19
ATOM	1098	CE2	TYR	73	-7.802	1.952	5.865	1.00	3.68
ATOM	1099	HE2	TYR	73	-8.789	2.093	5.452	1.00	4.49
ATOM	1100	CZ	TYR	73	-7.545	0.866	6.710	1.00	3.90
ATOM	1101	OH	TYR	73	-8.554	-0.024	7.013	1.00	5.00
ATOM	1102	HH	TYR	73	-8.689	-0.590	6.249	1.00	5.22
ATOM	1103	C	TYR	73	-2.755	5.273	6.609	1.00	0.95
ATOM	1104	O	TYR	73	-2.219	5.127	5.529	1.00	1.21
ATOM	1105	N	PRO	74	-2.273	6.106	7.495	1.00	0.74
ATOM	1106	CA	PRO	74	-1.054	6.895	7.197	1.00	0.82
ATOM	1107	HA	PRO	74	-1.254	7.648	6.453	1.00	1.05
ATOM	1108	CB	PRO	74	-0.746	7.558	8.543	1.00	1.18
ATOM	1109	HB1	PRO	74	-0.786	8.631	8.438	1.00	1.46
ATOM	1110	HB2	PRO	74	0.239	7.261	8.876	1.00	1.28
ATOM	1111	CG	PRO	74	-1.795	7.105	9.566	1.00	1.35
ATOM	1112	HG1	PRO	74	-2.229	7.967	10.049	1.00	1.70
ATOM	1113	HG2	PRO	74	-1.330	6.468	10.305	1.00	1.61
ATOM	1114	CD	PRO	74	-2.889	6.328	8.828	1.00	1.04
ATOM	1115	HD2	PRO	74	-3.098	5.393	9.328	1.00	1.24
ATOM	1116	HD1	PRO	74	-3.778	6.929	8.733	1.00	1.14
ATOM	1117	C	PRO	74	0.097	5.988	6.765	1.00	0.65

19/75

ATOM	1118	O	PRO	74	0.136	4.822	7.106	1.00	0.66
ATOM	1119	N	PHE	75	1.038	6.503	6.032	1.00	0.56
ATOM	1120	HN	PHE	75	1.000	7.447	5.770	1.00	0.61
ATOM	1121	CA	PHE	75	2.179	5.651	5.605	1.00	0.45
ATOM	1122	HA	PHE	75	1.816	4.659	5.360	1.00	0.48
ATOM	1123	CB	PHE	75	2.859	6.266	4.379	1.00	0.42
ATOM	1124	HB1	PHE	75	3.761	5.718	4.153	1.00	0.44
ATOM	1125	HB2	PHE	75	3.104	7.298	4.582	1.00	0.45
ATOM	1126	CG	PHE	75	1.915	6.190	3.200	1.00	0.48
ATOM	1127	CD1	PHE	75	1.764	4.986	2.501	1.00	0.41
ATOM	1128	HD1	PHE	75	2.329	4.115	2.797	1.00	0.45
ATOM	1129	CD2	PHE	75	1.184	7.320	2.812	1.00	0.74
ATOM	1130	HD2	PHE	75	1.300	8.249	3.349	1.00	0.90
ATOM	1131	CE1	PHE	75	0.882	4.911	1.415	1.00	0.50
ATOM	1132	HE1	PHE	75	0.767	3.982	0.877	1.00	0.53
ATOM	1133	CE2	PHE	75	0.304	7.245	1.724	1.00	0.85
ATOM	1134	HE2	PHE	75	-0.258	8.117	1.423	1.00	1.09
ATOM	1135	CZ	PHE	75	0.154	6.041	1.026	1.00	0.69
ATOM	1136	HZ	PHE	75	-0.526	5.983	0.188	1.00	0.80
ATOM	1137	C	PHE	75	3.159	5.561	6.776	1.00	0.43
ATOM	1138	O	PHE	75	3.111	6.360	7.690	1.00	0.50
ATOM	1139	N	ASP	76	4.020	4.582	6.782	1.00	0.37
ATOM	1140	HN	ASP	76	4.028	3.929	6.050	1.00	0.32
ATOM	1141	CA	ASP	76	4.967	4.432	7.927	1.00	0.43
ATOM	1142	HA	ASP	76	4.551	4.906	8.804	1.00	0.50
ATOM	1143	CB	ASP	76	5.180	2.946	8.215	1.00	0.46
ATOM	1144	HB1	ASP	76	4.224	2.467	8.365	1.00	0.49
ATOM	1145	HB2	ASP	76	5.784	2.834	9.104	1.00	0.54
ATOM	1146	CG	ASP	76	5.892	2.295	7.028	1.00	0.38
ATOM	1147	OD1	ASP	76	6.468	1.236	7.218	1.00	0.45
ATOM	1148	OD2	ASP	76	5.846	2.864	5.950	1.00	0.30
ATOM	1149	C	ASP	76	6.314	5.074	7.596	1.00	0.42
ATOM	1150	O	ASP	76	7.314	4.770	8.216	1.00	0.54
ATOM	1151	N	GLY	77	6.347	5.958	6.632	1.00	0.35
ATOM	1152	HN	GLY	77	5.525	6.187	6.151	1.00	0.36
ATOM	1153	CA	GLY	77	7.634	6.625	6.267	1.00	0.38
ATOM	1154	HA1	GLY	77	8.378	6.388	7.004	1.00	0.45
ATOM	1155	HA2	GLY	77	7.484	7.696	6.238	1.00	0.44
ATOM	1156	C	GLY	77	8.084	6.131	4.884	1.00	0.31
ATOM	1157	O	GLY	77	7.262	5.767	4.068	1.00	0.37
ATOM	1158	N	PRO	78	9.370	6.117	4.603	1.00	0.33
ATOM	1159	CA	PRO	78	9.856	5.651	3.274	1.00	0.36
ATOM	1160	HA	PRO	78	9.435	6.254	2.488	1.00	0.42
ATOM	1161	CB	PRO	78	11.364	5.903	3.359	1.00	0.46
ATOM	1162	HB1	PRO	78	11.671	6.542	2.545	1.00	0.56
ATOM	1163	HB2	PRO	78	11.892	4.962	3.303	1.00	0.48
ATOM	1164	CG	PRO	78	11.675	6.592	4.694	1.00	0.64
ATOM	1165	HG1	PRO	78	11.965	7.616	4.516	1.00	0.87
ATOM	1166	HG2	PRO	78	12.478	6.068	5.194	1.00	0.83
ATOM	1167	CD	PRO	78	10.418	6.562	5.563	1.00	0.45
ATOM	1168	HD2	PRO	78	10.535	5.848	6.369	1.00	0.48
ATOM	1169	HD1	PRO	78	10.187	7.544	5.944	1.00	0.49
ATOM	1170	C	PRO	78	9.564	4.165	3.027	1.00	0.30
ATOM	1171	O	PRO	78	8.860	3.808	2.105	1.00	0.28
ATOM	1172	N	SER	79	10.102	3.297	3.840	1.00	0.31
ATOM	1173	HN	SER	79	10.670	3.604	4.577	1.00	0.35
ATOM	1174	CA	SER	79	9.855	1.837	3.647	1.00	0.30
ATOM	1175	HA	SER	79	9.916	1.599	2.595	1.00	0.30
ATOM	1176	CB	SER	79	10.911	1.037	4.410	1.00	0.37
ATOM	1177	HB1	SER	79	11.888	1.465	4.225	1.00	0.42
ATOM	1178	HB2	SER	79	10.901	0.013	4.076	1.00	0.39
ATOM	1179	OG	SER	79	10.617	1.080	5.800	1.00	0.38
ATOM	1180	HG	SER	79	11.173	1.752	6.201	1.00	0.98
ATOM	1181	C	SER	79	8.463	1.470	4.173	1.00	0.27
ATOM	1182	O	SER	79	7.888	2.183	4.971	1.00	0.25
ATOM	1183	N	GLY	80	7.927	0.356	3.734	1.00	0.31
ATOM	1184	HN	GLY	80	8.420	-0.200	3.095	1.00	0.37
ATOM	1185	CA	GLY	80	6.576	-0.081	4.207	1.00	0.30
ATOM	1186	HA1	GLY	80	6.224	0.586	4.977	1.00	0.31
ATOM	1187	HA2	GLY	80	6.646	-1.083	4.607	1.00	0.36
ATOM	1188	C	GLY	80	5.584	-0.070	3.042	1.00	0.25
ATOM	1189	O	GLY	80	5.850	-0.601	1.981	1.00	0.25
ATOM	1190	N	LEU	81	4.440	0.531	3.232	1.00	0.23
ATOM	1191	HN	LEU	81	4.246	0.951	4.096	1.00	0.25
ATOM	1192	CA	LEU	81	3.428	0.577	2.138	1.00	0.21
ATOM	1193	HA	LEU	81	3.259	-0.417	1.761	1.00	0.22
ATOM	1194	CB	LEU	81	2.123	1.164	2.692	1.00	0.24

20/75

ATOM	1195	HB1	LEU	81	1.587	1.658	1.896	1.00	0.25
ATOM	1196	HB2	LEU	81	2.356	1.881	3.465	1.00	0.29
ATOM	1197	CG	LEU	81	1.240	0.058	3.283	1.00	0.28
ATOM	1198	HG	LEU	81	1.856	-0.678	3.779	1.00	0.31
ATOM	1199	CD1	LEU	81	0.265	0.680	4.285	1.00	0.33
ATOM	1200	HD11	LEU	81	0.071	1.706	4.009	1.00	1.05
ATOM	1201	HD12	LEU	81	0.696	0.649	5.274	1.00	1.10
ATOM	1202	HD13	LEU	81	-0.662	0.125	4.278	1.00	1.06
ATOM	1203	CD2	LEU	81	0.426	-0.606	2.168	1.00	0.31
ATOM	1204	HD21	LEU	81	1.087	-0.997	1.412	1.00	1.02
ATOM	1205	HD22	LEU	81	-0.233	0.126	1.724	1.00	1.09
ATOM	1206	HD23	LEU	81	-0.161	-1.411	2.584	1.00	1.06
ATOM	1207	C	LEU	81	3.953	1.475	1.017	1.00	0.20
ATOM	1208	O	LEU	81	3.988	2.679	1.141	1.00	0.22
ATOM	1209	N	LEU	82	4.366	0.899	-0.078	1.00	0.18
ATOM	1210	HN	LEU	82	4.334	-0.077	-0.162	1.00	0.18
ATOM	1211	CA	LEU	82	4.901	1.728	-1.195	1.00	0.18
ATOM	1212	HA	LEU	82	5.519	2.520	-0.799	1.00	0.19
ATOM	1213	CB	LEU	82	5.728	0.840	-2.128	1.00	0.18
ATOM	1214	HB1	LEU	82	6.235	1.457	-2.854	1.00	0.20
ATOM	1215	HB2	LEU	82	5.071	0.151	-2.640	1.00	0.20
ATOM	1216	CG	LEU	82	6.763	0.050	-1.323	1.00	0.18
ATOM	1217	HG	LEU	82	6.262	-0.523	-0.556	1.00	0.22
ATOM	1218	CD1	LEU	82	7.513	-0.898	-2.259	1.00	0.17
ATOM	1219	HD11	LEU	82	8.102	-0.321	-2.957	1.00	0.97
ATOM	1220	HD12	LEU	82	6.802	-1.503	-2.802	1.00	0.95
ATOM	1221	HD13	LEU	82	8.163	-1.537	-1.681	1.00	0.98
ATOM	1222	CD2	LEU	82	7.764	1.010	-0.675	1.00	0.23
ATOM	1223	HD21	LEU	82	8.019	1.790	-1.375	1.00	1.03
ATOM	1224	HD22	LEU	82	8.657	0.466	-0.403	1.00	1.07
ATOM	1225	HD23	LEU	82	7.326	1.447	0.209	1.00	1.02
ATOM	1226	C	LEU	82	3.740	2.329	-1.986	1.00	0.19
ATOM	1227	O	LEU	82	3.882	3.341	-2.646	1.00	0.21
ATOM	1228	N	ALA	83	2.594	1.711	-1.919	1.00	0.21
ATOM	1229	HN	ALA	83	2.512	0.899	-1.376	1.00	0.24
ATOM	1230	CA	ALA	83	1.410	2.225	-2.662	1.00	0.22
ATOM	1231	HA	ALA	83	1.217	3.251	-2.381	1.00	0.22
ATOM	1232	CB	ALA	83	1.668	2.140	-4.171	1.00	0.23
ATOM	1233	HB1	ALA	83	2.522	2.746	-4.429	1.00	0.98
ATOM	1234	HB2	ALA	83	0.801	2.497	-4.705	1.00	1.00
ATOM	1235	HB3	ALA	83	1.860	1.113	-4.445	1.00	1.05
ATOM	1236	C	ALA	83	0.204	1.350	-2.317	1.00	0.27
ATOM	1237	O	ALA	83	0.342	0.301	-1.720	1.00	0.36
ATOM	1238	N	HIS	84	-0.976	1.762	-2.686	1.00	0.24
ATOM	1239	HN	HIS	84	-1.075	2.609	-3.170	1.00	0.20
ATOM	1240	CA	HIS	84	-2.173	0.933	-2.370	1.00	0.30
ATOM	1241	HA	HIS	84	-1.940	-0.108	-2.542	1.00	0.36
ATOM	1242	CB	HIS	84	-2.562	1.127	-0.903	1.00	0.40
ATOM	1243	HB1	HIS	84	-1.695	0.965	-0.278	1.00	0.48
ATOM	1244	HB2	HIS	84	-3.332	0.419	-0.638	1.00	0.45
ATOM	1245	CG	HIS	84	-3.074	2.525	-0.692	1.00	0.44
ATOM	1246	ND1	HIS	84	-4.384	2.781	-0.321	1.00	1.32
ATOM	1247	HD1	HIS	84	-5.084	2.112	-0.169	1.00	2.02
ATOM	1248	CD2	HIS	84	-2.465	3.752	-0.788	1.00	0.74
ATOM	1249	HD2	HIS	84	-1.432	3.915	-1.060	1.00	1.58
ATOM	1250	CE1	HIS	84	-4.521	4.114	-0.208	1.00	1.21
ATOM	1251	HE1	HIS	84	-5.441	4.606	0.071	1.00	1.87
ATOM	1252	NE2	HIS	84	-3.381	4.754	-0.482	1.00	0.53
ATOM	1253	C	HIS	84	-3.337	1.343	-3.274	1.00	0.25
ATOM	1254	O	HIS	84	-3.347	2.417	-3.843	1.00	0.23
ATOM	1255	N	ALA	85	-4.313	0.489	-3.417	1.00	0.27
ATOM	1256	HN	ALA	85	-4.279	-0.374	-2.954	1.00	0.34
ATOM	1257	CA	ALA	85	-5.474	0.817	-4.291	1.00	0.24
ATOM	1258	HA	ALA	85	-5.582	1.890	-4.364	1.00	0.22
ATOM	1259	CB	ALA	85	-5.236	0.231	-5.685	1.00	0.25
ATOM	1260	HB1	ALA	85	-5.079	-0.835	-5.605	1.00	1.05
ATOM	1261	HB2	ALA	85	-4.364	0.690	-6.126	1.00	1.05
ATOM	1262	HB3	ALA	85	-6.097	0.420	-6.308	1.00	1.06
ATOM	1263	C	ALA	85	-6.748	0.210	-3.698	1.00	0.26
ATOM	1264	O	ALA	85	-6.694	-0.611	-2.804	1.00	0.33
ATOM	1265	N	PHE	86	-7.892	0.605	-4.198	1.00	0.28
ATOM	1266	HN	PHE	86	-7.905	1.264	-4.922	1.00	0.31
ATOM	1267	CA	PHE	86	-9.179	0.053	-3.677	1.00	0.34
ATOM	1268	HA	PHE	86	-9.000	-0.443	-2.737	1.00	0.39
ATOM	1269	CB	PHE	86	-10.170	1.205	-3.471	1.00	0.36
ATOM	1270	HB1	PHE	86	-11.177	0.821	-3.459	1.00	0.42
ATOM	1271	HB2	PHE	86	-10.068	1.913	-4.279	1.00	0.33

21/75

ATOM	1272	CG	PHE	86	-9.877	1.896	-2.159	1.00	0.39
ATOM	1273	CD1	PHE	86	-8.784	2.764	-2.050	1.00	0.46
ATOM	1274	HD1	PHE	86	-8.146	2.939	-2.903	1.00	0.67
ATOM	1275	CD2	PHE	86	-10.703	1.670	-1.051	1.00	0.67
ATOM	1276	HD2	PHE	86	-11.546	1.001	-1.133	1.00	0.91
ATOM	1277	CE1	PHE	86	-8.516	3.406	-0.835	1.00	0.50
ATOM	1278	HE1	PHE	86	-7.673	4.075	-0.751	1.00	0.69
ATOM	1279	CE2	PHE	86	-10.435	2.311	0.165	1.00	0.74
ATOM	1280	HE2	PHE	86	-11.071	2.136	1.020	1.00	1.02
ATOM	1281	CZ	PHE	86	-9.342	3.179	0.273	1.00	0.54
ATOM	1282	HZ	PHE	86	-9.135	3.674	1.211	1.00	0.62
ATOM	1283	C	PHE	86	-9.746	-0.940	-4.710	1.00	0.36
ATOM	1284	O	PHE	86	-9.480	-0.812	-5.889	1.00	0.34
ATOM	1285	N	PRO	87	-10.516	-1.926	-4.293	1.00	0.43
ATOM	1286	CA	PRO	87	-11.082	-2.914	-5.257	1.00	0.46
ATOM	1287	HA	PRO	87	-10.296	-3.524	-5.665	1.00	0.53
ATOM	1288	CB	PRO	87	-11.990	-3.770	-4.370	1.00	0.60
ATOM	1289	HB1	PRO	87	-11.644	-4.792	-4.377	1.00	0.69
ATOM	1290	HB2	PRO	87	-13.004	-3.727	-4.742	1.00	0.73
ATOM	1291	CG	PRO	87	-11.943	-3.225	-2.937	1.00	0.58
ATOM	1292	HG1	PRO	87	-11.694	-4.022	-2.253	1.00	0.61
ATOM	1293	HG2	PRO	87	-12.905	-2.808	-2.676	1.00	0.66
ATOM	1294	CD	PRO	87	-10.872	-2.135	-2.861	1.00	0.50
ATOM	1295	HD2	PRO	87	-11.277	-1.235	-2.421	1.00	0.50
ATOM	1296	HD1	PRO	87	-10.014	-2.484	-2.309	1.00	0.52
ATOM	1297	C	PRO	87	-11.895	-2.246	-6.379	1.00	0.40
ATOM	1298	O	PRO	87	-12.221	-1.078	-6.299	1.00	0.42
ATOM	1299	N	PRO	88	-12.221	-2.981	-7.419	1.00	0.44
ATOM	1300	CA	PRO	88	-13.007	-2.416	-8.554	1.00	0.48
ATOM	1301	HA	PRO	88	-12.443	-1.645	-9.053	1.00	0.52
ATOM	1302	CB	PRO	88	-13.163	-3.622	-9.488	1.00	0.61
ATOM	1303	HB1	PRO	88	-12.604	-3.449	-10.395	1.00	0.83
ATOM	1304	HB2	PRO	88	-14.204	-3.772	-9.728	1.00	0.74
ATOM	1305	CG	PRO	88	-12.609	-4.863	-8.781	1.00	0.57
ATOM	1306	HG1	PRO	88	-11.945	-5.395	-9.446	1.00	0.71
ATOM	1307	HG2	PRO	88	-13.425	-5.508	-8.488	1.00	0.64
ATOM	1308	CD	PRO	88	-11.835	-4.413	-7.540	1.00	0.56
ATOM	1309	HD2	PRO	88	-12.146	-4.977	-6.671	1.00	0.62
ATOM	1310	HD1	PRO	88	-10.773	-4.503	-7.702	1.00	0.65
ATOM	1311	C	PRO	88	-14.372	-1.873	-8.109	1.00	0.47
ATOM	1312	O	PRO	88	-15.380	-2.551	-8.172	1.00	0.88
ATOM	1313	N	GLY	89	-14.400	-0.647	-7.661	1.00	0.63
ATOM	1314	HN	GLY	89	-13.571	-0.129	-7.626	1.00	1.01
ATOM	1315	CA	GLY	89	-15.681	-0.026	-7.209	1.00	0.65
ATOM	1316	HA1	GLY	89	-15.536	0.422	-6.239	1.00	0.62
ATOM	1317	HA2	GLY	89	-16.455	-0.778	-7.148	1.00	0.78
ATOM	1318	C	GLY	89	-16.092	1.057	-8.210	1.00	0.74
ATOM	1319	O	GLY	89	-15.541	1.151	-9.289	1.00	0.84
ATOM	1320	N	PRO	90	-17.044	1.878	-7.852	1.00	0.95
ATOM	1321	CA	PRO	90	-17.499	2.973	-8.750	1.00	1.19
ATOM	1322	HA	PRO	90	-17.819	2.565	-9.697	1.00	1.37
ATOM	1323	CB	PRO	90	-18.720	3.532	-7.990	1.00	1.55
ATOM	1324	HB1	PRO	90	-19.602	3.432	-8.605	1.00	1.85
ATOM	1325	HB2	PRO	90	-18.572	4.567	-7.740	1.00	1.74
ATOM	1326	CG	PRO	90	-18.913	2.724	-6.702	1.00	1.46
ATOM	1327	HG1	PRO	90	-19.828	2.155	-6.763	1.00	1.60
ATOM	1328	HG2	PRO	90	-18.959	3.396	-5.857	1.00	1.57
ATOM	1329	CD	PRO	90	-17.729	1.769	-6.539	1.00	1.17
ATOM	1330	HD2	PRO	90	-17.083	2.099	-5.736	1.00	1.17
ATOM	1331	HD1	PRO	90	-18.067	0.759	-6.375	1.00	1.28
ATOM	1332	C	PRO	90	-16.375	4.011	-8.972	1.00	1.14
ATOM	1333	O	PRO	90	-15.269	3.649	-9.320	1.00	1.53
ATOM	1334	N	ASN	91	-16.624	5.282	-8.790	1.00	1.17
ATOM	1335	HN	ASN	91	-17.514	5.578	-8.517	1.00	1.40
ATOM	1336	CA	ASN	91	-15.541	6.286	-9.008	1.00	1.38
ATOM	1337	HA	ASN	91	-15.147	6.169	-10.005	1.00	1.58
ATOM	1338	CB	ASN	91	-16.116	7.700	-8.857	1.00	1.87
ATOM	1339	HB1	ASN	91	-15.336	8.372	-8.532	1.00	2.33
ATOM	1340	HB2	ASN	91	-16.908	7.686	-8.122	1.00	1.96
ATOM	1341	CG	ASN	91	-16.678	8.184	-10.197	1.00	2.69
ATOM	1342	OD1	ASN	91	-16.132	7.890	-11.242	1.00	3.20
ATOM	1343	ND2	ASN	91	-17.748	8.931	-10.212	1.00	3.47
ATOM	1344	HD21	ASN	91	-18.186	9.176	-9.370	1.00	3.59
ATOM	1345	HD22	ASN	91	-18.112	9.249	-11.064	1.00	4.20
ATOM	1346	C	ASN	91	-14.404	6.098	-7.992	1.00	1.15
ATOM	1347	O	ASN	91	-13.242	6.135	-8.344	1.00	1.26
ATOM	1348	N	TYR	92	-14.719	5.924	-6.736	1.00	1.01

22/75

ATOM	1349	HN	TYR	92	-15.660	5.916	-6.462	1.00	1.08
ATOM	1350	CA	TYR	92	-13.639	5.768	-5.711	1.00	0.97
ATOM	1351	HA	TYR	92	-12.994	6.632	-5.739	1.00	1.14
ATOM	1352	CB	TYR	92	-14.262	5.652	-4.319	1.00	1.09
ATOM	1353	HB1	TYR	92	-13.543	5.214	-3.643	1.00	1.62
ATOM	1354	HB2	TYR	92	-15.135	5.020	-4.369	1.00	1.45
ATOM	1355	CG	TYR	92	-14.656	7.018	-3.810	1.00	1.52
ATOM	1356	CD1	TYR	92	-13.672	7.979	-3.549	1.00	2.14
ATOM	1357	HD1	TYR	92	-12.631	7.747	-3.719	1.00	2.46
ATOM	1358	CD2	TYR	92	-16.006	7.320	-3.588	1.00	2.44
ATOM	1359	HD2	TYR	92	-16.766	6.580	-3.789	1.00	2.86
ATOM	1360	CE1	TYR	92	-14.037	9.241	-3.066	1.00	3.06
ATOM	1361	HE1	TYR	92	-13.278	9.982	-2.865	1.00	3.78
ATOM	1362	CE2	TYR	92	-16.370	8.582	-3.107	1.00	3.33
ATOM	1363	HE2	TYR	92	-17.411	8.815	-2.936	1.00	4.19
ATOM	1364	CZ	TYR	92	-15.386	9.542	-2.846	1.00	3.50
ATOM	1365	OH	TYR	92	-15.746	10.786	-2.368	1.00	4.57
ATOM	1366	HH	TYR	92	-15.602	10.791	-1.419	1.00	4.91
ATOM	1367	C	TYR	92	-12.808	4.508	-5.966	1.00	0.78
ATOM	1368	O	TYR	92	-11.605	4.506	-5.798	1.00	0.81
ATOM	1369	N	GLY	93	-13.436	3.430	-6.337	1.00	0.64
ATOM	1370	HN	GLY	93	-14.410	3.441	-6.445	1.00	0.70
ATOM	1371	CA	GLY	93	-12.674	2.170	-6.560	1.00	0.51
ATOM	1372	HA1	GLY	93	-13.366	1.366	-6.740	1.00	0.51
ATOM	1373	HA2	GLY	93	-12.090	1.947	-5.678	1.00	0.51
ATOM	1374	C	GLY	93	-11.739	2.310	-7.761	1.00	0.49
ATOM	1375	O	GLY	93	-11.832	3.242	-8.534	1.00	0.61
ATOM	1376	N	GLY	94	-10.844	1.373	-7.923	1.00	0.45
ATOM	1377	HN	GLY	94	-10.799	0.627	-7.288	1.00	0.44
ATOM	1378	CA	GLY	94	-9.902	1.420	-9.075	1.00	0.55
ATOM	1379	HA1	GLY	94	-10.459	1.569	-9.988	1.00	0.63
ATOM	1380	HA2	GLY	94	-9.363	0.485	-9.133	1.00	0.58
ATOM	1381	C	GLY	94	-8.905	2.569	-8.901	1.00	0.60
ATOM	1382	O	GLY	94	-8.109	2.838	-9.772	1.00	1.14
ATOM	1383	N	ASP	95	-8.933	3.252	-7.790	1.00	0.24
ATOM	1384	HN	ASP	95	-9.581	3.028	-7.089	1.00	0.52
ATOM	1385	CA	ASP	95	-7.976	4.382	-7.597	1.00	0.24
ATOM	1386	HA	ASP	95	-7.888	4.939	-8.518	1.00	0.28
ATOM	1387	CB	ASP	95	-8.493	5.303	-6.491	1.00	0.26
ATOM	1388	HB1	ASP	95	-9.500	5.617	-6.724	1.00	0.28
ATOM	1389	HB2	ASP	95	-7.853	6.170	-6.415	1.00	0.30
ATOM	1390	CG	ASP	95	-8.494	4.549	-5.162	1.00	0.28
ATOM	1391	OD1	ASP	95	-8.543	5.200	-4.132	1.00	1.08
ATOM	1392	OD2	ASP	95	-8.440	3.331	-5.198	1.00	1.14
ATOM	1393	C	ASP	95	-6.605	3.827	-7.202	1.00	0.23
ATOM	1394	O	ASP	95	-6.479	2.683	-6.815	1.00	0.24
ATOM	1395	N	ALA	96	-5.573	4.626	-7.297	1.00	0.23
ATOM	1396	HN	ALA	96	-5.692	5.546	-7.614	1.00	0.23
ATOM	1397	CA	ALA	96	-4.215	4.131	-6.926	1.00	0.25
ATOM	1398	HA	ALA	96	-4.307	3.360	-6.175	1.00	0.25
ATOM	1399	CB	ALA	96	-3.527	3.553	-8.164	1.00	0.30
ATOM	1400	HB1	ALA	96	-2.528	3.236	-7.905	1.00	1.08
ATOM	1401	HB2	ALA	96	-3.476	4.309	-8.934	1.00	1.08
ATOM	1402	HB3	ALA	96	-4.090	2.706	-8.528	1.00	1.03
ATOM	1403	C	ALA	96	-3.375	5.284	-6.372	1.00	0.25
ATOM	1404	O	ALA	96	-3.222	6.313	-7.005	1.00	0.29
ATOM	1405	N	HXS	97	-2.831	5.113	-5.192	1.00	0.25
ATOM	1406	HN	HXS	97	-2.976	4.271	-4.710	1.00	0.28
ATOM	1407	CA	HXS	97	-1.996	6.187	-4.574	1.00	0.27
ATOM	1408	HA	HXS	97	-2.010	7.068	-5.198	1.00	0.28
ATOM	1409	CB	HXS	97	-2.564	6.537	-3.197	1.00	0.33
ATOM	1410	HB1	HXS	97	-1.969	7.319	-2.750	1.00	0.44
ATOM	1411	HB2	HXS	97	-2.540	5.661	-2.566	1.00	0.39
ATOM	1412	CG	HXS	97	-3.983	7.009	-3.349	1.00	0.37
ATOM	1413	ND1	HXS	97	-4.697	7.052	-2.163	1.00	0.80
ATOM	1414	CD2	HXS	97	-4.783	7.420	-4.384	1.00	0.55
ATOM	1415	HD2	HXS	97	-4.517	7.497	-5.428	1.00	0.94
ATOM	1416	CE1	HXS	97	-5.918	7.487	-2.498	1.00	0.86
ATOM	1417	HE1	HXS	97	-6.724	7.632	-1.795	1.00	1.24
ATOM	1418	NE2	HXS	97	-6.018	7.722	-3.819	1.00	0.59
ATOM	1419	HE2	HXS	97	-6.812	8.044	-4.294	1.00	0.72
ATOM	1420	C	HXS	97	-0.552	5.700	-4.420	1.00	0.26
ATOM	1421	O	HXS	97	-0.299	4.525	-4.237	1.00	0.39
ATOM	1422	N	PHE	98	0.391	6.604	-4.496	1.00	0.18
ATOM	1423	HN	PHE	98	0.147	7.540	-4.648	1.00	0.23
ATOM	1424	CA	PHE	98	1.832	6.230	-4.360	1.00	0.17
ATOM	1425	HA	PHE	98	1.921	5.190	-4.085	1.00	0.18

23/75

ATOM	1426	CB	PHE	98	2.543	6.472	-5.691	1.00	0.18
ATOM	1427	HB1	PHE	98	3.611	6.464	-5.536	1.00	0.21
ATOM	1428	HB2	PHE	98	2.243	7.431	-6.085	1.00	0.20
ATOM	1429	CG	PHE	98	2.169	5.391	-6.674	1.00	0.19
ATOM	1430	CD1	PHE	98	3.114	4.428	-7.048	1.00	0.22
ATOM	1431	HD1	PHE	98	4.110	4.456	-6.631	1.00	0.25
ATOM	1432	CD2	PHE	98	0.880	5.355	-7.214	1.00	0.22
ATOM	1433	HD2	PHE	98	0.151	6.098	-6.924	1.00	0.24
ATOM	1434	CE1	PHE	98	2.768	3.429	-7.963	1.00	0.25
ATOM	1435	HE1	PHE	98	3.496	2.685	-8.252	1.00	0.29
ATOM	1436	CE2	PHE	98	0.533	4.355	-8.127	1.00	0.26
ATOM	1437	HE2	PHE	98	-0.462	4.327	-8.542	1.00	0.31
ATOM	1438	CZ	PHE	98	1.478	3.392	-8.503	1.00	0.26
ATOM	1439	HZ	PHE	98	1.214	2.622	-9.211	1.00	0.30
ATOM	1440	C	PHE	98	2.487	7.104	-3.286	1.00	0.17
ATOM	1441	O	PHE	98	2.081	8.226	-3.058	1.00	0.19
ATOM	1442	N	ASP	99	3.498	6.604	-2.625	1.00	0.19
ATOM	1443	HN	ASP	99	3.813	5.693	-2.820	1.00	0.22
ATOM	1444	CA	ASP	99	4.167	7.424	-1.570	1.00	0.20
ATOM	1445	HA	ASP	99	3.421	7.956	-0.998	1.00	0.20
ATOM	1446	CB	ASP	99	4.973	6.516	-0.638	1.00	0.25
ATOM	1447	HB1	ASP	99	5.567	7.122	0.029	1.00	0.28
ATOM	1448	HB2	ASP	99	5.624	5.884	-1.226	1.00	0.30
ATOM	1449	CG	ASP	99	4.023	5.646	0.180	1.00	0.41
ATOM	1450	OD1	ASP	99	2.838	5.680	-0.100	1.00	0.89
ATOM	1451	OD2	ASP	99	4.497	4.968	1.079	1.00	0.27
ATOM	1452	C	ASP	99	5.123	8.426	-2.224	1.00	0.21
ATOM	1453	O	ASP	99	6.020	8.054	-2.954	1.00	0.25
ATOM	1454	N	ASP	100	4.946	9.694	-1.962	1.00	0.23
ATOM	1455	HN	ASP	100	4.222	9.976	-1.365	1.00	0.23
ATOM	1456	CA	ASP	100	5.857	10.710	-2.565	1.00	0.29
ATOM	1457	HA	ASP	100	6.169	10.379	-3.545	1.00	0.31
ATOM	1458	CB	ASP	100	5.127	12.049	-2.684	1.00	0.34
ATOM	1459	HB1	ASP	100	5.130	12.544	-1.727	1.00	0.34
ATOM	1460	HB2	ASP	100	4.109	11.879	-2.999	1.00	0.34
ATOM	1461	CG	ASP	100	5.844	12.929	-3.710	1.00	0.43
ATOM	1462	OD1	ASP	100	5.240	13.887	-4.164	1.00	1.21
ATOM	1463	OD2	ASP	100	6.984	12.630	-4.025	1.00	1.12
ATOM	1464	C	ASP	100	7.085	10.885	-1.667	1.00	0.30
ATOM	1465	O	ASP	100	8.032	11.559	-2.018	1.00	0.32
ATOM	1466	N	ASP	101	7.074	10.280	-0.510	1.00	0.31
ATOM	1467	HN	ASP	101	6.298	9.741	-0.249	1.00	0.32
ATOM	1468	CA	ASP	101	8.236	10.407	0.415	1.00	0.33
ATOM	1469	HA	ASP	101	8.647	11.403	0.345	1.00	0.36
ATOM	1470	CB	ASP	101	7.778	10.142	1.851	1.00	0.39
ATOM	1471	HB1	ASP	101	8.641	10.060	2.495	1.00	0.41
ATOM	1472	HB2	ASP	101	7.216	9.220	1.884	1.00	0.39
ATOM	1473	CG	ASP	101	6.896	11.296	2.330	1.00	0.45
ATOM	1474	OD1	ASP	101	7.027	12.380	1.786	1.00	1.25
ATOM	1475	OD2	ASP	101	6.104	11.076	3.231	1.00	1.09
ATOM	1476	C	ASP	101	9.304	9.385	0.028	1.00	0.30
ATOM	1477	O	ASP	101	10.411	9.405	0.529	1.00	0.29
ATOM	1478	N	GLU	102	8.971	8.484	-0.849	1.00	0.30
ATOM	1479	HN	GLU	102	8.068	8.484	-1.230	1.00	0.31
ATOM	1480	CA	GLU	102	9.950	7.444	-1.266	1.00	0.29
ATOM	1481	HA	GLU	102	10.649	7.263	-0.463	1.00	0.30
ATOM	1482	CB	GLU	102	9.195	6.155	-1.585	1.00	0.35
ATOM	1483	HB1	GLU	102	9.873	5.437	-2.020	1.00	0.36
ATOM	1484	HB2	GLU	102	8.397	6.368	-2.282	1.00	0.40
ATOM	1485	CG	GLU	102	8.611	5.584	-0.293	1.00	0.46
ATOM	1486	HG1	GLU	102	8.020	6.342	0.200	1.00	1.18
ATOM	1487	HG2	GLU	102	9.415	5.276	0.356	1.00	1.03
ATOM	1488	CD	GLU	102	7.724	4.381	-0.616	1.00	0.83
ATOM	1489	OE1	GLU	102	7.601	4.060	-1.786	1.00	1.63
ATOM	1490	OE2	GLU	102	7.184	3.801	0.314	1.00	0.87
ATOM	1491	C	GLU	102	10.707	7.917	-2.508	1.00	0.25
ATOM	1492	O	GLU	102	10.359	8.910	-3.115	1.00	0.25
ATOM	1493	N	THR	103	11.741	7.213	-2.886	1.00	0.25
ATOM	1494	HN	THR	103	12.003	6.416	-2.379	1.00	0.28
ATOM	1495	CA	THR	103	12.525	7.620	-4.088	1.00	0.23
ATOM	1496	HA	THR	103	12.356	8.665	-4.301	1.00	0.23
ATOM	1497	CB	THR	103	14.016	7.383	-3.824	1.00	0.27
ATOM	1498	HB	THR	103	14.169	6.359	-3.521	1.00	0.30
ATOM	1499	OG1	THR	103	14.455	8.252	-2.789	1.00	0.29
ATOM	1500	HG1	THR	103	15.334	8.564	-3.016	1.00	0.86
ATOM	1501	CG2	THR	103	14.820	7.656	-5.098	1.00	0.29
ATOM	1502	HG21	THR	103	15.864	7.777	-4.846	1.00	1.00

24/75

ATOM	1503	HG22	THR	103	14.457	8.557	-5.569	1.00	1.08
ATOM	1504	HG23	THR	103	14.710	6.824	-5.779	1.00	1.01
ATOM	1505	C	THR	103	12.083	6.777	-5.281	1.00	0.22
ATOM	1506	O	THR	103	12.417	5.614	-5.394	1.00	0.23
ATOM	1507	N	TRP	104	11.332	7.358	-6.175	1.00	0.21
ATOM	1508	HN	TRP	104	11.076	8.297	-6.063	1.00	0.23
ATOM	1509	CA	TRP	104	10.867	6.598	-7.364	1.00	0.21
ATOM	1510	HA	TRP	104	10.750	5.556	-7.104	1.00	0.20
ATOM	1511	CB	TRP	104	9.525	7.165	-7.831	1.00	0.23
ATOM	1512	HB1	TRP	104	9.188	6.623	-8.702	1.00	0.24
ATOM	1513	HB2	TRP	104	9.641	8.210	-8.078	1.00	0.25
ATOM	1514	CG	TRP	104	8.520	7.018	-6.731	1.00	0.24
ATOM	1515	CD1	TRP	104	8.098	8.019	-5.924	1.00	0.31
ATOM	1516	HD1	TRP	104	8.427	9.045	-5.972	1.00	0.36
ATOM	1517	CD2	TRP	104	7.811	5.821	-6.300	1.00	0.21
ATOM	1518	NE1	TRP	104	7.176	7.512	-5.026	1.00	0.31
ATOM	1519	HE1	TRP	104	6.718	8.030	-4.331	1.00	0.36
ATOM	1520	CE2	TRP	104	6.963	6.162	-5.220	1.00	0.24
ATOM	1521	CE3	TRP	104	7.819	4.486	-6.739	1.00	0.18
ATOM	1522	HE3	TRP	104	8.458	4.198	-7.559	1.00	0.19
ATOM	1523	CZ2	TRP	104	6.153	5.213	-4.596	1.00	0.23
ATOM	1524	HZ2	TRP	104	5.515	5.499	-3.774	1.00	0.27
ATOM	1525	CZ3	TRP	104	7.005	3.527	-6.114	1.00	0.20
ATOM	1526	HZ3	TRP	104	7.019	2.504	-6.460	1.00	0.23
ATOM	1527	CH2	TRP	104	6.173	3.891	-5.045	1.00	0.21
ATOM	1528	HH2	TRP	104	5.548	3.150	-4.568	1.00	0.23
ATOM	1529	C	TRP	104	11.911	6.732	-8.474	1.00	0.21
ATOM	1530	O	TRP	104	12.276	7.824	-8.864	1.00	0.24
ATOM	1531	N	THR	105	12.403	5.630	-8.973	1.00	0.20
ATOM	1532	HN	THR	105	12.098	4.763	-8.633	1.00	0.19
ATOM	1533	CA	THR	105	13.437	5.685	-10.048	1.00	0.21
ATOM	1534	HA	THR	105	13.415	6.652	-10.525	1.00	0.24
ATOM	1535	CB	THR	105	14.817	5.459	-9.428	1.00	0.21
ATOM	1536	HB	THR	105	15.018	6.233	-8.704	1.00	0.21
ATOM	1537	OG1	THR	105	15.806	5.497	-10.447	1.00	0.24
ATOM	1538	HG1	THR	105	15.882	6.404	-10.752	1.00	0.86
ATOM	1539	CG2	THR	105	14.846	4.101	-8.729	1.00	0.21
ATOM	1540	HG21	THR	105	15.178	4.233	-7.711	1.00	1.04
ATOM	1541	HG22	THR	105	15.524	3.442	-9.249	1.00	1.07
ATOM	1542	HG23	THR	105	13.854	3.674	-8.731	1.00	0.99
ATOM	1543	C	THR	105	13.166	4.597	-11.087	1.00	0.23
ATOM	1544	O	THR	105	12.521	3.606	-10.808	1.00	0.23
ATOM	1545	N	SER	106	13.668	4.769	-12.282	1.00	0.26
ATOM	1546	HN	SER	106	14.194	5.572	-12.480	1.00	0.29
ATOM	1547	CA	SER	106	13.454	3.739	-13.337	1.00	0.29
ATOM	1548	HA	SER	106	12.570	3.163	-13.111	1.00	0.30
ATOM	1549	CB	SER	106	13.290	4.423	-14.695	1.00	0.35
ATOM	1550	HB1	SER	106	14.249	4.467	-15.193	1.00	1.09
ATOM	1551	HB2	SER	106	12.916	5.424	-14.554	1.00	0.96
ATOM	1552	OG	SER	106	12.365	3.685	-15.483	1.00	1.44
ATOM	1553	HG	SER	106	11.671	4.285	-15.766	1.00	1.97
ATOM	1554	C	SER	106	14.674	2.817	-13.372	1.00	0.28
ATOM	1555	O	SER	106	14.669	1.781	-14.006	1.00	0.31
ATOM	1556	N	SER	107	15.715	3.187	-12.677	1.00	0.26
ATOM	1557	HN	SER	107	15.687	4.023	-12.166	1.00	0.25
ATOM	1558	CA	SER	107	16.940	2.340	-12.641	1.00	0.27
ATOM	1559	HA	SER	107	17.018	1.778	-13.560	1.00	0.29
ATOM	1560	CB	SER	107	18.175	3.226	-12.474	1.00	0.28
ATOM	1561	HB1	SER	107	18.292	3.847	-13.353	1.00	1.12
ATOM	1562	HB2	SER	107	19.049	2.609	-12.355	1.00	1.04
ATOM	1563	OG	SER	107	18.017	4.040	-11.320	1.00	1.29
ATOM	1564	HG	SER	107	18.556	4.827	-11.436	1.00	1.82
ATOM	1565	C	SER	107	16.836	1.376	-11.460	1.00	0.26
ATOM	1566	O	SER	107	15.829	1.324	-10.781	1.00	0.26
ATOM	1567	N	SER	108	17.859	0.609	-11.203	1.00	0.28
ATOM	1568	HN	SER	108	18.666	0.658	-11.757	1.00	0.31
ATOM	1569	CA	SER	108	17.788	-0.342	-10.061	1.00	0.30
ATOM	1570	HA	SER	108	16.775	-0.706	-9.967	1.00	0.30
ATOM	1571	CB	SER	108	18.728	-1.527	-10.330	1.00	0.36
ATOM	1572	HB1	SER	108	19.561	-1.505	-9.642	1.00	1.09
ATOM	1573	HB2	SER	108	19.103	-1.468	-11.338	1.00	0.95
ATOM	1574	OG	SER	108	18.005	-2.741	-10.176	1.00	1.47
ATOM	1575	HG	SER	108	18.550	-3.456	-10.513	1.00	2.00
ATOM	1576	C	SER	108	18.181	0.390	-8.767	1.00	0.28
ATOM	1577	O	SER	108	19.279	0.265	-8.261	1.00	0.33
ATOM	1578	N	LYS	109	17.272	1.157	-8.224	1.00	0.24
ATOM	1579	HN	LYS	109	16.392	1.244	-8.646	1.00	0.22

25/75

ATOM	1580	CA	LYS	109	17.561	1.897	-6.960	1.00	0.23
ATOM	1581	HA	LYS	109	18.275	1.341	-6.370	1.00	0.25
ATOM	1582	CB	LYS	109	18.123	3.293	-7.268	1.00	0.24
ATOM	1583	HB1	LYS	109	18.172	3.868	-6.355	1.00	0.27
ATOM	1584	HB2	LYS	109	17.472	3.793	-7.970	1.00	0.25
ATOM	1585	CG	LYS	109	19.525	3.177	-7.868	1.00	0.30
ATOM	1586	HG1	LYS	109	19.476	2.615	-8.785	1.00	0.54
ATOM	1587	HG2	LYS	109	20.177	2.675	-7.170	1.00	0.70
ATOM	1588	CD	LYS	109	20.072	4.574	-8.169	1.00	0.75
ATOM	1589	HD1	LYS	109	20.124	5.144	-7.254	1.00	1.27
ATOM	1590	HD2	LYS	109	19.420	5.074	-8.870	1.00	1.27
ATOM	1591	CE	LYS	109	21.475	4.453	-8.770	1.00	1.13
ATOM	1592	HE1	LYS	109	21.396	4.264	-9.830	1.00	1.68
ATOM	1593	HE2	LYS	109	22.000	3.636	-8.297	1.00	1.68
ATOM	1594	NZ	LYS	109	22.224	5.721	-8.545	1.00	1.79
ATOM	1595	HZ1	LYS	109	21.689	6.516	-8.948	1.00	2.22
ATOM	1596	HZ2	LYS	109	23.155	5.660	-9.006	1.00	2.17
ATOM	1597	HZ3	LYS	109	22.351	5.873	-7.525	1.00	2.34
ATOM	1598	C	LYS	109	16.259	2.052	-6.175	1.00	0.21
ATOM	1599	O	LYS	109	15.190	2.110	-6.747	1.00	0.20
ATOM	1600	N	GLY	110	16.338	2.124	-4.873	1.00	0.23
ATOM	1601	HN	GLY	110	17.212	2.079	-4.432	1.00	0.26
ATOM	1602	CA	GLY	110	15.099	2.283	-4.056	1.00	0.22
ATOM	1603	HA1	GLY	110	14.751	3.302	-4.124	1.00	0.23
ATOM	1604	HA2	GLY	110	15.316	2.044	-3.024	1.00	0.25
ATOM	1605	C	GLY	110	14.013	1.342	-4.581	1.00	0.19
ATOM	1606	O	GLY	110	14.281	0.216	-4.949	1.00	0.20
ATOM	1607	N	TYR	111	12.789	1.801	-4.626	1.00	0.17
ATOM	1608	HN	TYR	111	12.599	2.716	-4.330	1.00	0.18
ATOM	1609	CA	TYR	111	11.683	0.941	-5.136	1.00	0.15
ATOM	1610	HA	TYR	111	11.975	-0.098	-5.088	1.00	0.16
ATOM	1611	CB	TYR	111	10.437	1.162	-4.277	1.00	0.15
ATOM	1612	HB1	TYR	111	9.633	0.540	-4.641	1.00	0.15
ATOM	1613	HB2	TYR	111	10.143	2.200	-4.330	1.00	0.16
ATOM	1614	CG	TYR	111	10.745	0.798	-2.844	1.00	0.17
ATOM	1615	CD1	TYR	111	10.648	-0.533	-2.422	1.00	0.17
ATOM	1616	HD1	TYR	111	10.354	-1.301	-3.121	1.00	0.17
ATOM	1617	CD2	TYR	111	11.127	1.794	-1.936	1.00	0.20
ATOM	1618	HD2	TYR	111	11.201	2.821	-2.261	1.00	0.23
ATOM	1619	CE1	TYR	111	10.933	-0.868	-1.093	1.00	0.19
ATOM	1620	HE1	TYR	111	10.858	-1.895	-0.767	1.00	0.20
ATOM	1621	CE2	TYR	111	11.412	1.459	-0.607	1.00	0.22
ATOM	1622	HE2	TYR	111	11.706	2.227	0.093	1.00	0.26
ATOM	1623	CZ	TYR	111	11.315	0.127	-0.185	1.00	0.21
ATOM	1624	OH	TYR	111	11.595	-0.204	1.125	1.00	0.23
ATOM	1625	HH	TYR	111	12.543	-0.121	1.255	1.00	0.95
ATOM	1626	C	TYR	111	11.374	1.321	-6.588	1.00	0.14
ATOM	1627	O	TYR	111	10.949	2.424	-6.871	1.00	0.15
ATOM	1628	N	ASN	112	11.581	0.421	-7.511	1.00	0.15
ATOM	1629	HN	ASN	112	11.924	-0.464	-7.264	1.00	0.17
ATOM	1630	CA	ASN	112	11.295	0.739	-8.939	1.00	0.16
ATOM	1631	HA	ASN	112	11.870	1.605	-9.235	1.00	0.16
ATOM	1632	CB	ASN	112	11.677	-0.450	-9.822	1.00	0.19
ATOM	1633	HB1	ASN	112	11.025	-1.276	-9.607	1.00	0.22
ATOM	1634	HB2	ASN	112	12.698	-0.739	-9.622	1.00	0.19
ATOM	1635	CG	ASN	112	11.531	-0.060	-11.295	1.00	0.24
ATOM	1636	OD1	ASN	112	10.446	0.248	-11.748	1.00	0.96
ATOM	1637	ND2	ASN	112	12.583	-0.059	-12.067	1.00	1.06
ATOM	1638	HD21	ASN	112	13.458	-0.308	-11.704	1.00	1.80
ATOM	1639	HD22	ASN	112	12.497	0.189	-13.012	1.00	1.08
ATOM	1640	C	ASN	112	9.803	1.040	-9.108	1.00	0.15
ATOM	1641	O	ASN	112	8.953	0.310	-8.637	1.00	0.14
ATOM	1642	N	LEU	113	9.482	2.112	-9.777	1.00	0.15
ATOM	1643	HN	LEU	113	10.187	2.684	-10.145	1.00	0.16
ATOM	1644	CA	LEU	113	8.049	2.475	-9.984	1.00	0.15
ATOM	1645	HA	LEU	113	7.582	2.620	-9.025	1.00	0.14
ATOM	1646	CB	LEU	113	7.981	3.781	-10.791	1.00	0.16
ATOM	1647	HB1	LEU	113	8.513	3.646	-11.721	1.00	0.17
ATOM	1648	HB2	LEU	113	8.452	4.571	-10.226	1.00	0.16
ATOM	1649	CG	LEU	113	6.523	4.177	-11.095	1.00	0.17
ATOM	1650	HG	LEU	113	6.041	3.387	-11.652	1.00	0.18
ATOM	1651	CD1	LEU	113	5.748	4.421	-9.793	1.00	0.18
ATOM	1652	HD11	LEU	113	4.841	4.969	-10.007	1.00	0.99
ATOM	1653	HD12	LEU	113	6.359	4.991	-9.110	1.00	1.00
ATOM	1654	HD13	LEU	113	5.490	3.474	-9.343	1.00	0.97
ATOM	1655	CD2	LEU	113	6.526	5.457	-11.943	1.00	0.20
ATOM	1656	HD21	LEU	113	6.115	6.277	-11.374	1.00	1.06

26/75

ATOM	1657	HD22	LEU	113	5.930	5.302	-12.830	1.00	1.03
ATOM	1658	HD23	LEU	113	7.539	5.696	-12.231	1.00	1.00
ATOM	1659	C	LEU	113	7.320	1.361	-10.743	1.00	0.15
ATOM	1660	O	LEU	113	6.203	1.014	-10.419	1.00	0.15
ATOM	1661	N	PHE	114	7.928	0.817	-11.762	1.00	0.16
ATOM	1662	HN	PHE	114	8.822	1.123	-12.020	1.00	0.17
ATOM	1663	CA	PHE	114	7.245	-0.250	-12.555	1.00	0.17
ATOM	1664	HA	PHE	114	6.338	0.151	-12.980	1.00	0.18
ATOM	1665	CB	PHE	114	8.159	-0.720	-13.685	1.00	0.21
ATOM	1666	HB1	PHE	114	9.077	-1.108	-13.271	1.00	0.22
ATOM	1667	HB2	PHE	114	8.380	0.111	-14.340	1.00	0.22
ATOM	1668	CG	PHE	114	7.457	-1.807	-14.464	1.00	0.24
ATOM	1669	CD1	PHE	114	7.545	-3.135	-14.031	1.00	0.35
ATOM	1670	HD1	PHE	114	8.105	-3.376	-13.147	1.00	0.43
ATOM	1671	CD2	PHE	114	6.724	-1.494	-15.613	1.00	0.24
ATOM	1672	HD2	PHE	114	6.655	-0.470	-15.950	1.00	0.28
ATOM	1673	CE1	PHE	114	6.902	-4.149	-14.741	1.00	0.39
ATOM	1674	HE1	PHE	114	6.975	-5.171	-14.402	1.00	0.50
ATOM	1675	CE2	PHE	114	6.078	-2.512	-16.327	1.00	0.26
ATOM	1676	HE2	PHE	114	5.511	-2.273	-17.214	1.00	0.30
ATOM	1677	CZ	PHE	114	6.168	-3.839	-15.890	1.00	0.32
ATOM	1678	HZ	PHE	114	5.670	-4.623	-16.438	1.00	0.35
ATOM	1679	C	PHE	114	6.900	-1.452	-11.676	1.00	0.17
ATOM	1680	O	PHE	114	5.842	-2.034	-11.806	1.00	0.17
ATOM	1681	N	LEU	115	7.774	-1.846	-10.797	1.00	0.18
ATOM	1682	HN	LEU	115	8.631	-1.380	-10.706	1.00	0.18
ATOM	1683	CA	LEU	115	7.463	-3.028	-9.946	1.00	0.20
ATOM	1684	HA	LEU	115	7.297	-3.882	-10.579	1.00	0.21
ATOM	1685	CB	LEU	115	8.634	-3.304	-8.984	1.00	0.23
ATOM	1686	HB1	LEU	115	8.237	-3.650	-8.041	1.00	0.26
ATOM	1687	HB2	LEU	115	9.172	-2.387	-8.821	1.00	0.22
ATOM	1688	CG	LEU	115	9.612	-4.369	-9.539	1.00	0.28
ATOM	1689	HG	LEU	115	10.397	-4.525	-8.812	1.00	0.33
ATOM	1690	CD1	LEU	115	8.886	-5.702	-9.749	1.00	0.36
ATOM	1691	HD11	LEU	115	9.551	-6.514	-9.498	1.00	0.99
ATOM	1692	HD12	LEU	115	8.578	-5.795	-10.779	1.00	1.11
ATOM	1693	HD13	LEU	115	8.017	-5.740	-9.109	1.00	1.13
ATOM	1694	CD2	LEU	115	10.249	-3.903	-10.859	1.00	0.30
ATOM	1695	HD21	LEU	115	10.497	-4.761	-11.466	1.00	1.10
ATOM	1696	HD22	LEU	115	11.149	-3.351	-10.645	1.00	1.06
ATOM	1697	HD23	LEU	115	9.567	-3.272	-11.395	1.00	1.01
ATOM	1698	C	LEU	115	6.194	-2.748	-9.136	1.00	0.19
ATOM	1699	O	LEU	115	5.280	-3.548	-9.106	1.00	0.20
ATOM	1700	N	VAL	116	6.130	-1.624	-8.475	1.00	0.18
ATOM	1701	HN	VAL	116	6.879	-0.993	-8.508	1.00	0.18
ATOM	1702	CA	VAL	116	4.919	-1.305	-7.664	1.00	0.19
ATOM	1703	HA	VAL	116	4.686	-2.146	-7.028	1.00	0.21
ATOM	1704	CB	VAL	116	5.203	-0.078	-6.794	1.00	0.20
ATOM	1705	HB	VAL	116	5.581	0.722	-7.414	1.00	0.19
ATOM	1706	CG1	VAL	116	3.914	0.381	-6.103	1.00	0.22
ATOM	1707	HG11	VAL	116	3.253	0.832	-6.828	1.00	1.05
ATOM	1708	HG12	VAL	116	4.155	1.105	-5.339	1.00	1.05
ATOM	1709	HG13	VAL	116	3.426	-0.470	-5.650	1.00	1.03
ATOM	1710	CG2	VAL	116	6.246	-0.443	-5.737	1.00	0.21
ATOM	1711	HG21	VAL	116	7.188	-0.654	-6.221	1.00	1.02
ATOM	1712	HG22	VAL	116	5.917	-1.317	-5.194	1.00	0.98
ATOM	1713	HG23	VAL	116	6.370	0.382	-5.052	1.00	1.03
ATOM	1714	C	VAL	116	3.724	-1.020	-8.582	1.00	0.18
ATOM	1715	O	VAL	116	2.615	-1.433	-8.312	1.00	0.19
ATOM	1716	N	ALA	117	3.934	-0.307	-9.659	1.00	0.17
ATOM	1717	HN	ALA	117	4.833	0.028	-9.859	1.00	0.16
ATOM	1718	CA	ALA	117	2.796	0.007	-10.572	1.00	0.17
ATOM	1719	HA	ALA	117	2.064	0.598	-10.044	1.00	0.19
ATOM	1720	CB	ALA	117	3.306	0.795	-11.780	1.00	0.18
ATOM	1721	HB1	ALA	117	4.378	0.709	-11.840	1.00	1.05
ATOM	1722	HB2	ALA	117	3.033	1.834	-11.674	1.00	1.01
ATOM	1723	HB3	ALA	117	2.863	0.397	-12.682	1.00	0.98
ATOM	1724	C	ALA	117	2.150	-1.291	-11.058	1.00	0.17
ATOM	1725	O	ALA	117	0.956	-1.480	-10.951	1.00	0.19
ATOM	1726	N	ALA	118	2.931	-2.187	-11.588	1.00	0.16
ATOM	1727	HN	ALA	118	3.893	-2.015	-11.663	1.00	0.16
ATOM	1728	CA	ALA	118	2.366	-3.472	-12.083	1.00	0.17
ATOM	1729	HA	ALA	118	1.643	-3.273	-12.859	1.00	0.19
ATOM	1730	CB	ALA	118	3.491	-4.335	-12.653	1.00	0.17
ATOM	1731	HB1	ALA	118	3.125	-5.338	-12.812	1.00	1.05
ATOM	1732	HB2	ALA	118	4.316	-4.358	-11.956	1.00	1.02
ATOM	1733	HB3	ALA	118	3.824	-3.824	-12.824	1.00	1.02

27/75

ATOM	1734	C	ALA	118	1.687	-4.220	-10.935	1.00	0.17
ATOM	1735	O	ALA	118	0.699	-4.901	-11.124	1.00	0.18
ATOM	1736	N	HIS	119	2.225	-4.123	-9.751	1.00	0.16
ATOM	1737	HN	HIS	119	3.035	-3.585	-9.623	1.00	0.16
ATOM	1738	CA	HIS	119	1.627	-4.855	-8.599	1.00	0.17
ATOM	1739	HA	HIS	119	1.576	-5.907	-8.833	1.00	0.18
ATOM	1740	CB	HIS	119	2.513	-4.655	-7.368	1.00	0.19
ATOM	1741	HB1	HIS	119	2.547	-3.605	-7.116	1.00	0.19
ATOM	1742	HB2	HIS	119	3.512	-5.005	-7.584	1.00	0.20
ATOM	1743	CG	HIS	119	1.950	-5.431	-6.210	1.00	0.21
ATOM	1744	ND1	HIS	119	2.228	-6.775	-6.020	1.00	0.26
ATOM	1745	HD1	HIS	119	2.791	-7.336	-6.593	1.00	0.30
ATOM	1746	CD2	HIS	119	1.128	-5.067	-5.172	1.00	0.20
ATOM	1747	HD2	HIS	119	0.719	-4.079	-5.019	1.00	0.21
ATOM	1748	CE1	HIS	119	1.585	-7.168	-4.906	1.00	0.27
ATOM	1749	HE1	HIS	119	1.622	-8.171	-4.509	1.00	0.33
ATOM	1750	NE2	HIS	119	0.899	-6.166	-4.350	1.00	0.23
ATOM	1751	C	HIS	119	0.215	-4.333	-8.299	1.00	0.17
ATOM	1752	O	HIS	119	-0.721	-5.101	-8.185	1.00	0.18
ATOM	1753	N	GLU	120	0.043	-3.044	-8.160	1.00	0.18
ATOM	1754	HN	GLU	120	0.801	-2.430	-8.248	1.00	0.18
ATOM	1755	CA	GLU	120	-1.322	-2.520	-7.860	1.00	0.20
ATOM	1756	HA	GLU	120	-1.666	-2.977	-6.943	1.00	0.21
ATOM	1757	CB	GLU	120	-1.294	-0.999	-7.668	1.00	0.22
ATOM	1758	HB1	GLU	120	-0.719	-0.763	-6.785	1.00	0.37
ATOM	1759	HB2	GLU	120	-2.302	-0.635	-7.542	1.00	0.33
ATOM	1760	CG	GLU	120	-0.663	-0.314	-8.875	1.00	0.41
ATOM	1761	HG1	GLU	120	-1.125	-0.668	-9.781	1.00	0.63
ATOM	1762	HG2	GLU	120	0.393	-0.531	-8.895	1.00	0.87
ATOM	1763	CD	GLU	120	-0.875	1.194	-8.757	1.00	0.94
ATOM	1764	OE1	GLU	120	-0.757	1.703	-7.654	1.00	1.67
ATOM	1765	OE2	GLU	120	-1.151	1.816	-9.769	1.00	1.56
ATOM	1766	C	GLU	120	-2.291	-2.903	-8.984	1.00	0.20
ATOM	1767	O	GLU	120	-3.432	-3.238	-8.737	1.00	0.21
ATOM	1768	N	PHE	121	-1.853	-2.872	-10.217	1.00	0.19
ATOM	1769	HN	PHE	121	-0.928	-2.608	-10.405	1.00	0.19
ATOM	1770	CA	PHE	121	-2.767	-3.251	-11.331	1.00	0.21
ATOM	1771	HA	PHE	121	-3.628	-2.600	-11.317	1.00	0.23
ATOM	1772	CB	PHE	121	-2.053	-3.130	-12.685	1.00	0.22
ATOM	1773	HB1	PHE	121	-2.576	-3.726	-13.419	1.00	0.24
ATOM	1774	HB2	PHE	121	-1.041	-3.493	-12.587	1.00	0.21
ATOM	1775	CG	PHE	121	-2.026	-1.684	-13.141	1.00	0.25
ATOM	1776	CD1	PHE	121	-0.804	-1.019	-13.308	1.00	0.27
ATOM	1777	HD1	PHE	121	0.121	-1.535	-13.113	1.00	0.40
ATOM	1778	CD2	PHE	121	-3.227	-1.007	-13.403	1.00	0.45
ATOM	1779	HD2	PHE	121	-4.173	-1.513	-13.281	1.00	0.60
ATOM	1780	CE1	PHE	121	-0.781	0.314	-13.733	1.00	0.29
ATOM	1781	HE1	PHE	121	0.163	0.824	-13.862	1.00	0.39
ATOM	1782	CE2	PHE	121	-3.202	0.327	-13.828	1.00	0.49
ATOM	1783	HE2	PHE	121	-4.127	0.847	-14.029	1.00	0.68
ATOM	1784	CZ	PHE	121	-1.979	0.988	-13.993	1.00	0.34
ATOM	1785	HZ	PHE	121	-1.961	2.017	-14.321	1.00	0.38
ATOM	1786	C	PHE	121	-3.228	-4.693	-11.120	1.00	0.20
ATOM	1787	O	PHE	121	-4.374	-5.027	-11.344	1.00	0.21
ATOM	1788	N	GLY	122	-2.344	-5.551	-10.690	1.00	0.18
ATOM	1789	HN	GLY	122	-1.424	-5.262	-10.514	1.00	0.17
ATOM	1790	CA	GLY	122	-2.737	-6.970	-10.464	1.00	0.20
ATOM	1791	HA1	GLY	122	-1.890	-7.523	-10.092	1.00	0.21
ATOM	1792	HA2	GLY	122	-3.072	-7.404	-11.394	1.00	0.21
ATOM	1793	C	GLY	122	-3.867	-7.022	-9.435	1.00	0.20
ATOM	1794	O	GLY	122	-4.823	-7.756	-9.589	1.00	0.22
ATOM	1795	N	HIS	123	-3.778	-6.240	-8.392	1.00	0.20
ATOM	1796	HN	HIS	123	-3.005	-5.644	-8.287	1.00	0.20
ATOM	1797	CA	HIS	123	-4.864	-6.243	-7.371	1.00	0.22
ATOM	1798	HA	HIS	123	-5.047	-7.255	-7.042	1.00	0.23
ATOM	1799	CB	HIS	123	-4.456	-5.382	-6.174	1.00	0.25
ATOM	1800	HB1	HIS	123	-5.324	-5.180	-5.564	1.00	0.30
ATOM	1801	HB2	HIS	123	-4.041	-4.449	-6.527	1.00	0.25
ATOM	1802	CG	HIS	123	-3.427	-6.108	-5.354	1.00	0.27
ATOM	1803	ND1	HIS	123	-3.736	-7.247	-4.628	1.00	0.37
ATOM	1804	HD1	HIS	123	-4.611	-7.685	-4.581	1.00	0.45
ATOM	1805	CD2	HIS	123	-2.096	-5.866	-5.125	1.00	0.25
ATOM	1806	HD2	HIS	123	-1.532	-5.046	-5.545	1.00	0.27
ATOM	1807	CE1	HIS	123	-2.614	-7.644	-4.001	1.00	0.38
ATOM	1808	HE1	HIS	123	-2.553	-8.514	-3.367	1.00	0.47
ATOM	1809	NE2	HIS	123	-1.584	-6.837	-4.269	1.00	0.29
ATOM	1810	C	HIS	123	-6.137	-6.671	-7.882	1.00	0.33

28/75

ATOM	1811	O	HIS	123	-7.229	-6.148	-7.755	1.00	0.25
ATOM	1812	N	SER	124	-6.002	-4.646	-8.788	1.00	0.23
ATOM	1813	HN	SER	124	-5.110	-4.278	-8.962	1.00	0.22
ATOM	1814	CA	SER	124	-7.196	-4.030	-9.429	1.00	0.25
ATOM	1815	HA	SER	124	-7.928	-3.790	-8.672	1.00	0.27
ATOM	1816	CB	SER	124	-6.778	-2.751	-10.156	1.00	0.27
ATOM	1817	HB1	SER	124	-6.219	-2.119	-9.478	1.00	0.29
ATOM	1818	HB2	SER	124	-7.654	-2.224	-10.494	1.00	0.29
ATOM	1819	OG	SER	124	-5.975	-3.091	-11.279	1.00	0.25
ATOM	1820	HG	SER	124	-6.545	-3.131	-12.050	1.00	0.88
ATOM	1821	C	SER	124	-7.805	-5.006	-10.437	1.00	0.24
ATOM	1822	O	SER	124	-8.975	-4.932	-10.755	1.00	0.26
ATOM	1823	N	LEU	125	-7.022	-5.913	-10.952	1.00	0.22
ATOM	1824	HN	LEU	125	-6.078	-5.953	-10.690	1.00	0.21
ATOM	1825	CA	LEU	125	-7.562	-6.879	-11.949	1.00	0.23
ATOM	1826	HA	LEU	125	-8.285	-6.374	-12.568	1.00	0.24
ATOM	1827	CB	LEU	125	-6.420	-7.398	-12.827	1.00	0.22
ATOM	1828	HB1	LEU	125	-6.759	-8.247	-13.398	1.00	0.24
ATOM	1829	HB2	LEU	125	-5.594	-7.698	-12.197	1.00	0.22
ATOM	1830	CG	LEU	125	-5.956	-6.280	-13.779	1.00	0.22
ATOM	1831	HG	LEU	125	-5.928	-5.343	-13.241	1.00	0.24
ATOM	1832	CD1	LEU	125	-4.556	-6.601	-14.302	1.00	0.25
ATOM	1833	HD11	LEU	125	-4.588	-7.515	-14.874	1.00	0.99
ATOM	1834	HD12	LEU	125	-3.879	-6.719	-13.471	1.00	1.00
ATOM	1835	HD13	LEU	125	-4.215	-5.794	-14.933	1.00	1.05
ATOM	1836	CD2	LEU	125	-6.913	-6.155	-14.976	1.00	0.24
ATOM	1837	HD21	LEU	125	-7.793	-5.604	-14.682	1.00	1.05
ATOM	1838	HD22	LEU	125	-7.201	-7.135	-15.324	1.00	1.00
ATOM	1839	HD23	LEU	125	-6.415	-5.627	-15.775	1.00	1.03
ATOM	1840	C	LEU	125	-8.256	-8.044	-11.234	1.00	0.24
ATOM	1841	O	LEU	125	-8.790	-8.935	-11.864	1.00	0.33
ATOM	1842	N	GLY	126	-8.277	-8.035	-9.927	1.00	0.24
ATOM	1843	HN	GLY	126	-7.858	-7.298	-9.435	1.00	0.29
ATOM	1844	CA	GLY	126	-8.968	-9.132	-9.185	1.00	0.27
ATOM	1845	HA1	GLY	126	-9.748	-9.545	-9.807	1.00	0.29
ATOM	1846	HA2	GLY	126	-9.408	-8.727	-8.285	1.00	0.29
ATOM	1847	C	GLY	126	-7.985	-10.245	-8.809	1.00	0.26
ATOM	1848	O	GLY	126	-8.377	-11.268	-8.283	1.00	0.30
ATOM	1849	N	LEU	127	-6.719	-10.068	-9.063	1.00	0.23
ATOM	1850	HN	LEU	127	-6.410	-9.239	-9.484	1.00	0.22
ATOM	1851	CA	LEU	127	-5.744	-11.138	-8.700	1.00	0.25
ATOM	1852	HA	LEU	127	-6.212	-12.099	-8.815	1.00	0.28
ATOM	1853	CB	LEU	127	-4.507	-11.052	-9.602	1.00	0.23
ATOM	1854	HB1	LEU	127	-3.733	-11.696	-9.211	1.00	0.25
ATOM	1855	HB2	LEU	127	-4.156	-10.033	-9.602	1.00	0.22
ATOM	1856	CG	LEU	127	-4.844	-11.471	-11.045	1.00	0.24
ATOM	1857	HG	LEU	127	-5.707	-10.915	-11.384	1.00	0.23
ATOM	1858	CD1	LEU	127	-3.646	-11.159	-11.962	1.00	0.24
ATOM	1859	HD11	LEU	127	-4.001	-10.692	-12.868	1.00	1.00
ATOM	1860	HD12	LEU	127	-3.126	-12.073	-12.208	1.00	1.02
ATOM	1861	HD13	LEU	127	-2.962	-10.491	-11.460	1.00	1.03
ATOM	1862	CD2	LEU	127	-5.150	-12.980	-11.109	1.00	0.30
ATOM	1863	HD21	LEU	127	-5.021	-13.334	-12.121	1.00	1.04
ATOM	1864	HD22	LEU	127	-6.169	-13.159	-10.805	1.00	1.11
ATOM	1865	HD23	LEU	127	-4.478	-13.515	-10.454	1.00	1.03
ATOM	1866	C	LEU	127	-5.315	-10.969	-7.241	1.00	0.28
ATOM	1867	O	LEU	127	-5.245	-9.872	-6.723	1.00	0.32
ATOM	1868	N	ASP	128	-5.027	-12.059	-6.581	1.00	0.32
ATOM	1869	HN	ASP	128	-5.093	-12.928	-7.029	1.00	0.34
ATOM	1870	CA	ASP	128	-4.598	-11.997	-5.154	1.00	0.39
ATOM	1871	HA	ASP	128	-4.882	-11.046	-4.728	1.00	0.40
ATOM	1872	CB	ASP	128	-5.271	-13.130	-4.375	1.00	0.48
ATOM	1873	HB1	ASP	128	-4.779	-14.064	-4.600	1.00	0.48
ATOM	1874	HB2	ASP	128	-6.311	-13.193	-4.661	1.00	0.50
ATOM	1875	CG	ASP	128	-5.171	-12.854	-2.873	1.00	0.55
ATOM	1876	OD1	ASP	128	-4.082	-12.980	-2.339	1.00	1.23
ATOM	1877	OD2	ASP	128	-6.185	-12.521	-2.283	1.00	1.22
ATOM	1878	C	ASP	128	-3.078	-12.159	-5.082	1.00	0.37
ATOM	1879	O	ASP	128	-2.424	-12.387	-6.080	1.00	0.59
ATOM	1880	N	HIS	129	-2.507	-12.042	-3.914	1.00	0.23
ATOM	1881	HN	HIS	129	-3.048	-11.856	-3.118	1.00	0.32
ATOM	1882	CA	HIS	129	-1.029	-12.189	-3.797	1.00	0.22
ATOM	1883	HA	HIS	129	-0.543	-11.439	-4.401	1.00	0.21
ATOM	1884	CB	HIS	129	-0.606	-12.019	-2.335	1.00	0.23
ATOM	1885	HB1	HIS	129	0.430	-12.302	-2.227	1.00	0.24
ATOM	1886	HB2	HIS	129	-1.217	-12.653	-1.710	1.00	0.25
ATOM	1887	CG	HIS	129	-0.779	-10.585	-1.912	1.00	0.22

29/75

ATOM	1888	ND1	HIS	129	-1.862	-10.161	-1.156	1.00	0.35
ATOM	1889	HD1	HIS	129	-2.602	-10.720	-0.841	1.00	0.53
ATOM	1890	CD2	HIS	129	-0.007	-9.468	-2.116	1.00	0.34
ATOM	1891	HD2	HIS	129	0.918	-9.447	-2.673	1.00	0.54
ATOM	1892	CE1	HIS	129	-1.711	-8.842	-0.936	1.00	0.31
ATOM	1893	HE1	HIS	129	-2.406	-8.239	-0.370	1.00	0.44
ATOM	1894	NE2	HIS	129	-0.597	-8.369	-1.501	1.00	0.28
ATOM	1895	C	HIS	129	-0.614	-13.584	-4.277	1.00	0.24
ATOM	1896	O	HIS	129	-1.267	-14.568	-3.991	1.00	0.28
ATOM	1897	N	SER	130	0.474	-13.671	-4.999	1.00	0.24
ATOM	1898	HN	SER	130	0.984	-12.862	-5.210	1.00	0.23
ATOM	1899	CA	SER	130	0.949	-14.996	-5.498	1.00	0.29
ATOM	1900	HA	SER	130	0.139	-15.710	-5.464	1.00	0.33
ATOM	1901	CB	SER	130	1.442	-14.852	-6.938	1.00	0.32
ATOM	1902	HB1	SER	130	2.201	-14.082	-6.982	1.00	0.31
ATOM	1903	HB2	SER	130	0.618	-14.577	-7.576	1.00	0.35
ATOM	1904	OG	SER	130	1.980	-16.092	-7.378	1.00	0.40
ATOM	1905	HG	SER	130	1.254	-16.714	-7.469	1.00	0.97
ATOM	1906	C	SER	130	2.096	-15.484	-4.609	1.00	0.28
ATOM	1907	O	SER	130	2.801	-14.696	-4.009	1.00	0.29
ATOM	1908	N	LYS	131	2.287	-16.775	-4.514	1.00	0.30
ATOM	1909	HN	LYS	131	1.705	-17.393	-5.003	1.00	0.32
ATOM	1910	CA	LYS	131	3.386	-17.310	-3.656	1.00	0.32
ATOM	1911	HA	LYS	131	3.665	-16.567	-2.923	1.00	0.34
ATOM	1912	CB	LYS	131	2.903	-18.572	-2.936	1.00	0.39
ATOM	1913	HB1	LYS	131	3.714	-18.988	-2.355	1.00	0.42
ATOM	1914	HB2	LYS	131	2.572	-19.298	-3.664	1.00	0.40
ATOM	1915	CG	LYS	131	1.743	-18.214	-2.003	1.00	0.45
ATOM	1916	HG1	LYS	131	0.932	-17.798	-2.581	1.00	0.79
ATOM	1917	HG2	LYS	131	2.077	-17.488	-1.276	1.00	1.01
ATOM	1918	CD	LYS	131	1.255	-19.472	-1.280	1.00	1.18
ATOM	1919	HD1	LYS	131	2.064	-19.890	-0.698	1.00	1.86
ATOM	1920	HD2	LYS	131	0.921	-20.199	-2.006	1.00	1.66
ATOM	1921	CE	LYS	131	0.096	-19.108	-0.349	1.00	1.52
ATOM	1922	HE1	LYS	131	-0.788	-18.908	-0.937	1.00	1.92
ATOM	1923	HE2	LYS	131	0.355	-18.229	0.222	1.00	1.93
ATOM	1924	NZ	LYS	131	-0.174	-20.242	0.581	1.00	2.23
ATOM	1925	HZ1	LYS	131	-1.103	-20.109	1.030	1.00	2.72
ATOM	1926	HZ2	LYS	131	0.565	-20.272	1.313	1.00	2.53
ATOM	1927	HZ3	LYS	131	-0.174	-21.135	0.050	1.00	2.72
ATOM	1928	C	LYS	131	4.604	-17.649	-4.521	1.00	0.31
ATOM	1929	O	LYS	131	5.612	-18.116	-4.027	1.00	0.34
ATOM	1930	N	ASP	132	4.532	-17.411	-5.804	1.00	0.29
ATOM	1931	HN	ASP	132	3.717	-17.028	-6.190	1.00	0.28
ATOM	1932	CA	ASP	132	5.703	-17.719	-6.674	1.00	0.30
ATOM	1933	HA	ASP	132	6.187	-18.601	-6.302	1.00	0.32
ATOM	1934	CB	ASP	132	5.225	-17.970	-8.108	1.00	0.32
ATOM	1935	HB1	ASP	132	4.727	-17.090	-8.483	1.00	0.31
ATOM	1936	HB2	ASP	132	4.539	-18.804	-8.118	1.00	0.34
ATOM	1937	CG	ASP	132	6.430	-18.289	-8.996	1.00	0.35
ATOM	1938	OD1	ASP	132	6.457	-19.371	-9.558	1.00	1.10
ATOM	1939	OD2	ASP	132	7.306	-17.446	-9.097	1.00	1.15
ATOM	1940	C	ASP	132	6.656	-16.501	-6.659	1.00	0.28
ATOM	1941	O	ASP	132	6.226	-15.399	-6.939	1.00	0.28
ATOM	1942	N	PRO	133	7.930	-16.658	-6.328	1.00	0.30
ATOM	1943	CA	PRO	133	8.852	-15.484	-6.296	1.00	0.31
ATOM	1944	HA	PRO	133	8.517	-14.766	-5.566	1.00	0.32
ATOM	1945	CB	PRO	133	10.173	-16.097	-5.832	1.00	0.36
ATOM	1946	HB1	PRO	133	10.441	-15.694	-4.867	1.00	0.36
ATOM	1947	HB2	PRO	133	10.949	-15.869	-6.549	1.00	0.41
ATOM	1948	CG	PRO	133	10.007	-17.615	-5.721	1.00	0.42
ATOM	1949	HG1	PRO	133	10.293	-17.940	-4.732	1.00	0.51
ATOM	1950	HG2	PRO	133	10.630	-18.103	-6.457	1.00	0.51
ATOM	1951	CD	PRO	133	8.540	-17.972	-5.969	1.00	0.35
ATOM	1952	HD2	PRO	133	8.456	-18.679	-6.785	1.00	0.34
ATOM	1953	HD1	PRO	133	8.091	-18.362	-5.069	1.00	0.38
ATOM	1954	C	PRO	133	9.032	-14.810	-7.662	1.00	0.31
ATOM	1955	O	PRO	133	9.496	-13.691	-7.749	1.00	0.34
ATOM	1956	N	GLY	134	8.684	-15.477	-8.729	1.00	0.32
ATOM	1957	HN	GLY	134	8.320	-16.382	-8.647	1.00	0.35
ATOM	1958	CA	GLY	134	8.860	-14.856	-10.074	1.00	0.34
ATOM	1959	HA1	GLY	134	9.048	-15.630	-10.803	1.00	0.37
ATOM	1960	HA2	GLY	134	9.701	-14.177	-10.047	1.00	0.36
ATOM	1961	C	GLY	134	7.598	-14.087	-10.471	1.00	0.29
ATOM	1962	O	GLY	134	7.563	-13.420	-11.486	1.00	0.29
ATOM	1963	N	ALA	135	6.563	-14.168	-9.683	1.00	0.27
ATOM	1964	HN	ALA	135	6.607	-14.709	-8.867	1.00	0.28

30/75

ATOM	1965	CA	ALA	135	5.312	-13.434	-10.026	1.00	0.24
ATOM	1966	HA	ALA	135	5.199	-13.401	-11.099	1.00	0.25
ATOM	1967	CB	ALA	135	4.109	-14.151	-9.410	1.00	0.25
ATOM	1968	HB1	ALA	135	3.633	-14.765	-10.160	1.00	1.07
ATOM	1969	HB2	ALA	135	3.405	-13.421	-9.041	1.00	1.01
ATOM	1970	HB3	ALA	135	4.442	-14.774	-8.593	1.00	1.04
ATOM	1971	C	ALA	135	5.388	-12.007	-9.479	1.00	0.21
ATOM	1972	O	ALA	135	5.968	-11.760	-8.440	1.00	0.23
ATOM	1973	N	LEU	136	4.799	-11.067	-10.164	1.00	0.22
ATOM	1974	HN	LEU	136	4.330	-11.286	-10.996	1.00	0.24
ATOM	1975	CA	LEU	136	4.830	-9.660	-9.676	1.00	0.23
ATOM	1976	HA	LEU	136	5.842	-9.382	-9.427	1.00	0.25
ATOM	1977	CB	LEU	136	4.279	-8.724	-10.761	1.00	0.25
ATOM	1978	HB1	LEU	136	4.193	-7.724	-10.365	1.00	0.27
ATOM	1979	HB2	LEU	136	3.302	-9.072	-11.064	1.00	0.26
ATOM	1980	CG	LEU	136	5.213	-8.709	-11.980	1.00	0.26
ATOM	1981	HG	LEU	136	5.312	-9.713	-12.368	1.00	0.29
ATOM	1982	CD1	LEU	136	4.624	-7.801	-13.063	1.00	0.29
ATOM	1983	HD11	LEU	136	3.546	-7.848	-13.030	1.00	1.06
ATOM	1984	HD12	LEU	136	4.967	-8.126	-14.033	1.00	1.05
ATOM	1985	HD13	LEU	136	4.944	-6.784	-12.893	1.00	1.06
ATOM	1986	CD2	LEU	136	6.592	-8.176	-11.578	1.00	0.32
ATOM	1987	HD21	LEU	136	6.485	-7.477	-10.762	1.00	1.05
ATOM	1988	HD22	LEU	136	7.046	-7.677	-12.422	1.00	1.09
ATOM	1989	HD23	LEU	136	7.220	-8.998	-11.269	1.00	0.97
ATOM	1990	C	LEU	136	3.954	-9.556	-8.427	1.00	0.25
ATOM	1991	O	LEU	136	4.201	-8.761	-7.542	1.00	0.30
ATOM	1992	N	MET	137	2.924	-10.353	-8.357	1.00	0.28
ATOM	1993	HN	MET	137	2.744	-10.981	-9.087	1.00	0.31
ATOM	1994	CA	MET	137	2.016	-10.309	-7.177	1.00	0.33
ATOM	1995	HA	MET	137	1.768	-9.283	-6.959	1.00	0.38
ATOM	1996	CB	MET	137	0.734	-11.087	-7.494	1.00	0.42
ATOM	1997	HB1	MET	137	0.118	-11.136	-6.615	1.00	0.57
ATOM	1998	HB2	MET	137	0.995	-12.089	-7.803	1.00	0.50
ATOM	1999	CG	MET	137	-0.035	-10.391	-8.625	1.00	0.58
ATOM	2000	HG1	MET	137	-0.909	-10.975	-8.875	1.00	1.13
ATOM	2001	HG2	MET	137	0.601	-10.311	-9.494	1.00	1.22
ATOM	2002	SD	MET	137	-0.551	-8.729	-8.108	1.00	0.83
ATOM	2003	CE	MET	137	-2.048	-9.184	-7.194	1.00	0.39
ATOM	2004	HE1	MET	137	-2.231	-8.450	-6.426	1.00	1.14
ATOM	2005	HE2	MET	137	-1.927	-10.151	-6.741	1.00	1.07
ATOM	2006	HE3	MET	137	-2.885	-9.212	-7.872	1.00	1.06
ATOM	2007	C	MET	137	2.700	-10.925	-5.951	1.00	0.27
ATOM	2008	O	MET	137	2.050	-11.287	-4.990	1.00	0.28
ATOM	2009	N	PHE	138	4.000	-11.042	-5.964	1.00	0.25
ATOM	2010	HN	PHE	138	4.514	-10.741	-6.743	1.00	0.28
ATOM	2011	CA	PHE	138	4.699	-11.628	-4.785	1.00	0.23
ATOM	2012	HA	PHE	138	4.225	-12.557	-4.534	1.00	0.26
ATOM	2013	CB	PHE	138	6.167	-11.877	-5.152	1.00	0.25
ATOM	2014	HB1	PHE	138	6.710	-10.945	-5.104	1.00	0.24
ATOM	2015	HB2	PHE	138	6.221	-12.270	-6.156	1.00	0.27
ATOM	2016	CG	PHE	138	6.790	-12.873	-4.194	1.00	0.28
ATOM	2017	CD1	PHE	138	6.295	-14.184	-4.113	1.00	0.32
ATOM	2018	HD1	PHE	138	5.465	-14.490	-4.731	1.00	0.33
ATOM	2019	CD2	PHE	138	7.871	-12.486	-3.392	1.00	0.30
ATOM	2020	HD2	PHE	138	8.256	-11.481	-3.455	1.00	0.30
ATOM	2021	CE1	PHE	138	6.881	-15.100	-3.230	1.00	0.38
ATOM	2022	HE1	PHE	138	6.500	-16.109	-3.168	1.00	0.42
ATOM	2023	CE2	PHE	138	8.455	-13.404	-2.511	1.00	0.36
ATOM	2024	HE2	PHE	138	9.288	-13.104	-1.894	1.00	0.39
ATOM	2025	CZ	PHE	138	7.960	-14.710	-2.430	1.00	0.39
ATOM	2026	HZ	PHE	138	8.411	-15.417	-1.749	1.00	0.44
ATOM	2027	C	PHE	138	4.601	-10.615	-3.615	1.00	0.20
ATOM	2028	O	PHE	138	4.874	-9.447	-3.808	1.00	0.22
ATOM	2029	N	PRO	139	4.185	-11.019	-2.421	1.00	0.22
ATOM	2030	CA	PRO	139	4.044	-10.048	-1.291	1.00	0.25
ATOM	2031	HA	PRO	139	3.262	-9.340	-1.509	1.00	0.27
ATOM	2032	CB	PRO	139	3.600	-10.936	-0.127	1.00	0.31
ATOM	2033	HB1	PRO	139	2.615	-10.638	0.199	1.00	0.38
ATOM	2034	HB2	PRO	139	4.299	-10.835	0.691	1.00	0.42
ATOM	2035	CG	PRO	139	3.562	-12.392	-0.597	1.00	0.33
ATOM	2036	HG1	PRO	139	2.588	-12.812	-0.396	1.00	0.41
ATOM	2037	HG2	PRO	139	4.317	-12.961	-0.074	1.00	0.42
ATOM	2038	CD	PRO	139	3.834	-12.435	-2.102	1.00	0.27
ATOM	2039	HD2	PRO	139	4.661	-13.100	-2.318	1.00	0.28
ATOM	2040	HD1	PRO	139	2.946	-12.732	-2.637	1.00	0.30
ATOM	2041	C	PRO	139	5.312	-13.434	-10.026	1.00	0.24

31/75

ATOM	2042	O	PRO	139	5.302	-8.351	-0.173	1.00	0.44
ATOM	2043	N	ILE	140	6.467	-9.726	-1.437	1.00	0.24
ATOM	2044	HN	ILE	140	6.474	-10.500	-2.038	1.00	0.37
ATOM	2045	CA	ILE	140	7.749	-9.031	-1.094	1.00	0.23
ATOM	2046	HA	ILE	140	7.572	-8.308	-0.312	1.00	0.24
ATOM	2047	CB	ILE	140	8.775	-10.054	-0.600	1.00	0.25
ATOM	2048	HB	ILE	140	8.978	-10.770	-1.379	1.00	0.25
ATOM	2049	CG1	ILE	140	8.207	-10.768	0.632	1.00	0.29
ATOM	2050	HG11	ILE	140	7.246	-11.196	0.384	1.00	0.32
ATOM	2051	HG12	ILE	140	8.084	-10.055	1.434	1.00	0.33
ATOM	2052	CG2	ILE	140	10.070	-9.332	-0.214	1.00	0.26
ATOM	2053	HG21	ILE	140	9.850	-8.567	0.517	1.00	1.04
ATOM	2054	HG22	ILE	140	10.505	-8.876	-1.090	1.00	1.06
ATOM	2055	HG23	ILE	140	10.768	-10.040	0.207	1.00	1.04
ATOM	2056	CD1	ILE	140	9.156	-11.883	1.082	1.00	0.30
ATOM	2057	HD11	ILE	140	9.716	-12.250	0.236	1.00	1.08
ATOM	2058	HD12	ILE	140	8.582	-12.691	1.511	1.00	0.98
ATOM	2059	HD13	ILE	140	9.838	-11.495	1.824	1.00	1.08
ATOM	2060	C	ILE	140	8.284	-8.301	-2.329	1.00	0.22
ATOM	2061	O	ILE	140	8.265	-8.817	-3.429	1.00	0.22
ATOM	2062	N	TYR	141	8.745	-7.092	-2.150	1.00	0.21
ATOM	2063	HN	TYR	141	8.736	-6.696	-1.254	1.00	0.22
ATOM	2064	CA	TYR	141	9.265	-6.303	-3.304	1.00	0.21
ATOM	2065	HA	TYR	141	8.560	-6.348	-4.120	1.00	0.20
ATOM	2066	CB	TYR	141	9.444	-4.847	-2.865	1.00	0.21
ATOM	2067	HB1	TYR	141	10.050	-4.810	-1.972	1.00	0.22
ATOM	2068	HB2	TYR	141	8.476	-4.413	-2.661	1.00	0.22
ATOM	2069	CG	TYR	141	10.122	-4.066	-3.962	1.00	0.23
ATOM	2070	CD1	TYR	141	11.515	-4.104	-4.089	1.00	0.25
ATOM	2071	HD1	TYR	141	12.104	-4.697	-3.404	1.00	0.26
ATOM	2072	CD2	TYR	141	9.359	-3.298	-4.848	1.00	0.24
ATOM	2073	HD2	TYR	141	8.284	-3.268	-4.750	1.00	0.25
ATOM	2074	CE1	TYR	141	12.146	-3.376	-5.103	1.00	0.28
ATOM	2075	HE1	TYR	141	13.221	-3.405	-5.201	1.00	0.32
ATOM	2076	CE2	TYR	141	9.989	-2.569	-5.862	1.00	0.27
ATOM	2077	HE2	TYR	141	9.401	-1.975	-6.544	1.00	0.30
ATOM	2078	CZ	TYR	141	11.383	-2.608	-5.990	1.00	0.29
ATOM	2079	OH	TYR	141	12.005	-1.892	-6.991	1.00	0.33
ATOM	2080	HH	TYR	141	12.781	-2.385	-7.269	1.00	0.90
ATOM	2081	C	TYR	141	10.615	-6.864	-3.761	1.00	0.22
ATOM	2082	O	TYR	141	11.522	-7.050	-2.973	1.00	0.23
ATOM	2083	N	THR	142	10.750	-7.130	-5.035	1.00	0.22
ATOM	2084	HN	THR	142	10.002	-6.968	-5.648	1.00	0.22
ATOM	2085	CA	THR	142	12.035	-7.675	-5.563	1.00	0.24
ATOM	2086	HA	THR	142	12.835	-7.447	-4.874	1.00	0.25
ATOM	2087	CB	THR	142	11.917	-9.193	-5.723	1.00	0.25
ATOM	2088	HB	THR	142	11.645	-9.635	-4.777	1.00	0.26
ATOM	2089	OG1	THR	142	13.165	-9.720	-6.152	1.00	0.29
ATOM	2090	HG1	THR	142	13.274	-9.505	-7.081	1.00	0.97
ATOM	2091	CG2	THR	142	10.840	-9.517	-6.760	1.00	0.25
ATOM	2092	HG21	THR	142	10.577	-10.562	-6.691	1.00	1.04
ATOM	2093	HG22	THR	142	11.217	-9.304	-7.749	1.00	1.05
ATOM	2094	HG23	THR	142	9.965	-8.913	-6.570	1.00	1.06
ATOM	2095	C	THR	142	12.339	-7.040	-6.924	1.00	0.23
ATOM	2096	O	THR	142	11.454	-6.810	-7.724	1.00	0.23
ATOM	2097	N	TYR	143	13.586	-6.758	-7.195	1.00	0.25
ATOM	2098	HN	TYR	143	14.285	-6.955	-6.538	1.00	0.27
ATOM	2099	CA	TYR	143	13.948	-6.144	-8.506	1.00	0.26
ATOM	2100	HA	TYR	143	13.174	-5.452	-8.804	1.00	0.25
ATOM	2101	CB	TYR	143	15.277	-5.395	-8.370	1.00	0.29
ATOM	2102	HB1	TYR	143	16.072	-6.104	-8.190	1.00	0.33
ATOM	2103	HB2	TYR	143	15.217	-4.704	-7.542	1.00	0.30
ATOM	2104	CG	TYR	143	15.563	-4.633	-9.642	1.00	0.27
ATOM	2105	CD1	TYR	143	14.931	-3.406	-9.880	1.00	0.25
ATOM	2106	HD1	TYR	143	14.234	-3.008	-9.156	1.00	0.26
ATOM	2107	CD2	TYR	143	16.466	-5.148	-10.581	1.00	0.31
ATOM	2108	HD2	TYR	143	16.954	-6.094	-10.398	1.00	0.35
ATOM	2109	CE1	TYR	143	15.201	-2.695	-11.055	1.00	0.26
ATOM	2110	HE1	TYR	143	14.713	-1.749	-11.238	1.00	0.28
ATOM	2111	CE2	TYR	143	16.735	-4.436	-11.756	1.00	0.31
ATOM	2112	HE2	TYR	143	17.432	-4.833	-12.480	1.00	0.36
ATOM	2113	CZ	TYR	143	16.103	-3.210	-11.994	1.00	0.28
ATOM	2114	OH	TYR	143	16.369	-2.509	-13.152	1.00	0.30
ATOM	2115	HH	TYR	143	17.068	-2.969	-13.624	1.00	0.95
ATOM	2116	C	TYR	143	14.080	-7.244	-9.563	1.00	0.27
ATOM	2117	O	TYR	143	14.552	-8.328	-9.283	1.00	0.31
ATOM	2118	N	THR	144	13.650	-6.876	-8.770	1.00	0.22

32/75

ATOM	2119	HN	THR	144	13.277	-6.096	-10.972	1.00	0.32
ATOM	2120	CA	THR	144	13.753	-8.008	-11.847	1.00	0.32
ATOM	2121	HA	THR	144	14.479	-8.758	-11.573	1.00	0.35
ATOM	2122	CB	THR	144	12.385	-8.666	-12.031	1.00	0.37
ATOM	2123	HB	THR	144	11.922	-8.814	-11.067	1.00	0.84
ATOM	2124	OG1	THR	144	12.549	-9.918	-12.683	1.00	1.00
ATOM	2125	HG1	THR	144	13.280	-9.836	-13.301	1.00	1.42
ATOM	2126	CG2	THR	144	11.499	-7.757	-12.882	1.00	0.82
ATOM	2127	HG21	THR	144	10.461	-7.991	-12.699	1.00	1.51
ATOM	2128	HG22	THR	144	11.724	-7.911	-13.927	1.00	1.24
ATOM	2129	HG23	THR	144	11.687	-6.726	-12.622	1.00	1.49
ATOM	2130	C	THR	144	14.169	-7.351	-13.165	1.00	0.34
ATOM	2131	O	THR	144	13.922	-6.183	-13.392	1.00	0.32
ATOM	2132	N	GLY	145	14.789	-8.094	-14.043	1.00	0.43
ATOM	2133	HN	GLY	145	14.971	-9.037	-13.846	1.00	0.49
ATOM	2134	CA	GLY	145	15.205	-7.510	-15.350	1.00	0.49
ATOM	2135	HA1	GLY	145	15.842	-8.207	-15.872	1.00	0.57
ATOM	2136	HA2	GLY	145	15.742	-6.587	-15.178	1.00	0.50
ATOM	2137	C	GLY	145	13.957	-7.233	-16.191	1.00	0.47
ATOM	2138	O	GLY	145	13.331	-8.138	-16.706	1.00	0.53
ATOM	2139	N	LYS	146	13.583	-5.990	-16.322	1.00	0.46
ATOM	2140	HN	LYS	146	14.097	-5.277	-15.889	1.00	0.48
ATOM	2141	CA	LYS	146	12.367	-5.653	-17.116	1.00	0.49
ATOM	2142	HA	LYS	146	11.578	-6.350	-16.876	1.00	0.51
ATOM	2143	CB	LYS	146	11.911	-4.235	-16.764	1.00	0.52
ATOM	2144	HB1	LYS	146	10.973	-4.032	-17.254	1.00	0.58
ATOM	2145	HB2	LYS	146	12.657	-3.533	-17.103	1.00	0.57
ATOM	2146	CG	LYS	146	11.744	-4.128	-15.238	1.00	0.55
ATOM	2147	HG1	LYS	146	12.690	-3.853	-14.798	1.00	0.83
ATOM	2148	HG2	LYS	146	11.442	-5.089	-14.849	1.00	1.14
ATOM	2149	CD	LYS	146	10.684	-3.077	-14.854	1.00	1.23
ATOM	2150	HD1	LYS	146	10.308	-3.309	-13.871	1.00	1.78
ATOM	2151	HD2	LYS	146	9.865	-3.098	-15.556	1.00	1.79
ATOM	2152	CE	LYS	146	11.298	-1.671	-14.828	1.00	2.01
ATOM	2153	HE1	LYS	146	11.615	-1.439	-13.822	1.00	2.47
ATOM	2154	HE2	LYS	146	10.556	-0.952	-15.143	1.00	2.39
ATOM	2155	NZ	LYS	146	12.468	-1.601	-15.745	1.00	2.91
ATOM	2156	HZ1	LYS	146	12.847	-0.633	-15.750	1.00	3.39
ATOM	2157	HZ2	LYS	146	12.170	-1.861	-16.707	1.00	3.28
ATOM	2158	HZ3	LYS	146	13.205	-2.257	-15.420	1.00	3.27
ATOM	2159	C	LYS	146	12.677	-5.732	-18.613	1.00	0.59
ATOM	2160	O	LYS	146	11.845	-5.426	-19.444	1.00	1.16
ATOM	2161	N	SER	147	13.868	-6.131	-18.967	1.00	0.78
ATOM	2162	HN	SER	147	14.530	-6.366	-18.283	1.00	1.26
ATOM	2163	CA	SER	147	14.226	-6.214	-20.413	1.00	0.87
ATOM	2164	HA	SER	147	14.141	-5.234	-20.859	1.00	1.03
ATOM	2165	CB	SER	147	15.667	-6.709	-20.554	1.00	0.95
ATOM	2166	HB1	SER	147	15.798	-7.158	-21.530	1.00	1.42
ATOM	2167	HB2	SER	147	15.871	-7.445	-19.794	1.00	1.34
ATOM	2168	OG	SER	147	16.561	-5.616	-20.395	1.00	1.71
ATOM	2169	HG	SER	147	17.097	-5.555	-21.190	1.00	2.16
ATOM	2170	C	SER	147	13.288	-7.185	-21.138	1.00	0.79
ATOM	2171	O	SER	147	12.747	-6.865	-22.178	1.00	1.40
ATOM	2172	N	HIS	148	13.098	-8.366	-20.605	1.00	0.66
ATOM	2173	HN	HIS	148	13.551	-8.602	-19.768	1.00	1.10
ATOM	2174	CA	HIS	148	12.199	-9.360	-21.272	1.00	0.65
ATOM	2175	HA	HIS	148	11.629	-8.874	-22.048	1.00	0.74
ATOM	2176	CB	HIS	148	13.041	-10.479	-21.887	1.00	0.79
ATOM	2177	HB1	HIS	148	12.401	-11.312	-22.138	1.00	1.14
ATOM	2178	HB2	HIS	148	13.786	-10.801	-21.174	1.00	1.30
ATOM	2179	CG	HIS	148	13.723	-9.980	-23.130	1.00	1.66
ATOM	2180	ND1	HIS	148	13.104	-9.116	-24.019	1.00	2.52
ATOM	2181	HD1	HIS	148	12.200	-8.747	-23.934	1.00	2.81
ATOM	2182	CD2	HIS	148	14.969	-10.226	-23.652	1.00	2.62
ATOM	2183	HD2	HIS	148	15.715	-10.867	-23.206	1.00	3.00
ATOM	2184	CE1	HIS	148	13.970	-8.875	-25.020	1.00	3.46
ATOM	2185	HE1	HIS	148	13.759	-8.233	-25.863	1.00	4.33
ATOM	2186	NE2	HIS	148	15.123	-9.528	-24.846	1.00	3.55
ATOM	2187	C	HIS	148	11.238	-9.971	-20.249	1.00	0.55
ATOM	2188	O	HIS	148	10.743	-11.064	-20.435	1.00	0.60
ATOM	2189	N	PHE	149	10.978	-9.293	-19.167	1.00	0.57
ATOM	2190	HN	PHE	149	11.392	-8.417	-19.021	1.00	0.73
ATOM	2191	CA	PHE	149	10.060	-9.871	-18.145	1.00	0.48
ATOM	2192	HA	PHE	149	10.416	-10.849	-17.857	1.00	0.51
ATOM	2193	CB	PHE	149	10.022	-8.967	-16.911	1.00	0.44
ATOM	2194	HB1	PHE	149	9.603	-8.008	-17.177	1.00	0.44
ATOM	2195	HB2	PHE	149	11.023	-8.831	-16.530	1.00	0.46

33/75

ATOM	2196	CG	PHE	149	9.161	-9.615	-15.851	1.00	0.40
ATOM	2197	CD1	PHE	149	7.766	-9.507	-15.919	1.00	0.36
ATOM	2198	HD1	PHE	149	7.305	-8.956	-16.726	1.00	0.38
ATOM	2199	CD2	PHE	149	9.757	-10.328	-14.804	1.00	0.42
ATOM	2200	HD2	PHE	149	10.832	-10.412	-14.750	1.00	0.48
ATOM	2201	CE1	PHE	149	6.969	-10.112	-14.941	1.00	0.35
ATOM	2202	HE1	PHE	149	5.894	-10.031	-14.996	1.00	0.37
ATOM	2203	CE2	PHE	149	8.958	-10.932	-13.825	1.00	0.40
ATOM	2204	HE2	PHE	149	9.417	-11.482	-13.016	1.00	0.45
ATOM	2205	CZ	PHE	149	7.564	-10.825	-13.894	1.00	0.37
ATOM	2206	HZ	PHE	149	6.948	-11.291	-13.140	1.00	0.38
ATOM	2207	C	PHE	149	8.641	-9.993	-18.706	1.00	0.43
ATOM	2208	O	PHE	149	8.080	-9.044	-19.217	1.00	0.45
ATOM	2209	N	MET	150	8.050	-11.153	-18.575	1.00	0.43
ATOM	2210	HN	MET	150	8.523	-11.888	-18.133	1.00	0.50
ATOM	2211	CA	MET	150	6.651	-11.357	-19.051	1.00	0.39
ATOM	2212	HA	MET	150	6.189	-10.400	-19.245	1.00	0.38
ATOM	2213	CB	MET	150	6.632	-12.207	-20.328	1.00	0.44
ATOM	2214	HB1	MET	150	5.610	-12.374	-20.632	1.00	0.45
ATOM	2215	HB2	MET	150	7.109	-13.157	-20.134	1.00	0.47
ATOM	2216	CG	MET	150	7.381	-11.477	-21.446	1.00	0.50
ATOM	2217	HG1	MET	150	8.401	-11.831	-21.485	1.00	0.98
ATOM	2218	HG2	MET	150	7.376	-10.415	-21.253	1.00	0.86
ATOM	2219	SD	MET	150	6.571	-11.806	-23.033	1.00	1.32
ATOM	2220	CE	MET	150	7.378	-13.384	-23.393	1.00	2.23
ATOM	2221	HE1	MET	150	7.326	-14.022	-22.521	1.00	2.66
ATOM	2222	HE2	MET	150	8.411	-13.211	-23.647	1.00	2.74
ATOM	2223	HE3	MET	150	6.879	-13.861	-24.225	1.00	2.74
ATOM	2224	C	MET	150	5.877	-12.071	-17.943	1.00	0.32
ATOM	2225	O	MET	150	6.435	-12.837	-17.183	1.00	0.32
ATOM	2226	N	LEU	151	4.605	-11.819	-17.827	1.00	0.28
ATOM	2227	HN	LEU	151	4.169	-11.188	-18.437	1.00	0.30
ATOM	2228	CA	LEU	151	3.821	-12.478	-16.746	1.00	0.24
ATOM	2229	HA	LEU	151	4.120	-12.064	-15.803	1.00	0.24
ATOM	2230	CB	LEU	151	2.327	-12.212	-16.966	1.00	0.24
ATOM	2231	HB1	LEU	151	1.765	-12.626	-16.145	1.00	0.25
ATOM	2232	HB2	LEU	151	2.012	-12.680	-17.887	1.00	0.28
ATOM	2233	CG	LEU	151	2.061	-10.703	-17.047	1.00	0.28
ATOM	2234	HG	LEU	151	2.900	-10.208	-17.512	1.00	0.52
ATOM	2235	CD1	LEU	151	0.804	-10.457	-17.881	1.00	0.35
ATOM	2236	HD11	LEU	151	0.506	-9.424	-17.788	1.00	1.07
ATOM	2237	HD12	LEU	151	0.007	-11.095	-17.526	1.00	1.02
ATOM	2238	HD13	LEU	151	1.009	-10.682	-18.917	1.00	1.17
ATOM	2239	CD2	LEU	151	1.848	-10.140	-15.638	1.00	0.46
ATOM	2240	HD21	LEU	151	2.078	-9.084	-15.635	1.00	1.14
ATOM	2241	HD22	LEU	151	2.495	-10.650	-14.941	1.00	1.16
ATOM	2242	HD23	LEU	151	0.820	-10.284	-15.345	1.00	1.11
ATOM	2243	C	LEU	151	4.076	-14.004	-16.794	1.00	0.24
ATOM	2244	O	LEU	151	3.879	-14.613	-17.826	1.00	0.28
ATOM	2245	N	PRO	152	4.504	-14.641	-15.711	1.00	0.22
ATOM	2246	CA	PRO	152	4.748	-16.112	-15.751	1.00	0.23
ATOM	2247	HA	PRO	152	5.480	-16.354	-16.503	1.00	0.24
ATOM	2248	CB	PRO	152	5.323	-16.404	-14.364	1.00	0.24
ATOM	2249	HB1	PRO	152	6.361	-16.686	-14.453	1.00	0.29
ATOM	2250	HB2	PRO	152	4.766	-17.208	-13.903	1.00	0.26
ATOM	2251	CG	PRO	152	5.209	-15.141	-13.507	1.00	0.32
ATOM	2252	HG1	PRO	152	6.166	-14.917	-13.061	1.00	0.44
ATOM	2253	HG2	PRO	152	4.473	-15.295	-12.730	1.00	0.41
ATOM	2254	CD	PRO	152	4.778	-13.976	-14.402	1.00	0.25
ATOM	2255	HD2	PRO	152	3.886	-13.507	-14.008	1.00	0.25
ATOM	2256	HD1	PRO	152	5.581	-13.263	-14.503	1.00	0.27
ATOM	2257	C	PRO	152	3.462	-16.915	-15.974	1.00	0.21
ATOM	2258	O	PRO	152	2.378	-16.371	-16.038	1.00	0.20
ATOM	2259	N	ASP	153	3.582	-18.209	-16.090	1.00	0.23
ATOM	2260	HN	ASP	153	4.468	-18.622	-16.031	1.00	0.25
ATOM	2261	CA	ASP	153	2.380	-19.063	-16.304	1.00	0.23
ATOM	2262	HA	ASP	153	1.890	-18.772	-17.221	1.00	0.23
ATOM	2263	CB	ASP	153	2.813	-20.526	-16.401	1.00	0.25
ATOM	2264	HB1	ASP	153	1.943	-21.163	-16.363	1.00	0.26
ATOM	2265	HB2	ASP	153	3.470	-20.762	-15.576	1.00	0.26
ATOM	2266	CG	ASP	153	3.550	-20.752	-17.722	1.00	0.27
ATOM	2267	OD1	ASP	153	4.768	-20.687	-17.717	1.00	1.08
ATOM	2268	OD2	ASP	153	2.884	-20.994	-18.715	1.00	1.14
ATOM	2269	C	ASP	153	1.409	-18.899	-15.133	1.00	0.21
ATOM	2270	O	ASP	153	0.208	-18.858	-15.310	1.00	0.21
ATOM	2271	N	ASP	154	1.919	-18.820	-13.935	1.00	0.21
ATOM	2272	HN	ASP	154	2.881	-18.866	-13.811	1.00	0.22

34/75

ATOM	2273	CA	ASP	154	1.025	-18.678	-12.752	1.00	0.21
ATOM	2274	HA	ASP	154	0.431	-19.572	-12.641	1.00	0.22
ATOM	2275	CB	ASP	154	1.880	-18.474	-11.496	1.00	0.23
ATOM	2276	HB1	ASP	154	2.466	-17.572	-11.602	1.00	0.22
ATOM	2277	HB2	ASP	154	2.541	-19.319	-11.370	1.00	0.25
ATOM	2278	CG	ASP	154	0.975	-18.347	-10.267	1.00	0.25
ATOM	2279	OD1	ASP	154	1.276	-18.982	-9.269	1.00	1.13
ATOM	2280	OD2	ASP	154	0.004	-17.613	-10.340	1.00	1.07
ATOM	2281	C	ASP	154	0.102	-17.473	-12.943	1.00	0.19
ATOM	2282	O	ASP	154	-1.095	-17.564	-12.759	1.00	0.19
ATOM	2283	N	ASP	155	0.645	-16.345	-13.303	1.00	0.19
ATOM	2284	HN	ASP	155	1.613	-16.288	-13.443	1.00	0.21
ATOM	2285	CA	ASP	155	-0.210	-15.140	-13.496	1.00	0.19
ATOM	2286	HA	ASP	155	-0.843	-15.011	-12.631	1.00	0.20
ATOM	2287	CB	ASP	155	0.683	-13.909	-13.653	1.00	0.21
ATOM	2288	HB1	ASP	155	0.087	-13.067	-13.969	1.00	0.22
ATOM	2289	HB2	ASP	155	1.443	-14.113	-14.393	1.00	0.22
ATOM	2290	CG	ASP	155	1.351	-13.588	-12.315	1.00	0.24
ATOM	2291	OD1	ASP	155	2.355	-12.896	-12.327	1.00	1.07
ATOM	2292	OD2	ASP	155	0.845	-14.038	-11.300	1.00	1.14
ATOM	2293	C	ASP	155	-1.087	-15.300	-14.744	1.00	0.19
ATOM	2294	O	ASP	155	-2.240	-14.918	-14.750	1.00	0.19
ATOM	2295	N	VAL	156	-0.555	-15.850	-15.802	1.00	0.19
ATOM	2296	HN	VAL	156	0.379	-16.147	-15.787	1.00	0.19
ATOM	2297	CA	VAL	156	-1.372	-16.013	-17.041	1.00	0.21
ATOM	2298	HA	VAL	156	-1.726	-15.044	-17.362	1.00	0.22
ATOM	2299	CB	VAL	156	-0.519	-16.630	-18.148	1.00	0.23
ATOM	2300	HB	VAL	156	-0.034	-17.521	-17.776	1.00	0.23
ATOM	2301	CG1	VAL	156	-1.416	-16.995	-19.333	1.00	0.27
ATOM	2302	HG11	VAL	156	-2.273	-16.338	-19.348	1.00	1.00
ATOM	2303	HG12	VAL	156	-1.747	-18.018	-19.235	1.00	1.05
ATOM	2304	HG13	VAL	156	-0.861	-16.882	-20.253	1.00	1.05
ATOM	2305	CG2	VAL	156	0.535	-15.618	-18.600	1.00	0.26
ATOM	2306	HG21	VAL	156	0.990	-15.162	-17.733	1.00	1.07
ATOM	2307	HG22	VAL	156	0.067	-14.856	-19.204	1.00	1.05
ATOM	2308	HG23	VAL	156	1.293	-16.123	-19.180	1.00	1.00
ATOM	2309	C	VAL	156	-2.574	-16.919	-16.754	1.00	0.20
ATOM	2310	O	VAL	156	-3.694	-16.615	-17.107	1.00	0.21
ATOM	2311	N	GLN	157	-2.356	-18.035	-16.124	1.00	0.20
ATOM	2312	HN	GLN	157	-1.447	-18.277	-15.847	1.00	0.20
ATOM	2313	CA	GLN	157	-3.498	-18.941	-15.824	1.00	0.22
ATOM	2314	HA	GLN	157	-3.987	-19.214	-16.747	1.00	0.24
ATOM	2315	CB	GLN	157	-2.995	-20.204	-15.117	1.00	0.24
ATOM	2316	HB1	GLN	157	-3.838	-20.774	-14.756	1.00	0.26
ATOM	2317	HB2	GLN	157	-2.368	-19.922	-14.282	1.00	0.23
ATOM	2318	CG	GLN	157	-2.184	-21.064	-16.095	1.00	0.25
ATOM	2319	HG1	GLN	157	-1.174	-20.686	-16.152	1.00	0.94
ATOM	2320	HG2	GLN	157	-2.636	-21.032	-17.074	1.00	0.87
ATOM	2321	CD	GLN	157	-2.152	-22.510	-15.598	1.00	1.19
ATOM	2322	OE1	GLN	157	-2.594	-22.799	-14.504	1.00	1.89
ATOM	2323	NE2	GLN	157	-1.646	-23.437	-16.364	1.00	1.96
ATOM	2324	HE21	GLN	157	-1.291	-23.203	-17.247	1.00	2.18
ATOM	2325	HE22	GLN	157	-1.624	-24.368	-16.058	1.00	2.65
ATOM	2326	C	GLN	157	-4.505	-18.214	-14.925	1.00	0.22
ATOM	2327	O	GLN	157	-5.702	-18.356	-15.077	1.00	0.24
ATOM	2328	N	GLY	158	-4.027	-17.456	-13.974	1.00	0.21
ATOM	2329	HN	GLY	158	-3.057	-17.370	-13.859	1.00	0.20
ATOM	2330	CA	GLY	158	-4.952	-16.741	-13.045	1.00	0.22
ATOM	2331	HA1	GLY	158	-4.380	-16.319	-12.232	1.00	0.22
ATOM	2332	HA2	GLY	158	-5.667	-17.446	-12.646	1.00	0.25
ATOM	2333	C	GLY	158	-5.704	-15.615	-13.766	1.00	0.20
ATOM	2334	O	GLY	158	-6.918	-15.552	-13.730	1.00	0.21
ATOM	2335	N	ILE	159	-5.007	-14.713	-14.405	1.00	0.18
ATOM	2336	HN	ILE	159	-4.028	-14.763	-14.418	1.00	0.18
ATOM	2337	CA	ILE	159	-5.713	-13.593	-15.097	1.00	0.19
ATOM	2338	HA	ILE	159	-6.301	-13.054	-14.375	1.00	0.20
ATOM	2339	CB	ILE	159	-4.679	-12.648	-15.735	1.00	0.19
ATOM	2340	HB	ILE	159	-3.950	-12.367	-14.988	1.00	0.20
ATOM	2341	CG1	ILE	159	-5.355	-11.384	-16.284	1.00	0.24
ATOM	2342	HG11	ILE	159	-6.308	-11.645	-16.717	1.00	0.26
ATOM	2343	HG12	ILE	159	-4.725	-10.952	-17.045	1.00	0.28
ATOM	2344	CG2	ILE	159	-3.968	-13.361	-16.880	1.00	0.21
ATOM	2345	HG21	ILE	159	-2.998	-12.914	-17.036	1.00	1.01
ATOM	2346	HG22	ILE	159	-4.556	-13.274	-17.781	1.00	1.01
ATOM	2347	HG23	ILE	159	-3.848	-14.398	-16.628	1.00	1.04
ATOM	2348	CD1	ILE	159	-5.571	-10.356	-15.166	1.00	0.27
ATOM	2349	HD11	ILE	159	-6.322	-9.644	-15.476	1.00	1.05

35/75

ATOM	2350	HD12	ILE	159	-4.644	-9.838	-14.978	1.00	1.06
ATOM	2351	HD13	ILE	159	-5.893	-10.848	-14.265	1.00	1.02
ATOM	2352	C	ILE	159	-6.644	-14.162	-16.173	1.00	0.21
ATOM	2353	O	ILE	159	-7.754	-13.700	-16.347	1.00	0.23
ATOM	2354	N	GLN	160	-6.215	-15.168	-16.885	1.00	0.22
ATOM	2355	HN	GLN	160	-5.322	-15.538	-16.726	1.00	0.21
ATOM	2356	CA	GLN	160	-7.097	-15.763	-17.930	1.00	0.27
ATOM	2357	HA	GLN	160	-7.457	-14.979	-18.580	1.00	0.29
ATOM	2358	CB	GLN	160	-6.317	-16.786	-18.756	1.00	0.31
ATOM	2359	HB1	GLN	160	-6.999	-17.334	-19.389	1.00	0.35
ATOM	2360	HB2	GLN	160	-5.809	-17.472	-18.093	1.00	0.30
ATOM	2361	CG	GLN	160	-5.289	-16.062	-19.626	1.00	0.34
ATOM	2362	HG1	GLN	160	-4.606	-15.512	-18.997	1.00	0.92
ATOM	2363	HG2	GLN	160	-5.799	-15.378	-20.290	1.00	0.91
ATOM	2364	CD	GLN	160	-4.508	-17.087	-20.451	1.00	1.11
ATOM	2365	OE1	GLN	160	-4.451	-18.248	-20.100	1.00	1.88
ATOM	2366	NE2	GLN	160	-3.901	-16.704	-21.540	1.00	1.83
ATOM	2367	HE21	GLN	160	-3.947	-15.767	-21.824	1.00	2.13
ATOM	2368	HE22	GLN	160	-3.398	-17.353	-22.075	1.00	2.46
ATOM	2369	C	GLN	160	-8.290	-16.447	-17.261	1.00	0.28
ATOM	2370	O	GLN	160	-9.386	-16.449	-17.779	1.00	0.31
ATOM	2371	N	SER	161	-8.086	-17.035	-16.117	1.00	0.27
ATOM	2372	HN	SER	161	-7.193	-17.030	-15.714	1.00	0.25
ATOM	2373	CA	SER	161	-9.213	-17.718	-15.424	1.00	0.30
ATOM	2374	HA	SER	161	-9.658	-18.444	-16.089	1.00	0.34
ATOM	2375	CB	SER	161	-8.690	-18.427	-14.174	1.00	0.33
ATOM	2376	HB1	SER	161	-7.861	-19.067	-14.444	1.00	0.35
ATOM	2377	HB2	SER	161	-9.476	-19.024	-13.741	1.00	0.36
ATOM	2378	OG	SER	161	-8.267	-17.455	-13.227	1.00	0.33
ATOM	2379	HG	SER	161	-9.045	-16.986	-12.915	1.00	0.94
ATOM	2380	C	SER	161	-10.267	-16.684	-15.019	1.00	0.30
ATOM	2381	O	SER	161	-11.433	-16.997	-14.882	1.00	0.35
ATOM	2382	N	LEU	162	-9.867	-15.457	-14.815	1.00	0.27
ATOM	2383	HN	LEU	162	-8.920	-15.225	-14.921	1.00	0.26
ATOM	2384	CA	LEU	162	-10.852	-14.413	-14.405	1.00	0.29
ATOM	2385	HA	LEU	162	-11.637	-14.869	-13.821	1.00	0.33
ATOM	2386	CB	LEU	162	-10.141	-13.350	-13.563	1.00	0.28
ATOM	2387	HB1	LEU	162	-10.802	-12.509	-13.411	1.00	0.29
ATOM	2388	HB2	LEU	162	-9.256	-13.017	-14.086	1.00	0.27
ATOM	2389	CG	LEU	162	-9.736	-13.937	-12.206	1.00	0.30
ATOM	2390	HG	LEU	162	-9.157	-14.836	-12.367	1.00	0.30
ATOM	2391	CD1	LEU	162	-8.883	-12.918	-11.450	1.00	0.33
ATOM	2392	HD11	LEU	162	-8.496	-13.370	-10.549	1.00	1.03
ATOM	2393	HD12	LEU	162	-9.490	-12.063	-11.191	1.00	1.01
ATOM	2394	HD13	LEU	162	-8.062	-12.601	-12.075	1.00	1.12
ATOM	2395	CD2	LEU	162	-10.980	-14.272	-11.374	1.00	0.33
ATOM	2396	HD21	LEU	162	-11.227	-15.315	-11.502	1.00	1.05
ATOM	2397	HD22	LEU	162	-11.812	-13.664	-11.697	1.00	1.09
ATOM	2398	HD23	LEU	162	-10.776	-14.078	-10.332	1.00	1.01
ATOM	2399	C	LEU	162	-11.461	-13.742	-15.643	1.00	0.30
ATOM	2400	O	LEU	162	-12.664	-13.615	-15.757	1.00	0.36
ATOM	2401	N	TYR	163	-10.645	-13.300	-16.564	1.00	0.27
ATOM	2402	HN	TYR	163	-9.677	-13.404	-16.452	1.00	0.26
ATOM	2403	CA	TYR	163	-11.188	-12.626	-17.783	1.00	0.31
ATOM	2404	HA	TYR	163	-12.144	-12.182	-17.549	1.00	0.33
ATOM	2405	CB	TYR	163	-10.219	-11.531	-18.236	1.00	0.29
ATOM	2406	HB1	TYR	163	-10.562	-11.112	-19.170	1.00	0.32
ATOM	2407	HB2	TYR	163	-9.234	-11.952	-18.371	1.00	0.29
ATOM	2408	CG	TYR	163	-10.162	-10.444	-17.190	1.00	0.25
ATOM	2409	CD1	TYR	163	-9.223	-10.520	-16.155	1.00	0.23
ATOM	2410	HD1	TYR	163	-8.545	-11.359	-16.103	1.00	0.23
ATOM	2411	CD2	TYR	163	-11.042	-9.357	-17.258	1.00	0.27
ATOM	2412	HD2	TYR	163	-11.767	-9.298	-18.056	1.00	0.30
ATOM	2413	CE1	TYR	163	-9.164	-9.511	-15.187	1.00	0.24
ATOM	2414	HE1	TYR	163	-8.439	-9.571	-14.388	1.00	0.25
ATOM	2415	CE2	TYR	163	-10.984	-8.348	-16.289	1.00	0.27
ATOM	2416	HE2	TYR	163	-11.663	-7.510	-16.340	1.00	0.30
ATOM	2417	CZ	TYR	163	-10.044	-8.425	-15.253	1.00	0.27
ATOM	2418	OH	TYR	163	-9.985	-7.430	-14.299	1.00	0.31
ATOM	2419	HH	TYR	163	-10.344	-7.782	-13.481	1.00	0.99
ATOM	2420	C	TYR	163	-11.367	-13.647	-18.909	1.00	0.37
ATOM	2421	O	TYR	163	-11.953	-13.357	-19.933	1.00	0.43
ATOM	2422	N	GLY	164	-10.865	-14.836	-18.729	1.00	0.38
ATOM	2423	HN	GLY	164	-10.394	-15.046	-17.896	1.00	0.35
ATOM	2424	CA	GLY	164	-11.001	-15.877	-19.789	1.00	0.47
ATOM	2425	HA1	GLY	164	-11.851	-15.651	-20.413	1.00	0.53

36/75

ATOM	2427	C	GLY	164	-9.735	-15.902	-20.648	1.00	0.55
ATOM	2428	O	GLY	164	-9.761	-15.580	-21.819	1.00	1.01
TER	2429		GLY	164					
HETATM	2430	ZN	ZN	166	-0.218	-6.515	-2.613	1.00	0.24
HETATM	2431	ZN	ZN	167	-3.506	6.833	-0.714	1.00	0.97
HETATM	2432	CA	CA	168	6.060	3.350	3.030	1.00	0.23
HETATM	2433	C1	WAY	169	2.180	-4.315	1.627	0.00	0.30
HETATM	2434	C2	WAY	169	0.865	-4.629	1.215	0.00	0.33
HETATM	2435	1CE1	WAY	169	-0.170	-4.517	2.143	0.00	0.38
HETATM	2436	1CZ	WAY	169	0.074	-4.157	3.457	0.00	0.40
HETATM	2437	1CE2	WAY	169	1.355	-3.807	3.841	0.00	0.38
HETATM	2438	C6	WAY	169	2.395	-3.805	2.922	0.00	0.33
HETATM	2439	1HE1	WAY	169	-1.190	-4.713	1.839	0.00	0.42
HETATM	2440	1HZ	WAY	169	-0.734	-4.151	4.173	0.00	0.45
HETATM	2441	1HE2	WAY	169	1.535	-3.534	4.872	0.00	0.42
HETATM	2442	C10	WAY	169	0.444	-5.080	-0.136	0.00	0.36
HETATM	2443	O11	WAY	169	0.467	-6.264	-0.463	0.00	0.58
HETATM	2444	N12	WAY	169	-0.019	-4.195	-1.032	0.00	0.61
HETATM	2445	O13	WAY	169	-0.045	-4.608	-2.371	0.00	0.68
HETATM	2446	H14	WAY	169	-0.357	-3.297	-0.743	0.00	0.88
HETATM	2447	H15	WAY	169	-0.953	-4.727	-2.645	0.00	1.13
HETATM	2448	1CH1	WAY	169	3.728	-3.247	3.360	0.00	0.37
HETATM	2449	1HH1	WAY	169	3.702	-2.162	3.422	0.00	1.07
HETATM	2450	1HH2	WAY	169	4.519	-3.516	2.664	0.00	1.06
HETATM	2451	1HH3	WAY	169	4.013	-3.623	4.339	0.00	1.11
HETATM	2452	N20	WAY	169	3.274	-4.485	0.819	0.00	0.29
HETATM	2453	S21	WAY	169	3.865	-3.175	0.021	0.00	0.25
HETATM	2454	2CB	WAY	169	3.882	-5.812	0.684	0.00	0.32
HETATM	2455	2CE1	WAY	169	7.334	-6.241	2.178	0.00	1.09
HETATM	2456	2CZ	WAY	169	6.971	-6.520	3.488	0.00	0.53
HETATM	2457	N25	WAY	169	5.697	-6.659	3.876	0.00	1.47
HETATM	2458	2CD2	WAY	169	4.747	-6.451	2.954	0.00	1.37
HETATM	2459	C27	WAY	169	5.010	-6.084	1.640	0.00	0.36
HETATM	2460	2CD1	WAY	169	6.338	-5.982	1.250	0.00	1.14
HETATM	2461	2HE1	WAY	169	8.374	-6.224	1.881	0.00	1.94
HETATM	2462	2HZ	WAY	169	7.752	-6.630	4.227	0.00	0.61
HETATM	2463	2HD2	WAY	169	3.708	-6.570	3.227	0.00	2.23
HETATM	2464	2HD1	WAY	169	6.599	-5.706	0.239	0.00	2.05
HETATM	2465	2HB1	WAY	169	4.245	-5.905	-0.339	0.00	0.31
HETATM	2466	2HB2	WAY	169	3.095	-6.552	0.832	0.00	0.34
HETATM	2467	C35	WAY	169	4.187	-3.617	-1.665	0.00	0.23
HETATM	2468	3CD1	WAY	169	3.310	-3.216	-2.661	0.00	0.25
HETATM	2469	3CE1	WAY	169	3.622	-3.465	-3.992	0.00	0.27
HETATM	2470	C38	WAY	169	4.769	-4.183	-4.326	0.00	0.24
HETATM	2471	3CE2	WAY	169	5.602	-4.644	-3.308	0.00	0.23
HETATM	2472	3CD2	WAY	169	5.315	-4.359	-1.979	0.00	0.23
HETATM	2473	3HD1	WAY	169	2.392	-2.714	-2.389	0.00	0.29
HETATM	2474	3HE1	WAY	169	2.961	-3.091	-4.758	0.00	0.31
HETATM	2475	3HE2	WAY	169	6.481	-5.228	-3.535	0.00	0.26
HETATM	2476	3HD2	WAY	169	5.959	-4.707	-1.184	0.00	0.27
HETATM	2477	O45	WAY	169	5.078	-4.439	-5.664	0.00	0.27
HETATM	2478	3CH	WAY	169	6.245	-5.202	-5.904	0.00	0.28
HETATM	2479	3HH1	WAY	169	6.379	-5.372	-6.973	0.00	0.31
HETATM	2480	3HH2	WAY	169	6.178	-6.172	-5.407	0.00	0.28
HETATM	2481	3HH3	WAY	169	7.127	-4.683	-5.526	0.00	0.29
HETATM	2482	O50	WAY	169	5.123	-2.847	0.614	0.00	0.27
HETATM	2483	O51	WAY	169	2.834	-2.186	0.004	0.00	0.25
END									

37/75

	Atom Type	Res.		X	Y	Z	Occ.	B	MOL.
ATOM	1 CB	THR	7	73.468	27.410	6.079	1.00	42.70	A ₁₃
ATOM	2 OG1	THR	7	72.149	27.911	6.358	1.00	37.82	A ₁₃
ATOM	4 CG2	THR	7	73.843	26.297	7.068	1.00	25.79	A ₁₃
ATOM	5 C	THR	7	75.936	28.076	6.227	1.00	28.29	A ₁₃
ATOM	6 O	THR	7	76.497	28.090	7.332	1.00	22.94	A ₁₃
ATOM	9 N	THR	7	74.360	29.396	4.862	1.00	20.25	A ₁₃
ATOM	11 CA	THR	7	74.501	28.593	6.099	1.00	21.49	A ₁₃
ATOM	12 N	LEU	8	76.547	27.691	5.099	1.00	32.90	A ₁₃
ATOM	14 CA	LEU	8	77.915	27.150	5.105	1.00	31.85	A ₁₃
ATOM	15 GB	LEU	8	77.952	25.759	4.438	1.00	21.38	A ₁₃
ATOM	16 CG	LEU	8	78.016	25.576	2.910	1.00	29.31	A ₁₃
ATOM	17 CD1	LEU	8	79.463	25.509	2.425	1.00	16.78	A ₁₃
ATOM	18 CD2	LEU	8	77.334	24.292	2.527	1.00	23.37	A ₁₃
ATOM	19 C	LEU	8	78.956	28.070	4.465	1.00	24.01	A ₁₃
ATOM	20 O	LEU	8	78.835	28.415	3.293	1.00	26.18	A ₁₃
ATOM	21 N	LYS	9	79.980	28.424	5.251	1.00	36.26	A ₁₃
ATOM	23 CA	LYS	9	81.106	29.298	4.867	1.00	23.24	A ₁₃
ATOM	24 CB	LYS	9	82.438	28.521	4.977	1.00	25.52	A ₁₃
ATOM	25 CG	LYS	9	82.767	27.570	3.815	1.00	19.05	A ₁₃
ATOM	26 CD	LYS	9	83.661	28.243	2.753	1.00	31.69	A ₁₃
ATOM	27 CE	LYS	9	83.451	27.688	1.323	1.00	25.30	A ₁₃
ATOM	28 NZ	LYS	9	82.056	27.938	0.797	1.00	20.65	A ₁₃
ATOM	32 C	LYS	9	81.042	30.073	3.526	1.00	31.41	A ₁₃
ATOM	33 O	LYS	9	80.764	29.505	2.466	1.00	22.31	A ₁₃
ATOM	34 N	TRP	10	81.327	31.372	3.573	1.00	15.84	A ₁₃
ATOM	36 CA	TRP	10	81.312	32.172	2.361	1.00	10.58	A ₁₃
ATOM	37 CB	TRP	10	81.636	33.620	2.680	1.00	21.39	A ₁₃
ATOM	38 CG	TRP	10	80.529	34.337	3.343	1.00	22.84	A ₁₃
ATOM	39 CD2	TRP	10	79.479	35.074	2.697	1.00	20.41	A ₁₃
ATOM	40 CE2	TRP	10	78.676	35.631	3.718	1.00	24.50	A ₁₃
ATOM	41 CE3	TRP	10	79.142	35.320	1.357	1.00	13.29	A ₁₃
ATOM	42 CD1	TRP	10	80.327	34.469	4.682	1.00	13.40	A ₁₃
ATOM	43 NE1	TRP	10	79.220	35.253	4.919	1.00	18.40	A ₁₃
ATOM	45 CZ2	TRP	10	77.550	36.418	3.442	1.00	12.63	A ₁₃
ATOM	46 CZ3	TRP	10	78.021	36.105	1.083	1.00	19.89	A ₁₃
ATOM	47 CH2	TRP	10	77.242	36.641	2.120	1.00	13.62	A ₁₃
ATOM	48 C	TRP	10	82.377	31.594	1.455	1.00	22.95	A ₁₃
ATOM	49 O	TRP	10	83.450	31.221	1.920	1.00	16.28	A ₁₃
ATOM	50 N	SER	11	82.087	31.533	0.167	1.00	14.81	A ₁₃
ATOM	52 CA	SER	11	83.017	30.975	-0.801	1.00	19.50	A ₁₃
ATOM	53 CB	SER	11	82.282	30.596	-2.086	1.00	24.36	A ₁₃
ATOM	54 OG	SER	11	81.605	29.353	-1.958	1.00	40.49	A ₁₃
ATOM	56 C	SER	11	84.190	31.867	-1.134	1.00	16.53	A ₁₃
ATOM	57 O	SER	11	85.132	31.423	-1.779	1.00	23.48	A ₁₃
ATOM	58 N	LYS	12	84.153	33.113	-0.686	1.00	12.50	A ₁₃
ATOM	60 CA	LYS	12	85.232	34.057	-0.961	1.00	17.05	A ₁₃
ATOM	61 CB	LYS	12	84.741	35.168	-1.891	1.00	17.32	A ₁₃
ATOM	62 CG	LYS	12	83.526	35.898	-1.350	1.00	18.49	A ₁₃
ATOM	63 CD	LYS	12	82.788	36.644	-2.446	1.00	18.29	A ₁₃
ATOM	64 CE	LYS	12	81.534	37.282	-1.888	1.00	18.44	A ₁₃
ATOM	65 NZ	LYS	12	80.805	38.094	-2.895	1.00	16.65	A ₁₃
ATOM	69 C	LYS	12	85.687	34.662	0.344	1.00	11.16	A ₁₃
ATOM	70 O	LYS	12	84.946	34.637	1.319	1.00	12.63	A ₁₃
ATOM	71 N	MET	13	85.915	35.185	0.355	1.00	15.52	A ₁₃
ATOM	73 CA	MET	13	87.516	35.801	1.537	1.00	11.04	A ₁₃
ATOM	74 CB	MET	13	89.028	35.547	1.565	1.00	16.57	A ₁₃
ATOM	75 CG	MET	13	89.431	34.082	1.707	1.00	20.92	A ₁₃
ATOM	76 SD	MET	13	88.905	33.235	3.227	1.00	20.10	A ₁₃
ATOM	77 CE	MET	13	87.486	32.313	2.604	1.00	16.29	A ₁₃
ATOM	78 C	MET	13	87.258	37.296	1.572	1.00	13.23	A ₁₃
ATOM	79 O	MET	13	87.247	37.916	2.634	1.00	22.80	A ₁₃
ATOM	80 N	ASN	14	87.111	37.875	0.389	1.00	15.02	A ₁₃
ATOM	82 CA	ASN	14	86.853	39.294	0.241	1.00	33.02	A ₁₃
ATOM	83 CB	ASN	14	87.445	39.801	-1.082	1.00	19.42	A ₁₃
ATOM	84 CG	ASN	14	88.925	39.482	-1.217	1.00	30.32	A ₁₃
ATOM	85 OD1	ASN	14	89.343	38.341	-1.031	1.00	30.12	A ₁₃
ATOM	86 ND2	ASN	14	89.723	40.489	-1.549	1.00	28.22	A ₁₃
ATOM	89 C	ASN	14	85.337	39.482	0.277	1.00	27.58	A ₁₃
ATOM	90 O	ASN	14	84.606	38.935	-0.568	1.00	28.01	A ₁₃
ATOM	91 N	LEU	15	84.868	40.212	1.287	1.00	19.06	A ₁₃
ATOM	93 CA	LEU	15	83.444	40.450	1.459	1.00	20.03	A ₁₃
ATOM	94 CB	LEU	15	82.930	39.690	2.691	1.00	19.55	A ₁₃
ATOM	95 CG	LEU	15	83.027	38.166	2.593	1.00	19.02	A ₁₃
ATOM	96 CD1	LEU	15	83.216	37.555	3.962	1.00	17.48	A ₁₃
ATOM	97 CD2	LEU	15	81.799	37.604	1.903	1.00	23.43	A ₁₃
ATOM	98 C	LEU	15	83.161	41.928	1.609	1.00	19.52	A ₁₃
ATOM	99 O	LEU	15	83.980	42.676	2.130	1.00	15.98	A ₁₃

FIG. 5

38/75

ATOM	100	N	THR	16	81.983	42.343	1.162	1.00	21.22	A_13
ATOM	102	CA	THR	16	81.578	43.736	1.252	1.00	10.00	A_13
ATOM	103	CB	THR	16	81.194	44.257	-0.109	1.00	10.00	A_13
ATOM	104	OG1	THR	16	80.225	43.370	-0.681	1.00	22.43	A_13
ATOM	106	CG2	THR	16	82.427	44.383	-1.009	1.00	15.42	A_13
ATOM	107	C	THR	16	80.368	43.869	2.184	1.00	14.48	A_13
ATOM	108	O	THR	16	79.647	42.897	2.445	1.00	15.74	A_13
ATOM	109	N	TYR	17	80.176	45.065	2.716	1.00	15.89	A_13
ATOM	111	CA	TYR	17	79.064	45.340	3.604	1.00	13.19	A_13
ATOM	112	CB	TYR	17	79.480	45.195	5.067	1.00	21.42	A_13
ATOM	113	CG	TYR	17	80.448	46.236	5.580	1.00	26.23	A_13
ATOM	114	CD1	TYR	17	81.824	46.081	5.412	1.00	16.37	A_13
ATOM	115	CE1	TYR	17	82.724	46.981	5.988	1.00	12.90	A_13
ATOM	116	CD2	TYR	17	79.990	47.329	6.331	1.00	17.15	A_13
ATOM	117	CE2	TYR	17	80.880	48.235	6.912	1.00	24.15	A_13
ATOM	118	CZ	TYR	17	82.244	48.057	6.743	1.00	23.38	A_13
ATOM	119	OH	TYR	17	83.121	48.942	7.343	1.00	19.47	A_13
ATOM	121	C	TYR	17	78.573	46.740	3.343	1.00	10.00	A_13
ATOM	122	O	TYR	17	79.298	47.559	2.782	1.00	19.27	A_13
ATOM	123	N	ARG	18	77.349	47.019	3.762	1.00	18.52	A_13
ATOM	125	CA	ARG	18	76.762	48.332	3.577	1.00	10.00	A_13
ATOM	126	CB	ARG	18	75.970	48.363	2.274	1.00	10.00	A_13
ATOM	127	CG	ARG	18	75.134	49.619	2.094	1.00	14.01	A_13
ATOM	128	CD	ARG	18	74.266	49.524	0.846	1.00	13.91	A_13
ATOM	129	NE	ARG	18	73.298	50.615	0.782	1.00	13.55	A_13
ATOM	131	CZ	ARG	18	72.165	50.571	0.092	1.00	10.00	A_13
ATOM	132	NH1	ARG	18	71.855	49.488	-0.602	1.00	14.30	A_13
ATOM	135	NH2	ARG	18	71.331	51.604	0.125	1.00	28.79	A_13
ATOM	138	C	ARG	18	75.842	48.640	4.741	1.00	10.65	A_13
ATOM	139	O	ARG	18	75.037	47.796	5.141	1.00	12.86	A_13
ATOM	140	N	ILE	19	76.014	49.814	5.332	1.00	25.54	A_13
ATOM	142	CA	ILE	19	75.169	50.265	6.436	1.00	24.52	A_13
ATOM	143	CB	ILE	19	75.944	51.236	7.350	1.00	18.37	A_13
ATOM	144	CG2	ILE	19	75.034	51.765	8.485	1.00	13.87	A_13
ATOM	145	CG1	ILE	19	77.204	50.545	7.888	1.00	27.67	A_13
ATOM	146	CD1	ILE	19	78.203	51.501	8.557	1.00	22.81	A_13
ATOM	147	C	ILE	19	74.062	51.027	5.698	1.00	21.11	A_13
ATOM	148	O	ILE	19	74.261	52.179	5.300	1.00	10.00	A_13
ATOM	149	N	VAL	20	72.916	50.378	5.487	1.00	19.76	A_13
ATOM	151	CA	VAL	20	71.829	51.014	4.735	1.00	18.20	A_13
ATOM	152	CB	VAL	20	70.774	49.983	4.193	1.00	15.42	A_13
ATOM	153	CG1	VAL	20	71.384	48.570	4.088	1.00	10.00	A_13
ATOM	154	CG2	VAL	20	69.496	50.030	4.992	1.00	18.62	A_13
ATOM	155	C	VAL	20	71.175	52.206	5.443	1.00	11.67	A_13
ATOM	156	O	VAL	20	70.652	53.110	4.798	1.00	18.36	A_13
ATOM	157	N	ASN	21	71.153	52.187	6.773	1.00	10.94	A_13
ATOM	159	CA	ASN	21	70.609	53.316	7.544	1.00	11.99	A_13
ATOM	160	CB	ASN	21	69.078	53.307	7.675	1.00	10.00	A_13
ATOM	161	CG	ASN	21	68.533	51.978	8.107	1.00	14.93	A_13
ATOM	162	OD1	ASN	21	67.627	51.449	7.486	1.00	21.54	A_13
ATOM	163	ND2	ASN	21	69.105	51.408	9.148	1.00	10.00	A_13
ATOM	166	C	ASN	21	71.291	53.382	8.897	1.00	18.90	A_13
ATOM	167	O	ASN	21	72.006	52.447	9.283	1.00	12.49	A_13
ATOM	168	N	TYR	22	71.053	54.471	9.618	1.00	17.47	A_13
ATOM	170	CA	TYR	22	71.681	54.708	10.910	1.00	24.85	A_13
ATOM	171	CB	TYR	22	72.556	55.954	10.818	1.00	13.52	A_13
ATOM	172	CG	TYR	22	73.791	55.748	9.991	1.00	10.00	A_13
ATOM	173	CD1	TYR	22	75.033	55.600	10.598	1.00	14.05	A_13
ATOM	174	CE1	TYR	22	76.180	55.370	9.841	1.00	13.69	A_13
ATOM	175	CD2	TYR	22	73.717	55.663	8.608	1.00	10.00	A_13
ATOM	176	CE2	TYR	22	74.848	55.432	7.847	1.00	17.10	A_13
ATOM	177	CZ	TYR	22	76.077	55.288	8.476	1.00	14.43	A_13
ATOM	178	OH	TYR	22	77.204	55.072	7.737	1.00	10.00	A_13
ATOM	180	C	TYR	22	70.726	54.862	12.076	1.00	25.95	A_13
ATOM	181	O	TYR	22	69.593	55.311	11.916	1.00	10.00	A_13
ATOM	182	N	THR	23	71.187	54.483	13.259	1.00	20.30	A_13
ATOM	184	CA	THR	23	70.367	54.606	14.450	1.00	29.11	A_13
ATOM	185	CB	THR	23	70.821	53.635	15.584	1.00	10.90	A_13
ATOM	186	OG1	THR	23	70.136	53.968	16.792	1.00	10.00	A_13
ATOM	188	CG2	THR	23	72.328	53.752	15.852	1.00	16.51	A_13
ATOM	189	C	THR	23	70.459	56.038	14.959	1.00	18.14	A_13
ATOM	190	O	THR	23	71.360	56.785	14.575	1.00	10.00	A_13
ATOM	191	N	PRO	24	69.433	56.487	15.691	1.00	12.76	A_13
ATOM	192	CD	PRO	24	68.061	55.950	15.716	1.00	15.26	A_13
ATOM	193	CA	PRO	24	69.453	57.844	16.232	1.00	22.70	A_13
ATOM	194	CB	PRO	24	67.985	58.086	16.585	1.00	28.52	A_13
ATOM	195	CG	PRO	24	67.448	56.706	16.841	1.00	15.78	A_13

39/75

ATOM	196	C	PRO	24	70.346	57.945	17.475	1.00	24.52	A_13
ATOM	197	O	PRO	24	70.790	59.040	17.831	1.00	10.00	A_13
ATOM	198	N	ASP	25	70.614	56.797	18.105	1.00	11.82	A_13
ATOM	200	CA	ASP	25	71.416	56.721	19.336	1.00	12.31	A_13
ATOM	201	CB	ASP	25	71.339	55.317	19.917	1.00	25.26	A_13
ATOM	202	CG	ASP	25	69.927	54.782	19.977	1.00	10.00	A_13
ATOM	203	OD1	ASP	25	69.783	53.567	20.159	1.00	20.90	A_13
ATOM	204	OD2	ASP	25	68.960	55.558	19.841	1.00	18.45	A_13
ATOM	205	C	ASP	25	72.891	57.113	19.286	1.00	14.34	A_13
ATOM	206	O	ASP	25	73.449	57.511	20.301	1.00	11.77	A_13
ATOM	207	N	MET	26	73.546	56.873	18.157	1.00	20.78	A_13
ATOM	209	CA	MET	26	74.960	57.208	18.010	1.00	20.03	A_13
ATOM	210	CB	MET	26	75.791	55.928	17.916	1.00	13.86	A_13
ATOM	211	CG	MET	26	75.966	55.181	19.231	1.00	19.00	A_13
ATOM	212	SD	MET	26	76.043	53.404	18.941	1.00	14.67	A_13
ATOM	213	CE	MET	26	77.737	53.223	18.385	1.00	19.74	A_13
ATOM	214	C	MET	26	75.157	58.047	16.754	1.00	13.32	A_13
ATOM	215	O	MET	26	74.274	58.086	15.900	1.00	16.81	A_13
ATOM	216	N	THR	27	76.285	58.749	16.656	1.00	10.29	A_13
ATOM	218	CA	THR	27	76.568	59.564	15.470	1.00	17.00	A_13
ATOM	219	CB	THR	27	77.710	60.596	15.700	1.00	11.79	A_13
ATOM	220	OG1	THR	27	78.969	59.921	15.729	1.00	23.77	A_13
ATOM	222	CG2	THR	27	77.519	61.342	17.020	1.00	21.98	A_13
ATOM	223	C	THR	27	76.996	58.634	14.347	1.00	13.37	A_13
ATOM	224	O	THR	27	77.411	57.500	14.608	1.00	11.05	A_13
ATOM	225	N	HIS	28	76.972	59.124	13.113	1.00	10.00	A_13
ATOM	227	CA	HIS	28	77.362	58.300	11.980	1.00	10.96	A_13
ATOM	228	CB	HIS	28	77.240	59.071	10.657	1.00	16.07	A_13
ATOM	229	CG	HIS	28	75.829	59.382	10.264	1.00	15.53	A_13
ATOM	230	CD2	HIS	28	74.707	59.531	11.016	1.00	21.47	A_13
ATOM	231	ND1	HIS	28	75.440	59.597	8.959	1.00	30.32	A_13
ATOM	233	CE1	HIS	28	74.149	59.868	8.920	1.00	19.38	A_13
ATOM	234	NE2	HIS	28	73.680	59.833	10.160	1.00	29.43	A_13
ATOM	236	C	HIS	28	78.769	57.735	12.151	1.00	14.80	A_13
ATOM	237	O	HIS	28	79.005	56.568	11.851	1.00	28.24	A_13
ATOM	238	N	SER	29	79.703	58.548	12.634	1.00	14.00	A_13
ATOM	240	CA	SER	29	81.068	58.070	12.854	1.00	19.57	A_13
ATOM	241	CB	SER	29	82.001	59.219	13.242	1.00	17.84	A_13
ATOM	242	OG	SER	29	82.383	59.936	12.084	1.00	28.25	A_13
ATOM	244	C	SER	29	81.134	56.983	13.917	1.00	15.23	A_13
ATOM	245	O	SER	29	81.818	55.973	13.733	1.00	13.73	A_13
ATOM	246	N	GLU	30	80.428	57.182	15.027	1.00	27.71	A_13
ATOM	248	CA	GLU	30	80.430	56.186	16.100	1.00	23.60	A_13
ATOM	249	CB	GLU	30	79.571	56.635	17.289	1.00	21.72	A_13
ATOM	250	CG	GLU	30	80.048	57.913	17.973	1.00	24.07	A_13
ATOM	251	CD	GLU	30	79.205	58.279	19.185	1.00	21.06	A_13
ATOM	252	OE1	GLU	30	79.784	58.660	20.218	1.00	46.95	A_13
ATOM	253	OE2	GLU	30	77.963	58.185	19.119	1.00	18.27	A_13
ATOM	254	C	GLU	30	79.895	54.877	15.553	1.00	18.75	A_13
ATOM	255	O	GLU	30	80.456	53.809	15.815	1.00	13.06	A_13
ATOM	256	N	VAL	31	78.839	54.970	14.746	1.00	16.23	A_13
ATOM	258	CA	VAL	31	78.225	53.781	14.146	1.00	22.33	A_13
ATOM	259	CB	VAL	31	76.899	54.135	13.390	1.00	23.53	A_13
ATOM	260	CG1	VAL	31	76.384	52.920	12.628	1.00	14.39	A_13
ATOM	261	CG2	VAL	31	75.829	54.587	14.377	1.00	10.00	A_13
ATOM	262	C	VAL	31	79.208	53.040	13.216	1.00	20.29	A_13
ATOM	263	O	VAL	31	79.330	51.814	13.282	1.00	14.02	A_13
ATOM	264	N	GLU	32	79.913	53.790	12.370	1.00	23.94	A_13
ATOM	266	CA	GLU	32	80.887	53.219	11.446	1.00	10.18	A_13
ATOM	267	CB	GLU	32	81.406	54.285	10.502	1.00	16.50	A_13
ATOM	268	CG	GLU	32	80.424	54.605	9.427	1.00	20.84	A_13
ATOM	269	CD	GLU	32	80.330	56.080	9.155	1.00	22.31	A_13
ATOM	270	OE1	GLU	32	79.285	56.509	8.639	1.00	29.39	A_13
ATOM	271	OE2	GLU	32	81.294	56.812	9.458	1.00	22.01	A_13
ATOM	272	C	GLU	32	82.056	52.565	12.137	1.00	18.93	A_13
ATOM	273	O	GLU	32	82.474	51.470	11.753	1.00	24.42	A_13
ATOM	274	N	LYS	33	82.610	53.241	13.139	1.00	19.78	A_13
ATOM	276	CA	LYS	33	83.726	52.661	13.873	1.00	28.68	A_13
ATOM	277	CB	LYS	33	84.340	53.681	14.837	1.00	18.54	A_13
ATOM	278	CG	LYS	33	85.016	54.855	14.135	1.00	31.19	A_13
ATOM	279	CD	LYS	33	86.135	54.425	13.148	1.00	40.31	A_13
ATOM	280	CE	LYS	33	85.600	53.972	11.785	1.00	21.99	A_13
ATOM	281	NZ	LYS	33	86.646	53.779	10.773	1.00	33.20	A_13
ATOM	285	C	LYS	33	83.242	51.407	14.594	1.00	12.66	A_13
ATOM	286	O	LYS	33	83.892	50.361	14.552	1.00	15.54	A_13
ATOM	287	N	ALA	34	82.036	51.481	15.148	1.00	20.70	A_13
ATOM	289	CA	ALA	34	81.453	50.344	15.843	1.00	10.00	A_13

40/75

ATOM	290	CB	ALA	34	80.040	50.651	16.279	1.00	18.59	A_13
ATOM	291	C	ALA	34	81.468	49.119	14.940	1.00	13.45	A_13
ATOM	292	O	ALA	34	82.067	48.095	15.284	1.00	15.90	A_13
ATOM	293	N	PHE	35	80.857	49.234	13.766	1.00	19.57	A_13
ATOM	295	CA	PHE	35	80.802	48.112	12.812	1.00	26.77	A_13
ATOM	296	CB	PHE	35	79.837	48.423	11.660	1.00	17.34	A_13
ATOM	297	CG	PHE	35	78.390	48.477	12.077	1.00	30.55	A_13
ATOM	298	CD1	PHE	35	77.838	47.464	12.863	1.00	26.58	A_13
ATOM	299	CD2	PHE	35	77.570	49.512	11.653	1.00	10.00	A_13
ATOM	300	CE1	PHE	35	76.494	47.485	13.212	1.00	12.45	A_13
ATOM	301	CE2	PHE	35	76.224	49.538	12.002	1.00	17.92	A_13
ATOM	302	CZ	PHE	35	75.684	48.525	12.777	1.00	13.29	A_13
ATOM	303	C	PHE	35	82.170	47.754	12.236	1.00	11.31	A_13
ATOM	304	O	PHE	35	82.493	46.573	12.034	1.00	11.37	A_13
ATOM	305	N	LYS	36	82.962	48.778	11.945	1.00	17.06	A_13
ATOM	307	CA	LYS	36	84.293	48.573	11.400	1.00	17.41	A_13
ATOM	308	CB	LYS	36	84.991	49.922	11.208	1.00	11.20	A_13
ATOM	309	CG	LYS	36	86.282	49.792	10.439	1.00	28.84	A_13
ATOM	310	CD	LYS	36	87.246	50.917	10.738	1.00	24.52	A_13
ATOM	311	CE	LYS	36	88.542	50.703	9.978	1.00	12.87	A_13
ATOM	312	NZ	LYS	36	88.264	50.536	8.514	1.00	23.69	A_13
ATOM	316	C	LYS	36	85.122	47.685	12.345	1.00	16.09	A_13
ATOM	317	O	LYS	36	85.701	46.686	11.938	1.00	21.50	A_13
ATOM	318	N	LYS	37	85.173	48.057	13.613	1.00	12.42	A_13
ATOM	320	CA	LYS	37	85.926	47.303	14.591	1.00	12.36	A_13
ATOM	321	CB	LYS	37	85.953	48.066	15.917	1.00	13.65	A_13
ATOM	322	CG	LYS	37	86.744	47.374	17.028	1.00	13.38	A_13
ATOM	323	CD	LYS	37	88.192	47.125	16.616	1.00	38.32	A_13
ATOM	324	CE	LYS	37	88.750	45.825	17.205	1.00	34.46	A_13
ATOM	325	NZ	LYS	37	88.234	44.576	16.557	1.00	12.49	A_13
ATOM	329	C	LYS	37	85.372	45.887	14.786	1.00	17.04	A_13
ATOM	330	O	LYS	37	86.131	44.958	15.053	1.00	18.14	A_13
ATOM	331	N	ALA	38	84.061	45.711	14.649	1.00	24.47	A_13
ATOM	333	CA	ALA	38	83.452	44.392	14.822	1.00	11.03	A_13
ATOM	334	CB	ALA	38	81.941	44.504	14.890	1.00	14.71	A_13
ATOM	335	C	ALA	38	83.900	43.451	13.697	1.00	20.27	A_13
ATOM	336	O	ALA	38	84.143	42.266	13.936	1.00	18.80	A_13
ATOM	337	N	PHE	39	84.021	43.971	12.477	1.00	22.58	A_13
ATOM	339	CA	PHE	39	84.492	43.158	11.355	1.00	18.87	A_13
ATOM	340	CB	PHE	39	84.350	43.899	10.027	1.00	19.91	A_13
ATOM	341	CG	PHE	39	82.993	43.783	9.414	1.00	10.00	A_13
ATOM	342	CD1	PHE	39	82.266	44.915	9.097	1.00	17.54	A_13
ATOM	343	CD2	PHE	39	82.438	42.533	9.143	1.00	15.92	A_13
ATOM	344	CE1	PHE	39	81.008	44.808	8.520	1.00	20.75	A_13
ATOM	345	CE2	PHE	39	81.186	42.418	8.569	1.00	10.00	A_13
ATOM	346	CZ	PHE	39	80.467	43.555	8.252	1.00	10.00	A_13
ATOM	347	C	PHE	39	85.955	42.827	11.589	1.00	16.52	A_13
ATOM	348	O	PHE	39	86.382	41.689	11.387	1.00	19.70	A_13
ATOM	349	N	LYS	40	86.699	43.822	12.072	1.00	21.31	A_13
ATOM	351	CA	LYS	40	88.117	43.673	12.369	1.00	20.07	A_13
ATOM	352	CB	LYS	40	88.703	44.967	12.927	1.00	13.77	A_13
ATOM	353	CG	LYS	40	90.192	44.885	13.171	1.00	11.54	A_13
ATOM	354	CD	LYS	40	90.757	46.242	13.507	1.00	10.34	A_13
ATOM	355	CE	LYS	40	92.236	46.142	13.838	1.00	11.24	A_13
ATOM	356	NZ	LYS	40	92.468	45.518	15.179	1.00	27.33	A_13
ATOM	360	C	LYS	40	88.352	42.534	13.337	1.00	12.06	A_13
ATOM	361	O	LYS	40	89.252	41.719	13.124	1.00	25.09	A_13
ATOM	362	N	VAL	41	87.495	42.418	14.349	1.00	12.26	A_13
ATOM	364	CA	VAL	41	87.630	41.331	15.325	1.00	17.89	A_13
ATOM	365	CB	VAL	41	86.351	41.205	16.216	1.00	10.00	A_13
ATOM	366	CG1	VAL	41	86.298	39.865	16.894	1.00	23.82	A_13
ATOM	367	CG2	VAL	41	86.329	42.274	17.259	1.00	17.65	A_13
ATOM	368	C	VAL	41	87.822	40.009	14.560	1.00	23.06	A_13
ATOM	369	O	VAL	41	88.664	39.168	14.912	1.00	11.82	A_13
ATOM	370	N	TRP	42	87.069	39.871	13.471	1.00	21.42	A_13
ATOM	372	CA	TRP	42	87.085	38.666	12.661	1.00	21.32	A_13
ATOM	373	CB	TRP	42	85.713	38.476	12.009	1.00	18.84	A_13
ATOM	374	CG	TRP	42	84.605	38.387	13.025	1.00	25.92	A_13
ATOM	375	CD2	TRP	42	84.437	37.369	14.024	1.00	16.65	A_13
ATOM	376	CE2	TRP	42	83.260	37.680	14.737	1.00	17.58	A_13
ATOM	377	CE3	TRP	42	85.165	36.223	14.380	1.00	11.14	A_13
ATOM	378	CD1	TRP	42	83.563	39.249	13.179	1.00	10.00	A_13
ATOM	379	NE1	TRP	42	82.755	38.832	14.200	1.00	10.91	A_13
ATOM	381	CZ2	TRP	42	82.785	36.879	15.793	1.00	14.81	A_13
ATOM	382	CZ3	TRP	42	84.691	35.425	15.436	1.00	23.68	A_13
ATOM	383	CH2	TRP	42	83.513	35.759	16.125	1.00	12.75	A_13
ATOM	384	C	TRP	42	88.190	38.600	11.623	1.00	27.45	A_13

41/75

ATOM	385	O	TRP	42	88.834	37.556	11.472	1.00	11.84	A_13
ATOM	386	N	SER	43	88.413	39.702	10.909	1.00	25.46	A_13
ATOM	388	CA	SER	43	89.449	39.740	9.881	1.00	19.61	A_13
ATOM	389	CB	SER	43	89.342	40.993	8.991	1.00	16.16	A_13
ATOM	390	OG	SER	43	89.495	42.199	9.709	1.00	26.34	A_13
ATOM	392	C	SER	43	90.837	39.615	10.491	1.00	11.53	A_13
ATOM	393	O	SER	43	91.758	39.119	9.834	1.00	17.99	A_13
ATOM	394	N	ASP	44	90.949	39.973	11.771	1.00	10.00	A_13
ATOM	396	CA	ASP	44	92.206	39.908	12.505	1.00	16.90	A_13
ATOM	397	CB	ASP	44	92.057	40.588	13.857	1.00	17.79	A_13
ATOM	398	CG	ASP	44	92.544	42.013	13.839	1.00	15.93	A_13
ATOM	399	OD1	ASP	44	92.605	42.618	14.920	1.00	17.21	A_13
ATOM	400	OD2	ASP	44	92.874	42.533	12.754	1.00	19.50	A_13
ATOM	401	C	ASP	44	92.781	38.523	12.729	1.00	26.12	A_13
ATOM	402	O	ASP	44	93.996	38.362	12.897	1.00	21.21	A_13
ATOM	403	N	VAL	45	91.911	37.523	12.745	1.00	20.89	A_13
ATOM	405	CA	VAL	45	92.353	36.161	12.996	1.00	27.53	A_13
ATOM	406	CB	VAL	45	91.853	35.678	14.381	1.00	16.30	A_13
ATOM	407	CG1	VAL	45	92.557	36.472	15.504	1.00	10.00	A_13
ATOM	408	CG2	VAL	45	90.348	35.857	14.495	1.00	10.86	A_13
ATOM	409	C	VAL	45	91.928	35.187	11.911	1.00	24.33	A_13
ATOM	410	O	VAL	45	91.864	33.978	12.157	1.00	18.84	A_13
ATOM	411	N	THR	46	91.750	35.705	10.694	1.00	16.30	A_13
ATOM	413	CA	THR	46	91.293	34.893	9.574	1.00	14.48	A_13
ATOM	414	CB	THR	46	89.750	34.796	9.662	1.00	22.05	A_13
ATOM	415	OG1	THR	46	89.279	33.609	9.028	1.00	31.53	A_13
ATOM	417	CG2	THR	46	89.112	36.014	9.040	1.00	10.99	A_13
ATOM	418	C	THR	46	91.716	35.575	8.257	1.00	25.10	A_13
ATOM	419	O	THR	46	92.022	36.764	8.256	1.00	17.64	A_13
ATOM	420	N	PRO	47	91.688	34.845	7.114	1.00	15.31	A_13
ATOM	421	CD	PRO	47	91.459	33.398	6.985	1.00	17.94	A_13
ATOM	422	CA	PRO	47	92.069	35.416	5.815	1.00	21.50	A_13
ATOM	423	CB	PRO	47	92.199	34.182	4.911	1.00	17.57	A_13
ATOM	424	CG	PRO	47	92.369	33.041	5.848	1.00	27.45	A_13
ATOM	425	C	PRO	47	90.991	36.348	5.256	1.00	21.44	A_13
ATOM	426	O	PRO	47	91.095	36.788	4.116	1.00	11.08	A_13
ATOM	427	N	LEU	48	89.918	36.567	6.018	1.00	10.00	A_13
ATOM	429	CA	LEU	48	88.826	37.434	5.581	1.00	22.09	A_13
ATOM	430	CB	LEU	48	87.575	37.212	6.432	1.00	15.92	A_13
ATOM	431	CG	LEU	48	86.848	35.867	6.435	1.00	13.58	A_13
ATOM	432	CD1	LEU	48	85.931	35.811	7.654	1.00	25.90	A_13
ATOM	433	CD2	LEU	48	86.073	35.666	5.157	1.00	16.47	A_13
ATOM	434	C	LEU	48	89.156	38.916	5.641	1.00	21.20	A_13
ATOM	435	O	LEU	48	89.936	39.366	6.480	1.00	17.28	A_13
ATOM	436	N	ASN	49	88.569	39.670	4.723	1.00	26.12	A_13
ATOM	438	CA	ASN	49	88.738	41.112	4.717	1.00	26.84	A_13
ATOM	439	CB	ASN	49	89.936	41.569	3.885	1.00	18.29	A_13
ATOM	440	CG	ASN	49	90.010	40.912	2.568	1.00	22.55	A_13
ATOM	441	OD1	ASN	49	90.928	40.131	2.305	1.00	24.41	A_13
ATOM	442	ND2	ASN	49	89.068	41.235	1.693	1.00	46.51	A_13
ATOM	445	C	ASN	49	87.416	41.705	4.259	1.00	12.18	A_13
ATOM	446	O	ASN	49	86.732	41.128	3.400	1.00	20.77	A_13
ATOM	447	N	PHE	50	87.025	42.802	4.900	1.00	21.39	A_13
ATOM	449	CA	PHE	50	85.738	43.439	4.642	1.00	10.00	A_13
ATOM	450	CB	PHE	50	84.914	43.440	5.932	1.00	11.45	A_13
ATOM	451	CG	PHE	50	84.863	42.098	6.629	1.00	10.63	A_13
ATOM	452	CD1	PHE	50	85.886	41.705	7.490	1.00	10.00	A_13
ATOM	453	CD2	PHE	50	83.809	41.216	6.395	1.00	14.63	A_13
ATOM	454	CE1	PHE	50	85.858	40.457	8.097	1.00	26.88	A_13
ATOM	455	CE2	PHE	50	83.773	39.963	7.000	1.00	21.13	A_13
ATOM	456	CZ	PHE	50	84.801	39.581	7.852	1.00	10.30	A_13
ATOM	457	C	PHE	50	85.867	44.842	4.093	1.00	22.56	A_13
ATOM	458	O	PHE	50	86.638	45.644	4.612	1.00	19.33	A_13
ATOM	459	N	THR	51	85.099	45.129	3.044	1.00	21.47	A_13
ATOM	461	CA	THR	51	85.125	46.433	2.371	1.00	24.21	A_13
ATOM	462	CB	THR	51	85.602	46.306	0.895	1.00	15.39	A_13
ATOM	463	OG1	THR	51	86.950	45.811	0.853	1.00	24.33	A_13
ATOM	465	CG2	THR	51	85.551	47.654	0.192	1.00	25.47	A_13
ATOM	466	C	THR	51	83.735	47.048	2.359	1.00	22.17	A_13
ATOM	467	O	THR	51	82.766	46.421	1.912	1.00	20.53	A_13
ATOM	468	N	ARG	52	83.653	48.294	2.797	1.00	16.53	A_13
ATOM	470	CA	ARG	52	82.393	49.004	2.871	1.00	10.00	A_13
ATOM	471	CB	ARG	52	82.490	50.085	3.939	1.00	10.00	A_13
ATOM	472	CG	ARG	52	81.201	50.778	4.259	1.00	12.47	A_13
ATOM	473	CD	ARG	52	81.462	51.879	5.278	1.00	19.61	A_13
ATOM	474	NE	ARG	52	80.371	52.836	5.333	1.00	30.55	A_13
ATOM	476	CZ	ARG	52	80.489	54.074	5.795	1.00	24.06	A_13

42/75

ATOM	477	NH1	ARG	52	81.661	54.508	6.257	1.00	21.24	A_13
ATOM	480	NH2	ARG	52	79.421	54.862	5.829	1.00	27.78	A_13
ATOM	483	C	ARG	52	81.980	49.620	1.540	1.00	30.22	A_13
ATOM	484	O	ARG	52	82.782	50.269	0.859	1.00	16.27	A_13
ATOM	485	N	LEU	53	80.730	49.372	1.161	1.00	21.07	A_13
ATOM	487	CA	LEU	53	80.159	49.914	-0.062	1.00	15.73	A_13
ATOM	488	CB	LEU	53	79.435	48.831	-0.868	1.00	11.53	A_13
ATOM	489	CG	LEU	53	80.304	47.770	-1.530	1.00	10.00	A_13
ATOM	490	CD1	LEU	53	79.429	46.790	-2.296	1.00	13.21	A_13
ATOM	491	CD2	LEU	53	81.280	48.443	-2.448	1.00	12.78	A_13
ATOM	492	C	LEU	53	79.149	50.932	0.421	1.00	10.00	A_13
ATOM	493	O	LEU	53	78.463	50.713	1.411	1.00	13.62	A_13
ATOM	494	N	HIS	54	79.043	52.041	-0.283	1.00	15.73	A_13
ATOM	496	CA	HIS	54	78.102	53.065	0.126	1.00	12.47	A_13
ATOM	497	CB	HIS	54	78.765	54.435	0.011	1.00	15.18	A_13
ATOM	498	CG	HIS	54	79.967	54.589	0.884	1.00	21.27	A_13
ATOM	499	CD2	HIS	54	81.207	54.056	0.798	1.00	25.30	A_13
ATOM	500	ND1	HIS	54	79.951	55.338	2.043	1.00	16.48	A_13
ATOM	502	CE1	HIS	54	81.127	55.255	2.633	1.00	21.62	A_13
ATOM	503	NE2	HIS	54	81.910	54.482	1.899	1.00	29.91	A_13
ATOM	505	C	HIS	54	76.796	53.044	-0.664	1.00	15.50	A_13
ATOM	506	O	HIS	54	75.914	53.849	-0.403	1.00	21.80	A_13
ATOM	507	N	ASP	55	76.707	52.178	-1.671	1.00	18.31	A_13
ATOM	509	CA	ASP	55	75.509	52.077	-2.502	1.00	17.23	A_13
ATOM	510	CB	ASP	55	75.645	52.928	-3.773	1.00	19.94	A_13
ATOM	511	CG	ASP	55	75.864	54.393	-3.495	1.00	26.81	A_13
ATOM	512	OD1	ASP	55	75.059	54.991	-2.741	1.00	35.97	A_13
ATOM	513	OD2	ASP	55	76.839	54.948	-4.058	1.00	25.09	A_13
ATOM	514	C	ASP	55	75.343	50.645	-2.970	1.00	21.50	A_13
ATOM	515	O	ASP	55	76.286	49.862	-2.929	1.00	17.45	A_13
ATOM	516	N	GLY	56	74.160	50.337	-3.489	1.00	10.31	A_13
ATOM	518	CA	GLY	56	73.897	49.014	-4.014	1.00	13.67	A_13
ATOM	519	C	GLY	56	73.842	47.869	-3.030	1.00	17.61	A_13
ATOM	520	O	GLY	56	73.683	48.065	-1.825	1.00	12.57	A_13
ATOM	521	N	ILE	57	73.943	46.653	-3.560	1.00	22.27	A_13
ATOM	523	CA	ILE	57	73.895	45.460	-2.737	1.00	11.39	A_13
ATOM	524	CB	ILE	57	72.941	44.391	-3.347	1.00	22.87	A_13
ATOM	525	CG2	ILE	57	73.365	42.995	-2.955	1.00	22.98	A_13
ATOM	526	CG1	ILE	57	71.522	44.582	-2.787	1.00	30.87	A_13
ATOM	527	CD1	ILE	57	71.002	46.022	-2.796	1.00	28.15	A_13
ATOM	528	C	ILE	57	75.289	44.919	-2.446	1.00	22.32	A_13
ATOM	529	O	ILE	57	76.140	44.849	-3.332	1.00	25.00	A_13
ATOM	530	N	ALA	58	75.517	44.631	-1.168	1.00	25.02	A_13
ATOM	532	CA	ALA	58	76.773	44.105	-0.669	1.00	15.45	A_13
ATOM	533	CB	ALA	58	77.366	45.060	0.358	1.00	11.62	A_13
ATOM	534	C	ALA	58	76.438	42.780	-0.006	1.00	12.08	A_13
ATOM	535	O	ALA	58	75.289	42.521	0.307	1.00	13.30	A_13
ATOM	536	N	ASP	59	77.449	41.968	0.247	1.00	14.79	A_13
ATOM	538	CA	ASP	59	77.245	40.675	0.880	1.00	18.50	A_13
ATOM	539	CB	ASP	59	78.608	39.974	1.093	1.00	10.83	A_13
ATOM	540	CG	ASP	59	79.425	39.858	-0.210	1.00	23.35	A_13
ATOM	541	OD1	ASP	59	80.598	40.266	-0.236	1.00	17.98	A_13
ATOM	542	OD2	ASP	59	78.896	39.379	-1.230	1.00	16.89	A_13
ATOM	543	C	ASP	59	76.480	40.806	2.200	1.00	13.69	A_13
ATOM	544	O	ASP	59	75.402	40.227	2.380	1.00	15.93	A_13
ATOM	545	N	ILE	60	77.025	41.596	3.109	1.00	13.15	A_13
ATOM	547	CA	ILE	60	76.422	41.800	4.412	1.00	12.20	A_13
ATOM	548	CB	ILE	60	77.500	41.695	5.508	1.00	12.12	A_13
ATOM	549	CG2	ILE	60	76.921	42.060	6.864	1.00	19.27	A_13
ATOM	550	CG1	ILE	60	78.118	40.287	5.481	1.00	10.00	A_13
ATOM	551	CD1	ILE	60	79.330	40.120	6.360	1.00	10.00	A_13
ATOM	552	C	ILE	60	75.743	43.164	4.456	1.00	17.78	A_13
ATOM	553	O	ILE	60	76.410	44.193	4.478	1.00	18.65	A_13
ATOM	554	N	MET	61	74.416	43.168	4.431	1.00	12.54	A_13
ATOM	556	CA	MET	61	73.640	44.416	4.476	1.00	12.86	A_13
ATOM	557	CB	MET	61	72.385	44.314	3.604	1.00	18.16	A_13
ATOM	558	CG	MET	61	72.634	43.979	2.141	1.00	10.00	A_13
ATOM	559	SD	MET	61	73.374	45.314	1.251	1.00	10.69	A_13
ATOM	560	CE	MET	61	71.836	46.299	0.764	1.00	10.00	A_13
ATOM	561	C	MET	61	73.239	44.666	5.921	1.00	10.15	A_13
ATOM	562	O	MET	61	72.584	43.838	6.547	1.00	18.13	A_13
ATOM	563	N	ILE	62	73.706	45.784	6.456	1.00	15.60	A_13
ATOM	565	CA	ILE	62	73.452	46.170	7.837	1.00	18.55	A_13
ATOM	566	CB	ILE	62	74.723	46.828	8.437	1.00	10.00	A_13
ATOM	567	CG2	ILE	62	74.498	47.163	9.900	1.00	26.36	A_13
ATOM	568	CG1	ILE	62	75.936	45.897	8.302	1.00	11.04	A_13
ATOM	569	CD1	ILE	62	77.228	46.481	8.891	1.00	10.00	A_13

43/75

ATOM	570	C	ILE	62	72.289	47.172	7.920	1.00	17.99	A_13
ATOM	571	O	ILE	62	72.335	48.208	7.264	1.00	12.72	A_13
ATOM	572	N	SER	63	71.285	46.896	8.751	1.00	10.00	A_13
ATOM	574	CA	SER	63	70.149	47.803	8.882	1.00	12.52	A_13
ATOM	575	CB	SER	63	69.016	47.364	7.956	1.00	13.06	A_13
ATOM	576	OG	SER	63	68.448	46.146	8.415	1.00	27.90	A_13
ATOM	578	C	SER	63	69.625	47.854	10.314	1.00	13.14	A_13
ATOM	579	O	SER	63	69.869	46.951	11.101	1.00	22.10	A_13
ATOM	580	N	PHE	64	68.919	48.932	10.640	1.00	21.17	A_13
ATOM	582	CA	PHE	64	68.317	49.139	11.954	1.00	22.01	A_13
ATOM	583	CB	PHE	64	68.777	50.468	12.574	1.00	10.98	A_13
ATOM	584	CG	PHE	64	70.189	50.448	13.092	1.00	10.00	A_13
ATOM	585	CD1	PHE	64	70.473	49.885	14.322	1.00	10.00	A_13
ATOM	586	CD2	PHE	64	71.229	51.016	12.357	1.00	16.56	A_13
ATOM	587	CE1	PHE	64	71.777	49.885	14.825	1.00	10.00	A_13
ATOM	588	CE2	PHE	64	72.540	51.025	12.846	1.00	10.00	A_13
ATOM	589	CZ	PHE	64	72.812	50.459	14.081	1.00	18.83	A_13
ATOM	590	C	PHE	64	66.825	49.207	11.675	1.00	22.55	A_13
ATOM	591	O	PHE	64	66.405	49.940	10.779	1.00	19.49	A_13
ATOM	592	N	GLY	65	66.031	48.485	12.453	1.00	13.69	A_13
ATOM	594	CA	GLY	65	64.593	48.491	12.238	1.00	10.70	A_13
ATOM	595	C	GLY	65	63.894	48.138	13.521	1.00	12.62	A_13
ATOM	596	O	GLY	65	64.559	47.777	14.491	1.00	18.29	A_13
ATOM	597	N	ILE	66	62.577	48.309	13.565	1.00	13.69	A_13
ATOM	599	CA	ILE	66	61.803	47.968	14.760	1.00	21.58	A_13
ATOM	600	CB	ILE	66	61.227	49.228	15.503	1.00	30.51	A_13
ATOM	601	CG2	ILE	66	62.351	50.110	16.025	1.00	10.43	A_13
ATOM	602	CG1	ILE	66	60.332	50.062	14.586	1.00	14.56	A_13
ATOM	603	CD1	ILE	66	59.587	51.149	15.333	1.00	16.94	A_13
ATOM	604	C	ILE	66	60.662	47.030	14.361	1.00	10.81	A_13
ATOM	605	O	ILE	66	60.311	46.962	13.188	1.00	10.00	A_13
ATOM	606	N	LYS	67	60.143	46.271	15.330	1.00	10.00	A_13
ATOM	608	CA	LYS	67	59.036	45.327	15.103	1.00	10.23	A_13
ATOM	609	CB	LYS	67	57.689	46.042	15.268	1.00	10.29	A_13
ATOM	610	CG	LYS	67	57.584	46.895	16.510	1.00	14.63	A_13
ATOM	611	CD	LYS	67	57.646	46.056	17.774	1.00	14.94	A_13
ATOM	612	CE	LYS	67	57.382	46.923	18.986	1.00	22.99	A_13
ATOM	613	NZ	LYS	67	57.480	46.174	20.258	1.00	28.27	A_13
ATOM	617	C	LYS	67	59.113	44.633	13.726	1.00	17.91	A_13
ATOM	618	O	LYS	67	60.167	44.106	13.366	1.00	24.16	A_13
ATOM	619	N	GLU	68	58.027	44.690	12.949	1.00	12.72	A_13
ATOM	621	CA	GLU	68	57.960	44.067	11.624	1.00	16.06	A_13
ATOM	622	CB	GLU	68	56.505	44.019	11.128	1.00	26.89	A_13
ATOM	623	CG	GLU	68	55.566	43.258	12.087	1.00	36.97	A_13
ATOM	624	CD	GLU	68	54.217	43.973	12.381	1.00	41.61	A_13
ATOM	625	OE1	GLU	68	53.289	43.921	11.537	1.00	17.31	A_13
ATOM	626	OE2	GLU	68	54.074	44.561	13.485	1.00	26.72	A_13
ATOM	627	C	GLU	68	58.823	44.911	10.705	1.00	22.50	A_13
ATOM	628	O	GLU	68	58.587	46.093	10.532	1.00	20.64	A_13
ATOM	629	N	HIS	69	59.848	44.315	10.120	1.00	16.43	A_13
ATOM	631	CA	HIS	69	60.732	45.102	9.283	1.00	13.69	A_13
ATOM	632	CB	HIS	69	61.930	45.603	10.103	1.00	10.97	A_13
ATOM	633	CG	HIS	69	62.786	44.502	10.643	1.00	24.02	A_13
ATOM	634	CD2	HIS	69	63.873	43.876	10.133	1.00	10.00	A_13
ATOM	635	ND1	HIS	69	62.512	43.876	11.839	1.00	17.68	A_13
ATOM	637	CE1	HIS	69	63.384	42.912	12.041	1.00	12.53	A_13
ATOM	638	NE2	HIS	69	64.228	42.888	11.020	1.00	10.00	A_13
ATOM	639	C	HIS	69	61.214	44.469	7.983	1.00	21.28	A_13
ATOM	640	O	HIS	69	62.314	44.780	7.529	1.00	18.74	A_13
ATOM	641	N	GLY	70	60.451	43.537	7.411	1.00	13.11	A_13
ATOM	643	CA	GLY	70	60.832	42.968	6.127	1.00	10.00	A_13
ATOM	644	C	GLY	70	61.262	41.533	5.936	1.00	10.00	A_13
ATOM	645	O	GLY	70	61.523	41.125	4.794	1.00	15.12	A_13
ATOM	646	N	ASP	71	61.412	40.768	7.012	1.00	19.99	A_13
ATOM	648	CA	ASP	71	61.842	39.381	6.862	1.00	19.99	A_13
ATOM	649	CB	ASP	71	63.332	39.223	7.218	1.00	10.00	A_13
ATOM	650	CG	ASP	71	63.672	39.752	8.592	1.00	23.52	A_13
ATOM	651	OD1	ASP	71	64.846	40.110	8.803	1.00	13.38	A_13
ATOM	652	OD2	ASP	71	62.774	39.812	9.464	1.00	12.94	A_13
ATOM	653	C	ASP	71	60.998	38.377	7.632	1.00	22.07	A_13
ATOM	654	O	ASP	71	61.319	37.190	7.649	1.00	24.45	A_13
ATOM	655	N	PHE	72	59.946	38.865	8.292	1.00	14.15	A_13
ATOM	657	CA	PHE	72	59.040	38.035	9.094	1.00	10.00	A_13
ATOM	658	CB	PHE	72	58.410	36.905	8.272	1.00	10.00	A_13
ATOM	659	CG	PHE	72	57.360	37.387	7.332	1.00	10.00	A_13
ATOM	660	CD1	PHE	72	56.115	37.773	7.815	1.00	23.01	A_13
ATOM	661	CD2	PHE	72	57.624	37.507	5.973	1.00	12.52	A_13

44/75

ATOM	662	CE1	PHE	72	55.144	38.290	6.950	1.00	18.99	A_13
ATOM	663	CE2	PHE	72	56.662	38.023	5.091	1.00	13.37	A_13
ATOM	664	CZ	PHE	72	55.420	38.413	5.576	1.00	22.50	A_13
ATOM	665	C	PHE	72	59.634	37.523	10.392	1.00	16.31	A_13
ATOM	666	O	PHE	72	59.111	36.596	11.021	1.00	15.64	A_13
ATOM	667	N	TYR	73	60.737	38.141	10.793	1.00	18.10	A_13
ATOM	669	CA	TYR	73	61.407	37.827	12.046	1.00	14.01	A_13
ATOM	670	CB	TYR	73	62.845	37.331	11.803	1.00	21.08	A_13
ATOM	671	CG	TYR	73	62.915	35.965	11.138	1.00	22.48	A_13
ATOM	672	CD1	TYR	73	63.579	35.788	9.923	1.00	30.23	A_13
ATOM	673	CE1	TYR	73	63.615	34.538	9.291	1.00	24.04	A_13
ATOM	674	CD2	TYR	73	62.288	34.856	11.710	1.00	19.23	A_13
ATOM	675	CE2	TYR	73	62.320	33.606	11.083	1.00	29.35	A_13
ATOM	676	CZ	TYR	73	62.984	33.460	9.875	1.00	12.50	A_13
ATOM	677	OH	TYR	73	63.018	32.246	9.241	1.00	17.89	A_13
ATOM	679	C	TYR	73	61.360	39.203	12.721	1.00	22.00	A_13
ATOM	680	O	TYR	73	62.365	39.919	12.819	1.00	10.93	A_13
ATOM	681	N	PRO	74	60.175	39.570	13.221	1.00	19.94	A_13
ATOM	682	CD	PRO	74	58.969	38.723	13.278	1.00	15.69	A_13
ATOM	683	CA	PRO	74	59.934	40.843	13.886	1.00	16.75	A_13
ATOM	684	CB	PRO	74	58.417	40.836	14.067	1.00	17.27	A_13
ATOM	685	CG	PRO	74	58.131	39.407	14.335	1.00	16.24	A_13
ATOM	686	C	PRO	74	60.640	41.037	15.216	1.00	17.39	A_13
ATOM	687	O	PRO	74	60.779	40.105	16.023	1.00	10.00	A_13
ATOM	688	N	PHE	75	61.098	42.264	15.431	1.00	10.00	A_13
ATOM	690	CA	PHE	75	61.743	42.618	16.675	1.00	16.45	A_13
ATOM	691	CB	PHE	75	62.613	43.865	16.512	1.00	20.71	A_13
ATOM	692	CG	PHE	75	63.931	43.590	15.841	1.00	23.32	A_13
ATOM	693	CD1	PHE	75	64.694	42.482	16.200	1.00	12.03	A_13
ATOM	694	CD2	PHE	75	64.405	44.420	14.842	1.00	22.30	A_13
ATOM	695	CE1	PHE	75	65.905	42.214	15.572	1.00	17.64	A_13
ATOM	696	CE2	PHE	75	65.622	44.148	14.208	1.00	15.43	A_13
ATOM	697	CZ	PHE	75	66.367	43.044	14.576	1.00	10.00	A_13
ATOM	698	C	PHE	75	60.632	42.784	17.707	1.00	25.73	A_13
ATOM	699	O	PHE	75	59.443	42.778	17.370	1.00	18.57	A_13
ATOM	700	N	ASP	76	61.009	43.002	18.952	1.00	20.50	A_13
ATOM	702	CA	ASP	76	60.023	43.049	20.006	1.00	13.89	A_13
ATOM	703	CB	ASP	76	60.241	41.805	20.873	1.00	20.69	A_13
ATOM	704	CG	ASP	76	61.672	41.685	21.378	1.00	22.52	A_13
ATOM	705	OD1	ASP	76	61.947	40.771	22.174	1.00	20.06	A_13
ATOM	706	OD2	ASP	76	62.525	42.506	20.998	1.00	10.69	A_13
ATOM	707	C	ASP	76	59.971	44.277	20.900	1.00	25.20	A_13
ATOM	708	O	ASP	76	59.397	44.207	21.986	1.00	29.52	A_13
ATOM	709	N	GLY	77	60.585	45.379	20.488	1.00	10.00	A_13
ATOM	711	CA	GLY	77	60.575	46.553	21.334	1.00	10.00	A_13
ATOM	712	C	GLY	77	61.769	46.514	22.266	1.00	10.00	A_13
ATOM	713	O	GLY	77	62.735	45.797	21.987	1.00	18.49	A_13
ATOM	714	N	PRO	78	61.785	47.344	23.322	1.00	16.07	A_13
ATOM	715	CD	PRO	78	60.790	48.426	23.505	1.00	15.88	A_13
ATOM	716	CA	PRO	78	62.855	47.439	24.330	1.00	16.23	A_13
ATOM	717	CB	PRO	78	62.261	48.391	25.363	1.00	22.96	A_13
ATOM	718	CG	PRO	78	61.470	49.349	24.501	1.00	22.37	A_13
ATOM	719	C	PRO	78	63.150	46.090	24.969	1.00	25.32	A_13
ATOM	720	O	PRO	78	62.227	45.356	25.272	1.00	20.04	A_13
ATOM	721	N	SER	79	64.432	45.750	25.099	1.00	20.93	A_13
ATOM	723	CA	SER	79	64.878	44.478	25.689	1.00	20.51	A_13
ATOM	724	CB	SER	79	64.364	44.311	27.131	1.00	23.69	A_13
ATOM	725	OG	SER	79	65.028	45.211	28.006	1.00	33.37	A_13
ATOM	727	C	SER	79	64.557	43.248	24.863	1.00	20.39	A_13
ATOM	728	O	SER	79	64.124	43.362	23.708	1.00	17.27	A_13
ATOM	729	N	GLY	80	64.825	42.071	25.415	1.00	13.38	A_13
ATOM	731	CA	GLY	80	64.564	40.850	24.678	1.00	10.11	A_13
ATOM	732	C	GLY	80	65.471	40.808	23.458	1.00	13.15	A_13
ATOM	733	O	GLY	80	66.614	41.251	23.538	1.00	31.80	A_13
ATOM	734	N	LEU	81	64.939	40.393	22.310	1.00	29.05	A_13
ATOM	736	CA	LEU	81	65.720	40.317	21.078	1.00	29.63	A_13
ATOM	737	CB	LEU	81	64.789	40.033	19.905	1.00	19.67	A_13
ATOM	738	CG	LEU	81	65.121	38.872	18.971	1.00	21.79	A_13
ATOM	739	CD1	LEU	81	64.215	38.980	17.773	1.00	23.87	A_13
ATOM	740	CD2	LEU	81	66.590	38.918	18.518	1.00	22.09	A_13
ATOM	741	C	LEU	81	66.442	41.649	20.835	1.00	19.25	A_13
ATOM	742	O	LEU	81	65.808	42.700	20.872	1.00	14.95	A_13
ATOM	743	N	LEU	82	67.760	41.599	20.657	1.00	25.03	A_13
ATOM	745	CA	LEU	82	68.573	42.795	20.421	1.00	27.35	A_13
ATOM	746	CB	LEU	82	69.868	42.747	21.244	1.00	12.74	A_13
ATOM	747	CG	LEU	82	69.802	42.748	22.773	1.00	16.50	A_13
ATOM	748	CD1	LEU	82	68.590	43.520	23.263	1.00	17.99	A_13

45/75

ATOM	749	CD2	LEU	82	69.744	41.343	23.279	1.00	13.28	A_13
ATOM	750	C	LEU	82	68.938	42.945	18.949	1.00	24.79	A_13
ATOM	751	O	LEU	82	68.812	44.039	18.363	1.00	14.36	A_13
ATOM	752	N	ALA	83	69.387	41.839	18.359	1.00	21.15	A_13
ATOM	754	CA	ALA	83	69.790	41.819	16.961	1.00	15.64	A_13
ATOM	755	CB	ALA	83	71.180	42.410	16.820	1.00	15.74	A_13
ATOM	756	C	ALA	83	69.806	40.400	16.444	1.00	19.37	A_13
ATOM	757	O	ALA	83	69.864	39.458	17.227	1.00	20.42	A_13
ATOM	758	N	HIS	84	69.746	40.252	15.126	1.00	10.72	A_13
ATOM	760	CA	HIS	84	69.808	38.939	14.502	1.00	20.51	A_13
ATOM	761	CB	HIS	84	68.454	38.185	14.476	1.00	12.34	A_13
ATOM	762	CG	HIS	84	67.361	38.849	13.679	1.00	24.79	A_13
ATOM	763	CD2	HIS	84	67.381	39.489	12.488	1.00	10.00	A_13
ATOM	764	ND1	HIS	84	66.052	38.869	14.104	1.00	13.50	A_13
ATOM	766	CE1	HIS	84	65.307	39.497	13.210	1.00	14.37	A_13
ATOM	767	NE2	HIS	84	66.087	39.886	12.220	1.00	15.00	A_13
ATOM	768	C	HIS	84	70.418	39.088	13.130	1.00	22.78	A_13
ATOM	769	O	HIS	84	70.338	40.162	12.532	1.00	10.00	A_13
ATOM	770	N	ALA	85	71.086	38.027	12.685	1.00	13.43	A_13
ATOM	772	CA	ALA	85	71.746	37.983	11.402	1.00	10.00	A_13
ATOM	773	CB	ALA	85	73.234	38.132	11.596	1.00	10.05	A_13
ATOM	774	C	ALA	85	71.426	36.661	10.721	1.00	17.89	A_13
ATOM	775	O	ALA	85	70.900	35.746	11.346	1.00	19.43	A_13
ATOM	776	N	PHE	86	71.697	36.585	9.425	1.00	13.49	A_13
ATOM	778	CA	PHE	86	71.459	35.372	8.651	1.00	12.49	A_13
ATOM	779	CB	PHE	86	70.739	35.728	7.344	1.00	10.00	A_13
ATOM	780	CG	PHE	86	69.348	36.240	7.529	1.00	19.96	A_13
ATOM	781	CD1	PHE	86	68.252	35.434	7.212	1.00	21.89	A_13
ATOM	782	CD2	PHE	86	69.119	37.530	8.003	1.00	10.63	A_13
ATOM	783	CE1	PHE	86	66.946	35.900	7.364	1.00	16.59	A_13
ATOM	784	CE2	PHE	86	67.829	38.009	8.158	1.00	19.06	A_13
ATOM	785	CZ	PHE	86	66.732	37.194	7.838	1.00	24.79	A_13
ATOM	786	C	PHE	86	72.802	34.721	8.298	1.00	11.05	A_13
ATOM	787	O	PHE	86	73.774	35.435	8.041	1.00	25.56	A_13
ATOM	788	N	PRO	87	72.892	33.375	8.304	1.00	19.41	A_13
ATOM	789	CD	PRO	87	71.876	32.383	8.717	1.00	17.25	A_13
ATOM	790	CA	PRO	87	74.149	32.686	7.956	1.00	29.29	A_13
ATOM	791	CB	PRO	87	73.800	31.198	8.135	1.00	18.88	A_13
ATOM	792	CG	PRO	87	72.329	31.160	7.939	1.00	20.17	A_13
ATOM	793	C	PRO	87	74.562	32.999	6.503	1.00	10.00	A_13
ATOM	794	O	PRO	87	73.728	33.448	5.703	1.00	20.68	A_13
ATOM	795	N	PRO	88	75.814	32.701	6.120	1.00	10.00	A_13
ATOM	796	CD	PRO	88	76.796	31.854	6.831	1.00	19.58	A_13
ATOM	797	CA	PRO	88	76.280	32.977	4.756	1.00	12.43	A_13
ATOM	798	CB	PRO	88	77.600	32.201	4.676	1.00	18.69	A_13
ATOM	799	CG	PRO	88	78.073	32.163	6.098	1.00	18.48	A_13
ATOM	800	C	PRO	88	75.304	32.510	3.672	1.00	24.39	A_13
ATOM	801	O	PRO	88	74.596	31.522	3.854	1.00	16.92	A_13
ATOM	802	N	GLY	89	75.266	33.230	2.560	1.00	10.73	A_13
ATOM	804	CA	GLY	89	74.386	32.868	1.471	1.00	10.00	A_13
ATOM	805	C	GLY	89	73.960	34.127	0.772	1.00	10.94	A_13
ATOM	806	O	GLY	89	74.143	35.218	1.307	1.00	19.86	A_13
ATOM	807	N	PRO	90	73.390	34.019	-0.432	1.00	26.31	A_13
ATOM	808	CD	PRO	90	73.090	32.792	-1.192	1.00	18.46	A_13
ATOM	809	CA	PRO	90	72.960	35.212	-1.163	1.00	25.07	A_13
ATOM	810	CB	PRO	90	72.670	34.651	-2.556	1.00	15.47	A_13
ATOM	811	CG	PRO	90	72.108	33.289	-2.236	1.00	24.63	A_13
ATOM	812	C	PRO	90	71.726	35.879	-0.543	1.00	20.41	A_13
ATOM	813	O	PRO	90	71.176	35.390	0.442	1.00	17.00	A_13
ATOM	814	N	ASN	91	71.303	37.000	-1.125	1.00	18.43	A_13
ATOM	816	CA	ASN	91	70.127	37.721	-0.653	1.00	14.03	A_13
ATOM	817	CB	ASN	91	68.863	36.932	-0.999	1.00	15.26	A_13
ATOM	818	CG	ASN	91	68.860	36.430	-2.439	1.00	36.74	A_13
ATOM	819	OD1	ASN	91	68.497	35.282	-2.701	1.00	29.56	A_13
ATOM	820	ND2	ASN	91	69.265	37.286	-3.376	1.00	27.03	A_13
ATOM	823	C	ASN	91	70.226	37.986	0.849	1.00	24.66	A_13
ATOM	824	O	ASN	91	71.257	38.479	1.313	1.00	17.43	A_13
ATOM	825	N	TYR	92	69.198	37.632	1.622	1.00	17.69	A_13
ATOM	827	CA	TYR	92	69.233	37.876	3.061	1.00	10.17	A_13
ATOM	828	CB	TYR	92	67.942	37.428	3.744	1.00	16.78	A_13
ATOM	829	CG	TYR	92	66.786	38.364	3.523	1.00	26.17	A_13
ATOM	830	CD1	TYR	92	66.015	38.803	4.581	1.00	17.79	A_13
ATOM	831	CE1	TYR	92	64.947	39.678	4.380	1.00	29.60	A_13
ATOM	832	CD2	TYR	92	66.467	38.818	2.250	1.00	25.90	A_13
ATOM	833	CE2	TYR	92	65.406	39.691	2.040	1.00	30.60	A_13
ATOM	834	CZ	TYR	92	64.647	40.117	3.107	1.00	12.31	A_13
ATOM	835	OH	TYR	92	63.575	40.967	2.886	1.00	26.07	A_13

46/75

ATOM	837	C	TYR	92	70.427	37.245	3.763	1.00	11.94	A_13
ATOM	838	O	TYR	92	70.752	37.617	4.882	1.00	17.58	A_13
ATOM	839	N	GLY	93	71.095	36.311	3.097	1.00	24.67	A_13
ATOM	841	CA	GLY	93	72.250	35.666	3.691	1.00	18.05	A_13
ATOM	842	C	GLY	93	73.295	36.681	4.116	1.00	10.00	A_13
ATOM	843	O	GLY	93	73.573	37.656	3.391	1.00	10.13	A_13
ATOM	844	N	GLY	94	73.812	36.495	5.328	1.00	12.44	A_13
ATOM	846	CA	GLY	94	74.827	37.372	5.872	1.00	10.00	A_13
ATOM	847	C	GLY	94	74.358	38.694	6.456	1.00	17.29	A_13
ATOM	848	O	GLY	94	75.052	39.271	7.284	1.00	14.53	A_13
ATOM	849	N	ASP	95	73.221	39.206	5.993	1.00	10.00	A_13
ATOM	851	CA	ASP	95	72.689	40.485	6.472	1.00	16.35	A_13
ATOM	852	CB	ASP	95	71.332	40.777	5.814	1.00	10.00	A_13
ATOM	853	CG	ASP	95	71.421	40.904	4.309	1.00	14.54	A_13
ATOM	854	OD1	ASP	95	70.406	41.256	3.673	1.00	11.86	A_13
ATOM	855	OD2	ASP	95	72.502	40.647	3.753	1.00	15.39	A_13
ATOM	856	C	ASP	95	72.548	40.523	7.994	1.00	22.31	A_13
ATOM	857	O	ASP	95	72.279	39.497	8.635	1.00	10.88	A_13
ATOM	858	N	ALA	96	72.703	41.711	8.566	1.00	18.45	A_13
ATOM	860	CA	ALA	96	72.609	41.877	10.011	1.00	15.08	A_13
ATOM	861	CB	ALA	96	73.982	42.244	10.587	1.00	19.20	A_13
ATOM	862	C	ALA	96	71.587	42.961	10.345	1.00	14.91	A_13
ATOM	863	O	ALA	96	71.702	44.092	9.876	1.00	10.00	A_13
ATOM	864	N	HIS	97	70.635	42.646	11.215	1.00	14.01	A_13
ATOM	866	CA	HIS	97	69.599	43.620	11.581	1.00	11.35	A_13
ATOM	867	CB	HIS	97	68.207	43.083	11.203	1.00	20.32	A_13
ATOM	868	CG	HIS	97	68.027	42.786	9.742	1.00	15.00	A_13
ATOM	869	CD2	HIS	97	68.734	43.186	8.654	1.00	10.00	A_13
ATOM	870	ND1	HIS	97	67.014	41.978	9.257	1.00	14.03	A_13
ATOM	871	CE1	HIS	97	67.108	41.895	7.936	1.00	10.00	A_13
ATOM	872	NE2	HIS	97	68.142	42.618	7.552	1.00	17.10	A_13
ATOM	874	C	HIS	97	69.650	43.952	13.078	1.00	13.37	A_13
ATOM	875	O	HIS	97	69.736	43.055	13.908	1.00	13.48	A_13
ATOM	876	N	PHE	98	69.596	45.237	13.423	1.00	21.01	A_13
ATOM	878	CA	PHE	98	69.634	45.668	14.823	1.00	11.27	A_13
ATOM	879	CB	PHE	98	70.817	46.615	15.055	1.00	10.00	A_13
ATOM	880	CG	PHE	98	72.138	46.011	14.703	1.00	20.49	A_13
ATOM	881	CD1	PHE	98	72.984	45.524	15.707	1.00	17.49	A_13
ATOM	882	CD2	PHE	98	72.506	45.853	13.365	1.00	13.51	A_13
ATOM	883	CE1	PHE	98	74.171	44.888	15.382	1.00	20.00	A_13
ATOM	884	CE2	PHE	98	73.693	45.215	13.024	1.00	10.00	A_13
ATOM	885	CZ	PHE	98	74.527	44.728	14.029	1.00	10.00	A_13
ATOM	886	C	PHE	98	68.336	46.336	15.245	1.00	25.38	A_13
ATOM	887	O	PHE	98	67.815	47.218	14.552	1.00	10.00	A_13
ATOM	888	N	ASP	99	67.817	45.924	16.394	1.00	21.68	A_13
ATOM	890	CA	ASP	99	66.567	46.476	16.886	1.00	10.00	A_13
ATOM	891	CB	ASP	99	66.039	45.604	18.010	1.00	10.00	A_13
ATOM	892	CG	ASP	99	64.648	45.998	18.473	1.00	14.00	A_13
ATOM	893	OD1	ASP	99	64.104	45.272	19.329	1.00	15.19	A_13
ATOM	894	OD2	ASP	99	64.089	47.011	18.001	1.00	17.01	A_13
ATOM	895	C	ASP	99	66.817	47.871	17.391	1.00	13.06	A_13
ATOM	896	O	ASP	99	67.528	48.056	18.374	1.00	10.00	A_13
ATOM	897	N	ASP	100	66.203	48.856	16.746	1.00	15.56	A_13
ATOM	899	CA	ASP	100	66.397	50.232	17.177	1.00	18.23	A_13
ATOM	900	CB	ASP	100	66.121	51.228	16.041	1.00	15.05	A_13
ATOM	901	CG	ASP	100	67.275	52.180	15.838	1.00	11.67	A_13
ATOM	902	OD1	ASP	100	67.602	52.516	14.683	1.00	21.07	A_13
ATOM	903	OD2	ASP	100	67.879	52.569	16.860	1.00	14.72	A_13
ATOM	904	C	ASP	100	65.610	50.572	18.445	1.00	10.00	A_13
ATOM	905	O	ASP	100	65.767	51.635	19.009	1.00	17.18	A_13
ATOM	906	N	ASP	101	64.755	49.669	18.895	1.00	14.57	A_13
ATOM	908	CA	ASP	101	64.031	49.924	20.123	1.00	17.59	A_13
ATOM	909	CB	ASP	101	62.769	49.051	20.236	1.00	12.50	A_13
ATOM	910	CG	ASP	101	61.532	49.721	19.606	1.00	17.12	A_13
ATOM	911	OD1	ASP	101	60.599	49.023	19.179	1.00	10.39	A_13
ATOM	912	OD2	ASP	101	61.480	50.962	19.536	1.00	18.09	A_13
ATOM	913	C	ASP	101	64.994	49.766	21.306	1.00	19.33	A_13
ATOM	914	O	ASP	101	64.610	49.972	22.456	1.00	10.00	A_13
ATOM	915	N	GLU	102	66.213	49.301	21.019	1.00	16.15	A_13
ATOM	917	CA	GLU	102	67.267	49.194	22.044	1.00	13.43	A_13
ATOM	918	CB	GLU	102	68.264	48.085	21.720	1.00	18.25	A_13
ATOM	919	CG	GLU	102	67.697	46.704	21.636	1.00	10.00	A_13
ATOM	920	CD	GLU	102	66.650	46.467	22.672	1.00	11.18	A_13
ATOM	921	OE1	GLU	102	66.872	46.746	23.870	1.00	16.09	A_13
ATOM	922	OE2	GLU	102	65.572	46.033	22.271	1.00	26.76	A_13
ATOM	923	C	GLU	102	68.070	50.495	22.007	1.00	11.07	A_13
ATOM	924	O	GLU	102	68.103	51.161	20.971	1.00	13.97	A_13

47/75

ATOM	925	N	THR	103	68.774	50.823	23.091	1.00	22.82	A_13
ATOM	927	CA	THR	103	69.606	52.034	23.102	1.00	13.45	A_13
ATOM	928	CB	THR	103	69.571	52.793	24.459	1.00	20.78	A_13
ATOM	929	OG1	THR	103	68.236	53.228	24.745	1.00	10.69	A_13
ATOM	931	CG2	THR	103	70.445	54.046	24.378	1.00	19.45	A_13
ATOM	932	C	THR	103	71.030	51.571	22.822	1.00	12.42	A_13
ATOM	933	O	THR	103	71.639	50.896	23.642	1.00	19.81	A_13
ATOM	934	N	TRP	104	71.525	51.854	21.626	1.00	10.00	A_13
ATOM	936	CA	TRP	104	72.873	51.448	21.248	1.00	13.61	A_13
ATOM	937	CB	TRP	104	72.943	51.221	19.739	1.00	29.21	A_13
ATOM	938	CG	TRP	104	71.970	50.174	19.313	1.00	21.39	A_13
ATOM	939	CD2	TRP	104	72.101	48.760	19.501	1.00	25.13	A_13
ATOM	940	CE2	TRP	104	70.937	48.156	18.964	1.00	28.84	A_13
ATOM	941	CE3	TRP	104	73.088	47.941	20.070	1.00	13.36	A_13
ATOM	942	CD1	TRP	104	70.765	50.372	18.694	1.00	21.59	A_13
ATOM	943	NE1	TRP	104	70.139	49.163	18.484	1.00	19.91	A_13
ATOM	945	CZ2	TRP	104	70.738	46.768	18.977	1.00	10.00	A_13
ATOM	946	CZ3	TRP	104	72.888	46.568	20.084	1.00	14.54	A_13
ATOM	947	CH2	TRP	104	71.720	45.995	19.539	1.00	11.93	A_13
ATOM	948	C	TRP	104	73.912	52.453	21.725	1.00	16.59	A_13
ATOM	949	O	TRP	104	73.707	53.671	21.642	1.00	12.90	A_13
ATOM	950	N	THR	105	75.013	51.949	22.268	1.00	20.85	A_13
ATOM	952	CA	THR	105	76.040	52.831	22.794	1.00	12.38	A_13
ATOM	953	CB	THR	105	75.974	52.890	24.322	1.00	14.39	A_13
ATOM	954	OG1	THR	105	76.345	51.609	24.849	1.00	16.42	A_13
ATOM	956	CG2	THR	105	74.575	53.273	24.797	1.00	12.17	A_13
ATOM	957	C	THR	105	77.437	52.378	22.457	1.00	10.00	A_13
ATOM	958	O	THR	105	77.644	51.261	22.012	1.00	18.98	A_13
ATOM	959	N	SER	106	78.385	53.277	22.704	1.00	26.01	A_13
ATOM	961	CA	SER	106	79.809	53.043	22.502	1.00	17.80	A_13
ATOM	962	CB	SER	106	80.466	54.284	21.888	1.00	20.63	A_13
ATOM	963	OG	SER	106	79.744	54.756	20.763	1.00	38.89	A_13
ATOM	965	C	SER	106	80.435	52.779	23.880	1.00	34.75	A_13
ATOM	966	O	SER	106	81.652	52.884	24.042	1.00	33.01	A_13
ATOM	967	N	SER	107	79.590	52.494	24.875	1.00	25.87	A_13
ATOM	969	CA	SER	107	80.032	52.221	26.240	1.00	19.68	A_13
ATOM	970	CB	SER	107	80.082	53.510	27.061	1.00	23.47	A_13
ATOM	971	OG	SER	107	78.819	54.158	27.096	1.00	33.70	A_13
ATOM	973	C	SER	107	79.100	51.200	26.892	1.00	13.60	A_13
ATOM	974	O	SER	107	78.460	50.418	26.193	1.00	16.40	A_13
ATOM	975	N	SER	108	79.028	51.205	28.221	1.00	17.31	A_13
ATOM	977	CA	SER	108	78.188	50.259	28.949	1.00	20.12	A_13
ATOM	978	CB	SER	108	78.745	50.009	30.364	1.00	22.63	A_13
ATOM	979	OG	SER	108	78.444	51.061	31.271	1.00	27.69	A_13
ATOM	981	C	SER	108	76.702	50.606	29.076	1.00	19.98	A_13
ATOM	982	O	SER	108	75.921	49.785	29.562	1.00	35.96	A_13
ATOM	983	N	LYS	109	76.311	51.820	28.713	1.00	16.24	A_13
ATOM	985	CA	LYS	109	74.907	52.186	28.847	1.00	11.10	A_13
ATOM	986	CB	LYS	109	74.740	53.688	28.690	1.00	12.41	A_13
ATOM	987	CG	LYS	109	73.555	54.239	29.462	1.00	32.67	A_13
ATOM	988	CD	LYS	109	73.353	55.732	29.258	1.00	25.94	A_13
ATOM	989	CE	LYS	109	74.535	56.599	29.749	1.00	25.11	A_13
ATOM	990	NZ	LYS	109	74.225	58.070	29.636	1.00	22.70	A_13
ATOM	994	C	LYS	109	74.138	51.424	27.773	1.00	21.67	A_13
ATOM	995	O	LYS	109	74.667	51.210	26.694	1.00	32.76	A_13
ATOM	996	N	GLY	110	72.932	50.955	28.081	1.00	29.60	A_13
ATOM	998	CA	GLY	110	72.156	50.206	27.096	1.00	10.31	A_13
ATOM	999	C	GLY	110	72.965	49.043	26.542	1.00	20.08	A_13
ATOM	1000	O	GLY	110	73.672	48.362	27.285	1.00	11.17	A_13
ATOM	1001	N	TYR	111	72.924	48.859	25.227	1.00	12.05	A_13
ATOM	1003	CA	TYR	111	73.665	47.791	24.583	1.00	13.45	A_13
ATOM	1004	CB	TYR	111	72.713	46.871	23.806	1.00	21.16	A_13
ATOM	1005	CG	TYR	111	71.776	46.101	24.716	1.00	12.28	A_13
ATOM	1006	CD1	TYR	111	70.455	46.510	24.906	1.00	14.85	A_13
ATOM	1007	CE1	TYR	111	69.618	45.837	25.795	1.00	19.08	A_13
ATOM	1008	CD2	TYR	111	72.232	44.995	25.435	1.00	21.86	A_13
ATOM	1009	CE2	TYR	111	71.405	44.314	26.324	1.00	10.00	A_13
ATOM	1010	CZ	TYR	111	70.101	44.740	26.505	1.00	18.51	A_13
ATOM	1011	OH	TYR	111	69.282	44.077	27.398	1.00	14.32	A_13
ATOM	1013	C	TYR	111	74.779	48.335	23.695	1.00	16.73	A_13
ATOM	1014	O	TYR	111	74.540	49.105	22.764	1.00	11.98	A_13
ATOM	1015	N	ASN	112	76.008	47.930	23.999	1.00	11.80	A_13
ATOM	1017	CA	ASN	112	77.184	48.357	23.240	1.00	16.37	A_13
ATOM	1018	CB	ASN	112	78.453	47.867	23.927	1.00	27.52	A_13
ATOM	1019	CG	ASN	112	79.701	48.460	23.324	1.00	20.16	A_13
ATOM	1020	OD1	ASN	112	80.327	47.861	22.447	1.00	20.99	A_13
ATOM	1021	ND2	ASN	112	80.082	49.640	23.801	1.00	15.12	A_13

48/75

ATOM	1024	C	ASN	112	77.137	47.809	21.813	1.00	18.08	A_13
ATOM	1025	O	ASN	112	77.288	46.606	21.592	1.00	12.69	A_13
ATOM	1026	N	LEU	113	76.972	48.700	20.844	1.00	11.15	A_13
ATOM	1028	CA	LEU	113	76.878	48.296	19.461	1.00	10.00	A_13
ATOM	1029	CB	LEU	113	76.718	49.526	18.568	1.00	10.24	A_13
ATOM	1030	CG	LEU	113	76.325	49.262	17.106	1.00	15.67	A_13
ATOM	1031	CD1	LEU	113	75.155	48.296	17.050	1.00	26.54	A_13
ATOM	1032	CD2	LEU	113	75.967	50.555	16.415	1.00	15.60	A_13
ATOM	1033	C	LEU	113	78.037	47.403	18.986	1.00	25.17	A_13
ATOM	1034	O	LEU	113	77.799	46.380	18.336	1.00	17.24	A_13
ATOM	1035	N	PHE	114	79.274	47.759	19.327	1.00	28.89	A_13
ATOM	1037	CA	PHE	114	80.442	46.974	18.910	1.00	19.15	A_13
ATOM	1038	CB	PHE	114	81.753	47.579	19.434	1.00	14.60	A_13
ATOM	1039	CG	PHE	114	82.923	46.627	19.374	1.00	18.53	A_13
ATOM	1040	CD1	PHE	114	83.419	46.175	18.144	1.00	26.13	A_13
ATOM	1041	CD2	PHE	114	83.514	46.162	20.547	1.00	17.22	A_13
ATOM	1042	CE1	PHE	114	84.475	45.271	18.086	1.00	10.43	A_13
ATOM	1043	CE2	PHE	114	84.571	45.259	20.502	1.00	16.51	A_13
ATOM	1044	CZ	PHE	114	85.052	44.815	19.260	1.00	15.54	A_13
ATOM	1045	C	PHE	114	80.359	45.508	19.306	1.00	10.00	A_13
ATOM	1046	O	PHE	114	80.437	44.625	18.445	1.00	33.07	A_13
ATOM	1047	N	LEU	115	80.206	45.249	20.600	1.00	12.18	A_13
ATOM	1049	CA	LEU	115	80.113	43.877	21.103	1.00	10.59	A_13
ATOM	1050	CB	LEU	115	79.874	43.895	22.616	1.00	14.14	A_13
ATOM	1051	CG	LEU	115	81.082	43.937	23.578	1.00	34.39	A_13
ATOM	1052	CD1	LEU	115	82.337	44.354	22.863	1.00	14.93	A_13
ATOM	1053	CD2	LEU	115	80.815	44.836	24.793	1.00	13.42	A_13
ATOM	1054	C	LEU	115	79.019	43.080	20.379	1.00	12.06	A_13
ATOM	1055	O	LEU	115	79.298	42.109	19.675	1.00	13.35	A_13
ATOM	1056	N	VAL	116	77.786	43.558	20.459	1.00	13.11	A_13
ATOM	1058	CA	VAL	116	76.678	42.875	19.814	1.00	12.97	A_13
ATOM	1059	CB	VAL	116	75.343	43.569	20.129	1.00	28.07	A_13
ATOM	1060	CG1	VAL	116	74.200	42.926	19.340	1.00	17.32	A_13
ATOM	1061	CG2	VAL	116	75.074	43.491	21.617	1.00	22.14	A_13
ATOM	1062	C	VAL	116	76.862	42.724	18.313	1.00	10.00	A_13
ATOM	1063	O	VAL	116	76.473	41.716	17.755	1.00	14.68	A_13
ATOM	1064	N	ALA	117	77.481	43.706	17.667	1.00	10.80	A_13
ATOM	1066	CA	ALA	117	77.726	43.662	16.224	1.00	18.28	A_13
ATOM	1067	CB	ALA	117	78.223	45.014	15.727	1.00	14.94	A_13
ATOM	1068	C	ALA	117	78.735	42.579	15.863	1.00	25.24	A_13
ATOM	1069	O	ALA	117	78.562	41.872	14.861	1.00	18.50	A_13
ATOM	1070	N	ALA	118	79.795	42.458	16.665	1.00	24.40	A_13
ATOM	1072	CA	ALA	118	80.829	41.451	16.422	1.00	11.80	A_13
ATOM	1073	CB	ALA	118	81.945	41.590	17.447	1.00	19.28	A_13
ATOM	1074	C	ALA	118	80.178	40.056	16.496	1.00	10.00	A_13
ATOM	1075	O	ALA	118	80.426	39.183	15.660	1.00	10.00	A_13
ATOM	1076	N	HIS	119	79.309	39.875	17.487	1.00	19.01	A_13
ATOM	1078	CA	HIS	119	78.587	38.624	17.674	1.00	14.36	A_13
ATOM	1079	CB	HIS	119	77.725	38.751	18.924	1.00	10.00	A_13
ATOM	1080	CG	HIS	119	76.796	37.602	19.166	1.00	10.00	A_13
ATOM	1081	CD2	HIS	119	75.691	37.187	18.498	1.00	14.94	A_13
ATOM	1082	ND1	HIS	119	76.905	36.783	20.263	1.00	20.37	A_13
ATOM	1084	CE1	HIS	119	75.917	35.909	20.270	1.00	17.53	A_13
ATOM	1085	NE2	HIS	119	75.161	36.134	19.208	1.00	17.55	A_13
ATOM	1086	C	HIS	119	77.741	38.339	16.419	1.00	10.00	A_13
ATOM	1087	O	HIS	119	77.779	37.245	15.856	1.00	10.64	A_13
ATOM	1088	N	GLU	120	77.004	39.343	15.968	1.00	22.95	A_13
ATOM	1090	CA	GLU	120	76.174	39.224	14.775	1.00	23.96	A_13
ATOM	1091	CB	GLU	120	75.429	40.545	14.502	1.00	17.19	A_13
ATOM	1092	CG	GLU	120	74.373	40.889	15.555	1.00	16.14	A_13
ATOM	1093	CD	GLU	120	73.492	39.691	15.929	1.00	10.00	A_13
ATOM	1094	OE1	GLU	120	73.478	39.354	17.122	1.00	17.94	A_13
ATOM	1095	OE2	GLU	120	72.844	39.078	15.047	1.00	17.03	A_13
ATOM	1096	C	GLU	120	76.992	38.832	13.549	1.00	11.45	A_13
ATOM	1097	O	GLU	120	76.594	37.946	12.772	1.00	13.34	A_13
ATOM	1098	N	PHE	121	78.127	39.498	13.353	1.00	10.00	A_13
ATOM	1100	CA	PHE	121	78.959	39.187	12.216	1.00	14.70	A_13
ATOM	1101	CB	PHE	121	80.040	40.245	12.039	1.00	10.00	A_13
ATOM	1102	CG	PHE	121	79.481	41.623	11.792	1.00	21.57	A_13
ATOM	1103	CD1	PHE	121	80.235	42.764	12.069	1.00	16.73	A_13
ATOM	1104	CD2	PHE	121	78.164	41.788	11.331	1.00	13.91	A_13
ATOM	1105	CE1	PHE	121	79.682	44.054	11.891	1.00	11.69	A_13
ATOM	1106	CE2	PHE	121	77.615	43.066	11.152	1.00	18.93	A_13
ATOM	1107	CZ	PHE	121	78.373	44.192	11.436	1.00	10.00	A_13
ATOM	1108	C	PHE	121	79.505	37.756	12.283	1.00	17.14	A_13
ATOM	1109	O	PHE	121	79.642	37.104	11.256	1.00	13.04	A_13
ATOM	1110	N	GLY	122	79.738	37.245	13.490	1.00	16.60	A_13

49/75

ATOM	1112	CA	GLY	122	80.202	35.872	13.627	1.00	19.45	A_13
ATOM	1113	C	GLY	122	79.162	34.982	12.966	1.00	18.55	A_13
ATOM	1114	O	GLY	122	79.500	33.988	12.306	1.00	10.03	A_13
ATOM	1115	N	HIS	123	77.892	35.361	13.140	1.00	18.22	A_13
ATOM	1117	CA	HIS	123	76.753	34.665	12.525	1.00	16.31	A_13
ATOM	1118	CB	HIS	123	75.424	35.224	13.031	1.00	11.35	A_13
ATOM	1119	CG	HIS	123	75.049	34.768	14.403	1.00	10.33	A_13
ATOM	1120	CD2	HIS	123	74.552	35.454	15.457	1.00	16.64	A_13
ATOM	1121	ND1	HIS	123	75.097	33.450	14.782	1.00	18.04	A_13
ATOM	1123	CE1	HIS	123	74.638	33.332	16.017	1.00	16.66	A_13
ATOM	1124	NE2	HIS	123	74.301	34.533	16.450	1.00	25.32	A_13
ATOM	1125	C	HIS	123	76.771	34.853	10.997	1.00	13.66	A_13
ATOM	1126	O	HIS	123	76.565	33.901	10.246	1.00	10.82	A_13
ATOM	1127	N	SER	124	77.006	36.082	10.539	1.00	13.57	A_13
ATOM	1129	CA	SER	124	77.030	36.368	9.099	1.00	12.03	A_13
ATOM	1130	CB	SER	124	77.311	37.863	8.832	1.00	10.35	A_13
ATOM	1131	OG	SER	124	76.399	38.706	9.510	1.00	14.26	A_13
ATOM	1133	C	SER	124	78.117	35.548	8.422	1.00	21.45	A_13
ATOM	1134	O	SER	124	78.079	35.333	7.210	1.00	10.00	A_13
ATOM	1135	N	LEU	125	79.091	35.108	9.216	1.00	10.00	A_13
ATOM	1137	CA	LEU	125	80.222	34.340	8.707	1.00	19.28	A_13
ATOM	1138	CB	LEU	125	81.521	34.754	9.422	1.00	22.39	A_13
ATOM	1139	CG	LEU	125	81.849	36.258	9.340	1.00	10.00	A_13
ATOM	1140	CD1	LEU	125	83.063	36.622	10.190	1.00	10.00	A_13
ATOM	1141	CD2	LEU	125	82.029	36.651	7.873	1.00	10.00	A_13
ATOM	1142	C	LEU	125	79.986	32.851	8.843	1.00	10.00	A_13
ATOM	1143	O	LEU	125	80.759	32.056	8.329	1.00	23.27	A_13
ATOM	1144	N	GLY	126	78.932	32.477	9.563	1.00	22.87	A_13
ATOM	1146	CA	GLY	126	78.604	31.070	9.720	1.00	17.27	A_13
ATOM	1147	C	GLY	126	78.781	30.464	11.094	1.00	11.71	A_13
ATOM	1148	O	GLY	126	78.784	29.244	11.236	1.00	24.16	A_13
ATOM	1149	N	LEU	127	78.972	31.297	12.105	1.00	18.95	A_13
ATOM	1151	CA	LEU	127	79.152	30.790	13.457	1.00	22.84	A_13
ATOM	1152	CB	LEU	127	80.113	31.693	14.252	1.00	11.92	A_13
ATOM	1153	CG	LEU	127	81.244	30.969	14.983	1.00	18.83	A_13
ATOM	1154	CD1	LEU	127	82.096	30.197	13.979	1.00	16.63	A_13
ATOM	1155	CD2	LEU	127	82.104	31.970	15.760	1.00	22.15	A_13
ATOM	1156	C	LEU	127	77.802	30.699	14.163	1.00	21.02	A_13
ATOM	1157	O	LEU	127	76.996	31.629	14.098	1.00	14.68	A_13
ATOM	1158	N	ASP	128	77.563	29.572	14.828	1.00	18.87	A_13
ATOM	1160	CA	ASP	128	76.336	29.345	15.571	1.00	16.46	A_13
ATOM	1161	CB	ASP	128	75.996	27.855	15.540	1.00	17.60	A_13
ATOM	1162	CG	ASP	128	74.577	27.552	15.996	1.00	23.55	A_13
ATOM	1163	OD1	ASP	128	73.796	28.488	16.258	1.00	10.00	A_13
ATOM	1164	OD2	ASP	128	74.236	26.355	16.087	1.00	32.36	A_13
ATOM	1165	C	ASP	128	76.634	29.803	16.995	1.00	10.00	A_13
ATOM	1166	O	ASP	128	77.650	30.420	17.244	1.00	29.54	A_13
ATOM	1167	N	HIS	129	75.714	29.565	17.912	1.00	10.00	A_13
ATOM	1169	CA	HIS	129	75.910	29.955	19.289	1.00	10.00	A_13
ATOM	1170	CB	HIS	129	74.582	30.033	20.029	1.00	21.30	A_13
ATOM	1171	CG	HIS	129	73.798	31.282	19.761	1.00	24.16	A_13
ATOM	1172	CD2	HIS	129	74.180	32.585	19.725	1.00	10.00	A_13
ATOM	1173	ND1	HIS	129	72.460	31.263	19.476	1.00	21.70	A_13
ATOM	1175	CE1	HIS	129	72.031	32.501	19.271	1.00	10.27	A_13
ATOM	1176	NE2	HIS	129	73.057	33.319	19.407	1.00	14.37	A_13
ATOM	1177	C	HIS	129	76.780	28.947	19.992	1.00	30.04	A_13
ATOM	1178	O	HIS	129	76.624	27.730	19.822	1.00	22.13	A_13
ATOM	1179	N	SER	130	77.628	29.468	20.860	1.00	18.60	A_13
ATOM	1181	CA	SER	130	78.534	28.662	21.636	1.00	10.79	A_13
ATOM	1182	CB	SER	130	79.849	29.435	21.816	1.00	21.31	A_13
ATOM	1183	OG	SER	130	80.782	28.731	22.616	1.00	16.34	A_13
ATOM	1185	C	SER	130	77.898	28.368	22.987	1.00	31.13	A_13
ATOM	1186	O	SER	130	76.962	29.060	23.440	1.00	15.87	A_13
ATOM	1187	N	LYS	131	78.402	27.319	23.619	1.00	13.13	A_13
ATOM	1189	CA	LYS	131	77.924	26.925	24.928	1.00	13.21	A_13
ATOM	1190	CB	LYS	131	77.656	25.414	24.990	1.00	18.85	A_13
ATOM	1191	CG	LYS	131	78.689	24.541	24.303	1.00	32.55	A_13
ATOM	1192	CD	LYS	131	78.547	24.601	22.790	1.00	41.54	A_13
ATOM	1193	CE	LYS	131	79.909	24.672	22.117	1.00	19.64	A_13
ATOM	1194	NZ	LYS	131	80.747	25.799	22.617	1.00	13.47	A_13
ATOM	1198	C	LYS	131	78.922	27.379	25.982	1.00	10.00	A_13
ATOM	1199	O	LYS	131	78.666	27.260	27.185	1.00	13.35	A_13
ATOM	1200	N	ASP	132	80.025	27.968	25.519	1.00	13.47	A_13
ATOM	1202	CA	ASP	132	81.097	28.487	26.375	1.00	10.04	A_13
ATOM	1203	CB	ASP	132	82.376	28.617	25.522	1.00	18.14	A_13
ATOM	1204	CG	ASP	132	83.649	28.821	26.345	1.00	16.54	A_13
ATOM	1205	OD1	ASP	132	84.645	28.132	26.028	1.00	36.08	A_13

50/75

ATOM	1206	OD2	ASP	132	83.685	29.660	27.276	1.00	15.60	A_13
ATOM	1207	C	ASP	132	80.603	29.875	26.836	1.00	18.74	A_13
ATOM	1208	O	ASP	132	80.559	30.816	26.038	1.00	14.61	A_13
ATOM	1209	N	PRO	133	80.305	30.039	28.142	1.00	15.61	A_13
ATOM	1210	CD	PRO	133	80.617	29.127	29.251	1.00	21.19	A_13
ATOM	1211	CA	PRO	133	79.818	31.320	28.662	1.00	10.00	A_13
ATOM	1212	CB	PRO	133	79.542	31.007	30.135	1.00	10.00	A_13
ATOM	1213	CG	PRO	133	80.633	30.063	30.450	1.00	30.94	A_13
ATOM	1214	C	PRO	133	80.834	32.444	28.511	1.00	22.87	A_13
ATOM	1215	O	PRO	133	80.526	33.574	28.742	1.00	21.65	A_13
ATOM	1216	N	GLY	134	82.070	32.115	28.174	1.00	20.95	A_13
ATOM	1218	CA	GLY	134	83.055	33.167	28.028	1.00	15.22	A_13
ATOM	1219	C	GLY	134	83.182	33.578	26.581	1.00	34.54	A_13
ATOM	1220	O	GLY	134	83.962	34.488	26.252	1.00	18.06	A_13
ATOM	1221	N	ALA	135	82.490	32.846	25.706	1.00	21.09	A_13
ATOM	1223	CA	ALA	135	82.547	33.110	24.263	1.00	27.50	A_13
ATOM	1224	CB	ALA	135	82.131	31.858	23.453	1.00	10.00	A_13
ATOM	1225	C	ALA	135	81.722	34.308	23.814	1.00	21.74	A_13
ATOM	1226	O	ALA	135	80.641	34.556	24.328	1.00	13.84	A_13
ATOM	1227	N	LEU	136	82.220	34.990	22.787	1.00	19.10	A_13
ATOM	1229	CA	LEU	136	81.540	36.140	22.203	1.00	21.65	A_13
ATOM	1230	CB	LEU	136	82.448	36.803	21.161	1.00	10.00	A_13
ATOM	1231	CG	LEU	136	81.964	37.898	20.201	1.00	17.22	A_13
ATOM	1232	CD1	LEU	136	81.250	37.296	19.024	1.00	24.18	A_13
ATOM	1233	CD2	LEU	136	81.113	38.896	20.905	1.00	10.00	A_13
ATOM	1234	C	LEU	136	80.250	35.632	21.558	1.00	19.32	A_13
ATOM	1235	O	LEU	136	79.266	36.359	21.458	1.00	26.20	A_13
ATOM	1236	N	MET	137	80.297	34.409	21.029	1.00	10.00	A_13
ATOM	1238	CA	MET	137	79.123	33.791	20.423	1.00	10.02	A_13
ATOM	1239	CB	MET	137	79.507	32.691	19.428	1.00	15.14	A_13
ATOM	1240	CG	MET	137	80.181	33.223	18.169	1.00	16.42	A_13
ATOM	1241	SD	MET	137	79.366	34.665	17.397	1.00	10.65	A_13
ATOM	1242	CE	MET	137	77.848	34.005	16.975	1.00	10.87	A_13
ATOM	1243	C	MET	137	78.122	33.256	21.447	1.00	12.70	A_13
ATOM	1244	O	MET	137	77.187	32.539	21.087	1.00	10.00	A_13
ATOM	1245	N	PHE	138	78.295	33.627	22.713	1.00	18.70	A_13
ATOM	1247	CA	PHE	138	77.370	33.196	23.759	1.00	24.08	A_13
ATOM	1248	CB	PHE	138	77.954	33.448	25.159	1.00	24.15	A_13
ATOM	1249	CG	PHE	138	77.306	32.617	26.240	1.00	29.38	A_13
ATOM	1250	CD1	PHE	138	76.694	33.222	27.336	1.00	27.07	A_13
ATOM	1251	CD2	PHE	138	77.253	31.226	26.123	1.00	21.37	A_13
ATOM	1252	CE1	PHE	138	76.033	32.455	28.289	1.00	30.35	A_13
ATOM	1253	CE2	PHE	138	76.599	30.458	27.065	1.00	19.58	A_13
ATOM	1254	CZ	PHE	138	75.986	31.070	28.154	1.00	17.69	A_13
ATOM	1255	C	PHE	138	76.074	33.992	23.513	1.00	14.20	A_13
ATOM	1256	O	PHE	138	76.115	35.105	23.014	1.00	10.27	A_13
ATOM	1257	N	PRO	139	74.899	33.366	23.730	1.00	13.04	A_13
ATOM	1258	CD	PRO	139	74.664	31.975	24.131	1.00	11.17	A_13
ATOM	1259	CA	PRO	139	73.619	34.043	23.504	1.00	18.27	A_13
ATOM	1260	CB	PRO	139	72.625	32.875	23.384	1.00	14.33	A_13
ATOM	1261	CG	PRO	139	73.474	31.634	23.305	1.00	24.22	A_13
ATOM	1262	C	PRO	139	73.162	35.018	24.584	1.00	16.51	A_13
ATOM	1263	O	PRO	139	72.023	35.467	24.535	1.00	24.45	A_13
ATOM	1264	N	ILE	140	74.034	35.375	25.524	1.00	23.16	A_13
ATOM	1266	CA	ILE	140	73.652	36.290	26.604	1.00	25.00	A_13
ATOM	1267	CB	ILE	140	73.688	35.559	27.966	1.00	12.10	A_13
ATOM	1268	CG2	ILE	140	73.336	36.519	29.085	1.00	12.62	A_13
ATOM	1269	CG1	ILE	140	72.738	34.341	27.904	1.00	22.67	A_13
ATOM	1270	CD1	ILE	140	72.827	33.353	29.073	1.00	27.73	A_13
ATOM	1271	C	ILE	140	74.584	37.489	26.621	1.00	30.64	A_13
ATOM	1272	O	ILE	140	75.778	37.317	26.682	1.00	23.16	A_13
ATOM	1273	N	TYR	141	74.033	38.694	26.532	1.00	21.05	A_13
ATOM	1275	CA	TYR	141	74.851	39.901	26.528	1.00	20.10	A_13
ATOM	1276	CB	TYR	141	74.017	41.122	26.129	1.00	17.66	A_13
ATOM	1277	CG	TYR	141	74.784	42.433	26.103	1.00	22.24	A_13
ATOM	1278	CD1	TYR	141	74.711	43.318	27.171	1.00	18.07	A_13
ATOM	1279	CE1	TYR	141	75.386	44.527	27.144	1.00	19.84	A_13
ATOM	1280	CD2	TYR	141	75.563	42.798	24.999	1.00	18.08	A_13
ATOM	1281	CE2	TYR	141	76.244	44.008	24.961	1.00	10.00	A_13
ATOM	1282	CZ	TYR	141	76.149	44.867	26.038	1.00	25.17	A_13
ATOM	1283	OH	TYR	141	76.814	46.070	26.043	1.00	30.78	A_13
ATOM	1285	C	TYR	141	75.533	40.169	27.852	1.00	19.61	A_13
ATOM	1286	O	TYR	141	74.910	40.146	28.913	1.00	16.08	A_13
ATOM	1287	N	THR	142	76.817	40.476	27.772	1.00	26.26	A_13
ATOM	1289	CA	THR	142	77.612	40.788	28.944	1.00	24.52	A_13
ATOM	1290	CB	THR	142	78.498	39.568	29.362	1.00	10.00	A_13
ATOM	1291	OG1	THR	142	77.664	38.587	29.981	1.00	37.30	A_13

51/75

ATOM	1293	CG2	THR	142	79.543	39.961	30.390	1.00	14.88	A_13
ATOM	1294	C	THR	142	78.467	41.976	28.580	1.00	25.46	A_13
ATOM	1295	O	THR	142	78.980	42.058	27.464	1.00	10.00	A_13
ATOM	1296	N	TYR	143	78.575	42.947	29.476	1.00	20.23	A_13
ATOM	1298	CA	TYR	143	79.412	44.079	29.133	1.00	32.69	A_13
ATOM	1299	CB	TYR	143	79.024	45.363	29.854	1.00	35.01	A_13
ATOM	1300	CG	TYR	143	79.834	46.531	29.347	1.00	16.01	A_13
ATOM	1301	CD1	TYR	143	79.776	46.910	27.998	1.00	12.56	A_13
ATOM	1302	CE1	TYR	143	80.554	47.961	27.510	1.00	19.23	A_13
ATOM	1303	CD2	TYR	143	80.690	47.230	30.196	1.00	19.43	A_13
ATOM	1304	CE2	TYR	143	81.478	48.287	29.719	1.00	15.52	A_13
ATOM	1305	CZ	TYR	143	81.403	48.643	28.376	1.00	12.56	A_13
ATOM	1306	OH	TYR	143	82.193	49.654	27.892	1.00	18.85	A_13
ATOM	1308	C	TYR	143	80.871	43.754	29.382	1.00	25.10	A_13
ATOM	1309	O	TYR	143	81.373	43.846	30.503	1.00	28.90	A_13
ATOM	1310	N	THR	144	81.539	43.375	28.303	1.00	35.25	A_13
ATOM	1312	CA	THR	144	82.946	43.029	28.336	1.00	38.86	A_13
ATOM	1313	CB	THR	144	83.158	41.568	27.873	1.00	23.22	A_13
ATOM	1314	OG1	THR	144	82.129	41.219	26.934	1.00	35.22	A_13
ATOM	1316	CG2	THR	144	83.105	40.616	29.082	1.00	17.53	A_13
ATOM	1317	C	THR	144	83.720	44.017	27.488	1.00	21.63	A_13
ATOM	1318	O	THR	144	84.434	43.651	26.556	1.00	37.44	A_13
ATOM	1319	N	GLY	145	83.504	45.288	27.798	1.00	14.47	A_13
ATOM	1321	CA	GLY	145	84.200	46.375	27.131	1.00	24.39	A_13
ATOM	1322	C	GLY	145	84.119	46.536	25.628	1.00	41.65	A_13
ATOM	1323	O	GLY	145	84.053	45.565	24.877	1.00	42.39	A_13
ATOM	1324	N	LYS	146	84.122	47.792	25.195	1.00	33.04	A_13
ATOM	1326	CA	LYS	146	84.059	48.103	23.778	1.00	29.29	A_13
ATOM	1327	CB	LYS	146	83.260	49.392	23.539	1.00	26.47	A_13
ATOM	1328	CG	LYS	146	83.087	49.721	22.059	1.00	33.24	A_13
ATOM	1329	CD	LYS	146	82.812	51.194	21.833	1.00	13.70	A_13
ATOM	1330	CE	LYS	146	82.620	51.497	20.343	1.00	18.35	A_13
ATOM	1331	NZ	LYS	146	83.766	51.122	19.477	1.00	30.66	A_13
ATOM	1335	C	LYS	146	85.491	48.297	23.308	1.00	41.61	A_13
ATOM	1336	O	LYS	146	86.028	49.412	23.382	1.00	46.44	A_13
ATOM	1337	N	SER	147	86.130	47.206	22.898	1.00	34.67	A_13
ATOM	1339	CA	SER	147	87.509	47.258	22.416	1.00	30.76	A_13
ATOM	1340	CB	SER	147	87.624	48.258	21.249	1.00	24.56	A_13
ATOM	1341	OG	SER	147	86.638	48.002	20.257	1.00	31.81	A_13
ATOM	1343	C	SER	147	88.464	47.626	23.567	1.00	33.60	A_13
ATOM	1344	O	SER	147	88.789	48.806	23.789	1.00	39.96	A_13
ATOM	1345	N	HIS	148	88.862	46.611	24.331	1.00	36.71	A_13
ATOM	1347	CA	HIS	148	89.778	46.769	25.467	1.00	34.40	A_13
ATOM	1348	CB	HIS	148	89.307	47.862	26.438	1.00	26.40	A_13
ATOM	1349	CG	HIS	148	90.251	49.022	26.537	1.00	39.11	A_13
ATOM	1350	CD2	HIS	148	90.929	49.542	27.588	1.00	30.52	A_13
ATOM	1351	ND1	HIS	148	90.635	49.767	25.437	1.00	37.71	A_13
ATOM	1353	CE1	HIS	148	91.511	50.681	25.807	1.00	29.04	A_13
ATOM	1354	NE2	HIS	148	91.707	50.567	27.110	1.00	29.03	A_13
ATOM	1356	C	HIS	148	89.949	45.436	26.190	1.00	39.41	A_13
ATOM	1357	O	HIS	148	90.134	45.373	27.411	1.00	35.01	A_13
ATOM	1358	N	PHE	149	89.840	44.386	25.383	1.00	25.35	A_13
ATOM	1360	CA	PHE	149	89.996	42.966	25.721	1.00	30.54	A_13
ATOM	1361	CB	PHE	149	88.788	42.423	26.495	1.00	33.34	A_13
ATOM	1362	CG	PHE	149	88.951	42.440	27.996	1.00	31.37	A_13
ATOM	1363	CD1	PHE	149	89.387	41.302	28.673	1.00	30.46	A_13
ATOM	1364	CD2	PHE	149	88.624	43.575	28.740	1.00	40.67	A_13
ATOM	1365	CE1	PHE	149	89.492	41.293	30.075	1.00	18.92	A_13
ATOM	1366	CE2	PHE	149	88.728	43.574	30.136	1.00	23.23	A_13
ATOM	1367	CZ	PHE	149	89.161	42.430	30.803	1.00	17.03	A_13
ATOM	1368	C	PHE	149	90.026	42.366	24.295	1.00	41.76	A_13
ATOM	1369	O	PHE	149	89.967	43.119	23.307	1.00	40.43	A_13
ATOM	1370	N	MET	150	90.132	41.050	24.142	1.00	31.30	A_13
ATOM	1372	CA	MET	150	90.152	40.531	22.779	1.00	20.65	A_13
ATOM	1373	CB	MET	150	91.588	40.195	22.352	1.00	28.29	A_13
ATOM	1374	CG	MET	150	92.494	41.436	22.188	1.00	34.71	A_13
ATOM	1375	SD	MET	150	91.750	42.780	21.185	1.00	67.91	A_13
ATOM	1376	CE	MET	150	92.512	42.498	19.518	1.00	22.43	A_13
ATOM	1377	C	MET	150	89.201	39.370	22.497	1.00	21.51	A_13
ATOM	1378	O	MET	150	88.498	38.901	23.391	1.00	25.37	A_13
ATOM	1379	N	LEU	151	89.159	38.938	21.240	1.00	13.78	A_13
ATOM	1381	CA	LEU	151	88.313	37.825	20.834	1.00	14.73	A_13
ATOM	1382	CB	LEU	151	88.435	37.589	19.321	1.00	15.49	A_13
ATOM	1383	CG	LEU	151	87.535	36.511	18.691	1.00	27.05	A_13
ATOM	1384	CD1	LEU	151	86.070	36.915	18.847	1.00	10.98	A_13
ATOM	1385	CD2	LEU	151	87.879	36.310	17.208	1.00	15.73	A_13
ATOM	1386	C	LEU	151	88.732	36.563	21.600	1.00	25.01	A_13

52/75

ATOM	1337	O	LEU	151	89.912	36.178	21.589	1.00	17.37	A_13
ATOM	1388	N	PRO	152	87.777	35.927	22.306	1.00	10.37	A_13
ATOM	1389	CD	PRO	152	86.425	36.450	22.575	1.00	15.35	A_13
ATOM	1390	CA	PRO	152	88.030	34.712	23.087	1.00	11.49	A_13
ATOM	1391	CB	PRO	152	86.658	34.412	23.702	1.00	15.98	A_13
ATOM	1392	CG	PRO	152	86.083	35.789	23.898	1.00	27.60	A_13
ATOM	1393	C	PRO	152	88.533	33.553	22.230	1.00	18.06	A_13
ATOM	1394	O	PRO	152	88.160	33.430	21.063	1.00	16.21	A_13
ATOM	1395	N	ASP	153	89.350	32.696	22.836	1.00	15.86	A_13
ATOM	1397	CA	ASP	153	89.933	31.526	22.185	1.00	20.25	A_13
ATOM	1398	CB	ASP	153	90.632	30.630	23.227	1.00	18.17	A_13
ATOM	1399	CG	ASP	153	91.843	31.301	23.908	1.00	24.01	A_13
ATOM	1400	OD1	ASP	153	92.517	32.159	23.284	1.00	14.96	A_13
ATOM	1401	OD2	ASP	153	92.131	30.937	25.077	1.00	20.20	A_13
ATOM	1402	C	ASP	153	88.887	30.678	21.452	1.00	24.64	A_13
ATOM	1403	O	ASP	153	89.113	30.221	20.330	1.00	13.51	A_13
ATOM	1404	N	ASP	154	87.757	30.453	22.114	1.00	24.11	A_13
ATOM	1406	CA	ASP	154	86.664	29.657	21.577	1.00	19.19	A_13
ATOM	1407	CB	ASP	154	85.527	29.632	22.587	1.00	18.27	A_13
ATOM	1408	CG	ASP	154	84.406	28.751	22.161	1.00	24.26	A_13
ATOM	1409	OD1	ASP	154	83.314	29.291	21.950	1.00	20.97	A_13
ATOM	1410	OD2	ASP	154	84.609	27.530	22.031	1.00	20.32	A_13
ATOM	1411	C	ASP	154	86.162	30.170	20.229	1.00	18.99	A_13
ATOM	1412	O	ASP	154	86.043	29.408	19.277	1.00	22.56	A_13
ATOM	1413	N	ASP	155	85.873	31.465	20.158	1.00	16.11	A_13
ATOM	1415	CA	ASP	155	85.407	32.078	18.917	1.00	25.30	A_13
ATOM	1416	CB	ASP	155	85.011	33.527	19.158	1.00	13.32	A_13
ATOM	1417	CG	ASP	155	83.975	33.655	20.249	1.00	11.19	A_13
ATOM	1418	OD1	ASP	155	84.347	34.136	21.332	1.00	12.26	A_13
ATOM	1419	OD2	ASP	155	82.810	33.255	20.029	1.00	10.00	A_13
ATOM	1420	C	ASP	155	86.461	31.992	17.828	1.00	13.98	A_13
ATOM	1421	O	ASP	155	86.141	31.656	16.687	1.00	14.08	A_13
ATOM	1422	N	VAL	156	87.713	32.310	18.160	1.00	16.49	A_13
ATOM	1424	CA	VAL	156	88.771	32.201	17.159	1.00	27.34	A_13
ATOM	1425	CB	VAL	156	90.145	32.826	17.625	1.00	23.59	A_13
ATOM	1426	CG1	VAL	156	90.327	32.750	19.119	1.00	13.94	A_13
ATOM	1427	CG2	VAL	156	91.312	32.153	16.919	1.00	21.70	A_13
ATOM	1428	C	VAL	156	88.874	30.738	16.657	1.00	16.95	A_13
ATOM	1429	O	VAL	156	88.946	30.506	15.448	1.00	13.79	A_13
ATOM	1430	N	GLN	157	88.762	29.763	17.561	1.00	19.45	A_13
ATOM	1432	CA	GLN	157	88.796	28.352	17.154	1.00	30.53	A_13
ATOM	1433	CB	GLN	157	88.579	27.422	18.353	1.00	23.08	A_13
ATOM	1434	CG	GLN	157	89.633	27.521	19.452	1.00	24.83	A_13
ATOM	1435	CD	GLN	157	90.950	26.872	19.089	1.00	20.26	A_13
ATOM	1436	OE1	GLN	157	91.743	27.422	18.316	1.00	25.80	A_13
ATOM	1437	NE2	GLN	157	91.204	25.702	19.673	1.00	38.67	A_13
ATOM	1440	C	GLN	157	87.667	28.136	16.148	1.00	14.16	A_13
ATOM	1441	O	GLN	157	87.869	27.541	15.096	1.00	14.11	A_13
ATOM	1442	N	GLY	158	86.505	28.709	16.437	1.00	19.16	A_13
ATOM	1444	CA	GLY	158	85.361	28.584	15.551	1.00	12.79	A_13
ATOM	1445	C	GLY	158	85.510	29.144	14.143	1.00	24.46	A_13
ATOM	1446	O	GLY	158	85.181	28.449	13.177	1.00	18.77	A_13
ATOM	1447	N	ILE	159	85.936	30.403	13.989	1.00	22.41	A_13
ATOM	1449	CA	ILE	159	86.091	30.946	12.628	1.00	31.18	A_13
ATOM	1450	CB	ILE	159	86.300	32.508	12.532	1.00	23.53	A_13
ATOM	1451	CG2	ILE	159	84.991	33.203	12.177	1.00	17.28	A_13
ATOM	1452	CG1	ILE	159	87.022	33.063	13.758	1.00	15.28	A_13
ATOM	1453	CD1	ILE	159	88.507	32.949	13.707	1.00	14.71	A_13
ATOM	1454	C	ILE	159	87.226	30.280	11.875	1.00	10.56	A_13
ATOM	1455	O	ILE	159	87.167	30.139	10.653	1.00	18.79	A_13
ATOM	1456	N	GLN	160	88.287	29.927	12.590	1.00	20.71	A_13
ATOM	1458	CA	GLN	160	89.411	29.294	11.943	1.00	10.00	A_13
ATOM	1459	CB	GLN	160	90.640	29.274	12.855	1.00	10.00	A_13
ATOM	1460	CG	GLN	160	91.114	30.690	13.182	1.00	13.93	A_13
ATOM	1461	CD	GLN	160	92.402	30.754	13.981	1.00	25.61	A_13
ATOM	1462	OE1	GLN	160	92.814	29.786	14.629	1.00	19.40	A_13
ATOM	1463	NE2	GLN	160	93.042	31.915	13.950	1.00	24.78	A_13
ATOM	1466	C	GLN	160	89.000	27.917	11.477	1.00	10.00	A_13
ATOM	1467	O	GLN	160	89.458	27.481	10.432	1.00	21.73	A_13
ATOM	1468	N	SER	161	88.068	27.268	12.186	1.00	10.00	A_13
ATOM	1470	CA	SER	161	87.610	25.946	11.760	1.00	11.63	A_13
ATOM	1471	CB	SER	161	86.688	25.292	12.800	1.00	18.40	A_13
ATOM	1472	OG	SER	161	85.365	25.795	12.759	1.00	15.44	A_13
ATOM	1474	C	SER	161	86.913	26.048	10.396	1.00	26.18	A_13
ATOM	1475	O	SER	161	86.839	25.065	9.654	1.00	13.96	A_13
ATOM	1476	N	LEU	162	86.428	27.247	10.070	1.00	19.36	A_13
ATOM	1478	CA	LEU	162	85.749	27.493	8.808	1.00	17.21	A_13

53/75

ATOM	1479	CB	LEU	162	84.584	28.477	9.007	1.00	14.37	A_13
ATOM	1480	CG	LEU	162	83.489	28.144	10.021	1.00	31.09	A_13
ATOM	1481	CD1	LEU	162	82.596	29.351	10.217	1.00	14.96	A_13
ATOM	1482	CD2	LEU	162	82.672	26.949	9.548	1.00	23.87	A_13
ATOM	1483	C	LEU	162	86.654	28.080	7.744	1.00	11.98	A_13
ATOM	1484	O	LEU	162	86.596	27.680	6.584	1.00	15.25	A_13
ATOM	1485	N	TYR	163	87.459	29.063	8.135	1.00	26.64	A_13
ATOM	1487	CA	TYR	163	88.320	29.796	7.204	1.00	18.28	A_13
ATOM	1488	CB	TYR	163	87.977	31.289	7.277	1.00	26.89	A_13
ATOM	1489	CG	TYR	163	86.519	31.600	7.039	1.00	18.80	A_13
ATOM	1490	CD1	TYR	163	86.027	31.744	5.749	1.00	10.00	A_13
ATOM	1491	CE1	TYR	163	84.680	31.936	5.515	1.00	12.83	A_13
ATOM	1492	CD2	TYR	163	85.622	31.672	8.099	1.00	16.58	A_13
ATOM	1493	CE2	TYR	163	84.266	31.867	7.873	1.00	12.32	A_13
ATOM	1494	CZ	TYR	163	83.807	31.991	6.576	1.00	11.77	A_13
ATOM	1495	OH	TYR	163	82.472	32.141	6.331	1.00	21.93	A_13
ATOM	1497	C	TYR	163	89.818	29.669	7.397	1.00	15.67	A_13
ATOM	1498	O	TYR	163	90.590	30.089	6.526	1.00	18.92	A_13
ATOM	1499	N	GLY	164	90.225	29.096	8.525	1.00	18.34	A_13
ATOM	1501	CA	GLY	164	91.636	28.966	8.826	1.00	10.61	A_13
ATOM	1502	C	GLY	164	92.149	30.215	9.525	1.00	15.63	A_13
ATOM	1503	O	GLY	164	91.334	31.139	9.775	1.00	21.42	A_13
ATOM	1504	OT	GLY	164	93.353	30.250	9.858	1.00	21.99	A_13
ATOM	3009	ZN	ZN	166	73.275	35.223	18.371	1.00	27.40	AION
ATOM	3010	ZN	ZN	167	65.511	41.122	10.564	1.00	27.86	AION
ATOM	3011	CA	CA	168	64.285	44.152	21.635	1.00	11.76	AION
ATOM	3012	CA	CA	165	73.319	39.377	1.854	1.00	40.73	AION
ATOM	3017	C5	WAY	169	67.400	35.999	20.267	1.00	38.86	A693
ATOM	3018	CF1	WAY	169	66.626	35.606	19.161	1.00	30.96	A693
ATOM	3019	CH	WAY	169	67.199	35.400	17.901	1.00	41.17	A693
ATOM	3020	C2	WAY	169	68.561	35.623	17.728	1.00	36.26	A693
ATOM	3021	C3	WAY	169	69.339	36.039	18.811	1.00	35.73	A693
ATOM	3022	C4	WAY	169	68.807	36.216	20.078	1.00	33.71	A693
ATOM	3023	N20	WAY	169	69.699	36.617	21.141	1.00	33.16	A693
ATOM	3024	CD	WAY	169	70.137	35.640	22.189	1.00	29.78	A693
ATOM	3025	C23	WAY	169	68.986	34.739	22.685	1.00	25.69	A693
ATOM	3026	C28	WAY	169	68.187	35.088	23.798	1.00	31.72	A693
ATOM	3027	C27	WAY	169	67.141	34.238	24.205	1.00	33.61	A693
ATOM	3028	CM	WAY	169	66.921	33.061	23.490	1.00	32.16	A693
ATOM	3029	N25	WAY	169	67.703	32.748	22.426	1.00	42.39	A693
ATOM	3030	C24	WAY	169	68.709	33.546	22.016	1.00	27.88	A693
ATOM	3031	S21	WAY	169	69.757	38.213	21.577	1.00	24.43	A693
ATOM	3032	C16	WAY	169	71.513	38.570	21.438	1.00	29.69	A693
ATOM	3033	C21	WAY	169	72.032	39.163	20.269	1.00	19.32	A693
ATOM	3034	C20	WAY	169	73.400	39.453	20.169	1.00	11.82	A693
ATOM	3035	C19	WAY	169	74.267	39.156	21.241	1.00	19.50	A693
ATOM	3036	C18	WAY	169	73.748	38.564	22.402	1.00	11.88	A693
ATOM	3037	C17	WAY	169	72.382	38.272	22.507	1.00	26.57	A693
ATOM	3038	O33	WAY	169	75.623	39.445	21.141	1.00	16.99	A693
ATOM	3039	C36	WAY	169	76.504	39.509	22.271	1.00	12.69	A693
ATOM	3040	O15	WAY	169	69.030	39.032	20.657	1.00	13.98	A693
ATOM	3041	O14	WAY	169	69.419	38.338	22.942	1.00	22.94	A693
ATOM	3042	C7	WAY	169	70.780	36.256	18.621	1.00	30.48	A693
ATOM	3043	N9	WAY	169	71.192	36.946	17.553	1.00	10.00	A693
ATOM	3044	O10	WAY	169	72.581	37.127	17.426	1.00	38.25	A693
ATOM	3045	O8	WAY	169	71.614	35.847	19.414	1.00	39.46	A693
ATOM	3046	C29	WAY	169	66.584	36.175	21.566	1.00	46.13	A693
ATOM	1505	CB	THR	7	40.443	57.305	5.225	1.00	21.20	B_13
ATOM	1506	OG1	THR	7	39.149	56.999	5.762	1.00	25.31	B_13
ATOM	1508	CG2	THR	7	41.017	56.087	4.541	1.00	23.15	B_13
ATOM	1509	C	THR	7	40.920	59.113	6.901	1.00	32.45	B_13
ATOM	1510	O	THR	7	41.453	59.582	7.908	1.00	36.97	B_13
ATOM	1513	N	THR	7	41.386	56.786	7.488	1.00	34.12	B_13
ATOM	1515	CA	THR	7	41.371	57.761	6.365	1.00	26.16	B_13
ATOM	1516	N	LEU	8	39.907	59.694	6.265	1.00	23.60	B_13
ATOM	1518	CA	LEU	8	39.387	60.984	6.649	1.00	22.66	B_13
ATOM	1519	CB	LEU	8	38.113	60.848	7.503	1.00	21.78	B_13
ATOM	1520	CG	LEU	8	36.860	61.484	6.863	1.00	27.13	B_13
ATOM	1521	CD1	LEU	8	36.996	63.016	6.705	1.00	19.05	B_13
ATOM	1522	CD2	LEU	8	36.622	60.854	5.510	1.00	19.23	B_13
ATOM	1523	C	LEU	8	40.432	61.896	7.298	1.00	27.16	B_13
ATOM	1524	O	LEU	8	41.077	62.667	6.597	1.00	46.24	B_13
ATOM	1525	N	LYS	9	40.615	61.804	8.618	1.00	27.84	B_13
ATOM	1527	CA	LYS	9	41.572	62.674	9.306	1.00	15.20	B_13
ATOM	1528	CB	LYS	9	41.147	64.143	9.148	1.00	32.32	B_13
ATOM	1529	CG	LYS	9	39.663	64.342	8.853	1.00	29.47	B_13
ATOM	1530	CD	LYS	9	38.788	64.243	10.084	1.00	28.34	B_13

54/75

ATOM	1531	CE	LYS	9	38.830	65.556	10.842	1.00	18.48	B_13
ATOM	1532	NZ	LYS	9	38.732	66.725	9.888	1.00	33.19	B_13
ATOM	1536	C	LYS	9	41.809	62.384	10.780	1.00	20.69	B_13
ATOM	1537	O	LYS	9	41.268	61.428	11.334	1.00	25.62	B_13
ATOM	1538	N	TRP	10	42.654	63.208	11.390	1.00	12.09	B_13
ATOM	1540	CA	TRP	10	42.988	63.112	12.813	1.00	21.78	B_13
ATOM	1541	CB	TRP	10	44.403	63.660	13.048	1.00	23.03	B_13
ATOM	1542	CG	TRP	10	45.499	62.890	12.349	1.00	27.60	B_13
ATOM	1543	CD2	TRP	10	46.077	61.650	12.762	1.00	27.28	B_13
ATOM	1544	CE2	TRP	10	47.071	61.302	11.829	1.00	22.11	B_13
ATOM	1545	CE3	TRP	10	45.859	60.781	13.847	1.00	11.66	B_13
ATOM	1546	CD1	TRP	10	46.153	63.247	11.198	1.00	21.84	B_13
ATOM	1547	NE1	TRP	10	47.094	62.305	10.873	1.00	10.00	B_13
ATOM	1549	CZ2	TRP	10	47.847	60.143	11.929	1.00	25.24	B_13
ATOM	1550	CZ3	TRP	10	46.632	59.622	13.951	1.00	22.71	B_13
ATOM	1551	CH2	TRP	10	47.611	59.317	12.999	1.00	15.23	B_13
ATOM	1552	C	TRP	10	41.987	63.915	13.679	1.00	30.88	B_13
ATOM	1553	O	TRP	10	41.673	65.062	13.359	1.00	32.03	B_13
ATOM	1554	N	SER	11	41.495	63.316	14.765	1.00	35.64	B_13
ATOM	1556	CA	SER	11	40.548	63.981	15.665	1.00	30.37	B_13
ATOM	1557	CB	SER	11	39.498	62.995	16.176	1.00	31.03	B_13
ATOM	1558	OG	SER	11	38.485	62.815	15.202	1.00	41.11	B_13
ATOM	1560	C	SER	11	41.206	64.691	16.840	1.00	20.70	B_13
ATOM	1561	O	SER	11	40.558	65.002	17.838	1.00	36.52	B_13
ATOM	1562	N	LYS	12	42.504	64.910	16.731	1.00	23.56	B_13
ATOM	1564	CA	LYS	12	43.257	65.607	17.756	1.00	15.00	B_13
ATOM	1565	CB	LYS	12	43.991	64.631	18.688	1.00	18.58	B_13
ATOM	1566	CG	LYS	12	44.658	63.452	18.010	1.00	15.94	B_13
ATOM	1567	CD	LYS	12	45.456	62.589	19.007	1.00	23.03	B_13
ATOM	1568	CE	LYS	12	44.593	61.715	19.933	1.00	27.10	B_13
ATOM	1569	NZ	LYS	12	44.075	62.402	21.157	1.00	34.75	B_13
ATOM	1573	C	LYS	12	44.200	66.453	16.914	1.00	25.03	B_13
ATOM	1574	O	LYS	12	44.567	66.039	15.808	1.00	25.20	B_13
ATOM	1575	N	MET	13	44.536	67.647	17.401	1.00	18.44	B_13
ATOM	1577	CA	MET	13	45.377	68.582	16.663	1.00	24.63	B_13
ATOM	1578	CB	MET	13	44.864	70.015	16.880	1.00	13.15	B_13
ATOM	1579	CG	MET	13	43.421	70.253	16.419	1.00	21.56	B_13
ATOM	1580	SD	MET	13	43.167	70.131	14.616	1.00	31.39	B_13
ATOM	1581	CE	MET	13	41.433	69.678	14.474	1.00	24.70	B_13
ATOM	1582	C	MET	13	46.850	68.468	17.034	1.00	11.65	B_13
ATOM	1583	O	MET	13	47.728	68.815	16.247	1.00	14.33	B_13
ATOM	1584	N	ASN	14	47.103	67.985	18.242	1.00	16.99	B_13
ATOM	1586	CA	ASN	14	48.448	67.793	18.760	1.00	24.42	B_13
ATOM	1587	CB	ASN	14	48.437	68.006	20.268	1.00	17.84	B_13
ATOM	1588	CG	ASN	14	47.896	69.356	20.633	1.00	35.10	B_13
ATOM	1589	OD1	ASN	14	48.614	70.346	20.560	1.00	34.88	B_13
ATOM	1590	ND2	ASN	14	46.610	69.424	20.955	1.00	32.98	B_13
ATOM	1593	C	ASN	14	48.831	66.364	18.421	1.00	22.70	B_13
ATOM	1594	O	ASN	14	48.278	65.405	18.976	1.00	26.03	B_13
ATOM	1595	N	LEU	15	49.706	66.228	17.432	1.00	18.07	B_13
ATOM	1597	CA	LEU	15	50.144	64.912	16.969	1.00	29.36	B_13
ATOM	1598	CB	LEU	15	49.878	64.775	15.466	1.00	24.35	B_13
ATOM	1599	CG	LEU	15	48.380	64.762	15.162	1.00	19.51	B_13
ATOM	1600	CD1	LEU	15	48.079	65.469	13.852	1.00	27.59	B_13
ATOM	1601	CD2	LEU	15	47.902	63.326	15.163	1.00	19.66	B_13
ATOM	1602	C	LEU	15	51.613	64.704	17.257	1.00	28.48	B_13
ATOM	1603	O	LEU	15	52.341	65.657	17.552	1.00	22.28	B_13
ATOM	1604	N	THR	16	52.044	63.453	17.198	1.00	12.77	B_13
ATOM	1606	CA	THR	16	53.433	63.158	17.446	1.00	16.59	B_13
ATOM	1607	CB	THR	16	53.607	62.243	18.682	1.00	24.73	B_13
ATOM	1608	OG1	THR	16	52.912	61.005	18.481	1.00	12.79	B_13
ATOM	1610	CG2	THR	16	53.059	62.933	19.924	1.00	25.34	B_13
ATOM	1611	C	THR	16	54.038	62.515	16.214	1.00	21.94	B_13
ATOM	1612	O	THR	16	53.315	62.116	15.297	1.00	19.60	B_13
ATOM	1613	N	TYR	17	55.365	62.453	16.184	1.00	18.25	B_13
ATOM	1615	CA	TYR	17	56.092	61.810	15.097	1.00	19.54	B_13
ATOM	1616	CB	TYR	17	56.300	62.753	13.910	1.00	16.87	B_13
ATOM	1617	CG	TYR	17	57.277	63.892	14.116	1.00	27.90	B_13
ATOM	1618	CD1	TYR	17	56.839	65.135	14.587	1.00	13.93	B_13
ATOM	1619	CE1	TYR	17	57.700	66.221	14.652	1.00	17.08	B_13
ATOM	1620	CD2	TYR	17	58.613	63.764	13.723	1.00	14.99	B_13
ATOM	1621	CE2	TYR	17	59.479	64.841	13.777	1.00	25.98	B_13
ATOM	1622	CZ	TYR	17	59.017	66.075	14.242	1.00	33.12	B_13
ATOM	1623	OH	TYR	17	59.866	67.163	14.276	1.00	23.31	B_13
ATOM	1625	C	TYR	17	57.417	61.318	15.650	1.00	18.57	B_13
ATOM	1626	O	TYR	17	57.895	61.827	16.668	1.00	26.60	B_13
ATOM	1627	N	ARG	18	57.973	60.286	15.030	1.00	13.01	B_13

55/75

ATOM	1629	CA	ARG	18	59.245	59.750	15.492	1.00	18.74	B_13
ATOM	1630	CB	ARG	18	59.033	58.589	16.473	1.00	11.96	B_13
ATOM	1631	CG	ARG	18	60.320	57.911	16.970	1.00	15.06	B_13
ATOM	1632	CD	ARG	18	60.012	56.596	17.690	1.00	11.72	B_13
ATOM	1633	NE	ARG	18	61.165	55.689	17.752	1.00	10.00	B_13
ATOM	1635	CZ	ARG	18	61.134	54.428	18.181	1.00	24.87	B_13
ATOM	1636	NH1	ARG	18	60.004	53.882	18.614	1.00	13.34	B_13
ATOM	1639	NH2	ARG	18	62.247	53.703	18.169	1.00	20.03	B_13
ATOM	1642	C	ARG	18	60.076	59.309	14.307	1.00	13.14	B_13
ATOM	1643	O	ARG	18	59.598	58.588	13.434	1.00	14.10	B_13
ATOM	1644	N	ILE	19	61.304	59.813	14.252	1.00	15.55	B_13
ATOM	1646	CA	ILE	19	62.238	59.476	13.193	1.00	10.41	B_13
ATOM	1647	CB	ILE	19	63.307	60.603	13.054	1.00	17.20	B_13
ATOM	1648	CG2	ILE	19	64.273	60.307	11.903	1.00	16.57	B_13
ATOM	1649	CG1	ILE	19	62.613	61.952	12.836	1.00	15.47	B_13
ATOM	1650	CD1	ILE	19	63.543	63.110	12.783	1.00	14.99	B_13
ATOM	1651	C	ILE	19	62.870	58.166	13.673	1.00	10.00	B_13
ATOM	1652	O	ILE	19	63.829	58.179	14.434	1.00	10.00	B_13
ATOM	1653	N	VAL	20	62.289	57.037	13.276	1.00	17.84	B_13
ATOM	1655	CA	VAL	20	62.785	55.716	13.696	1.00	16.43	B_13
ATOM	1656	CB	VAL	20	61.911	54.570	13.138	1.00	13.17	B_13
ATOM	1657	CG1	VAL	20	62.519	53.208	13.493	1.00	10.00	B_13
ATOM	1658	CG2	VAL	20	60.521	54.673	13.698	1.00	10.00	B_13
ATOM	1659	C	VAL	20	64.268	55.449	13.387	1.00	16.02	B_13
ATOM	1660	O	VAL	20	65.001	54.909	14.218	1.00	21.07	B_13
ATOM	1661	N	ASN	21	64.698	55.762	12.177	1.00	10.00	B_13
ATOM	1663	CA	ASN	21	66.098	55.571	11.830	1.00	22.13	B_13
ATOM	1664	CB	ASN	21	66.392	54.128	11.386	1.00	19.75	B_13
ATOM	1665	CG	ASN	21	65.549	53.673	10.212	1.00	17.63	B_13
ATOM	1666	OD1	ASN	21	65.329	52.477	10.042	1.00	31.82	B_13
ATOM	1667	ND2	ASN	21	65.109	54.602	9.375	1.00	11.42	B_13
ATOM	1670	C	ASN	21	66.504	56.645	10.821	1.00	10.14	B_13
ATOM	1671	O	ASN	21	65.639	57.377	10.340	1.00	11.74	B_13
ATOM	1672	N	TYR	22	67.787	56.759	10.498	1.00	12.25	B_13
ATOM	1674	CA	TYR	22	68.233	57.829	9.602	1.00	12.46	B_13
ATOM	1675	CB	TYR	22	69.136	58.800	10.383	1.00	23.15	B_13
ATOM	1676	CG	TYR	22	68.461	59.584	11.492	1.00	21.95	B_13
ATOM	1677	CD1	TYR	22	68.221	60.945	11.348	1.00	22.29	B_13
ATOM	1678	CE1	TYR	22	67.625	61.678	12.347	1.00	10.00	B_13
ATOM	1679	CD2	TYR	22	68.077	58.974	12.687	1.00	13.42	B_13
ATOM	1680	CE2	TYR	22	67.471	59.710	13.693	1.00	14.69	B_13
ATOM	1681	CZ	TYR	22	67.254	61.064	13.505	1.00	12.89	B_13
ATOM	1682	OH	TYR	22	66.660	61.829	14.466	1.00	16.56	B_13
ATOM	1684	C	TYR	22	68.988	57.395	8.359	1.00	11.62	B_13
ATOM	1685	O	TYR	22	69.793	56.478	8.407	1.00	16.23	B_13
ATOM	1686	N	THR	23	68.792	58.111	7.261	1.00	10.39	B_13
ATOM	1688	CA	THR	23	69.503	57.800	6.024	1.00	20.36	B_13
ATOM	1689	CB	THR	23	68.909	58.582	4.829	1.00	16.21	B_13
ATOM	1690	OG1	THR	23	69.801	58.512	3.706	1.00	19.72	B_13
ATOM	1692	CG2	THR	23	68.663	60.039	5.206	1.00	16.62	B_13
ATOM	1693	C	THR	23	70.990	58.153	6.163	1.00	17.35	B_13
ATOM	1694	O	THR	23	71.377	58.958	7.024	1.00	13.88	B_13
ATOM	1695	N	PRO	24	71.852	57.503	5.364	1.00	15.86	B_13
ATOM	1696	CD	PRO	24	71.625	56.247	4.629	1.00	17.29	B_13
ATOM	1697	CA	PRO	24	73.287	57.796	5.436	1.00	15.96	B_13
ATOM	1698	CB	PRO	24	73.920	56.570	4.763	1.00	10.00	B_13
ATOM	1699	CG	PRO	24	72.891	55.504	4.905	1.00	15.15	B_13
ATOM	1700	C	PRO	24	73.635	59.069	4.668	1.00	27.08	B_13
ATOM	1701	O	PRO	24	74.698	59.656	4.869	1.00	19.47	B_13
ATOM	1702	N	ASP	25	72.728	59.489	3.794	1.00	16.99	B_13
ATOM	1704	CA	ASP	25	72.927	60.663	2.958	1.00	10.00	B_13
ATOM	1705	CB	ASP	25	71.792	60.758	1.953	1.00	11.53	B_13
ATOM	1706	CG	ASP	25	71.665	59.521	1.105	1.00	33.88	B_13
ATOM	1707	OD1	ASP	25	70.570	59.311	0.556	1.00	22.66	B_13
ATOM	1708	OD2	ASP	25	72.653	58.762	0.980	1.00	29.59	B_13
ATOM	1709	C	ASP	25	73.068	62.011	3.642	1.00	23.36	B_13
ATOM	1710	O	ASP	25	73.694	62.916	3.093	1.00	20.32	B_13
ATOM	1711	N	MET	26	72.480	62.158	4.826	1.00	18.44	B_13
ATOM	1713	CA	MET	26	72.510	63.432	5.537	1.00	13.83	B_13
ATOM	1714	CB	MET	26	71.154	64.151	5.368	1.00	10.00	B_13
ATOM	1715	CG	MET	26	70.782	64.491	3.913	1.00	28.32	B_13
ATOM	1716	SD	MET	26	69.016	64.786	3.599	1.00	12.18	B_13
ATOM	1717	CE	MET	26	68.395	63.255	3.887	1.00	37.25	B_13
ATOM	1718	C	MET	26	72.827	63.238	7.024	1.00	28.80	B_13
ATOM	1719	O	MET	26	72.839	62.107	7.533	1.00	20.90	B_13
ATOM	1720	N	THR	27	73.157	64.333	7.696	1.00	11.47	B_13
ATOM	1722	CA	THR	27	73.456	64.292	9.121	1.00	13.94	B_13

56/75

ATOM	1723	CB	THR	27	74.117	65.605	9.602	1.00	33.46	B_13
ATOM	1724	OG1	THR	27	73.209	66.702	9.415	1.00	10.00	B_13
ATOM	1726	CG2	THR	27	75.405	65.863	8.818	1.00	16.30	B_13
ATOM	1727	C	THR	27	72.135	64.113	9.861	1.00	10.67	B_13
ATOM	1728	O	THR	27	71.072	64.343	9.281	1.00	16.26	B_13
ATOM	1729	N	HIS	28	72.193	63.691	11.124	1.00	18.13	B_13
ATOM	1731	CA	HIS	28	70.986	63.514	11.915	1.00	10.00	B_13
ATOM	1732	CB	HIS	28	71.322	63.033	13.333	1.00	10.00	B_13
ATOM	1733	CG	HIS	28	71.793	61.608	13.401	1.00	22.65	B_13
ATOM	1734	CD2	HIS	28	72.893	61.003	12.889	1.00	22.73	B_13
ATOM	1735	ND1	HIS	28	71.103	60.627	14.080	1.00	19.90	B_13
ATOM	1737	CE1	HIS	28	71.755	59.481	13.985	1.00	16.52	B_13
ATOM	1738	NE2	HIS	28	72.843	59.681	13.268	1.00	20.38	B_13
ATOM	1740	C	HIS	28	70.281	64.870	11.957	1.00	29.38	B_13
ATOM	1741	O	HIS	28	69.074	64.941	11.742	1.00	17.20	B_13
ATOM	1742	N	SER	29	71.056	65.944	12.153	1.00	23.96	B_13
ATOM	1744	CA	SER	29	70.533	67.322	12.192	1.00	15.01	B_13
ATOM	1745	CB	SER	29	71.661	68.334	12.438	1.00	14.05	B_13
ATOM	1746	OG	SER	29	72.117	68.303	13.770	1.00	18.32	B_13
ATOM	1748	C	SER	29	69.808	67.729	10.909	1.00	10.95	B_13
ATOM	1749	O	SER	29	68.732	68.314	10.971	1.00	24.24	B_13
ATOM	1750	N	GLU	30	70.415	67.449	9.757	1.00	10.96	B_13
ATOM	1752	CA	GLU	30	69.820	67.786	8.470	1.00	10.00	B_13
ATOM	1753	CB	GLU	30	70.715	67.330	7.309	1.00	10.12	B_13
ATOM	1754	CG	GLU	30	71.967	68.143	7.042	1.00	22.31	B_13
ATOM	1755	CD	GLU	30	72.823	67.529	5.930	1.00	10.15	B_13
ATOM	1756	OE1	GLU	30	72.533	67.753	4.749	1.00	31.98	B_13
ATOM	1757	OE2	GLU	30	73.796	66.817	6.223	1.00	29.59	B_13
ATOM	1758	C	GLU	30	68.481	67.073	8.336	1.00	20.17	B_13
ATOM	1759	O	GLU	30	67.493	67.685	7.943	1.00	14.31	B_13
ATOM	1760	N	VAL	31	68.451	65.777	8.665	1.00	19.26	B_13
ATOM	1762	CA	VAL	31	67.228	64.989	8.536	1.00	14.22	B_13
ATOM	1763	CB	VAL	31	67.472	63.487	8.716	1.00	17.05	B_13
ATOM	1764	CG1	VAL	31	66.144	62.749	8.791	1.00	28.55	B_13
ATOM	1765	CG2	VAL	31	68.269	62.935	7.548	1.00	10.54	B_13
ATOM	1766	C	VAL	31	66.138	65.458	9.477	1.00	12.36	B_13
ATOM	1767	O	VAL	31	64.963	65.488	9.093	1.00	12.83	B_13
ATOM	1768	N	GLU	32	66.530	65.805	10.703	1.00	20.46	B_13
ATOM	1770	CA	GLU	32	65.596	66.306	11.710	1.00	16.04	B_13
ATOM	1771	CB	GLU	32	66.269	66.365	13.094	1.00	14.71	B_13
ATOM	1772	CG	GLU	32	66.512	64.985	13.741	1.00	23.30	B_13
ATOM	1773	CD	GLU	32	67.724	64.930	14.700	1.00	21.41	B_13
ATOM	1774	OE1	GLU	32	68.229	63.823	15.003	1.00	15.79	B_13
ATOM	1775	OE2	GLU	32	68.183	65.985	15.157	1.00	13.71	B_13
ATOM	1776	C	GLU	32	65.125	67.697	11.257	1.00	27.19	B_13
ATOM	1777	O	GLU	32	63.951	68.042	11.383	1.00	19.82	B_13
ATOM	1778	N	LYS	33	66.021	68.461	10.636	1.00	12.52	B_13
ATOM	1780	CA	LYS	33	65.663	69.786	10.171	1.00	13.00	B_13
ATOM	1781	CB	LYS	33	66.889	70.592	9.762	1.00	22.63	B_13
ATOM	1782	CG	LYS	33	66.581	72.054	9.560	1.00	18.24	B_13
ATOM	1783	CD	LYS	33	65.604	72.545	10.630	1.00	29.21	B_13
ATOM	1784	CE	LYS	33	66.185	72.429	12.048	1.00	41.79	B_13
ATOM	1785	NZ	LYS	33	65.181	71.939	13.054	1.00	20.17	B_13
ATOM	1789	C	LYS	33	64.698	69.686	9.023	1.00	10.62	B_13
ATOM	1790	O	LYS	33	63.734	70.437	8.971	1.00	22.94	B_13
ATOM	1791	N	ALA	34	64.915	68.707	8.150	1.00	10.00	B_13
ATOM	1793	CA	ALA	34	64.050	68.475	7.000	1.00	11.94	B_13
ATOM	1794	CB	ALA	34	64.611	67.374	6.100	1.00	10.00	B_13
ATOM	1795	C	ALA	34	62.640	68.115	7.423	1.00	10.00	B_13
ATOM	1796	O	ALA	34	61.675	68.650	6.878	1.00	15.32	B_13
ATOM	1797	N	PHE	35	62.510	67.208	8.387	1.00	21.32	B_13
ATOM	1799	CA	PHE	35	61.187	66.789	8.852	1.00	18.32	B_13
ATOM	1800	CB	PHE	35	61.267	65.451	9.614	1.00	25.48	B_13
ATOM	1801	CG	PHE	35	61.620	64.260	8.735	1.00	14.33	B_13
ATOM	1802	CD1	PHE	35	61.149	64.171	7.427	1.00	17.91	B_13
ATOM	1803	CD2	PHE	35	62.436	63.240	9.217	1.00	18.05	B_13
ATOM	1804	CE1	PHE	35	61.486	63.086	6.610	1.00	18.49	B_13
ATOM	1805	CE2	PHE	35	62.778	62.158	8.413	1.00	15.01	B_13
ATOM	1806	CZ	PHE	35	62.301	62.081	7.103	1.00	10.00	B_13
ATOM	1807	C	PHE	35	60.428	67.862	9.658	1.00	18.68	B_13
ATOM	1808	O	PHE	35	59.202	67.971	9.556	1.00	17.05	B_13
ATOM	1809	N	LYS	36	61.160	68.664	10.425	1.00	16.30	B_13
ATOM	1811	CA	LYS	36	60.579	69.749	11.229	1.00	19.34	B_13
ATOM	1812	CB	LYS	36	61.676	70.420	12.052	1.00	24.61	B_13
ATOM	1813	CG	LYS	36	61.200	71.293	13.191	1.00	18.38	B_13
ATOM	1814	CD	LYS	36	62.408	71.795	13.962	1.00	19.34	B_13
ATOM	1815	CE	LYS	36	62.067	72.267	15.356	1.00	21.80	B_13

57/75

ATOM	1816	NZ	LYS	36	63.299	72.615	16.118	1.00	27.76	B_13
ATOM	1820	C	LYS	36	59.924	70.770	10.301	1.00	10.19	B_13
ATOM	1821	O	LYS	36	58.788	71.183	10.528	1.00	14.95	B_13
ATOM	1822	N	LYS	37	60.630	71.134	9.233	1.00	15.89	B_13
ATOM	1824	CA	LYS	37	60.126	72.076	8.230	1.00	19.95	B_13
ATOM	1825	CB	LYS	37	61.202	72.386	7.189	1.00	10.00	B_13
ATOM	1826	CG	LYS	37	62.209	73.439	7.569	1.00	13.18	B_13
ATOM	1827	CD	LYS	37	62.869	73.966	6.311	1.00	28.86	B_13
ATOM	1828	CE	LYS	37	61.825	74.460	5.281	1.00	31.44	B_13
ATOM	1829	NZ	LYS	37	60.878	75.512	5.772	1.00	26.23	B_13
ATOM	1833	C	LYS	37	58.939	71.482	7.472	1.00	25.64	B_13
ATOM	1834	O	LYS	37	57.968	72.177	7.161	1.00	24.39	B_13
ATOM	1835	N	ALA	38	59.060	70.205	7.128	1.00	17.12	B_13
ATOM	1837	CA	ALA	38	58.031	69.493	6.381	1.00	16.06	B_13
ATOM	1838	CB	ALA	38	58.459	68.038	6.154	1.00	12.19	B_13
ATOM	1839	C	ALA	38	56.692	69.557	7.094	1.00	11.12	B_13
ATOM	1840	O	ALA	38	55.648	69.736	6.458	1.00	31.10	B_13
ATOM	1841	N	PHE	39	56.732	69.393	8.417	1.00	21.01	B_13
ATOM	1843	CA	PHE	39	55.540	69.446	9.257	1.00	10.85	B_13
ATOM	1844	CB	PHE	39	55.841	68.833	10.639	1.00	14.45	B_13
ATOM	1845	CG	PHE	39	55.851	67.325	10.659	1.00	21.88	B_13
ATOM	1846	CD1	PHE	39	57.016	66.625	10.954	1.00	16.88	B_13
ATOM	1847	CD2	PHE	39	54.675	66.599	10.442	1.00	22.14	B_13
ATOM	1848	CE1	PHE	39	57.010	65.223	11.037	1.00	17.95	B_13
ATOM	1849	CE2	PHE	39	54.655	65.190	10.522	1.00	17.22	B_13
ATOM	1850	CZ	PHE	39	55.823	64.503	10.823	1.00	13.51	B_13
ATOM	1851	C	PHE	39	55.044	70.898	9.426	1.00	19.98	B_13
ATOM	1852	O	PHE	39	53.839	71.160	9.393	1.00	14.30	B_13
ATOM	1853	N	LYS	40	55.981	71.826	9.611	1.00	20.03	B_13
ATOM	1855	CA	LYS	40	55.681	73.245	9.795	1.00	18.64	B_13
ATOM	1856	CB	LYS	40	56.989	74.011	10.020	1.00	19.28	B_13
ATOM	1857	CG	LYS	40	57.064	75.392	9.440	1.00	26.34	B_13
ATOM	1858	CD	LYS	40	58.288	76.093	9.974	1.00	18.46	B_13
ATOM	1859	CE	LYS	40	58.021	76.673	11.339	1.00	20.86	B_13
ATOM	1860	NZ	LYS	40	57.053	77.814	11.232	1.00	27.28	B_13
ATOM	1864	C	LYS	40	54.899	73.790	8.612	1.00	20.57	B_13
ATOM	1865	O	LYS	40	54.034	74.654	8.756	1.00	22.54	B_13
ATOM	1866	N	VAL	41	55.216	73.251	7.445	1.00	17.15	B_13
ATOM	1868	CA	VAL	41	54.565	73.576	6.184	1.00	19.19	B_13
ATOM	1869	CB	VAL	41	55.095	72.566	5.086	1.00	17.28	B_13
ATOM	1870	CG1	VAL	41	53.987	72.064	4.160	1.00	10.00	B_13
ATOM	1871	CG2	VAL	41	56.224	73.191	4.293	1.00	19.38	B_13
ATOM	1872	C	VAL	41	53.026	73.472	6.354	1.00	20.38	B_13
ATOM	1873	O	VAL	41	52.268	74.280	5.810	1.00	28.57	B_13
ATOM	1874	N	TRP	42	52.587	72.511	7.163	1.00	23.10	B_13
ATOM	1876	CA	TRP	42	51.166	72.265	7.403	1.00	19.29	B_13
ATOM	1877	CB	TRP	42	50.912	70.757	7.487	1.00	22.19	B_13
ATOM	1878	CG	TRP	42	51.437	70.007	6.313	1.00	19.32	B_13
ATOM	1879	CD2	TRP	42	50.836	69.909	5.015	1.00	31.02	B_13
ATOM	1880	CE2	TRP	42	51.659	69.067	4.238	1.00	22.49	B_13
ATOM	1881	CE3	TRP	42	49.677	70.448	4.434	1.00	15.54	B_13
ATOM	1882	CD1	TRP	42	52.571	69.251	6.269	1.00	14.04	B_13
ATOM	1883	NE1	TRP	42	52.710	68.681	5.027	1.00	13.55	B_13
ATOM	1885	CZ2	TRP	42	51.360	68.752	2.912	1.00	18.87	B_13
ATOM	1886	CZ3	TRP	42	49.383	70.132	3.116	1.00	13.33	B_13
ATOM	1887	CH2	TRP	42	50.219	69.294	2.370	1.00	20.30	B_13
ATOM	1888	C	TRP	42	50.617	72.926	8.660	1.00	24.68	B_13
ATOM	1889	O	TRP	42	49.455	73.339	8.688	1.00	20.93	B_13
ATOM	1890	N	SER	43	51.432	72.987	9.710	1.00	20.63	B_13
ATOM	1892	CA	SER	43	51.007	73.601	10.968	1.00	22.47	B_13
ATOM	1893	CB	SER	43	51.955	73.231	12.116	1.00	10.00	B_13
ATOM	1894	OG	SER	43	53.265	73.716	11.891	1.00	33.50	B_13
ATOM	1896	C	SER	43	50.913	75.122	10.829	1.00	14.99	B_13
ATOM	1897	O	SER	43	50.224	75.784	11.595	1.00	11.58	B_13
ATOM	1898	N	ASP	44	51.613	75.667	9.843	1.00	26.20	B_13
ATOM	1900	CA	ASP	44	51.595	77.100	9.617	1.00	22.11	B_13
ATOM	1901	CB	ASP	44	52.620	77.485	8.549	1.00	11.09	B_13
ATOM	1902	CG	ASP	44	54.000	77.751	9.125	1.00	18.45	B_13
ATOM	1903	OD1	ASP	44	54.903	78.114	8.347	1.00	17.67	B_13
ATOM	1904	OD2	ASP	44	54.195	77.602	10.345	1.00	21.44	B_13
ATOM	1905	C	ASP	44	50.216	77.575	9.190	1.00	32.83	B_13
ATOM	1906	O	ASP	44	49.795	78.677	9.549	1.00	34.78	B_13
ATOM	1907	N	VAL	45	49.508	76.735	8.439	1.00	31.40	B_13
ATOM	1909	CA	VAL	45	48.191	77.094	7.932	1.00	14.00	B_13
ATOM	1910	CB	VAL	45	48.121	76.872	6.401	1.00	15.73	B_13
ATOM	1911	CG1	VAL	45	49.123	77.755	5.707	1.00	19.37	B_13
ATOM	1912	CG2	VAL	45	48.407	75.409	6.055	1.00	10.00	B_13

58/75

ATOM	1913	C	VAL	45	47.054	76.333	8.575	1.00	18.43	B_13
ATOM	1914	O	VAL	45	45.954	76.304	8.026	1.00	26.09	B_13
ATOM	1915	N	THR	46	47.295	75.754	9.747	1.00	18.49	B_13
ATOM	1917	CA	THR	46	46.262	74.963	10.408	1.00	21.92	B_13
ATOM	1918	CB	THR	46	46.222	73.529	9.751	1.00	27.61	B_13
ATOM	1919	OG1	THR	46	44.876	73.047	9.661	1.00	28.78	B_13
ATOM	1921	CG2	THR	46	47.054	72.550	10.522	1.00	10.65	B_13
ATOM	1922	C	THR	46	46.505	74.931	11.932	1.00	18.41	B_13
ATOM	1923	O	THR	46	47.554	75.363	12.411	1.00	18.63	B_13
ATOM	1924	N	PRO	47	45.519	74.467	12.717	1.00	16.81	B_13
ATOM	1925	CD	PRO	47	44.113	74.209	12.348	1.00	32.80	B_13
ATOM	1926	CA	PRO	47	45.691	74.407	14.169	1.00	13.66	B_13
ATOM	1927	CB	PRO	47	44.256	74.489	14.675	1.00	30.52	B_13
ATOM	1928	CG	PRO	47	43.519	73.692	13.638	1.00	29.25	B_13
ATOM	1929	C	PRO	47	46.346	73.105	14.622	1.00	28.40	B_13
ATOM	1930	O	PRO	47	46.037	72.597	15.705	1.00	29.19	B_13
ATOM	1931	N	LEU	48	47.220	72.547	13.784	1.00	27.10	B_13
ATOM	1933	CA	LEU	48	47.915	71.302	14.124	1.00	21.49	B_13
ATOM	1934	CB	LEU	48	48.087	70.418	12.885	1.00	16.21	B_13
ATOM	1935	CG	LEU	48	46.924	69.476	12.538	1.00	15.14	B_13
ATOM	1936	CD1	LEU	48	45.618	70.049	13.000	1.00	26.83	B_13
ATOM	1937	CD2	LEU	48	46.894	69.206	11.035	1.00	32.93	B_13
ATOM	1938	C	LEU	48	49.262	71.611	14.771	1.00	16.35	B_13
ATOM	1939	O	LEU	48	49.885	72.648	14.498	1.00	26.65	B_13
ATOM	1940	N	ASN	49	49.691	70.744	15.669	1.00	18.84	B_13
ATOM	1942	CA	ASN	49	50.956	70.940	16.354	1.00	25.67	B_13
ATOM	1943	CB	ASN	49	50.741	71.205	17.846	1.00	23.64	B_13
ATOM	1944	CG	ASN	49	49.734	72.301	18.100	1.00	23.64	B_13
ATOM	1945	OD1	ASN	49	48.895	72.192	18.989	1.00	33.47	B_13
ATOM	1946	ND2	ASN	49	49.796	73.359	17.305	1.00	37.40	B_13
ATOM	1949	C	ASN	49	51.695	69.643	16.195	1.00	22.08	B_13
ATOM	1950	O	ASN	49	51.087	68.577	16.252	1.00	23.48	B_13
ATOM	1951	N	PHE	50	52.994	69.723	15.951	1.00	25.59	B_13
ATOM	1953	CA	PHE	50	53.762	68.510	15.806	1.00	19.57	B_13
ATOM	1954	CB	PHE	50	54.258	68.343	14.380	1.00	12.47	B_13
ATOM	1955	CG	PHE	50	53.161	68.024	13.432	1.00	14.47	B_13
ATOM	1956	CD1	PHE	50	52.665	68.989	12.581	1.00	17.81	B_13
ATOM	1957	CD2	PHE	50	52.566	66.770	13.445	1.00	14.44	B_13
ATOM	1958	CE1	PHE	50	51.585	68.705	11.754	1.00	23.43	B_13
ATOM	1959	CE2	PHE	50	51.488	66.482	12.624	1.00	20.62	B_13
ATOM	1960	CZ	PHE	50	50.999	67.447	11.781	1.00	13.34	B_13
ATOM	1961	C	PHE	50	54.858	68.419	16.826	1.00	23.56	B_13
ATOM	1962	O	PHE	50	55.720	69.299	16.922	1.00	20.28	B_13
ATOM	1963	N	THR	51	54.728	67.387	17.651	1.00	26.45	B_13
ATOM	1965	CA	THR	51	55.650	67.090	18.725	1.00	29.37	B_13
ATOM	1966	CB	THR	51	54.851	66.834	20.024	1.00	28.17	B_13
ATOM	1967	OG1	THR	51	53.946	65.738	19.824	1.00	40.86	B_13
ATOM	1969	CG2	THR	51	54.032	68.078	20.393	1.00	25.37	B_13
ATOM	1970	C	THR	51	56.435	65.838	18.331	1.00	21.26	B_13
ATOM	1971	O	THR	51	55.849	64.849	17.882	1.00	17.45	B_13
ATOM	1972	N	ARG	52	57.755	65.889	18.477	1.00	15.17	B_13
ATOM	1974	CA	ARG	52	58.604	64.752	18.126	1.00	20.79	B_13
ATOM	1975	CB	ARG	52	59.868	65.241	17.429	1.00	20.81	B_13
ATOM	1976	CG	ARG	52	60.871	64.160	17.110	1.00	19.06	B_13
ATOM	1977	CD	ARG	52	62.208	64.808	16.880	1.00	22.17	B_13
ATOM	1978	NE	ARG	52	63.293	63.848	16.904	1.00	18.57	B_13
ATOM	1980	CZ	ARG	52	64.563	64.160	17.108	1.00	10.00	B_13
ATOM	1981	NH1	ARG	52	64.915	65.414	17.315	1.00	19.35	B_13
ATOM	1984	NH2	ARG	52	65.488	63.214	17.039	1.00	35.90	B_13
ATOM	1987	C	ARG	52	58.995	63.903	19.328	1.00	22.29	B_13
ATOM	1988	O	ARG	52	59.326	64.433	20.387	1.00	24.98	B_13
ATOM	1989	N	LEU	53	59.013	62.586	19.140	1.00	19.90	B_13
ATOM	1991	CA	LEU	53	59.378	61.660	20.203	1.00	27.02	B_13
ATOM	1992	CB	LEU	53	58.279	60.625	20.434	1.00	16.80	B_13
ATOM	1993	CG	LEU	53	56.859	61.138	20.639	1.00	23.45	B_13
ATOM	1994	CD1	LEU	53	55.943	59.943	20.884	1.00	24.07	B_13
ATOM	1995	CD2	LEU	53	56.801	62.143	21.785	1.00	21.02	B_13
ATOM	1996	C	LEU	53	60.657	60.944	19.813	1.00	15.08	B_13
ATOM	1997	O	LEU	53	60.822	60.539	18.671	1.00	13.89	B_13
ATOM	1998	N	HIS	54	61.532	60.750	20.792	1.00	19.96	B_13
ATOM	2000	CA	HIS	54	62.812	60.079	20.568	1.00	28.80	B_13
ATOM	2001	CB	HIS	54	63.848	60.604	21.569	1.00	19.40	B_13
ATOM	2002	CG	HIS	54	64.113	62.075	21.431	1.00	31.96	B_13
ATOM	2003	CD2	HIS	54	63.365	63.060	20.883	1.00	21.32	B_13
ATOM	2004	ND1	HIS	54	65.292	62.662	21.835	1.00	33.94	B_13
ATOM	2006	CE1	HIS	54	65.260	63.949	21.539	1.00	18.64	B_13
ATOM	2007	NE2	HIS	54	64.103	64.218	20.960	1.00	19.56	B_13

59/75

ATOM	2009	C	HIS	54	62.695	58.555	20.647	1.00	13.04	B_13
ATOM	2010	O	HIS	54	63.620	57.850	20.282	1.00	19.90	B_13
ATOM	2011	N	ASP	55	61.586	58.076	21.219	1.00	17.27	B_13
ATOM	2013	CA	ASP	55	61.303	56.648	21.366	1.00	25.79	B_13
ATOM	2014	CB	ASP	55	62.099	56.038	22.533	1.00	29.40	B_13
ATOM	2015	CG	ASP	55	63.443	55.428	22.076	1.00	29.64	B_13
ATOM	2016	OD1	ASP	55	63.517	54.906	20.942	1.00	33.28	B_13
ATOM	2017	OD2	ASP	55	64.437	55.469	22.831	1.00	31.99	B_13
ATOM	2018	C	ASP	55	59.807	56.460	21.567	1.00	24.99	B_13
ATOM	2019	O	ASP	55	59.079	57.445	21.677	1.00	21.06	B_13
ATOM	2020	N	GLY	56	59.358	55.207	21.559	1.00	22.90	B_13
ATOM	2022	CA	GLY	56	57.954	54.877	21.737	1.00	21.80	B_13
ATOM	2023	C	GLY	56	57.155	54.926	20.447	1.00	14.48	B_13
ATOM	2024	O	GLY	56	57.720	55.108	19.379	1.00	19.38	B_13
ATOM	2025	N	ILE	57	55.841	54.742	20.545	1.00	11.78	B_13
ATOM	2027	CA	ILE	57	54.944	54.809	19.389	1.00	16.25	B_13
ATOM	2028	CB	ILE	57	53.737	53.804	19.510	1.00	22.94	B_13
ATOM	2029	CG2	ILE	57	52.442	54.417	18.955	1.00	24.79	B_13
ATOM	2030	CG1	ILE	57	54.025	52.505	18.744	1.00	25.63	B_13
ATOM	2031	CD1	ILE	57	53.586	52.520	17.240	1.00	17.48	B_13
ATOM	2032	C	ILE	57	54.410	56.238	19.301	1.00	18.78	B_13
ATOM	2033	O	ILE	57	53.866	56.777	20.270	1.00	11.40	B_13
ATOM	2034	N	ALA	58	54.598	56.842	18.140	1.00	14.67	B_13
ATOM	2036	CA	ALA	58	54.139	58.200	17.857	1.00	17.04	B_13
ATOM	2037	CB	ALA	58	55.270	59.015	17.245	1.00	10.00	B_13
ATOM	2038	C	ALA	58	53.048	58.009	16.825	1.00	25.41	B_13
ATOM	2039	O	ALA	58	52.956	56.940	16.243	1.00	22.59	B_13
ATOM	2040	N	ASP	59	52.211	59.020	16.609	1.00	13.36	B_13
ATOM	2042	CA	ASP	59	51.156	58.927	15.606	1.00	24.67	B_13
ATOM	2043	CB	ASP	59	50.348	60.237	15.545	1.00	10.00	B_13
ATOM	2044	CG	ASP	59	49.743	60.631	16.899	1.00	12.93	B_13
ATOM	2045	OD1	ASP	59	49.922	61.788	17.327	1.00	32.89	B_13
ATOM	2046	OD2	ASP	59	49.076	59.793	17.541	1.00	21.52	B_13
ATOM	2047	C	ASP	59	51.784	58.653	14.242	1.00	11.46	B_13
ATOM	2048	O	ASP	59	51.378	57.736	13.531	1.00	16.58	B_13
ATOM	2049	N	ILE	60	52.791	59.445	13.899	1.00	24.90	B_13
ATOM	2051	CA	ILE	60	53.494	59.346	12.624	1.00	12.17	B_13
ATOM	2052	CB	ILE	60	53.620	60.738	11.975	1.00	10.91	B_13
ATOM	2053	CG2	ILE	60	54.289	60.641	10.588	1.00	10.70	B_13
ATOM	2054	CG1	ILE	60	52.228	61.367	11.851	1.00	18.58	B_13
ATOM	2055	CD1	ILE	60	52.219	62.870	11.726	1.00	12.00	B_13
ATOM	2056	C	ILE	60	54.881	58.750	12.841	1.00	12.93	B_13
ATOM	2057	O	ILE	60	55.788	59.392	13.365	1.00	16.39	B_13
ATOM	2058	N	MET	61	55.015	57.485	12.483	1.00	19.08	B_13
ATOM	2060	CA	MET	61	56.275	56.784	12.617	1.00	16.97	B_13
ATOM	2061	CB	MET	61	56.011	55.328	13.035	1.00	23.79	B_13
ATOM	2062	CG	MET	61	55.313	55.172	14.422	1.00	12.37	B_13
ATOM	2063	SD	MET	61	56.389	55.360	15.913	1.00	31.01	B_13
ATOM	2064	CE	MET	61	57.204	53.749	15.861	1.00	14.93	B_13
ATOM	2065	C	MET	61	56.995	56.888	11.265	1.00	12.72	B_13
ATOM	2066	O	MET	61	56.438	56.538	10.216	1.00	15.31	B_13
ATOM	2067	N	ILE	62	58.170	57.518	11.294	1.00	16.64	B_13
ATOM	2069	CA	ILE	62	58.978	57.739	10.097	1.00	27.48	B_13
ATOM	2070	CB	ILE	62	59.557	59.181	10.060	1.00	10.00	B_13
ATOM	2071	CG2	ILE	62	60.191	59.462	8.717	1.00	18.65	B_13
ATOM	2072	CG1	ILE	62	58.460	60.203	10.342	1.00	18.51	B_13
ATOM	2073	CD1	ILE	62	58.983	61.499	10.931	1.00	16.23	B_13
ATOM	2074	C	ILE	62	60.155	56.787	10.046	1.00	15.06	B_13
ATOM	2075	O	ILE	62	60.873	56.606	11.033	1.00	10.73	B_13
ATOM	2076	N	SER	63	60.398	56.230	8.873	1.00	19.40	B_13
ATOM	2078	CA	SER	63	61.513	55.321	8.722	1.00	13.31	B_13
ATOM	2079	CB	SER	63	61.111	53.888	9.123	1.00	17.28	B_13
ATOM	2080	OG	SER	63	59.985	53.435	8.391	1.00	13.66	B_13
ATOM	2082	C	SER	63	62.086	55.339	7.315	1.00	19.86	B_13
ATOM	2083	O	SER	63	61.441	55.766	6.347	1.00	20.93	B_13
ATOM	2084	N	PHE	64	63.338	54.914	7.237	1.00	17.78	B_13
ATOM	2086	CA	PHE	64	64.072	54.823	5.989	1.00	18.81	B_13
ATOM	2087	CB	PHE	64	65.409	55.553	6.105	1.00	16.50	B_13
ATOM	2088	CG	PHE	64	65.278	57.054	6.171	1.00	22.54	B_13
ATOM	2089	CD1	PHE	64	65.321	57.817	5.013	1.00	20.48	B_13
ATOM	2090	CD2	PHE	64	65.155	57.708	7.395	1.00	24.76	B_13
ATOM	2091	CE1	PHE	64	65.246	59.207	5.071	1.00	13.94	B_13
ATOM	2092	CE2	PHE	64	65.079	59.105	7.461	1.00	14.29	B_13
ATOM	2093	CZ	PHE	64	65.128	59.847	6.298	1.00	10.16	B_13
ATOM	2094	C	PHE	64	64.293	53.336	5.823	1.00	10.30	B_13
ATOM	2095	O	PHE	64	64.571	52.637	6.799	1.00	14.11	B_13
ATOM	2096	N	GLY	65	64.121	52.842	4.610	1.00	13.58	B_13

60/75

ATOM	2098	CA	GLY	65	64.306	51.426	4.392	1.00	14.88	B_13
ATOM	2099	C	GLY	65	64.400	51.117	2.922	1.00	14.95	B_13
ATOM	2100	O	GLY	65	64.047	51.947	2.088	1.00	12.61	B_13
ATOM	2101	N	ILE	66	64.860	49.922	2.587	1.00	10.00	B_13
ATOM	2103	CA	ILE	66	64.995	49.555	1.187	1.00	19.70	B_13
ATOM	2104	CB	ILE	66	66.483	49.344	0.791	1.00	18.92	B_13
ATOM	2105	CG2	ILE	66	67.301	50.628	1.073	1.00	10.00	B_13
ATOM	2106	CG1	ILE	66	67.078	48.178	1.582	1.00	14.64	B_13
ATOM	2107	CD1	ILE	66	68.381	47.662	1.004	1.00	17.53	B_13
ATOM	2108	C	ILE	66	64.195	48.296	0.900	1.00	15.98	B_13
ATOM	2109	O	ILE	66	63.877	47.543	1.806	1.00	20.10	B_13
ATOM	2110	N	LYS	67	63.773	48.148	-0.349	1.00	18.78	B_13
ATOM	2112	CA	LYS	67	63.019	46.980	-0.787	1.00	14.73	B_13
ATOM	2113	CB	LYS	67	63.986	45.827	-1.073	1.00	22.08	B_13
ATOM	2114	CG	LYS	67	65.107	46.142	-2.066	1.00	15.53	B_13
ATOM	2115	CD	LYS	67	64.591	46.325	-3.487	1.00	16.76	B_13
ATOM	2116	CE	LYS	67	65.573	45.763	-4.523	1.00	21.90	B_13
ATOM	2117	NZ	LYS	67	66.975	46.257	-4.394	1.00	28.03	B_13
ATOM	2121	C	LYS	67	61.945	46.548	0.218	1.00	16.24	B_13
ATOM	2122	O	LYS	67	61.136	47.360	0.649	1.00	10.25	B_13
ATOM	2123	N	GLU	68	61.968	45.293	0.630	1.00	10.00	B_13
ATOM	2125	CA	GLU	68	60.986	44.787	1.570	1.00	10.00	B_13
ATOM	2126	CB	GLU	68	61.004	43.257	1.505	1.00	31.44	B_13
ATOM	2127	CG	GLU	68	59.733	42.550	1.696	1.00	27.13	B_13
ATOM	2128	CD	GLU	68	58.723	42.720	0.524	1.00	12.88	B_13
ATOM	2129	OE1	GLU	68	59.106	42.180	-0.613	1.00	14.05	B_13
ATOM	2130	OE2	GLU	68	57.681	43.274	0.753	1.00	38.61	B_13
ATOM	2131	C	GLU	68	61.402	45.292	2.954	1.00	32.89	B_13
ATOM	2132	O	GLU	68	62.541	45.099	3.390	1.00	19.77	B_13
ATOM	2133	N	HIS	69	60.467	45.918	3.659	1.00	15.43	B_13
ATOM	2135	CA	HIS	69	60.777	46.473	4.964	1.00	10.00	B_13
ATOM	2136	CB	HIS	69	61.173	47.928	4.802	1.00	15.60	B_13
ATOM	2137	CG	HIS	69	60.151	48.731	4.063	1.00	18.06	B_13
ATOM	2138	CD2	HIS	69	59.131	49.509	4.498	1.00	25.01	B_13
ATOM	2139	ND1	HIS	69	60.055	48.709	2.689	1.00	21.79	B_13
ATOM	2141	CE1	HIS	69	59.023	49.430	2.308	1.00	19.43	B_13
ATOM	2142	NE2	HIS	69	58.438	49.932	3.384	1.00	19.23	B_13
ATOM	2143	C	HIS	69	59.655	46.396	5.978	1.00	16.27	B_13
ATOM	2144	O	HIS	69	59.689	47.099	6.969	1.00	13.47	B_13
ATOM	2145	N	GLY	70	58.610	45.629	5.719	1.00	21.21	B_13
ATOM	2147	CA	GLY	70	57.567	45.520	6.720	1.00	15.93	B_13
ATOM	2148	C	GLY	70	56.147	45.784	6.287	1.00	13.13	B_13
ATOM	2149	O	GLY	70	55.283	45.986	7.147	1.00	12.19	B_13
ATOM	2150	N	ASP	71	55.891	45.805	4.983	1.00	10.00	B_13
ATOM	2152	CA	ASP	71	54.540	46.030	4.480	1.00	17.84	B_13
ATOM	2153	CB	ASP	71	54.086	47.490	4.636	1.00	21.86	B_13
ATOM	2154	CG	ASP	71	54.946	48.480	3.881	1.00	13.38	B_13
ATOM	2155	OD1	ASP	71	54.896	49.644	4.291	1.00	10.00	B_13
ATOM	2156	OD2	ASP	71	55.633	48.135	2.897	1.00	10.00	B_13
ATOM	2157	C	ASP	71	54.313	45.557	3.064	1.00	27.18	B_13
ATOM	2158	O	ASP	71	55.221	45.068	2.416	1.00	16.61	B_13
ATOM	2159	N	PHE	72	53.103	45.759	2.564	1.00	10.00	B_13
ATOM	2161	CA	PHE	72	52.788	45.317	1.213	1.00	19.60	B_13
ATOM	2162	CB	PHE	72	51.292	45.017	1.099	1.00	16.43	B_13
ATOM	2163	CG	PHE	72	50.849	43.779	1.851	1.00	27.69	B_13
ATOM	2164	CD1	PHE	72	51.399	42.532	1.561	1.00	22.33	B_13
ATOM	2165	CD2	PHE	72	49.848	43.855	2.823	1.00	27.58	B_13
ATOM	2166	CE1	PHE	72	50.955	41.383	2.225	1.00	22.03	B_13
ATOM	2167	CE2	PHE	72	49.403	42.709	3.486	1.00	21.82	B_13
ATOM	2168	CZ	PHE	72	49.957	41.473	3.184	1.00	10.00	B_13
ATOM	2169	C	PHE	72	53.225	46.313	0.130	1.00	18.56	B_13
ATOM	2170	O	PHE	72	52.840	46.190	-1.048	1.00	14.78	B_13
ATOM	2171	N	TYR	73	54.079	47.260	0.513	1.00	10.93	B_13
ATOM	2173	CA	TYR	73	54.558	48.295	-0.416	1.00	13.87	B_13
ATOM	2174	CB	TYR	73	53.943	49.649	-0.048	1.00	22.69	B_13
ATOM	2175	CG	TYR	73	52.439	49.581	0.007	1.00	16.43	B_13
ATOM	2176	CD1	TYR	73	51.774	49.385	1.219	1.00	18.21	B_13
ATOM	2177	CE1	TYR	73	50.386	49.219	1.257	1.00	35.13	B_13
ATOM	2178	CD2	TYR	73	51.683	49.618	-1.165	1.00	15.77	B_13
ATOM	2179	CE2	TYR	73	50.300	49.456	-1.133	1.00	39.16	B_13
ATOM	2180	CZ	TYR	73	49.663	49.258	0.080	1.00	28.27	B_13
ATOM	2181	OH	TYR	73	48.301	49.122	0.106	1.00	33.06	B_13
ATOM	2183	C	TYR	73	56.088	48.349	-0.425	1.00	18.05	B_13
ATOM	2184	O	TYR	73	56.721	49.339	0.003	1.00	10.00	B_13
ATOM	2185	N	PRO	74	56.702	47.287	-0.953	1.00	13.76	B_13
ATOM	2186	CD	PRO	74	56.063	46.221	-1.740	1.00	14.21	B_13
ATOM	2187	CA	PRO	74	58.158	47.183	-1.024	1.00	21.66	B_13

61/75

ATOM	2188	CB	PRO	74	58.353	45.768	-1.569	1.00	15.88	B_13
ATOM	2189	CG	PRO	74	57.225	45.653	-2.540	1.00	13.95	B_13
ATOM	2190	C	PRO	74	58.747	48.226	-1.959	1.00	27.68	B_13
ATOM	2191	O	PRO	74	58.173	48.526	-3.012	1.00	21.90	B_13
ATOM	2192	N	PHE	75	59.883	48.794	-1.562	1.00	20.91	B_13
ATOM	2194	CA	PHE	75	60.554	49.773	-2.395	1.00	15.84	B_13
ATOM	2195	CB	PHE	75	61.498	50.637	-1.548	1.00	11.67	B_13
ATOM	2196	CG	PHE	75	60.765	51.589	-0.641	1.00	14.42	B_13
ATOM	2197	CD1	PHE	75	59.831	52.484	-1.162	1.00	16.56	B_13
ATOM	2198	CD2	PHE	75	60.976	51.574	0.726	1.00	10.00	B_13
ATOM	2199	CE1	PHE	75	59.119	53.345	-0.327	1.00	11.14	B_13
ATOM	2200	CE2	PHE	75	60.274	52.423	1.558	1.00	10.28	B_13
ATOM	2201	CZ	PHE	75	59.340	53.316	1.027	1.00	10.00	B_13
ATOM	2202	C	PHE	75	61.236	49.068	-3.573	1.00	14.23	B_13
ATOM	2203	O	PHE	75	61.357	47.837	-3.582	1.00	18.64	B_13
ATOM	2204	N	ASP	76	61.742	49.845	-4.526	1.00	12.83	B_13
ATOM	2206	CA	ASP	76	62.330	49.287	-5.740	1.00	20.69	B_13
ATOM	2207	CB	ASP	76	61.394	49.644	-6.911	1.00	14.28	B_13
ATOM	2208	CG	ASP	76	61.212	51.144	-7.080	1.00	14.37	B_13
ATOM	2209	OD1	ASP	76	61.361	51.882	-6.095	1.00	22.32	B_13
ATOM	2210	OD2	ASP	76	60.941	51.597	-8.202	1.00	15.92	B_13
ATOM	2211	C	ASP	76	63.764	49.698	-6.104	1.00	19.31	B_13
ATOM	2212	O	ASP	76	64.056	49.864	-7.278	1.00	18.67	B_13
ATOM	2213	N	GLY	77	64.653	49.902	-5.132	1.00	10.00	B_13
ATOM	2215	CA	GLY	77	65.997	50.326	-5.501	1.00	10.00	B_13
ATOM	2216	C	GLY	77	65.989	51.790	-5.970	1.00	16.22	B_13
ATOM	2217	O	GLY	77	64.967	52.487	-5.752	1.00	17.04	B_13
ATOM	2218	N	PRO	78	67.080	52.305	-6.589	1.00	12.53	B_13
ATOM	2219	CD	PRO	78	68.319	51.564	-6.856	1.00	12.24	B_13
ATOM	2220	CA	PRO	78	67.207	53.691	-7.086	1.00	11.81	B_13
ATOM	2221	CB	PRO	78	68.546	53.678	-7.816	1.00	10.00	B_13
ATOM	2222	CG	PRO	78	69.316	52.693	-7.066	1.00	12.78	B_13
ATOM	2223	C	PRO	78	66.093	54.146	-8.027	1.00	10.00	B_13
ATOM	2224	O	PRO	78	65.621	53.381	-8.853	1.00	27.46	B_13
ATOM	2225	N	SER	79	65.641	55.386	-7.852	1.00	19.14	B_13
ATOM	2227	CA	SER	79	64.568	55.963	-8.669	1.00	10.00	B_13
ATOM	2228	CB	SER	79	64.970	56.033	-10.148	1.00	20.11	B_13
ATOM	2229	OG	SER	79	63.982	56.723	-10.901	1.00	23.87	B_13
ATOM	2231	C	SER	79	63.231	55.215	-8.507	1.00	31.68	B_13
ATOM	2232	O	SER	79	63.074	54.356	-7.606	1.00	26.48	B_13
ATOM	2233	N	GLY	80	62.250	55.589	-9.327	1.00	10.00	B_13
ATOM	2235	CA	GLY	80	60.940	54.969	-9.260	1.00	10.07	B_13
ATOM	2236	C	GLY	80	60.293	55.412	-7.968	1.00	30.72	B_13
ATOM	2237	O	GLY	80	60.347	56.600	-7.643	1.00	20.65	B_13
ATOM	2238	N	LEU	81	59.779	54.452	-7.193	1.00	23.74	B_13
ATOM	2240	CA	LEU	81	59.135	54.752	-5.917	1.00	13.14	B_13
ATOM	2241	CB	LEU	81	58.661	53.481	-5.213	1.00	16.20	B_13
ATOM	2242	CG	LEU	81	57.393	52.775	-5.687	1.00	17.33	B_13
ATOM	2243	CD1	LEU	81	57.554	52.277	-7.096	1.00	28.67	B_13
ATOM	2244	CD2	LEU	81	57.103	51.617	-4.745	1.00	27.02	B_13
ATOM	2245	C	LEU	81	60.122	55.466	-5.019	1.00	14.51	B_13
ATOM	2246	O	LEU	81	61.264	55.016	-4.846	1.00	16.24	B_13
ATOM	2247	N	LEU	82	59.692	56.590	-4.470	1.00	11.33	B_13
ATOM	2249	CA	LEU	82	60.540	57.381	-3.594	1.00	17.52	B_13
ATOM	2250	CB	LEU	82	60.442	58.861	-3.986	1.00	18.51	B_13
ATOM	2251	CG	LEU	82	61.355	59.499	-5.044	1.00	15.37	B_13
ATOM	2252	CD1	LEU	82	61.800	58.504	-6.104	1.00	17.05	B_13
ATOM	2253	CD2	LEU	82	60.639	60.744	-5.659	1.00	16.87	B_13
ATOM	2254	C	LEU	82	60.172	57.203	-2.127	1.00	10.00	B_13
ATOM	2255	O	LEU	82	61.045	57.056	-1.275	1.00	19.90	B_13
ATOM	2256	N	ALA	83	58.876	57.201	-1.840	1.00	18.16	B_13
ATOM	2258	CA	ALA	83	58.378	57.077	-0.472	1.00	13.17	B_13
ATOM	2259	CB	ALA	83	58.762	58.322	0.327	1.00	10.00	B_13
ATOM	2260	C	ALA	83	56.846	56.925	-0.500	1.00	10.00	B_13
ATOM	2261	O	ALA	83	56.209	57.155	-1.541	1.00	10.73	B_13
ATOM	2262	N	HIS	84	56.268	56.619	0.662	1.00	10.00	B_13
ATOM	2264	CA	HIS	84	54.811	56.472	0.810	1.00	23.81	B_13
ATOM	2265	CB	HIS	84	54.270	55.188	0.157	1.00	30.45	B_13
ATOM	2266	CG	HIS	84	54.848	53.925	0.711	1.00	17.68	B_13
ATOM	2267	CD2	HIS	84	54.856	53.415	1.964	1.00	10.00	B_13
ATOM	2268	ND1	HIS	84	55.525	53.025	-0.076	1.00	14.94	B_13
ATOM	2270	CE1	HIS	84	55.933	52.015	0.666	1.00	29.72	B_13
ATOM	2271	NE2	HIS	84	55.543	52.224	1.912	1.00	13.81	B_13
ATOM	2272	C	HIS	84	54.363	56.547	2.258	1.00	12.82	B_13
ATOM	2273	O	HIS	84	55.099	56.148	3.166	1.00	20.02	B_13
ATOM	2274	N	ALA	85	53.161	57.076	2.464	1.00	28.38	B_13
ATOM	2276	CA	ALA	85	52.584	57.230	3.796	1.00	18.64	B_13

62/75

ATOM	2277	CB	ALA	85	52.638	58.705	4.223	1.00	13.89	B_13
ATOM	2278	C	ALA	85	51.138	56.716	3.837	1.00	10.00	B_13
ATOM	2279	O	ALA	85	50.434	56.728	2.828	1.00	10.00	B_13
ATOM	2280	N	PHE	86	50.676	56.322	5.016	1.00	14.76	B_13
ATOM	2282	CA	PHE	86	49.316	55.811	5.143	1.00	17.96	B_13
ATOM	2283	CB	PHE	86	49.286	54.592	6.084	1.00	15.86	B_13
ATOM	2284	CG	PHE	86	50.320	53.542	5.748	1.00	26.30	B_13
ATOM	2285	CD1	PHE	86	49.973	52.398	5.042	1.00	22.30	B_13
ATOM	2286	CD2	PHE	86	51.654	53.730	6.090	1.00	27.63	B_13
ATOM	2287	CE1	PHE	86	50.938	51.472	4.681	1.00	27.85	B_13
ATOM	2288	CE2	PHE	86	52.620	52.810	5.731	1.00	13.97	B_13
ATOM	2289	CZ	PHE	86	52.266	51.683	5.027	1.00	23.08	B_13
ATOM	2290	C	PHE	86	48.427	56.924	5.669	1.00	13.02	B_13
ATOM	2291	O	PHE	86	48.870	57.747	6.466	1.00	15.02	B_13
ATOM	2292	N	PRO	87	47.174	57.006	5.186	1.00	17.55	B_13
ATOM	2293	CD	PRO	87	46.565	56.165	4.146	1.00	10.17	B_13
ATOM	2294	CA	PRO	87	46.228	58.041	5.628	1.00	32.09	B_13
ATOM	2295	CB	PRO	87	44.961	57.720	4.819	1.00	18.55	B_13
ATOM	2296	CG	PRO	87	45.115	56.277	4.481	1.00	18.86	B_13
ATOM	2297	C	PRO	87	45.995	57.955	7.139	1.00	25.18	B_13
ATOM	2298	O	PRO	87	46.284	56.919	7.752	1.00	18.18	B_13
ATOM	2299	N	PRO	88	45.462	59.032	7.760	1.00	11.49	B_13
ATOM	2300	CD	PRO	88	45.015	60.303	7.164	1.00	10.00	B_13
ATOM	2301	CA	PRO	88	45.217	59.034	9.202	1.00	19.03	B_13
ATOM	2302	CB	PRO	88	44.399	60.302	9.402	1.00	14.16	B_13
ATOM	2303	CG	PRO	88	44.939	61.196	8.357	1.00	16.39	B_13
ATOM	2304	C	PRO	88	44.500	57.787	9.733	1.00	25.43	B_13
ATOM	2305	O	PRO	88	43.670	57.165	9.044	1.00	15.90	B_13
ATOM	2306	N	GLY	89	44.865	57.422	10.955	1.00	26.28	B_13
ATOM	2308	CA	GLY	89	44.299	56.264	11.606	1.00	25.32	B_13
ATOM	2309	C	GLY	89	45.343	55.713	12.546	1.00	34.38	B_13
ATOM	2310	O	GLY	89	46.485	56.164	12.498	1.00	23.28	B_13
ATOM	2311	N	PRO	90	44.977	54.774	13.437	1.00	13.87	B_13
ATOM	2312	CD	PRO	90	43.613	54.259	13.631	1.00	16.36	B_13
ATOM	2313	CA	PRO	90	45.898	54.164	14.398	1.00	10.34	B_13
ATOM	2314	CB	PRO	90	44.963	53.360	15.300	1.00	15.93	B_13
ATOM	2315	CG	PRO	90	43.870	52.975	14.373	1.00	23.25	B_13
ATOM	2316	C	PRO	90	46.942	53.299	13.711	1.00	18.38	B_13
ATOM	2317	O	PRO	90	46.875	53.064	12.505	1.00	26.81	B_13
ATOM	2318	N	ASN	91	47.903	52.831	14.502	1.00	26.63	B_13
ATOM	2320	CA	ASN	91	49.022	52.010	14.033	1.00	21.91	B_13
ATOM	2321	CB	ASN	91	48.740	50.500	14.081	1.00	18.89	B_13
ATOM	2322	CG	ASN	91	47.437	50.117	13.448	1.00	22.49	B_13
ATOM	2323	OD1	ASN	91	47.335	50.017	12.237	1.00	29.37	B_13
ATOM	2324	ND2	ASN	91	46.438	49.858	14.273	1.00	28.01	B_13
ATOM	2327	C	ASN	91	49.656	52.438	12.721	1.00	20.07	B_13
ATOM	2328	O	ASN	91	50.301	53.479	12.681	1.00	21.24	B_13
ATOM	2329	N	TYR	92	49.423	51.716	11.633	1.00	20.15	B_13
ATOM	2331	CA	TYR	92	50.052	52.081	10.367	1.00	18.70	B_13
ATOM	2332	CB	TYR	92	49.905	50.953	9.344	1.00	14.48	B_13
ATOM	2333	CG	TYR	92	50.906	49.821	9.567	1.00	24.41	B_13
ATOM	2334	CD1	TYR	92	52.266	50.003	9.287	1.00	27.39	B_13
ATOM	2335	CE1	TYR	92	53.198	48.979	9.471	1.00	18.14	B_13
ATOM	2336	CD2	TYR	92	50.499	48.571	10.044	1.00	28.07	B_13
ATOM	2337	CE2	TYR	92	51.427	47.529	10.230	1.00	36.50	B_13
ATOM	2338	CZ	TYR	92	52.778	47.741	9.940	1.00	43.64	B_13
ATOM	2339	OH	TYR	92	53.694	46.710	10.105	1.00	32.21	B_13
ATOM	2341	C	TYR	92	49.633	53.431	9.797	1.00	21.78	B_13
ATOM	2342	O	TYR	92	50.384	54.049	9.040	1.00	12.55	B_13
ATOM	2343	N	GLY	93	48.464	53.916	10.198	1.00	15.83	B_13
ATOM	2345	CA	GLY	93	48.015	55.216	9.732	1.00	11.69	B_13
ATOM	2346	C	GLY	93	48.971	56.326	10.134	1.00	18.60	B_13
ATOM	2347	O	GLY	93	49.561	56.300	11.227	1.00	22.00	B_13
ATOM	2348	N	GLY	94	49.205	57.258	9.216	1.00	10.27	B_13
ATOM	2350	CA	GLY	94	50.099	58.365	9.492	1.00	18.36	B_13
ATOM	2351	C	GLY	94	51.567	58.061	9.234	1.00	15.54	B_13
ATOM	2352	O	GLY	94	52.334	58.967	8.938	1.00	17.55	B_13
ATOM	2353	N	ASP	95	51.977	56.801	9.351	1.00	17.69	B_13
ATOM	2355	CA	ASP	95	53.386	56.457	9.134	1.00	19.67	B_13
ATOM	2356	CB	ASP	95	53.637	54.986	9.444	1.00	15.96	B_13
ATOM	2357	CG	ASP	95	53.346	54.634	10.900	1.00	25.37	B_13
ATOM	2358	OD1	ASP	95	53.627	53.484	11.297	1.00	16.05	B_13
ATOM	2359	OD2	ASP	95	52.835	55.488	11.656	1.00	14.66	B_13
ATOM	2360	C	ASP	95	53.896	56.808	7.733	1.00	17.15	B_13
ATOM	2361	O	ASP	95	53.162	56.711	6.746	1.00	19.09	B_13
ATOM	2362	N	ALA	96	55.166	57.198	7.662	1.00	18.71	B_13
ATOM	2364	CA	ALA	96	55.803	57.581	6.400	1.00	19.97	B_13

63/75

ATOM	2365	CB	ALA	96	56.098	59.095	6.379	1.00	22.61	B_13
ATOM	2366	C	ALA	96	57.088	56.784	6.204	1.00	25.63	B_13
ATOM	2367	O	ALA	96	57.948	56.724	7.095	1.00	12.54	B_13
ATOM	2368	N	HIS	97	57.211	56.166	5.035	1.00	13.27	B_13
ATOM	2370	CA	HIS	97	58.375	55.357	4.730	1.00	25.28	B_13
ATOM	2371	CB	HIS	97	57.955	53.905	4.464	1.00	10.00	B_13
ATOM	2372	CG	HIS	97	57.264	53.257	5.624	1.00	12.02	B_13
ATOM	2373	CD2	HIS	97	57.214	53.603	6.929	1.00	10.00	B_13
ATOM	2374	ND1	HIS	97	56.516	52.104	5.499	1.00	12.91	B_13
ATOM	2375	CE1	HIS	97	56.038	51.770	6.688	1.00	10.00	B_13
ATOM	2376	NE2	HIS	97	56.445	52.664	7.571	1.00	10.64	B_13
ATOM	2378	C	HIS	97	59.069	55.959	3.520	1.00	13.82	B_13
ATOM	2379	O	HIS	97	58.415	56.273	2.517	1.00	12.27	B_13
ATOM	2380	N	PHE	98	60.379	56.154	3.647	1.00	10.67	B_13
ATOM	2382	CA	PHE	98	61.224	56.718	2.595	1.00	15.67	B_13
ATOM	2383	CB	PHE	98	61.970	57.938	3.156	1.00	10.76	B_13
ATOM	2384	CG	PHE	98	61.055	59.025	3.627	1.00	17.93	B_13
ATOM	2385	CD1	PHE	98	60.730	60.082	2.786	1.00	18.92	B_13
ATOM	2386	CD2	PHE	98	60.476	58.974	4.893	1.00	14.14	B_13
ATOM	2387	CE1	PHE	98	59.833	61.066	3.201	1.00	22.42	B_13
ATOM	2388	CE2	PHE	98	59.574	59.962	5.315	1.00	10.00	B_13
ATOM	2389	CZ	PHE	98	59.257	61.002	4.469	1.00	10.00	B_13
ATOM	2390	C	PHE	98	62.218	55.669	2.064	1.00	26.64	B_13
ATOM	2391	O	PHE	98	62.882	54.969	2.851	1.00	13.27	B_13
ATOM	2392	N	ASP	99	62.331	55.577	0.738	1.00	12.24	B_13
ATOM	2394	CA	ASP	99	63.229	54.612	0.102	1.00	10.00	B_13
ATOM	2395	CB	ASP	99	62.884	54.471	-1.385	1.00	10.00	B_13
ATOM	2396	CG	ASP	99	63.615	53.311	-2.067	1.00	22.86	B_13
ATOM	2397	OD1	ASP	99	63.170	52.890	-3.160	1.00	11.60	B_13
ATOM	2398	OD2	ASP	99	64.624	52.806	-1.528	1.00	21.20	B_13
ATOM	2399	C	ASP	99	64.677	55.046	0.264	1.00	12.66	B_13
ATOM	2400	O	ASP	99	65.121	56.010	-0.366	1.00	18.37	B_13
ATOM	2401	N	ASP	100	65.439	54.289	1.046	1.00	12.86	B_13
ATOM	2403	CA	ASP	100	66.833	54.642	1.260	1.00	14.46	B_13
ATOM	2404	CB	ASP	100	67.308	54.271	2.660	1.00	17.70	B_13
ATOM	2405	CG	ASP	100	68.006	55.437	3.358	1.00	16.15	B_13
ATOM	2406	OD1	ASP	100	68.091	55.447	4.602	1.00	15.74	B_13
ATOM	2407	OD2	ASP	100	68.470	56.354	2.655	1.00	27.08	B_13
ATOM	2408	C	ASP	100	67.793	54.171	0.179	1.00	13.66	B_13
ATOM	2409	O	ASP	100	68.961	53.932	0.416	1.00	19.54	B_13
ATOM	2410	N	ASP	101	67.254	53.954	-1.010	1.00	12.83	B_13
ATOM	2412	CA	ASP	101	68.074	53.590	-2.164	1.00	10.00	B_13
ATOM	2413	CB	ASP	101	67.471	52.413	-2.933	1.00	10.00	B_13
ATOM	2414	CG	ASP	101	67.997	51.065	-2.449	1.00	16.87	B_13
ATOM	2415	OD1	ASP	101	67.232	50.089	-2.458	1.00	19.89	B_13
ATOM	2416	OD2	ASP	101	69.184	50.968	-2.066	1.00	18.51	B_13
ATOM	2417	C	ASP	101	68.108	54.858	-3.029	1.00	26.72	B_13
ATOM	2418	O	ASP	101	68.602	54.853	-4.172	1.00	12.11	B_13
ATOM	2419	N	GLU	102	67.500	55.922	-2.496	1.00	13.76	B_13
ATOM	2421	CA	GLU	102	67.462	57.217	-3.161	1.00	12.54	B_13
ATOM	2422	CB	GLU	102	66.135	57.958	-2.916	1.00	13.01	B_13
ATOM	2423	CG	GLU	102	64.873	57.257	-3.381	1.00	15.50	B_13
ATOM	2424	CD	GLU	102	64.973	56.707	-4.791	1.00	29.02	B_13
ATOM	2425	OE1	GLU	102	65.640	57.307	-5.665	1.00	12.78	B_13
ATOM	2426	OE2	GLU	102	64.399	55.635	-5.021	1.00	12.36	B_13
ATOM	2427	C	GLU	102	68.544	58.040	-2.505	1.00	14.96	B_13
ATOM	2428	O	GLU	102	68.939	57.760	-1.371	1.00	10.00	B_13
ATOM	2429	N	THR	103	69.030	59.039	-3.228	1.00	19.38	B_13
ATOM	2431	CA	THR	103	70.021	59.957	-2.693	1.00	16.49	B_13
ATOM	2432	CB	THR	103	70.973	60.490	-3.801	1.00	19.31	B_13
ATOM	2433	OG1	THR	103	71.661	59.384	-4.399	1.00	25.44	B_13
ATOM	2435	CG2	THR	103	72.006	61.462	-3.212	1.00	10.75	B_13
ATOM	2436	C	THR	103	69.180	61.104	-2.141	1.00	12.91	B_13
ATOM	2437	O	THR	103	68.414	61.727	-2.867	1.00	13.59	B_13
ATOM	2438	N	TRP	104	69.252	61.322	-0.842	1.00	20.60	B_13
ATOM	2440	CA	TRP	104	68.497	62.388	-0.237	1.00	13.62	B_13
ATOM	2441	CB	TRP	104	67.852	61.902	1.063	1.00	22.66	B_13
ATOM	2442	CG	TRP	104	66.837	60.808	0.870	1.00	22.99	B_13
ATOM	2443	CD2	TRP	104	65.505	60.953	0.347	1.00	27.35	B_13
ATOM	2444	CE2	TRP	104	64.936	59.654	0.287	1.00	12.61	B_13
ATOM	2445	CE3	TRP	104	64.741	62.054	-0.079	1.00	11.89	B_13
ATOM	2446	CD1	TRP	104	67.013	59.473	1.108	1.00	17.89	B_13
ATOM	2447	NE1	TRP	104	65.876	58.775	0.755	1.00	14.24	B_13
ATOM	2449	CZ2	TRP	104	63.632	59.429	-0.186	1.00	10.00	B_13
ATOM	2450	CZ3	TRP	104	63.445	61.832	-0.549	1.00	22.21	B_13
ATOM	2451	CH2	TRP	104	62.904	60.527	-0.598	1.00	23.31	B_13
ATOM	2452	C	TRP	104	69.416	63.570	0.033	1.00	16.43	B_13

64/75

ATOM	2453	O	TRP	104	70.520	63.380	0.526	1.00	11.13	B_13
ATOM	2454	N	THR	105	68.960	64.775	-0.322	1.00	19.48	B_13
ATOM	2456	CA	THR	105	69.716	66.015	-0.097	1.00	10.40	B_13
ATOM	2457	CB	THR	105	70.153	66.749	-1.398	1.00	10.00	B_13
ATOM	2458	OG1	THR	105	69.305	66.401	-2.501	1.00	18.53	B_13
ATOM	2460	CG2	THR	105	71.596	66.484	-1.709	1.00	34.62	B_13
ATOM	2461	C	THR	105	68.904	67.062	0.641	1.00	20.82	B_13
ATOM	2462	O	THR	105	67.686	66.952	0.768	1.00	15.93	B_13
ATOM	2463	N	SER	106	69.621	68.073	1.125	1.00	38.37	B_13
ATOM	2465	CA	SER	106	69.029	69.222	1.791	1.00	20.77	B_13
ATOM	2466	CB	SER	106	69.979	69.778	2.862	1.00	17.95	B_13
ATOM	2467	OG	SER	106	70.281	68.825	3.864	1.00	29.88	B_13
ATOM	2469	C	SER	106	68.889	70.245	0.657	1.00	19.23	B_13
ATOM	2470	O	SER	106	68.202	71.260	0.782	1.00	21.34	B_13
ATOM	2471	N	SER	107	69.577	69.981	-0.450	1.00	18.73	B_13
ATOM	2473	CA	SER	107	69.533	70.884	-1.592	1.00	20.92	B_13
ATOM	2474	CB	SER	107	70.945	71.380	-1.927	1.00	19.84	B_13
ATOM	2475	OG	SER	107	71.556	71.957	-0.788	1.00	27.31	B_13
ATOM	2477	C	SER	107	68.848	70.284	-2.828	1.00	18.68	B_13
ATOM	2478	O	SER	107	67.660	69.953	-2.771	1.00	21.51	B_13
ATOM	2479	N	SER	108	69.623	70.038	-3.888	1.00	18.53	B_13
ATOM	2481	CA	SER	108	69.091	69.544	-5.152	1.00	16.21	B_13
ATOM	2482	CB	SER	108	69.285	70.632	-6.205	1.00	29.10	B_13
ATOM	2483	OG	SER	108	70.665	70.969	-6.271	1.00	21.47	B_13
ATOM	2485	C	SER	108	69.645	68.260	-5.745	1.00	17.68	B_13
ATOM	2486	O	SER	108	68.964	67.618	-6.541	1.00	19.67	B_13
ATOM	2487	N	LYS	109	70.895	67.919	-5.448	1.00	11.70	B_13
ATOM	2489	CA	LYS	109	71.468	66.721	-6.047	1.00	10.00	B_13
ATOM	2490	CB	LYS	109	72.994	66.748	-5.989	1.00	18.86	B_13
ATOM	2491	CG	LYS	109	73.657	65.833	-7.013	1.00	16.33	B_13
ATOM	2492	CD	LYS	109	75.143	65.726	-6.740	1.00	11.58	B_13
ATOM	2493	CE	LYS	109	75.787	64.655	-7.606	1.00	27.43	B_13
ATOM	2494	NZ	LYS	109	77.218	64.492	-7.251	1.00	35.03	B_13
ATOM	2498	C	LYS	109	70.916	65.428	-5.444	1.00	29.39	B_13
ATOM	2499	O	LYS	109	71.432	64.905	-4.449	1.00	29.95	B_13
ATOM	2500	N	GLY	110	69.852	64.922	-6.055	1.00	14.77	B_13
ATOM	2502	CA	GLY	110	69.227	63.705	-5.576	1.00	24.08	B_13
ATOM	2503	C	GLY	110	67.793	64.105	-5.342	1.00	20.25	B_13
ATOM	2504	O	GLY	110	67.203	64.737	-6.198	1.00	16.21	B_13
ATOM	2505	N	TYR	111	67.248	63.772	-4.182	1.00	10.00	B_13
ATOM	2507	CA	TYR	111	65.879	64.130	-3.845	1.00	24.52	B_13
ATOM	2508	CB	TYR	111	65.030	62.868	-3.688	1.00	22.46	B_13
ATOM	2509	CG	TYR	111	64.676	62.244	-4.999	1.00	10.83	B_13
ATOM	2510	CD1	TYR	111	65.380	61.155	-5.483	1.00	25.38	B_13
ATOM	2511	CE1	TYR	111	65.068	60.592	-6.720	1.00	18.68	B_13
ATOM	2512	CD2	TYR	111	63.646	62.769	-5.776	1.00	16.02	B_13
ATOM	2513	CE2	TYR	111	63.328	62.223	-7.013	1.00	31.72	B_13
ATOM	2514	CZ	TYR	111	64.041	61.131	-7.473	1.00	23.68	B_13
ATOM	2515	OH	TYR	111	63.711	60.550	-8.666	1.00	20.96	B_13
ATOM	2517	C	TYR	111	65.856	64.944	-2.553	1.00	22.83	B_13
ATOM	2518	O	TYR	111	66.410	64.518	-1.538	1.00	11.66	B_13
ATOM	2519	N	ASN	112	65.278	66.140	-2.611	1.00	17.47	B_13
ATOM	2521	CA	ASN	112	65.180	67.006	-1.431	1.00	15.77	B_13
ATOM	2522	CB	ASN	112	64.658	68.401	-1.817	1.00	15.93	B_13
ATOM	2523	CG	ASN	112	64.694	69.384	-0.657	1.00	10.00	B_13
ATOM	2524	OD1	ASN	112	63.757	69.465	0.132	1.00	15.33	B_13
ATOM	2525	ND2	ASN	112	65.754	70.180	-0.586	1.00	13.70	B_13
ATOM	2528	C	ASN	112	64.214	66.329	-0.472	1.00	17.73	B_13
ATOM	2529	O	ASN	112	63.007	66.243	-0.737	1.00	12.61	B_13
ATOM	2530	N	LEU	113	64.755	65.830	0.630	1.00	16.28	B_13
ATOM	2532	CA	LEU	113	63.962	65.121	1.619	1.00	15.93	B_13
ATOM	2533	CB	LEU	113	64.841	64.703	2.804	1.00	11.93	B_13
ATOM	2534	CG	LEU	113	64.719	63.352	3.521	1.00	17.15	B_13
ATOM	2535	CD1	LEU	113	65.002	63.640	4.987	1.00	10.00	B_13
ATOM	2536	CD2	LEU	113	63.370	62.667	3.362	1.00	16.08	B_13
ATOM	2537	C	LEU	113	62.802	65.994	2.085	1.00	14.61	B_13
ATOM	2538	O	LEU	113	61.673	65.528	2.161	1.00	17.98	B_13
ATOM	2539	N	PHE	114	63.073	67.267	2.346	1.00	16.81	B_13
ATOM	2541	CA	PHE	114	62.056	68.212	2.791	1.00	15.65	B_13
ATOM	2542	CB	PHE	114	62.638	69.630	2.888	1.00	22.16	B_13
ATOM	2543	CG	PHE	114	61.596	70.714	2.882	1.00	12.27	B_13
ATOM	2544	CD1	PHE	114	60.804	70.952	4.004	1.00	19.93	B_13
ATOM	2545	CD2	PHE	114	61.378	71.470	1.746	1.00	13.56	B_13
ATOM	2546	CE1	PHE	114	59.813	71.932	3.984	1.00	17.08	B_13
ATOM	2547	CE2	PHE	114	60.398	72.441	1.726	1.00	13.79	B_13
ATOM	2548	CZ	PHE	114	59.615	72.666	2.848	1.00	10.70	B_13
ATOM	2549	C	PHE	114	60.860	68.220	1.842	1.00	19.55	B_13

65/75

ATOM	2550	O	PHE	114	59.714	68.156	2.285	1.00	15.97	B_13
ATOM	2551	N	LEU	115	61.135	68.309	0.543	1.00	13.35	B_13
ATOM	2553	CA	LEU	115	60.096	68.323	-0.485	1.00	17.91	B_13
ATOM	2554	CB	LEU	115	60.741	68.462	-1.868	1.00	24.65	B_13
ATOM	2555	CG	LEU	115	60.501	69.739	-2.679	1.00	22.70	B_13
ATOM	2556	CD1	LEU	115	61.033	70.939	-1.943	1.00	17.98	B_13
ATOM	2557	CD2	LEU	115	61.148	69.624	-4.048	1.00	28.50	B_13
ATOM	2558	C	LEU	115	59.235	67.042	-0.443	1.00	21.61	B_13
ATOM	2559	O	LEU	115	58.002	67.093	-0.344	1.00	13.99	B_13
ATOM	2560	N	VAL	116	59.898	65.895	-0.511	1.00	11.14	B_13
ATOM	2562	CA	VAL	116	59.199	64.616	-0.482	1.00	22.27	B_13
ATOM	2563	CB	VAL	116	60.163	63.421	-0.772	1.00	17.40	B_13
ATOM	2564	CG1	VAL	116	59.437	62.086	-0.629	1.00	23.09	B_13
ATOM	2565	CG2	VAL	116	60.741	63.534	-2.169	1.00	12.16	B_13
ATOM	2566	C	VAL	116	58.502	64.414	0.864	1.00	10.00	B_13
ATOM	2567	O	VAL	116	57.368	63.950	0.911	1.00	16.18	B_13
ATOM	2568	N	ALA	117	59.153	64.803	1.954	1.00	10.00	B_13
ATOM	2570	CA	ALA	117	58.585	64.640	3.297	1.00	19.50	B_13
ATOM	2571	CB	ALA	117	59.608	64.995	4.352	1.00	11.81	B_13
ATOM	2572	C	ALA	117	57.309	65.455	3.505	1.00	30.87	B_13
ATOM	2573	O	ALA	117	56.327	64.955	4.053	1.00	10.00	B_13
ATOM	2574	N	ALA	118	57.322	66.714	3.087	1.00	24.62	B_13
ATOM	2576	CA	ALA	118	56.140	67.553	3.222	1.00	20.76	B_13
ATOM	2577	CB	ALA	118	56.407	68.917	2.654	1.00	16.19	B_13
ATOM	2578	C	ALA	118	54.968	66.894	2.485	1.00	20.54	B_13
ATOM	2579	O	ALA	118	53.843	66.889	2.981	1.00	22.12	B_13
ATOM	2580	N	HIS	119	55.255	66.315	1.321	1.00	10.00	B_13
ATOM	2582	CA	HIS	119	54.259	65.647	0.489	1.00	17.27	B_13
ATOM	2583	CB	HIS	119	54.909	65.263	-0.860	1.00	11.16	B_13
ATOM	2584	CG	HIS	119	54.006	64.530	-1.813	1.00	26.59	B_13
ATOM	2585	CD2	HIS	119	53.377	63.335	-1.706	1.00	16.63	B_13
ATOM	2586	ND1	HIS	119	53.723	64.995	-3.085	1.00	12.44	B_13
ATOM	2588	CE1	HIS	119	52.961	64.124	-3.715	1.00	14.58	B_13
ATOM	2589	NE2	HIS	119	52.734	63.101	-2.901	1.00	26.44	B_13
ATOM	2590	C	HIS	119	53.722	64.419	1.227	1.00	17.00	B_13
ATOM	2591	O	HIS	119	52.510	64.218	1.331	1.00	17.01	B_13
ATOM	2592	N	GLU	120	54.626	63.607	1.751	1.00	10.31	B_13
ATOM	2594	CA	GLU	120	54.231	62.401	2.466	1.00	12.32	B_13
ATOM	2595	CB	GLU	120	55.463	61.627	2.961	1.00	15.34	B_13
ATOM	2596	CG	GLU	120	56.354	61.078	1.848	1.00	10.00	B_13
ATOM	2597	CD	GLU	120	55.574	60.260	0.867	1.00	18.64	B_13
ATOM	2598	OE1	GLU	120	55.598	60.565	-0.348	1.00	18.08	B_13
ATOM	2599	OE2	GLU	120	54.920	59.308	1.320	1.00	14.49	B_13
ATOM	2600	C	GLU	120	53.347	62.777	3.635	1.00	12.41	B_13
ATOM	2601	O	GLU	120	52.323	62.130	3.888	1.00	26.62	B_13
ATOM	2602	N	PHE	121	53.750	63.813	4.359	1.00	10.29	B_13
ATOM	2604	CA	PHE	121	52.993	64.286	5.506	1.00	14.37	B_13
ATOM	2605	CB	PHE	121	53.780	65.344	6.270	1.00	20.10	B_13
ATOM	2606	CG	PHE	121	55.057	64.827	6.852	1.00	24.55	B_13
ATOM	2607	CD1	PHE	121	56.037	65.700	7.292	1.00	10.00	B_13
ATOM	2608	CD2	PHE	121	55.292	63.454	6.936	1.00	23.62	B_13
ATOM	2609	CE1	PHE	121	57.247	65.212	7.813	1.00	18.59	B_13
ATOM	2610	CE2	PHE	121	56.488	62.954	7.448	1.00	15.21	B_13
ATOM	2611	CZ	PHE	121	57.472	63.834	7.888	1.00	25.40	B_13
ATOM	2612	C	PHE	121	51.607	64.791	5.110	1.00	16.63	B_13
ATOM	2613	O	PHE	121	50.676	64.760	5.921	1.00	26.80	B_13
ATOM	2614	N	GLY	122	51.471	65.238	3.864	1.00	11.98	B_13
ATOM	2616	CA	GLY	122	50.175	65.664	3.380	1.00	12.95	B_13
ATOM	2617	C	GLY	122	49.284	64.427	3.381	1.00	13.71	B_13
ATOM	2618	O	GLY	122	48.113	64.483	3.753	1.00	13.74	B_13
ATOM	2619	N	HIS	123	49.859	63.284	3.016	1.00	16.90	B_13
ATOM	2621	CA	HIS	123	49.126	62.009	3.008	1.00	24.90	B_13
ATOM	2622	CB	HIS	123	49.918	60.918	2.279	1.00	18.28	B_13
ATOM	2623	CG	HIS	123	49.945	61.084	0.794	1.00	21.62	B_13
ATOM	2624	CD2	HIS	123	50.889	60.764	-0.119	1.00	13.04	B_13
ATOM	2625	ND1	HIS	123	48.887	61.618	0.093	1.00	17.18	B_13
ATOM	2627	CE1	HIS	123	49.176	61.621	-1.195	1.00	16.02	B_13
ATOM	2628	NE2	HIS	123	50.386	61.108	-1.353	1.00	15.58	B_13
ATOM	2629	C	HIS	123	48.864	61.562	4.446	1.00	19.74	B_13
ATOM	2630	O	HIS	123	47.744	61.179	4.785	1.00	15.41	B_13
ATOM	2631	N	SER	124	49.904	61.627	5.284	1.00	13.32	B_13
ATOM	2633	CA	SER	124	49.813	61.270	6.695	1.00	27.50	B_13
ATOM	2634	CB	SER	124	51.131	61.582	7.425	1.00	18.63	B_13
ATOM	2635	OG	SER	124	52.221	60.837	6.925	1.00	13.32	B_13
ATOM	2637	C	SEA	124	48.703	62.102	7.335	1.00	13.76	B_13
ATOM	2638	O	SER	124	48.061	61.677	8.306	1.00	20.65	B_13
ATOM	2639	N	LEU	125	48.481	63.300	6.814	1.00	13.33	B_13

66/75

ATOM	2641	CA	LEU	125	47.439	64.133	7.387	1.00	24.62	B_13
ATOM	2642	CB	LEU	125	47.893	65.592	7.436	1.00	20.76	B_13
ATOM	2643	CG	LEU	125	49.076	65.849	8.383	1.00	14.66	B_13
ATOM	2644	CD1	LEU	125	49.739	67.159	8.064	1.00	16.16	B_13
ATOM	2645	CD2	LEU	125	48.610	65.811	9.822	1.00	16.44	B_13
ATOM	2646	C	LEU	125	46.058	63.966	6.724	1.00	24.77	B_13
ATOM	2647	O	LEU	125	45.066	64.528	7.195	1.00	15.63	B_13
ATOM	2648	N	GLY	126	45.988	63.192	5.644	1.00	17.38	B_13
ATOM	2650	CA	GLY	126	44.700	62.968	5.001	1.00	22.41	B_13
ATOM	2651	C	GLY	126	44.453	63.487	3.603	1.00	13.20	B_13
ATOM	2652	O	GLY	126	43.349	63.366	3.096	1.00	20.86	B_13
ATOM	2653	N	LEU	127	45.452	64.079	2.972	1.00	12.39	B_13
ATOM	2655	CA	LEU	127	45.267	64.592	1.617	1.00	11.56	B_13
ATOM	2656	CB	LEU	127	45.965	65.947	1.467	1.00	19.19	B_13
ATOM	2657	CG	LEU	127	45.300	67.206	2.039	1.00	14.42	B_13
ATOM	2658	CD1	LEU	127	44.875	67.030	3.496	1.00	32.31	B_13
ATOM	2659	CD2	LEU	127	46.288	68.374	1.912	1.00	25.45	B_13
ATOM	2660	C	LEU	127	45.770	63.619	0.550	1.00	26.54	B_13
ATOM	2661	O	LEU	127	46.920	63.156	0.601	1.00	18.76	B_13
ATOM	2662	N	ASP	128	44.908	63.285	-0.407	1.00	28.54	B_13
ATOM	2664	CA	ASP	128	45.292	62.376	-1.480	1.00	10.89	B_13
ATOM	2665	CB	ASP	128	44.059	61.762	-2.136	1.00	15.95	B_13
ATOM	2666	CG	ASP	128	44.351	60.430	-2.794	1.00	23.44	B_13
ATOM	2667	OD1	ASP	128	43.377	59.735	-3.164	1.00	41.43	B_13
ATOM	2668	OD2	ASP	128	45.541	60.059	-2.918	1.00	18.12	B_13
ATOM	2669	C	ASP	128	46.060	63.203	-2.502	1.00	25.34	B_13
ATOM	2670	O	ASP	128	46.489	64.308	-2.213	1.00	16.36	B_13
ATOM	2671	N	HIS	129	46.283	62.645	-3.682	1.00	17.53	B_13
ATOM	2673	CA	HIS	129	47.001	63.366	-4.718	1.00	26.87	B_13
ATOM	2674	CB	HIS	129	47.495	62.398	-5.794	1.00	10.00	B_13
ATOM	2675	CG	HIS	129	48.729	61.645	-5.400	1.00	19.64	B_13
ATOM	2676	CD2	HIS	129	49.769	61.996	-4.609	1.00	19.96	B_13
ATOM	2677	ND1	HIS	129	49.012	60.373	-5.859	1.00	23.97	B_13
ATOM	2679	CE1	HIS	129	50.170	59.977	-5.372	1.00	17.95	B_13
ATOM	2680	NE2	HIS	129	50.658	60.944	-4.605	1.00	13.79	B_13
ATOM	2681	C	HIS	129	46.153	64.457	-5.360	1.00	39.97	B_13
ATOM	2682	O	HIS	129	45.011	64.220	-5.757	1.00	25.97	B_13
ATOM	2683	N	SER	130	46.743	65.640	-5.481	1.00	21.04	B_13
ATOM	2685	CA	SER	130	46.090	66.776	-6.109	1.00	16.72	B_13
ATOM	2686	CB	SER	130	46.847	68.058	-5.757	1.00	20.97	B_13
ATOM	2687	OG	SER	130	46.358	69.154	-6.502	1.00	25.52	B_13
ATOM	2689	C	SER	130	46.098	66.582	-7.622	1.00	24.66	B_13
ATOM	2690	O	SER	130	46.779	65.694	-8.145	1.00	29.24	B_13
ATOM	2691	N	LYS	131	45.315	67.403	-8.315	1.00	26.96	B_13
ATOM	2693	CA	LYS	131	45.253	67.358	-9.769	1.00	20.25	B_13
ATOM	2694	CB	LYS	131	43.796	67.379	-10.247	1.00	33.22	B_13
ATOM	2695	CG	LYS	131	43.159	68.775	-10.302	1.00	32.85	B_13
ATOM	2696	CD	LYS	131	43.335	69.436	-11.675	1.00	15.99	B_13
ATOM	2697	CE	LYS	131	43.023	70.919	-11.601	1.00	30.34	B_13
ATOM	2698	NZ	LYS	131	43.879	71.647	-10.600	1.00	30.44	B_13
ATOM	2702	C	LYS	131	45.998	68.602	-10.249	1.00	15.31	B_13
ATOM	2703	O	LYS	131	46.414	68.698	-11.402	1.00	30.72	B_13
ATOM	2704	N	ASP	132	46.191	69.536	-9.323	1.00	23.41	B_13
ATOM	2706	CA	ASP	132	46.869	70.798	-9.581	1.00	22.69	B_13
ATOM	2707	CB	ASP	132	46.641	71.726	-8.379	1.00	24.86	B_13
ATOM	2708	CG	ASP	132	46.819	73.200	-8.712	1.00	24.93	B_13
ATOM	2709	OD1	ASP	132	46.007	74.009	-8.208	1.00	29.71	B_13
ATOM	2710	OD2	ASP	132	47.766	73.555	-9.448	1.00	28.82	B_13
ATOM	2711	C	ASP	132	48.358	70.497	-9.728	1.00	14.97	B_13
ATOM	2712	O	ASP	132	49.047	70.235	-8.742	1.00	19.64	B_13
ATOM	2713	N	PRO	133	48.874	70.538	-10.964	1.00	16.94	B_13
ATOM	2714	CD	PRO	133	48.209	70.971	-12.199	1.00	21.42	B_13
ATOM	2715	CA	PRO	133	50.293	70.264	-11.215	1.00	19.34	B_13
ATOM	2716	CB	PRO	133	50.457	70.636	-12.690	1.00	20.48	B_13
ATOM	2717	CG	PRO	133	49.347	71.636	-12.929	1.00	21.80	B_13
ATOM	2718	C	PRO	133	51.237	71.059	-10.322	1.00	17.45	B_13
ATOM	2719	O	PRO	133	52.319	70.590	-10.006	1.00	23.30	B_13
ATOM	2720	N	GLY	134	50.799	72.246	-9.904	1.00	32.46	B_13
ATOM	2722	CA	GLY	134	51.610	73.104	-9.051	1.00	19.44	B_13
ATOM	2723	C	GLY	134	51.306	72.958	-7.569	1.00	22.33	B_13
ATOM	2724	O	GLY	134	51.556	73.877	-6.795	1.00	21.92	B_13
ATOM	2725	N	ALA	135	50.698	71.836	-7.190	1.00	34.71	B_13
ATOM	2727	CA	ALA	135	50.355	71.580	-5.794	1.00	18.35	B_13
ATOM	2728	CB	ALA	135	48.948	70.987	-5.690	1.00	14.30	B_13
ATOM	2729	C	ALA	135	51.370	70.616	-5.210	1.00	10.00	B_13
ATOM	2730	O	ALA	135	51.739	69.647	-5.858	1.00	17.52	B_13
ATOM	2731	N	LEU	136	51.727	70.842	-3.952	1.00	21.29	B_13

67/75

ATOM	2733	CA	LEU	136	52.692	70.015	-3.230	1.00	14.62	B_13
ATOM	2734	CB	LEU	136	52.738	70.458	-1.763	1.00	18.76	B_13
ATOM	2735	CG	LEU	136	54.007	70.308	-0.921	1.00	34.21	B_13
ATOM	2736	CD1	LEU	136	53.587	69.907	0.485	1.00	14.76	B_13
ATOM	2737	CD2	LEU	136	54.969	69.296	-1.508	1.00	11.64	B_13
ATOM	2738	C	LEU	136	52.232	68.564	-3.287	1.00	13.50	B_13
ATOM	2739	O	LEU	136	53.033	67.640	-3.238	1.00	19.04	B_13
ATOM	2740	N	MET	137	50.921	68.364	-3.281	1.00	17.54	B_13
ATOM	2742	CA	MET	137	50.360	67.019	-3.324	1.00	25.11	B_13
ATOM	2743	CB	MET	137	49.010	66.981	-2.599	1.00	19.80	B_13
ATOM	2744	CG	MET	137	49.083	67.312	-1.117	1.00	15.35	B_13
ATOM	2745	SD	MET	137	50.354	66.361	-0.262	1.00	11.22	B_13
ATOM	2746	CE	MET	137	49.882	64.680	-0.764	1.00	13.90	B_13
ATOM	2747	C	MET	137	50.254	66.387	-4.721	1.00	28.08	B_13
ATOM	2748	O	MET	137	49.730	65.268	-4.863	1.00	12.18	B_13
ATOM	2749	N	PHE	138	50.771	67.070	-5.743	1.00	10.00	B_13
ATOM	2751	CA	PHE	138	50.751	66.528	-7.097	1.00	12.27	B_13
ATOM	2752	CB	PHE	138	51.327	67.523	-8.094	1.00	19.38	B_13
ATOM	2753	CG	PHE	138	51.051	67.175	-9.534	1.00	25.74	B_13
ATOM	2754	CD1	PHE	138	52.090	67.077	-10.448	1.00	19.74	B_13
ATOM	2755	CD2	PHE	138	49.747	67.007	-9.990	1.00	24.46	B_13
ATOM	2756	CE1	PHE	138	51.843	66.824	-11.786	1.00	19.54	B_13
ATOM	2757	CE2	PHE	138	49.495	66.750	-11.335	1.00	24.12	B_13
ATOM	2758	CZ	PHE	138	50.544	66.664	-12.230	1.00	18.15	B_13
ATOM	2759	C	PHE	138	51.619	65.269	-7.068	1.00	25.93	B_13
ATOM	2760	O	PHE	138	52.658	65.226	-6.414	1.00	12.50	B_13
ATOM	2761	N	PRO	139	51.166	64.194	-7.714	1.00	25.17	B_13
ATOM	2762	CD	PRO	139	49.870	64.004	-8.392	1.00	10.00	B_13
ATOM	2763	CA	PRO	139	51.950	62.956	-7.713	1.00	18.48	B_13
ATOM	2764	CB	PRO	139	50.981	61.946	-8.339	1.00	15.96	B_13
ATOM	2765	CG	PRO	139	50.140	62.798	-9.250	1.00	18.82	B_13
ATOM	2766	C	PRO	139	53.299	62.950	-8.430	1.00	17.22	B_13
ATOM	2767	O	PRO	139	53.849	61.876	-8.661	1.00	36.93	B_13
ATOM	2768	N	ILE	140	53.844	64.114	-8.767	1.00	24.48	B_13
ATOM	2770	CA	ILE	140	55.118	64.155	-9.477	1.00	20.03	B_13
ATOM	2771	CB	ILE	140	54.996	64.807	-10.892	1.00	18.71	B_13
ATOM	2772	CG2	ILE	140	56.334	64.709	-11.639	1.00	23.96	B_13
ATOM	2773	CG1	ILE	140	53.932	64.113	-11.724	1.00	24.68	B_13
ATOM	2774	CD1	ILE	140	53.861	64.669	-13.125	1.00	25.83	B_13
ATOM	2775	C	ILE	140	56.109	64.992	-8.700	1.00	27.87	B_13
ATOM	2776	O	ILE	140	55.758	66.043	-8.248	1.00	22.39	B_13
ATOM	2777	N	TYR	141	57.332	64.512	-8.535	1.00	12.36	B_13
ATOM	2779	CA	TYR	141	58.350	65.281	-7.834	1.00	21.85	B_13
ATOM	2780	CB	TYR	141	59.418	64.353	-7.266	1.00	15.16	B_13
ATOM	2781	CG	TYR	141	60.592	65.096	-6.672	1.00	15.65	B_13
ATOM	2782	CD1	TYR	141	61.755	65.306	-7.407	1.00	18.56	B_13
ATOM	2783	CE1	TYR	141	62.836	65.967	-6.859	1.00	10.00	B_13
ATOM	2784	CD2	TYR	141	60.546	65.576	-5.366	1.00	11.42	B_13
ATOM	2785	CE2	TYR	141	61.626	66.236	-4.814	1.00	13.45	B_13
ATOM	2786	CZ	TYR	141	62.770	66.429	-5.567	1.00	10.00	B_13
ATOM	2787	OH	TYR	141	63.841	67.109	-5.016	1.00	18.97	B_13
ATOM	2789	C	TYR	141	59.042	66.270	-8.776	1.00	19.52	B_13
ATOM	2790	O	TYR	141	59.709	65.859	-9.727	1.00	21.37	B_13
ATOM	2791	N	THR	142	58.932	67.556	-8.465	1.00	23.99	B_13
ATOM	2793	CA	THR	142	59.573	68.616	-9.238	1.00	19.53	B_13
ATOM	2794	CB	THR	142	58.515	69.578	-9.807	1.00	10.00	B_13
ATOM	2795	OG1	THR	142	57.704	68.880	-10.756	1.00	37.02	B_13
ATOM	2797	CG2	THR	142	59.151	70.757	-10.457	1.00	34.35	B_13
ATOM	2798	C	THR	142	60.483	69.332	-8.235	1.00	19.89	B_13
ATOM	2799	O	THR	142	60.120	69.513	-7.076	1.00	25.67	B_13
ATOM	2800	N	TYR	143	61.699	69.677	-8.643	1.00	30.64	B_13
ATOM	2802	CA	TYR	143	62.609	70.344	-7.707	1.00	32.54	B_13
ATOM	2803	CB	TYR	143	64.091	70.190	-8.108	1.00	26.34	B_13
ATOM	2804	CG	TYR	143	65.008	71.048	-7.244	1.00	10.69	B_13
ATOM	2805	CD1	TYR	143	65.066	70.866	-5.852	1.00	16.37	B_13
ATOM	2806	CE1	TYR	143	65.801	71.738	-5.035	1.00	26.03	B_13
ATOM	2807	CD2	TYR	143	65.714	72.114	-7.795	1.00	17.36	B_13
ATOM	2808	CE2	TYR	143	66.451	73.006	-6.981	1.00	15.32	B_13
ATOM	2809	CZ	TYR	143	66.489	72.810	-5.610	1.00	10.00	B_13
ATOM	2810	OH	TYR	143	67.184	73.665	-4.790	1.00	27.84	B_13
ATOM	2812	C	TYR	143	62.330	71.815	-7.456	1.00	24.77	B_13
ATOM	2813	O	TYR	143	62.201	72.611	-8.399	1.00	26.19	B_13
ATOM	2814	N	THR	144	62.292	72.160	-6.170	1.00	22.23	B_13
ATOM	2816	CA	THR	144	62.103	73.533	-5.727	1.00	33.68	B_13
ATOM	2817	CB	THR	144	60.668	73.814	-5.189	1.00	28.06	B_13
ATOM	2818	OG1	THR	144	60.277	72.812	-4.241	1.00	38.14	B_13
ATOM	2820	CG2	THR	144	59.681	73.857	-6.346	1.00	48.73	B_13

68/75

ATOM	2821	C	THR	144	63.178	73.893	-4.695	1.00	35.52	B_13
ATOM	2822	O	THR	144	64.207	74.465	-5.064	1.00	39.57	B_13
ATOM	2823	N	GLY	145	62.967	73.552	-3.422	1.00	35.95	B_13
ATOM	2825	CA	GLY	145	63.967	73.872	-2.407	1.00	35.01	B_13
ATOM	2826	C	GLY	145	63.509	74.025	-0.965	1.00	26.81	B_13
ATOM	2827	O	GLY	145	62.566	74.773	-0.670	1.00	40.81	B_13
ATOM	2828	N	LYS	146	64.302	73.439	-0.066	1.00	27.13	B_13
ATOM	2830	CA	LYS	146	64.071	73.423	1.389	1.00	23.89	B_13
ATOM	2831	CB	LYS	146	65.163	72.548	2.049	1.00	29.08	B_13
ATOM	2832	CG	LYS	146	64.992	72.209	3.524	1.00	19.99	B_13
ATOM	2833	CD	LYS	146	66.079	71.224	3.913	1.00	20.44	B_13
ATOM	2834	CE	LYS	146	66.181	71.010	5.402	1.00	24.16	B_13
ATOM	2835	NZ	LYS	146	67.250	69.987	5.727	1.00	23.37	B_13
ATOM	2839	C	LYS	146	63.926	74.778	2.124	1.00	18.98	B_13
ATOM	2840	O	LYS	146	63.900	74.831	3.353	1.00	28.15	B_13
ATOM	2841	N	SER	147	63.826	75.871	1.382	1.00	35.50	B_13
ATOM	2843	CA	SER	147	63.661	77.185	1.992	1.00	31.59	B_13
ATOM	2844	CB	SER	147	64.988	77.673	2.594	1.00	27.05	B_13
ATOM	2845	OG	SER	147	65.996	77.756	1.586	1.00	48.28	B_13
ATOM	2847	C	SER	147	63.203	78.131	0.902	1.00	27.12	B_13
ATOM	2848	O	SER	147	62.743	79.251	1.168	1.00	33.75	B_13
ATOM	2849	N	HIS	148	63.248	77.644	-0.332	1.00	25.13	B_13
ATOM	2851	CA	HIS	148	62.872	78.465	-1.463	1.00	23.42	B_13
ATOM	2852	CB	HIS	148	63.704	78.076	-2.678	1.00	17.40	B_13
ATOM	2853	CG	HIS	148	65.174	78.020	-2.398	1.00	45.97	B_13
ATOM	2854	CD2	HIS	148	66.204	77.524	-3.121	1.00	27.24	B_13
ATOM	2855	ND1	HIS	148	65.724	78.476	-1.213	1.00	43.49	B_13
ATOM	2857	CE1	HIS	148	67.024	78.253	-1.218	1.00	30.28	B_13
ATOM	2858	NE2	HIS	148	67.342	77.676	-2.366	1.00	45.28	B_13
ATOM	2860	C	HIS	148	61.381	78.433	-1.796	1.00	47.15	B_13
ATOM	2861	O	HIS	148	60.936	79.166	-2.704	1.00	40.97	B_13
ATOM	2862	N	PHE	149	60.601	77.636	-1.053	1.00	48.76	B_13
ATOM	2864	CA	PHE	149	59.170	77.557	-1.347	1.00	32.44	B_13
ATOM	2865	CB	PHE	149	58.856	76.364	-2.269	1.00	27.77	B_13
ATOM	2866	CG	PHE	149	58.415	76.781	-3.657	1.00	24.63	B_13
ATOM	2867	CD1	PHE	149	57.826	75.874	-4.520	1.00	25.66	B_13
ATOM	2868	CD2	PHE	149	58.550	78.106	-4.072	1.00	30.89	B_13
ATOM	2869	CE1	PHE	149	57.376	76.277	-5.767	1.00	17.10	B_13
ATOM	2870	CE2	PHE	149	58.104	78.520	-5.311	1.00	18.57	B_13
ATOM	2871	CZ	PHE	149	57.513	77.608	-6.166	1.00	30.20	B_13
ATOM	2872	C	PHE	149	58.061	77.791	-0.308	1.00	27.40	B_13
ATOM	2873	O	PHE	149	58.299	77.971	0.892	1.00	29.69	B_13
ATOM	2874	N	MET	150	56.836	77.729	-0.822	1.00	28.66	B_13
ATOM	2876	CA	MET	150	55.621	78.027	-0.094	1.00	20.63	B_13
ATOM	2877	CB	MET	150	55.251	79.431	-0.503	1.00	25.60	B_13
ATOM	2878	CG	MET	150	55.599	79.691	-1.989	1.00	23.95	B_13
ATOM	2879	SD	MET	150	57.336	80.086	-2.296	1.00	76.68	B_13
ATOM	2880	CE	MET	150	57.209	81.473	-3.385	1.00	21.07	B_13
ATOM	2881	C	MET	150	54.436	77.118	-0.450	1.00	30.58	B_13
ATOM	2882	O	MET	150	54.104	76.948	-1.628	1.00	16.91	B_13
ATOM	2883	N	LEU	151	53.727	76.664	0.581	1.00	36.94	B_13
ATOM	2885	CA	LEU	151	52.576	75.772	0.431	1.00	25.68	B_13
ATOM	2886	CB	LEU	151	51.968	75.474	1.807	1.00	23.46	B_13
ATOM	2887	CG	LEU	151	51.087	74.232	1.927	1.00	24.21	B_13
ATOM	2888	CD1	LEU	151	51.936	72.998	1.657	1.00	21.54	B_13
ATOM	2889	CD2	LEU	151	50.487	74.150	3.314	1.00	19.89	B_13
ATOM	2890	C	LEU	151	51.498	76.322	-0.491	1.00	17.09	B_13
ATOM	2891	O	LEU	151	50.795	77.267	-0.136	1.00	35.38	B_13
ATOM	2892	N	PRO	152	51.338	75.727	-1.686	1.00	16.90	B_13
ATOM	2893	CD	PRO	152	52.154	74.643	-2.255	1.00	25.80	B_13
ATOM	2894	CA	PRO	152	50.334	76.170	-2.653	1.00	29.65	B_13
ATOM	2895	CB	PRO	152	50.447	75.110	-3.749	1.00	24.68	B_13
ATOM	2896	CG	PRO	152	51.892	74.791	-3.722	1.00	14.34	B_13
ATOM	2897	C	PRO	152	48.910	76.261	-2.087	1.00	10.00	B_13
ATOM	2898	O	PRO	152	48.543	75.505	-1.184	1.00	20.25	B_13
ATOM	2899	N	ASP	153	48.117	77.180	-2.639	1.00	19.53	B_13
ATOM	2901	CA	ASP	153	46.723	77.387	-2.226	1.00	15.90	B_13
ATOM	2902	CB	ASP	153	45.986	78.304	-3.213	1.00	22.34	B_13
ATOM	2903	CG	ASP	153	46.418	79.741	-3.115	1.00	28.86	B_13
ATOM	2904	OD1	ASP	153	47.016	80.115	-2.074	1.00	35.34	B_13
ATOM	2905	OD2	ASP	153	46.142	80.494	-4.084	1.00	30.09	B_13
ATOM	2906	C	ASP	153	45.953	76.084	-2.169	1.00	27.31	B_13
ATOM	2907	O	ASP	153	45.309	75.783	-1.167	1.00	23.50	B_13
ATOM	2908	N	ASP	154	46.000	75.339	-3.276	1.00	25.51	B_13
ATOM	2910	CA	ASP	154	45.316	74.063	-3.392	1.00	20.91	B_13
ATOM	2911	CB	ASP	154	45.745	73.364	-4.682	1.00	14.23	B_13
ATOM	2912	CG	ASP	154	45.033	72.062	-4.885	1.00	22.95	B_13

69/75

ATOM	2913	OD1	ASP	154	45.590	71.026	-4.516	1.00	17.80	B_13
ATOM	2914	OD2	ASP	154	43.904	72.076	-5.388	1.00	19.14	B_13
ATOM	2915	C	ASP	154	45.551	73.155	-2.173	1.00	26.95	B_13
ATOM	2916	O	ASP	154	44.629	72.491	-1.696	1.00	22.92	B_13
ATOM	2917	N	ASP	155	46.776	73.155	-1.654	1.00	23.56	B_13
ATOM	2919	CA	ASP	155	47.110	72.338	-0.490	1.00	28.69	B_13
ATOM	2920	CB	ASP	155	48.618	72.118	-0.388	1.00	12.87	B_13
ATOM	2921	CG	ASP	155	49.208	71.566	-1.676	1.00	24.35	B_13
ATOM	2922	OD1	ASP	155	49.705	72.369	-2.500	1.00	27.89	B_13
ATOM	2923	OD2	ASP	155	49.152	70.335	-1.875	1.00	16.96	B_13
ATOM	2924	C	ASP	155	46.582	72.976	0.781	1.00	25.41	B_13
ATOM	2925	O	ASP	155	46.055	72.275	1.656	1.00	13.36	B_13
ATOM	2926	N	VAL	156	46.733	74.296	0.891	1.00	16.99	B_13
ATOM	2928	CA	VAL	156	46.222	75.021	2.053	1.00	22.26	B_13
ATOM	2929	CB	VAL	156	46.340	76.571	1.901	1.00	25.69	B_13
ATOM	2930	CG1	VAL	156	45.811	77.249	3.158	1.00	14.95	B_13
ATOM	2931	CG2	VAL	156	47.768	77.007	1.641	1.00	17.52	B_13
ATOM	2932	C	VAL	156	44.727	74.705	2.129	1.00	10.00	B_13
ATOM	2933	O	VAL	156	44.224	74.234	3.145	1.00	22.47	B_13
ATOM	2934	N	GLN	157	44.033	74.980	1.029	1.00	16.19	B_13
ATOM	2936	CA	GLN	157	42.604	74.758	0.930	1.00	17.97	B_13
ATOM	2937	CB	GLN	157	42.108	75.039	-0.497	1.00	17.10	B_13
ATOM	2938	CG	GLN	157	40.804	75.852	-0.547	1.00	26.00	B_13
ATOM	2939	CD	GLN	157	40.949	77.284	-0.005	1.00	25.84	B_13
ATOM	2940	OE1	GLN	157	41.218	77.505	1.177	1.00	39.61	B_13
ATOM	2941	NE2	GLN	157	40.744	78.255	-0.875	1.00	32.22	B_13
ATOM	2944	C	GLN	157	42.347	73.324	1.309	1.00	18.69	B_13
ATOM	2945	O	GLN	157	41.368	73.015	1.982	1.00	10.00	B_13
ATOM	2946	N	GLY	158	43.272	72.460	0.903	1.00	31.05	B_13
ATOM	2948	CA	GLY	158	43.156	71.053	1.205	1.00	21.69	B_13
ATOM	2949	C	GLY	158	43.129	70.738	2.684	1.00	13.51	B_13
ATOM	2950	O	GLY	158	42.108	70.263	3.182	1.00	14.91	B_13
ATOM	2951	N	ILE	159	44.224	71.006	3.398	1.00	19.34	B_13
ATOM	2953	CA	ILE	159	44.268	70.686	4.827	1.00	19.14	B_13
ATOM	2954	CB	ILE	159	45.669	70.880	5.503	1.00	12.57	B_13
ATOM	2955	CG2	ILE	159	46.268	69.542	5.960	1.00	19.22	B_13
ATOM	2956	CG1	ILE	159	46.603	71.702	4.633	1.00	31.62	B_13
ATOM	2957	CD1	ILE	159	46.426	73.177	4.824	1.00	25.87	B_13
ATOM	2958	C	ILE	159	43.235	71.461	5.610	1.00	21.87	B_13
ATOM	2959	O	ILE	159	42.691	70.952	6.592	1.00	21.02	B_13
ATOM	2960	N	GLN	160	42.959	72.689	5.186	1.00	12.08	B_13
ATOM	2962	CA	GLN	160	41.967	73.483	5.874	1.00	11.43	B_13
ATOM	2963	CB	GLN	160	41.949	74.916	5.346	1.00	29.25	B_13
ATOM	2964	CG	GLN	160	43.158	75.737	5.827	1.00	22.01	B_13
ATOM	2965	CD	GLN	160	43.098	77.199	5.416	1.00	18.77	B_13
ATOM	2966	OE1	GLN	160	42.260	77.593	4.607	1.00	36.02	B_13
ATOM	2967	NE2	GLN	160	43.997	78.004	5.965	1.00	28.49	B_13
ATOM	2970	C	GLN	160	40.596	72.820	5.772	1.00	22.28	B_13
ATOM	2971	O	GLN	160	39.855	72.786	6.754	1.00	14.16	B_13
ATOM	2972	N	SER	161	40.304	72.183	4.634	1.00	32.89	B_13
ATOM	2974	CA	SER	161	39.005	71.537	4.474	1.00	29.25	B_13
ATOM	2975	CB	SER	161	38.847	70.901	3.085	1.00	19.70	B_13
ATOM	2976	OG	SER	161	39.594	69.706	2.946	1.00	24.88	B_13
ATOM	2978	C	SER	161	38.831	70.503	5.566	1.00	22.08	B_13
ATOM	2979	O	SER	161	37.745	70.340	6.118	1.00	26.26	B_13
ATOM	2980	N	LEU	162	39.931	69.852	5.919	1.00	19.14	B_13
ATOM	2982	CA	LEU	162	39.913	68.829	6.953	1.00	29.17	B_13
ATOM	2983	CB	LEU	162	41.081	67.852	6.767	1.00	12.08	B_13
ATOM	2984	CG	LEU	162	40.982	66.666	5.812	1.00	20.09	B_13
ATOM	2985	CD1	LEU	162	40.661	67.184	4.478	1.00	24.51	B_13
ATOM	2986	CD2	LEU	162	42.299	65.884	5.794	1.00	27.00	B_13
ATOM	2987	C	LEU	162	39.965	69.392	8.364	1.00	24.75	B_13
ATOM	2988	O	LEU	162	39.047	69.191	9.162	1.00	22.04	B_13
ATOM	2989	N	TYR	163	41.015	70.151	8.652	1.00	20.72	B_13
ATOM	2991	CA	TYR	163	41.211	70.689	9.980	1.00	10.00	B_13
ATOM	2992	CB	TYR	163	42.695	70.595	10.343	1.00	10.95	B_13
ATOM	2993	CG	TYR	163	43.221	69.167	10.209	1.00	10.00	B_13
ATOM	2994	CD1	TYR	163	43.114	68.261	11.264	1.00	37.53	B_13
ATOM	2995	CE1	TYR	163	43.452	66.913	11.103	1.00	26.00	B_13
ATOM	2996	CD2	TYR	163	43.703	68.689	8.990	1.00	23.78	B_13
ATOM	2997	CE2	TYR	163	44.048	67.342	8.822	1.00	17.88	B_13
ATOM	2998	CZ	TYR	163	43.914	66.461	9.879	1.00	24.28	B_13
ATOM	2999	OH	TYR	163	44.210	65.121	9.711	1.00	13.27	B_13
ATOM	3001	C	TYR	163	40.634	72.085	10.187	1.00	26.45	B_13
ATOM	3002	O	TYR	163	39.975	72.327	11.190	1.00	31.25	B_13
ATOM	3003	N	GLY	164	40.819	72.975	9.219	1.00	29.43	B_13
ATOM	3005	CA	GLY	164	40.291	74.324	9.340	1.00	30.64	B_13

70/75

ATOM	3006	C	GLY	164	41.402	75.344	9.424	1.00	30.89	B_13
ATOM	3007	O	GLY	164	41.101	76.564	9.368	1.00	26.89	B_13
ATOM	3008	OT	GLY	164	42.570	74.911	9.560	1.00	27.71	B_13
ATOM	3013	ZN	ZN	166	51.961	60.891	-2.865	1.00	28.31	BION
ATOM	3014	ZN	ZN	167	56.468	50.981	3.458	1.00	26.20	BION
ATOM	3015	CA	CA	168	63.096	53.752	-5.445	1.00	14.89	BION
ATOM	3016	CA	CA	165	50.705	55.618	13.085	1.00	15.79	BION
ATOM	3047	C5	WAY	169	54.585	56.119	-6.288	1.00	40.09	B693
ATOM	3048	CF1	WAY	169	54.019	54.934	-5.802	1.00	21.52	B693
ATOM	3049	CH	WAY	169	53.271	54.923	-4.624	1.00	32.32	B693
ATOM	3050	C2	WAY	169	53.100	56.104	-3.898	1.00	21.39	B693
ATOM	3051	C3	WAY	169	53.667	57.286	-4.369	1.00	18.26	B693
ATOM	3052	C4	WAY	169	54.402	57.308	-5.540	1.00	20.63	B693
ATOM	3053	N20	WAY	169	54.933	58.531	-5.964	1.00	22.15	B693
ATOM	3054	CD	WAY	169	54.297	59.340	-7.031	1.00	30.92	B693
ATOM	3055	C23	WAY	169	53.576	58.491	-8.087	1.00	20.75	B693
ATOM	3056	C28	WAY	169	54.224	58.114	-9.279	1.00	34.14	B693
ATOM	3057	C27	WAY	169	53.539	57.335	-10.228	1.00	33.99	B693
ATOM	3058	CM	WAY	169	52.209	56.944	-9.968	1.00	23.49	B693
ATOM	3059	N25	WAY	169	51.602	57.318	-8.814	1.00	23.61	B693
ATOM	3060	C24	WAY	169	52.246	58.071	-7.880	1.00	20.52	B693
ATOM	3061	S21	WAY	169	56.531	58.783	-5.660	1.00	20.46	B693
ATOM	3062	C16	WAY	169	56.457	60.446	-5.010	1.00	39.00	B693
ATOM	3063	C21	WAY	169	56.700	60.669	-3.634	1.00	28.79	B693
ATOM	3064	C20	WAY	169	56.656	61.967	-3.109	1.00	12.65	B693
ATOM	3065	C19	WAY	169	56.373	63.058	-3.946	1.00	15.68	B693
ATOM	3066	C18	WAY	169	56.126	62.828	-5.319	1.00	12.08	B693
ATOM	3067	C17	WAY	169	56.169	61.538	-5.852	1.00	15.19	B693
ATOM	3068	O33	WAY	169	56.337	64.360	-3.424	1.00	16.79	B693
ATOM	3069	C36	WAY	169	56.982	65.456	-4.084	1.00	20.80	B693
ATOM	3070	O15	WAY	169	56.973	57.923	-4.580	1.00	21.90	B693
ATOM	3071	O14	WAY	169	57.259	58.799	-6.913	1.00	10.86	B693
ATOM	3072	C7	WAY	169	53.486	58.556	-3.613	1.00	10.00	B693
ATOM	3073	N9	WAY	169	53.741	58.606	-2.303	1.00	10.00	B693
ATOM	3074	O10	WAY	169	53.539	59.846	-1.659	1.00	23.73	B693
ATOM	3075	O8	WAY	169	53.107	59.569	-4.154	1.00	15.89	B693
ATOM	3076	C29	WAY	169	55.383	55.968	-7.606	1.00	28.30	B693
ATOM	1	OH2	WAT	301	67.399	53.332	19.612	1.00	10.00	SOLV
ATOM	2	OH2	WAT	302	61.288	46.506	17.898	1.00	10.00	SOLV
ATOM	3	OH2	WAT	303	79.538	50.433	20.115	1.00	10.00	SOLV
ATOM	4	OH2	WAT	304	80.982	25.236	19.076	1.00	26.37	SOLV
ATOM	5	OH2	WAT	305	82.461	30.767	19.346	1.00	13.02	SOLV
ATOM	6	OH2	WAT	306	67.759	41.912	4.887	1.00	17.30	SOLV
ATOM	7	OH2	WAT	307	60.785	41.727	10.585	1.00	20.42	SOLV
ATOM	8	OH2	WAT	308	89.638	33.523	25.640	1.00	33.45	SOLV
ATOM	9	OH2	WAT	309	77.721	51.975	4.391	1.00	13.91	SOLV
ATOM	10	OH2	WAT	310	96.022	34.702	6.692	1.00	25.50	SOLV
ATOM	11	OH2	WAT	311	71.292	38.746	26.741	1.00	13.06	SOLV
ATOM	12	OH2	WAT	312	85.939	49.781	3.498	1.00	12.04	SOLV
ATOM	13	OH2	WAT	313	58.101	41.127	10.261	1.00	40.97	SOLV
ATOM	14	OH2	WAT	314	86.373	42.692	0.747	1.00	17.24	SOLV
ATOM	15	OH2	WAT	315	78.257	39.885	24.626	1.00	18.57	SOLV
ATOM	16	OH2	WAT	316	68.341	48.572	25.558	1.00	18.33	SOLV
ATOM	17	OH2	WAT	317	79.806	29.147	18.371	1.00	10.00	SOLV
ATOM	18	OH2	WAT	318	87.119	44.480	23.137	1.00	46.31	SOLV
ATOM	19	OH2	WAT	319	55.885	39.688	11.459	1.00	21.26	SOLV
ATOM	20	OH2	WAT	320	73.250	41.084	0.386	1.00	18.49	SOLV
ATOM	21	OH2	WAT	321	72.079	46.488	-6.835	1.00	27.48	SOLV
ATOM	22	OH2	WAT	322	71.923	37.638	-3.750	1.00	29.19	SOLV
ATOM	23	OH2	WAT	323	74.998	28.451	2.684	1.00	34.60	SOLV
ATOM	24	OH2	WAT	324	87.769	44.123	9.214	1.00	15.60	SOLV
ATOM	25	OH2	WAT	325	86.113	24.382	16.709	1.00	25.17	SOLV
ATOM	26	OH2	WAT	326	81.205	57.603	0.529	1.00	34.27	SOLV
ATOM	27	OH2	WAT	327	75.163	62.739	12.391	1.00	16.47	SOLV
ATOM	28	OH2	WAT	328	65.604	44.690	2.830	1.00	26.64	SOLV
ATOM	29	OH2	WAT	329	61.899	45.512	29.269	1.00	15.82	SOLV
ATOM	30	OH2	WAT	330	58.763	41.730	8.338	1.00	27.95	SOLV
ATOM	31	OH2	WAT	331	69.823	44.729	6.258	1.00	13.37	SOLV
ATOM	32	OH2	WAT	332	79.220	61.263	12.781	1.00	28.84	SOLV
ATOM	33	OH2	WAT	333	78.105	37.095	27.911	1.00	34.48	SOLV
ATOM	34	OH2	WAT	334	75.939	25.608	12.364	1.00	35.21	SOLV
ATOM	35	OH2	WAT	335	90.256	42.668	16.539	1.00	45.05	SOLV
ATOM	36	OH2	WAT	336	86.761	51.457	13.881	1.00	25.26	SOLV
ATOM	37	OH2	WAT	337	67.479	42.004	-5.009	1.00	33.30	SOLV
ATOM	38	OH2	WAT	338	82.018	50.963	8.823	1.00	19.80	SOLV
ATOM	39	OH2	WAT	339	80.278	32.895	-1.126	1.00	30.16	SOLV
ATOM	40	OH2	WAT	340	71.683	50.944	31.567	1.00	29.62	SOLV

71/75

ATOM	41	OH2	WAT	341	61.633	49.360	10.951	1.00	15.47	SOLV
ATOM	42	OH2	WAT	342	89.589	43.811	5.959	1.00	18.08	SOLV
ATOM	43	OH2	WAT	343	70.742	35.952	14.932	1.00	34.03	SOLV
ATOM	44	OH2	WAT	344	89.836	28.590	26.657	1.00	18.11	SOLV
ATOM	45	OH2	WAT	345	70.822	32.764	1.461	1.00	22.35	SOLV
ATOM	46	OH2	WAT	346	63.056	34.653	0.491	1.00	29.51	SOLV
ATOM	47	OH2	WAT	347	58.054	46.282	2.363	1.00	10.00	SOLV
ATOM	48	OH2	WAT	348	67.914	58.660	-6.267	1.00	18.30	SOLV
ATOM	49	OH2	WAT	349	70.170	56.725	0.575	1.00	11.89	SOLV
ATOM	50	OH2	WAT	350	55.922	73.897	0.623	1.00	18.86	SOLV
ATOM	51	OH2	WAT	351	73.489	53.195	2.061	1.00	24.35	SOLV
ATOM	52	OH2	WAT	352	58.033	50.530	19.075	1.00	25.52	SOLV
ATOM	53	OH2	WAT	353	63.245	57.302	17.340	1.00	13.88	SOLV
ATOM	54	OH2	WAT	354	58.442	71.334	-5.670	1.00	17.51	SOLV
ATOM	55	OH2	WAT	355	62.535	61.154	16.706	1.00	12.38	SOLV
ATOM	56	OH2	WAT	356	66.949	51.163	-10.284	1.00	17.92	SOLV
ATOM	57	OH2	WAT	357	57.588	54.191	9.850	1.00	17.88	SOLV
ATOM	58	OH2	WAT	358	64.836	48.085	4.627	1.00	17.80	SOLV
ATOM	59	OH2	WAT	359	66.445	61.785	19.640	1.00	24.12	SOLV
ATOM	60	OH2	WAT	360	55.740	42.557	0.533	1.00	27.32	SOLV
ATOM	61	OH2	WAT	361	74.075	57.146	13.179	1.00	18.01	SOLV
ATOM	62	OH2	WAT	362	46.987	69.315	-2.545	1.00	11.87	SOLV
ATOM	63	OH2	WAT	363	53.842	52.266	-2.612	1.00	25.20	SOLV
ATOM	64	OH2	WAT	364	33.425	65.313	-4.686	1.00	28.97	SOLV
ATOM	65	OH2	WAT	365	45.633	51.173	10.502	1.00	31.97	SOLV
ATOM	66	OH2	WAT	366	39.040	71.050	-0.722	1.00	20.81	SOLV
ATOM	67	OH2	WAT	367	54.517	67.335	-6.251	1.00	46.24	SOLV
ATOM	68	OH2	WAT	368	45.083	67.138	20.314	1.00	29.47	SOLV
ATOM	69	OH2	WAT	369	65.758	67.669	-6.655	1.00	14.69	SOLV
ATOM	70	OH2	WAT	370	44.943	78.174	12.948	1.00	23.88	SOLV
ATOM	71	OH2	WAT	371	37.141	57.403	1.723	1.00	23.72	SOLV
ATOM	72	OH2	WAT	372	62.407	66.806	13.368	1.00	13.36	SOLV
ATOM	73	OH2	WAT	373	50.776	47.263	5.661	1.00	38.22	SOLV
ATOM	74	OH2	WAT	374	56.697	47.264	11.752	1.00	24.75	SOLV
ATOM	75	OH2	WAT	375	42.566	60.884	15.739	1.00	16.25	SOLV
ATOM	76	OH2	WAT	376	59.299	74.342	13.838	1.00	31.27	SOLV
ATOM	77	OH2	WAT	377	72.976	63.691	-0.667	1.00	20.36	SOLV
ATOM	78	OH2	WAT	378	72.876	60.516	-6.752	1.00	34.24	SOLV
ATOM	79	OH2	WAT	379	63.998	68.760	16.371	1.00	19.04	SOLV
ATOM	80	OH2	WAT	380	44.947	66.728	-2.566	1.00	29.51	SOLV
ATOM	81	OH2	WAT	381	57.690	61.926	-9.414	1.00	29.01	SOLV
ATOM	82	OH2	WAT	382	44.595	80.810	5.831	1.00	27.43	SOLV
ATOM	83	OH2	WAT	383	78.065	36.583	24.121	1.00	14.08	SOLV
ATOM	84	OH2	WAT	384	42.289	64.651	-0.868	1.00	25.57	SOLV
ATOM	85	OH2	WAT	385	59.851	68.458	-12.381	1.00	30.18	SOLV
ATOM	86	OH2	WAT	386	53.784	72.644	-4.782	1.00	22.35	SOLV
ATOM	87	OH2	WAT	387	72.793	27.922	8.925	1.00	32.13	SOLV
ATOM	88	OH2	WAT	388	57.224	68.062	-6.072	1.00	17.87	SOLV
ATOM	89	OH2	WAT	389	45.210	44.988	4.285	1.00	25.10	SOLV
ATOM	90	OH2	WAT	390	49.413	53.782	1.546	1.00	21.68	SOLV
ATOM	91	OH2	WAT	391	45.232	59.677	1.393	1.00	19.25	SOLV
ATOM	92	OH2	WAT	392	42.551	59.954	5.056	1.00	27.30	SOLV
ATOM	93	OH2	WAT	393	58.412	43.750	3.948	1.00	58.70	SOLV
ATOM	94	OH2	WAT	394	56.942	54.199	-2.588	1.00	31.14	SOLV
ATOM	95	OH2	WAT	395	55.216	51.994	9.824	1.00	13.25	SOLV
ATOM	96	OH2	WAT	396	51.642	54.651	14.874	1.00	10.00	SOLV
ATOM	97	OH2	WAT	397	48.690	56.156	13.991	1.00	28.59	SOLV
ATOM	98	OH2	WAT	398	74.412	37.913	0.396	1.00	12.55	SOLV
ATOM	99	OH2	WAT	399	81.920	53.968	18.267	1.00	14.05	SOLV
ATOM	100	OH2	WAT	400	70.413	41.780	1.170	1.00	16.68	SOLV
ATOM	101	OH2	WAT	401	71.098	53.544	2.407	1.00	27.63	SOLV
ATOM	102	OH2	WAT	402	94.383	32.979	9.497	1.00	27.97	SOLV
ATOM	103	OH2	WAT	403	70.765	66.069	16.389	1.00	38.09	SOLV
ATOM	104	OH2	WAT	404	78.651	34.890	29.495	1.00	48.60	SOLV
ATOM	105	OH2	WAT	405	80.289	39.811	24.727	1.00	20.74	SOLV
ATOM	106	OH2	WAT	406	63.627	47.414	7.301	1.00	40.21	SOLV
ATOM	107	OH2	WAT	407	74.679	30.772	11.524	1.00	37.03	SOLV
ATOM	108	OH2	WAT	408	80.240	36.041	26.681	1.00	27.42	SOLV
ATOM	109	OH2	WAT	409	84.971	25.909	18.426	1.00	24.96	SOLV
ATOM	110	OH2	WAT	410	57.832	41.294	5.792	1.00	71.90	SOLV
ATOM	111	OH2	WAT	411	55.484	68.139	-9.086	1.00	48.47	SOLV
ATOM	112	OH2	WAT	412	65.535	68.260	2.400	1.00	26.24	SOLV
ATOM	113	OH2	WAT	413	80.085	42.291	-3.144	1.00	26.49	SOLV
ATOM	114	OH2	WAT	414	82.088	37.456	27.733	1.00	42.54	SOLV
ATOM	115	OH2	WAT	415	61.020	53.195	21.566	1.00	38.16	SOLV
ATOM	116	OH2	WAT	416	55.968	70.365	-5.096	1.00	28.42	SOLV
ATOM	117	OH2	WAT	417	51.619	57.620	-0.487	1.00	41.81	SOLV

72 / 75

ATOM	118	OH2	WAT	418	40.651	66.108	2.086	1.00	40.11	SOLV
ATOM	119	OH2	WAT	419	58.453	49.818	7.926	1.00	38.96	SOLV
ATOM	120	OH2	WAT	420	53.768	51.716	13.623	1.00	43.62	SOLV
ATOM	121	OH2	WAT	421	76.068	60.373	21.292	1.00	39.30	SOLV
ATOM	122	OH2	WAT	422	56.186	50.034	17.422	1.00	37.47	SOLV
END										

73/75

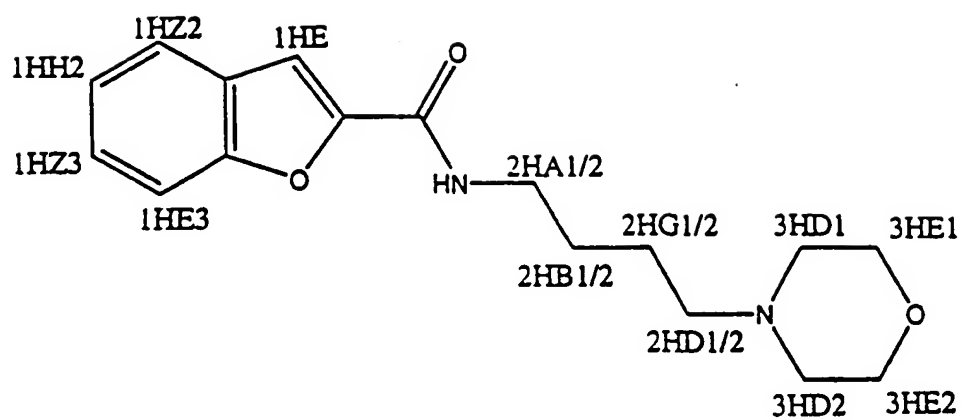


FIG. 6

74/75

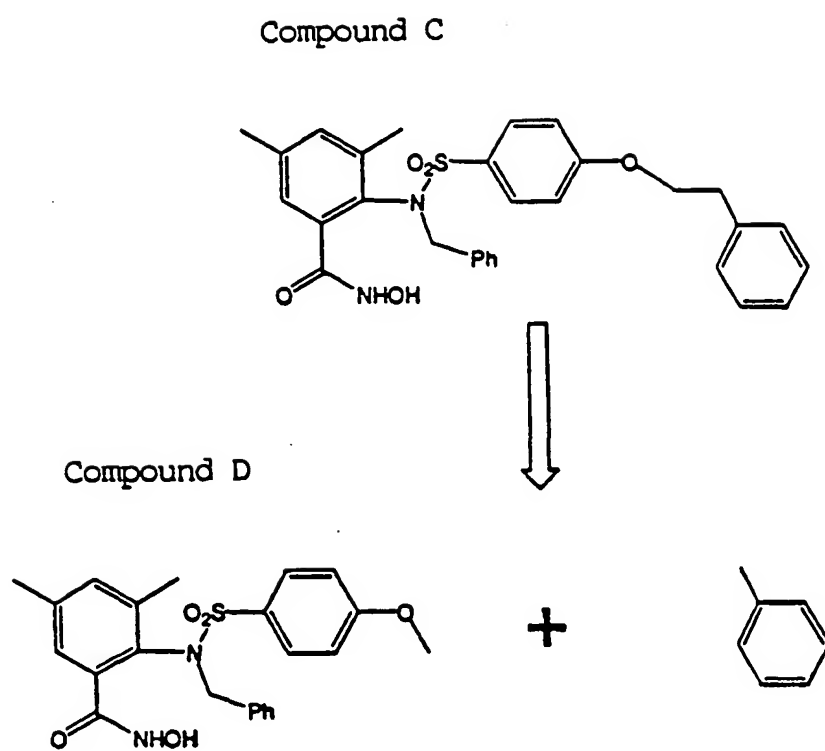


FIG. 7

75/75

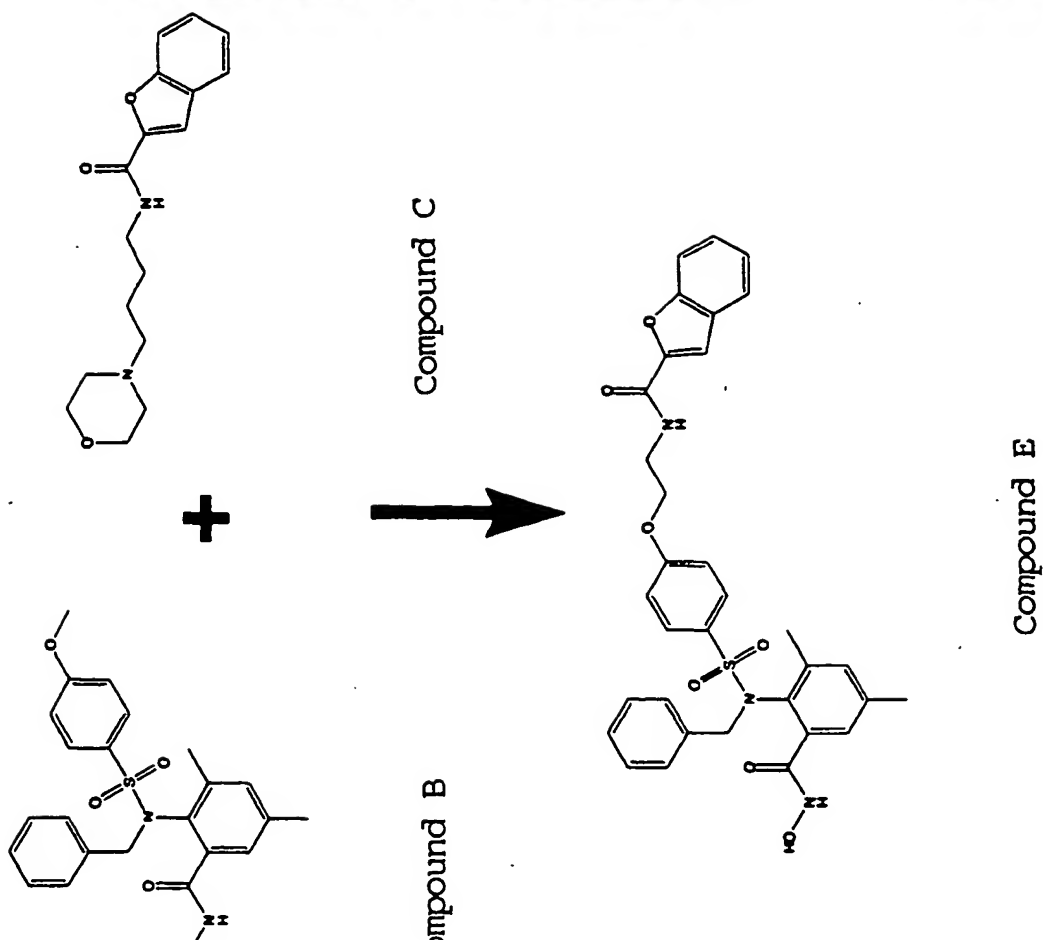
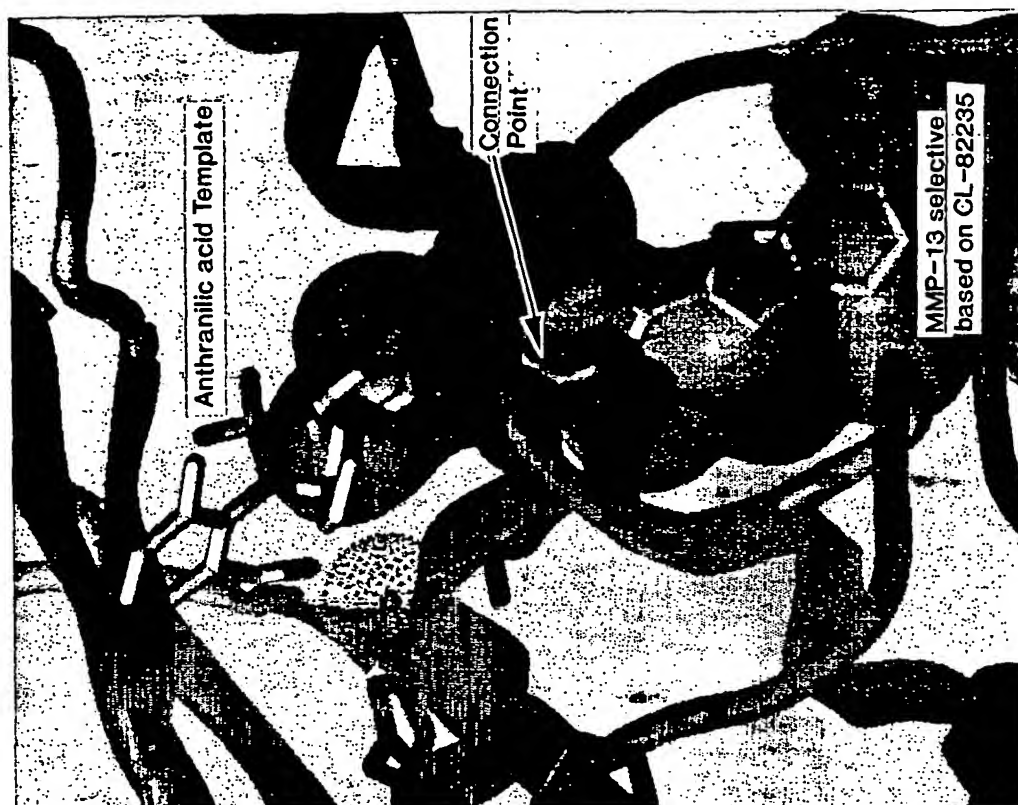


FIG. 8

INTERNATIONAL SEARCH REPORT

International application No.

PCT/US01/05150

A. CLASSIFICATION OF SUBJECT MATTER

IPC(7) : G01N 9/00, 33/48

US CL : 435/183; 702/22

According to International Patent Classification (IPC) or to both national classification and IPC

B. FIELDS SEARCHED

Minimum documentation searched (classification system followed by classification symbols)

U.S. : 435/183; 702/22

Documentation searched other than minimum documentation to the extent that such documents are included in the fields searched

NONE

Electronic data base consulted during the international search (name of data base and, where practicable, search terms used)

STN: WEST

C. DOCUMENTS CONSIDERED TO BE RELEVANT

Category*	Citation of document, with indication, where appropriate, of the relevant passages	Relevant to claim No.
X	GOMIS-RUTH, F.X. et al. The helping hand of collagenase-3 (MMP-13: 2.7, ANG> crystal structure of its C-terminal haemopexin-like domain. Journal Mol. Biol. 1996, Vol. 264, No. 3, pages 556-566, see entire document.	8-14
X	US 6,008,243 A (BENDER et al.) 28 December 1999(28.12.99), see entire document.	1-7, 15-20

☐ Further documents are listed in the continuation of Box C.
 ☐ See patent family annex.

* Special categories of cited documents:	*T* later document published after the international filing date or priority date and not in conflict with the application but cited to understand the principle or theory underlying the invention
A document defining the general state of the art which is not considered to be of particular relevance	*X* document of particular relevance; the claimed invention cannot be considered novel or cannot be considered to involve an inventive step when the document is taken alone
B earlier document published on or after the international filing date	*Y* document of particular relevance; the claimed invention cannot be considered to involve an inventive step when the document is combined with one or more other such documents, such combination being obvious to a person skilled in the art
L document which may throw doubts on priority claim(s) or which is cited to establish the publication date of another citation or other special reason (as specified)	*Z* document member of the same patent family
O document referring to an oral disclosure, use, exhibition or other means	
P document published prior to the international filing date but later than the priority date claimed	

Date of the actual completion of the international search

12 JULY 2001

Date of mailing of the international search report

30 JUL 2001

Name and mailing address of the ISA/US
Commissioner of Patents and Trademarks
Box PCT
Washington, D.C. 20231

Facsimile No. (703) 305-3230

Authorized official

AMY J. HARTTER

Telephone No. (703) 308-0196

INTERNATIONAL SEARCH REPORT

International application No.

PCT/US01/05150

Box I Observations where certain claims were found unsearchable (Continuation of item 1 of first sheet)

This international report has not been established in respect of certain claims under Article 17(2)(a) for the following reasons:

1. ☐ Claims Nos.:
because they relate to subject matter not required to be searched by this Authority, namely:
2. ☐ Claims Nos.:
because they relate to parts of the international application that do not comply with the prescribed requirements to such an extent that no meaningful international search can be carried out, specifically:
3. ☐ Claims Nos.:
because they are dependent claims and are not drafted in accordance with the second and third sentences of Rule 6.4(a).

Box II Observations where unity of invention is lacking (Continuation of item 2 of first sheet)

This International Searching Authority found multiple inventions in this international application, as follows:

Please See Extra Sheet.

1. ☒ As all required additional search fees were timely paid by the applicant, this international search report covers all searchable claims.
2. ☐ As all searchable claims could be searched without effort justifying an additional fee, this Authority did not invite payment of any additional fee.
3. ☐ As only some of the required additional search fees were timely paid by the applicant, this international search report covers only those claims for which fees were paid, specifically claims Nos.:
4. ☐ No required additional search fees were timely paid by the applicant. Consequently, this international search report is restricted to the invention first mentioned in the claims; it is covered by claims Nos.:

Remark on Protest

☐

The additional search fees were accompanied by the applicant's protest.

☐

No protest accompanied the payment of additional search fees.

BOX II. OBSERVATIONS WHERE UNITY OF INVENTION WAS LACKING

This ISA found multiple inventions as follows:

This application contains the following inventions or groups of inventions which are not so linked as to form a single inventive concept under PCT Rule 13.1. In order for these inventions to be searched the appropriate search fees must be paid.

Group I which consists of claims 1-7 is distinct as it addresses itself to the solution complex of the mixture of MMP-13 and the defined "Compound A." The solution is clearly distinct and different from the crystal complex, active site and methods that are claimed in succeeding groups and according claims.

Group II consists of claims 8-14. These claims pertain to the actual product of the crystal complex in its entirety. Thus it is distinct from Groups I and Groups 3-4. The group claims the whole crystal known as "Compound A" and the crystal is not in any other type of alternate environment or with any additional accoutrements.

Group III encompasses the claims of 15-20. These claims consist of the active site of the molecule of MMP-13. This chemical is a portion of the solution claimed in the first group and thus separate and distinct from the solution of Group I or the separate entity of "Compound A" that is claimed in Group 2. Thus these Groups are separate.

Group IV consists of claims 21-32 which claim a method of identifying an inhibitor or activator of the MMP-13 compound. The method that is embodied within this Group is clearly different from the preceding groups. Firstly the claims within Group 4 are directed toward a method of accomplishing the task of identifying different entities and not a product itself. Secondly its actions are addressed to entities outside the compound itself and not limited to "Compound A" of the MMP-13. Based on the aforementioned reasons and the distinct nature of the claims defined in each of the groups, the instant application has a lack of unity due to each group having a different Special Technical Feature as summarized above for each group.

This Page Blank (uspto)